




# COMPLETE PRODUCT CATALOGUE

Established in 1997. Proudly Canadian.   
Sold in over 85 countries worldwide.

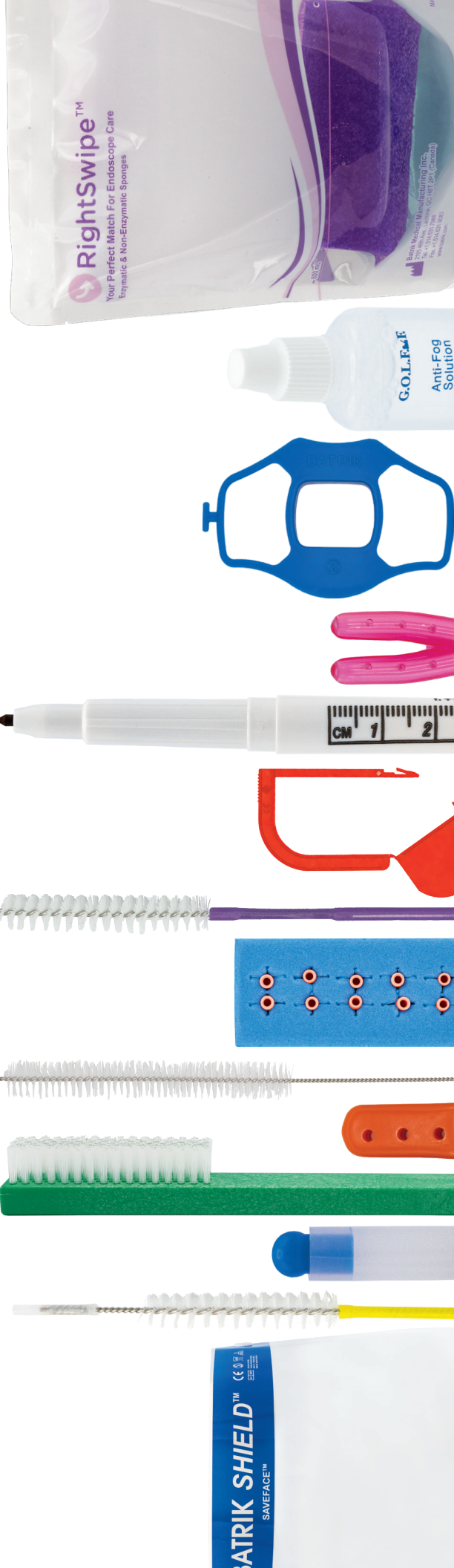

## ABOUT US

Proudly Canadian and established in 1997, we specialize in the manufacturing of enzymatic solutions offered in different formats as well as the fabrication of silicone and plastic products. Our products are sold in over 80 countries worldwide and used primarily in surgery, endoscopy, and central sterilization areas.

Our in-house controlled environment production facility in Montreal, Canada adheres to strictly enforced quality and regulatory control systems. We are registered as a device manufacturing establishment with the United States FDA and we are in full compliance with CE and Health Canada Regulations. We continually invest in upgrading our plant, machinery, and processes to ensure high quality products at high volumes. Most products are available off-the-shelf.

At Batrik, we are committed to innovation and providing our customers with the finest quality products backed by world-class customer care.

# CONTENTS

## SURGICAL PRODUCTS & GENERAL ACCESSORIES

04	SUTURE CLAMP TAGS	
05	CLAMP JAW COVERS	
06	VASCULAR SILICONE TIES	
07	CAUTERY TIP CLEANER	
08	SKINS™ <a href="#">Surgical MULTI-TIP Markers</a>	
09	PIN GUARDS	
10	SAFE-SOAK™ SPAC & MATS	
11	MAGNETIC MAT	
12	TOURNY™+ TOURNIQUETS	
14	SHOOZCOVERS™ SOFT LINE™	
16	SHOOZCOVERS™ TOUGH LINE™	
18	SHOOZCOVERS™ COATED LINE™	
20	G.O.L.F.F.™ ANTI-FOG SOLUTION	
21	WTFOG™ <a href="#">Reusable Anti-Fog Cloths</a>	
22	STACKABLE STEP STOOLS	
23	BLUE CLOUD-GEL™ ANTI-FATIGUE MATS	
24	BLUE-CLOUD NIMBUS™ ANTI-FATIGUE MATS	
25	PANDAPATCH™ EYE SAFETY TAPE	NEW!

## ENDOSCOPY PRODUCTS

27	G.O.L.F.F.™ ANTI-FOG SOLUTION	
28	BITE BLOCKS <a href="#">Endoscopy Bo2™ / Bitestrong™</a>	
29	ENDOSCOPY BITESTRONG™ BITE BLOCKS	
30	COUVRETTE™ VALVE COVERS	
31	DAYBRUSH™ CHANNEL BRUSHES	
32	ARC™ CHANNEL BRUSHES	
35	VALVE / CONTROL BRUSHES	
36	V-BRUSHES	

## INFECTION CONTROL & STERILE PROCESSING

37	TAG-IT™	NEW!
42	SILICONE PIN MAT	
43	CORNER POCKETS™ & SOFT POCKETS™	
44	LILA'S INSTRUMENT STAIN REMOVER	
45	INSTRUMENT PROTECTORS	
48	PrAmLi Power <sup>3</sup> LUBE™ INSTRUMENT LUBRICANT	
49	S.T.I.C.S.	

50	ID TAPE	
51	TAPE-B-GONE™ TAPE REMOVER	
52	PADLOCKS SECURITY TAGS	
56	SAVEFACE™ FACE SHIELDS	
57	SHIELD™ FACE SHIELD	
58	SCRUB PAD	
59	<b>Your First</b> <input checked="" type="checkbox"/> ™ ENZYMATIC & NON-ENZYMATIC CLEANERS	
60	<b>Your First</b> <input checked="" type="checkbox"/> ™ ENVIRO-SACHET™	
61	<b>Your First</b> <input checked="" type="checkbox"/> ™ FOAMING SPRAY	
62	<b>Your First</b> <input checked="" type="checkbox"/> ™ F1RST WIPES	
63	<b>Your First</b> <input checked="" type="checkbox"/> ™ DRY SPONGES	
64	<b>Your First</b> <input checked="" type="checkbox"/> ™ RIGHTSWIPE™ ENZYMATIC SPONGES BEDSIDE SPONGE KITS	NEW!
65	<b>Your First</b> <input checked="" type="checkbox"/> ™ RIGHTSWIPE™ NON-ENZYMATIC DETERGENT BEDSIDE SPONGE KITS	NEW!
66	<b>Your First</b> <input checked="" type="checkbox"/> ™ TRANSPORT BAGS FOR ENDOSCOPE AND MEDICAL INSTRUMENTS	

## BRUSHES


67	ARC™ WIRE STYLE BRUSHES	
71	DOUBLE ARC™ WIRE BRUSHES <a href="#">Asymmetrics / Symmetric</a>	
74	IN-BRUSH™ ANTIMICROBIAL TWISTED WIRE BRUSHES	
76	STAINLESS STEEL ARC™ MICRO BRUSHES	
77	INSTRUMENT CLEANING FLAT BRUSHES	
78	ORTHOPEDIC RASPS CLEANING BRUSH	
79	PRECISION BRUSH	
80	INSTRUMENT CLEANING BRUSHES	
81	ACETABULAR REAMER BRUSH	
82	ALL PURPOSE REUSABLE BRUSHES	
83	AUTOCLAVE CLEANING BRUSH KITS	

# BATRIK SUTURE CLAMP TAGS


Non-toxic, radiopaque, and latex free, our Suture Clamp Tags are designed to protect and safely hold sutures without any damage.

Additionally, its foam base, affixed to a sterile area with the clamps, allows for easy insertion.

- Packed in a convenient foam base with easy peel-off adhesive backing, even while wearing gloves
- Available in sterile and non-sterile
- Available in 2, 3, or 5 pairs
- Pre-counted

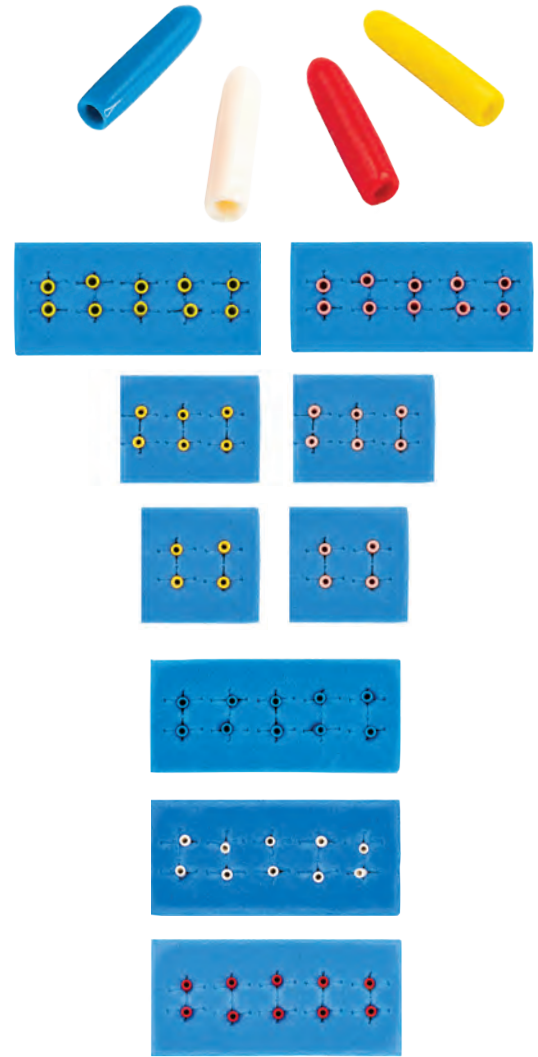

OUR "EASY PEEL CORNER" IS A FEATURE UNIQUE TO BATRIK

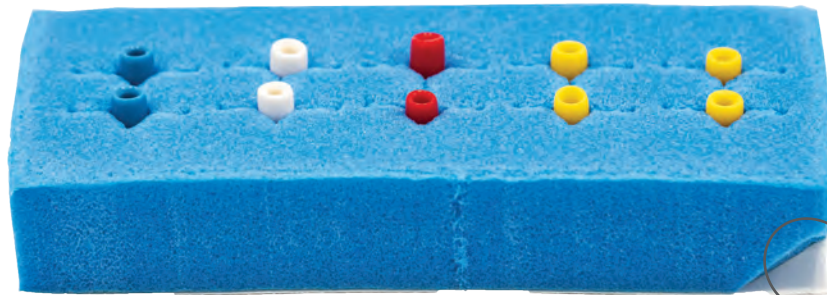

## PRODUCT SELECTION

ITEM #	STERILE	ITEM#	NON-STERILE	COLOR	DIMENSIONS	PACKAGING
11-1002		911-1002		Yellow	Standard ø (0.06" x 0.5" / 0.15 cm x 1.27 cm)	2 pairs / foam base
11-1003		911-1003		Yellow	Standard ø (0.06" x 0.5" L / 0.15 cm x 1.27 cm)	3 pairs / foam base
11-1005		911-1005		Yellow	Standard ø (0.06" x 0.5" L / 0.15 cm x 1.27 cm)	5 pairs / foam base
11-1025		911-1025		Assorted	Standard ø (0.06" x 0.5" L / 0.15 cm x 1.27 cm)	5 pairs / foam base
11-1045		911-1045		White	Standard ø (0.06" x 0.5" L / 0.15 cm x 1.27 cm)	5 pairs / foam base
11-1065		911-1065		Blue	Standard ø (0.06" x 0.5" L / 0.15 cm x 1.27 cm)	5 pairs / foam base
11-1085		911-1085		Red	Standard ø (0.06" x 0.5" L / 0.15 cm x 1.27 cm)	5 pairs / foam base
11-3002		911-3002		Pink	Large ø (0.09" x 0.75" L / 0.23 cm x 1.91 cm)	2 pairs / foam base
11-3003		911-3003		Pink	Large ø (0.09" x 0.75" L / 0.23 cm x 1.91 cm)	3 pairs / foam base
11-3005		911-3005		Pink	Large ø (0.09" x 0.75" L / 0.23 cm x 1.91 cm)	5 pairs / foam base

STERILE PACKAGING: 1 Foam Base / Pouch – 5 Pouches / Box

NON-STERILE PACKAGING: Sold per foam base

# BATRIK CLAMP JAW COVERS

Manufactured from medical grade silicone, our Clamp Jaw Covers reduce trauma to vessels requiring occlusion, particularly when a hemostatic clamp is utilized. The thread-like indents allow for gentle, constant and non-slip tension.

GOOD TO KNOW: OUR CLAMP JAW COVERS CONVERT ALMOST ANY INSTRUMENT TO ATTRAUMATIC



- Gamma sterilized
- Available in 3 colors and 3 widths to fit most clamp sizes
- Radiopaque and latex free

## PRODUCT SELECTION


ITEM # STERILE	ITEM # NON-STERILE	DIMENSIONS	COLOR
<b>12-1000</b>	<b>912-1000</b>	0.1" X 5" / 0.25 cm X 12.7 cm	Red
<b>12-1001</b>	<b>912-1001</b>	0.2" X 5" / 0.51 cm X 12.7 cm	White
<b>12-1002</b>	<b>912-1002</b>	0.3" X 5" / 0.76 cm X 12.7 cm	Blue

STERILE PACKAGING: 2 clamps per sterile pouch - 10 pouches per box

NON-STERILE PACKAGING: 20 clamps per bag (also available 100 units / bag)

# BATRIK VASCULAR SILICONE TIES

AVAILABLE IN 4 SIZES  
AND 4 COLORS

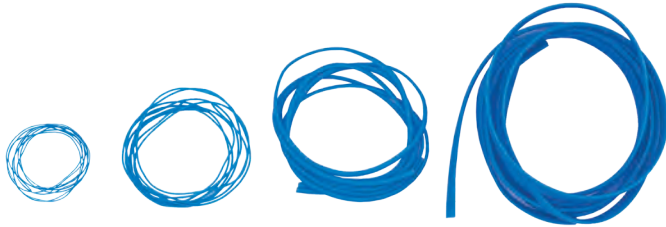


## PRODUCT SELECTION

ITEM #	DIMENSIONS	COLOR
10-1000	0.075" X 16" (0.2 cm x 40 cm)	Mini-Red
10-1001	0.075" X 16" (0.2 cm x 40 cm)	Mini-Blue
10-1002	0.075" X 16" (0.2 cm x 40 cm)	Mini-Yellow
10-1003	0.075" X 16" (0.2 cm x 40 cm)	Mini-White
10-1010	0.120" x 16" (0.3 cm x 40 cm)	Maxi-Red
10-1011	0.120" x 16" (0.3 cm x 40 cm)	Maxi-Blue
10-1012	0.120" x 16" (0.3 cm x 40 cm)	Maxi-Yellow
10-1013	0.120" x 16" (0.3 cm x 40 cm)	Maxi-White
10-1101	0.040" x 16" (0.1 cm x 40 cm)	Extra-Mini-Blue
10-1111	0.200" x 22" (0.5 cm x 56 cm)	Extra-Maxi-Blue
10-1113	0.200" x 22" (0.5 cm x 56 cm)	Extra-Maxi-White
10-1800	0.075" x 18" / 0.2 cm x 46 cm	Mini-Red
10-1801	0.075" x 18" / 0.2 cm x 46 cm	Mini-Blue
10-1802	0.075" x 18" / 0.2 cm x 46 cm	Mini-Yellow
10-1803	0.075" x 18" / 0.2 cm x 46 cm	Mini-White
10-1810	0.12" x 18" / 0.3 cm x 46 cm	Maxi-Red
10-1811	0.12" x 18" / 0.3 cm x 46 cm	Maxi-Blue
10-1812	0.120" x 18" / 0.3 cm x 46 cm	Maxi-Yellow
10-1813	0.120" x 18" / 0.3 cm x 46 cm	Maxi-White

PACKAGING: 2 ties per sterile pouch – 10 pouches per box

Our Vascular Silicone Ties permit safe and smooth retraction, occlusion and identification of vessels during surgical procedures.

Unlike clamps and retractors which can pinch the vessels flat, or even crack them, the BATRIK elliptical shape allows for the gentle handling of arteries, veins, nerves, tendons, spermatic cords, ureters, or whenever delicate retraction / occlusion is required.

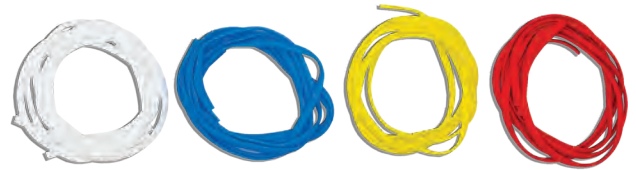

- Manufactured of medical grade silicone.
- Do not twist on themselves when unpacked.
- Easy to separate, eliminating the need for cardboard inserts.
- Sterile, radiopaque, and latex free.

## PRODUCT SELECTION

ITEM #	DIMENSIONS	COLOR
10-1831	0.040" x 18" / 0.1 cm x 46 cm	Mini-Blue
10-1841	0.200" x 18" / 0.5 cm x 46 cm	Maxi-Blue
10-1843	0.200" x 18" / 0.5 cm x 46 cm	Maxi-White
10-3010	0.120" x 30" (0.2 cm x 76 cm)	Maxi-Red
10-3200	0.075" x 32" (0.2 cm x 81 cm)	Mini-Red
10-3201	0.075" x 32" (0.2 cm x 81 cm)	Mini-Blue
10-3202	0.075" x 32" (0.2 cm x 81 cm)	Mini-Yellow
10-3203	0.075" x 32" (0.2 cm x 81 cm)	Mini-White
10-3210	0.120" x 32" (0.3 cm x 81 cm)	Maxi-Red
10-3211	0.120" x 32" (0.3 cm x 81 cm)	Maxi-Blue
10-3212	0.120" x 32" (0.3 cm x 81 cm)	Maxi-Yellow
10-3213	0.120" x 32" (0.3 cm x 81 cm)	Maxi-White

PACKAGING: 2 ties per sterile pouch – 10 pouches per box

# BATRIK CAUTERY TIP CLEANER

Our Cautery Tip Cleaner assists in removing material from cautery tips for both mono-polar and bi-polar electrodes.

- Foam pad with integrated abrasive material
- Adheres to disposable and non-disposable drapes
- Radiopaque
- Available exclusively in Non-Sterile
- Latex free



## PRODUCT SELECTION

ITEM #	DESCRIPTION	DIMENSIONS	PACKAGING
60-10000	Cautery Tip Cleaner Non-Sterile	1.8" x 1.8" / 4.5 cm x 4.5 cm	2,200 /Box

# BATRIK SKINS™ SURGICAL MULTI-TIP SKIN MARKERS

Our SKINS™ Surgical Multi-Tip Skin Markers are designed with a beveled tip to enable the creation of fine or bold tracing lines by holding the pen at different angles.

- Available with or without 6 labels and a flexible 6" calibrated (in/mm) ruler
- Available in non-sterile for pre-operative use, or sterile for use in the operating room
- Surgical prep resistant
- Made of lasting gentian violet ink that is non-toxic and non-tattooing
- Latex free



## PRODUCT SELECTION

ITEM #	TYPE	DESCRIPTION	PACKAGING
<b>S-3000</b>	Sterile	Multi-Tip fine marker with ruler and label set	50 / Box
<b>S-3001</b>	Sterile	Multi-Tip fine marker	50 / Box
<b>S-30005</b>	Sterile	Multi-Tip fine marker with ruler and label set (Volume Packaging)	500 / Box
<b>S-30015</b>	Sterile	Multi-Tip fine marker (Volume Packaging)	500 / Box
<b>S-93000</b>	Non-Sterile	Multi-Tip fine marker with ruler and label set	1000 / Box
<b>S-93001</b>	Non-Sterile	Multi-Tip fine marker	1000 / Box
<b>S-94000</b>	Non-Sterile	Ruler	100 / Bag
<b>S-94001</b>	Non-Sterile	Ruler (Bulk)	20,000 / Box
<b>S-95000</b>	Non-Sterile	Set of 3 labels	100 / Bag


SAVE by ordering our SKINS™ MARKERS in volume packaging. Cost effective and reduces waste.



# BATRIK PIN GUARDS



- Safely cover the sharp ends of your K-Wire and Steinmann Pins
- Color coded by wire size which eliminates the risk of selecting of the wrong size
- A removable sleeve is provided to avoid injury and facilitate application (sterile only)
- Non-sterile model is available with or without plastic sleeves
- Packaging designed to reduce waste
- Latex free

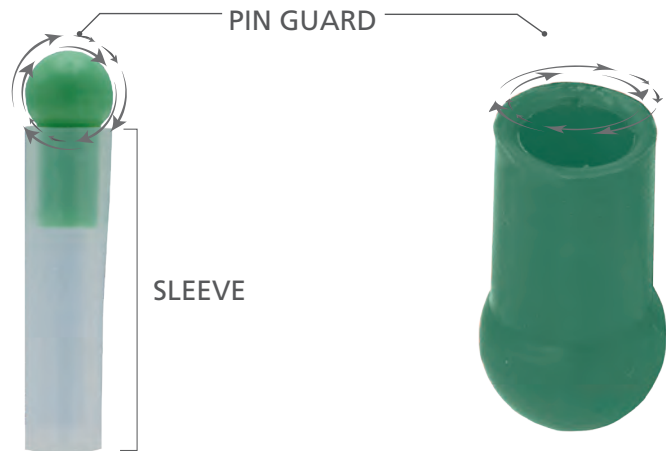

## PRODUCT SELECTION

ITEM # PIN GUARDS + SLEEVE	ITEM # PIN GUARDS + SLEEVE	ITEM # PIN GUARDS (NO SLEEVE)	COLOR	DIMENSIONS
<b>14-1028</b>	<b>914-1028 - SLV</b>	<b>914-1028</b>	Yellow	0.028" / 0.71 mm
<b>14-1035</b>	<b>914-1035 - SLV</b>	<b>914-1035</b>	Light Blue	0.035" / 0.89 mm
<b>14-1045</b>	<b>914-1045 - SLV</b>	<b>914-1045</b>	White	0.045" / 1.14 mm
<b>14-1054</b>	<b>914-1054 - SLV</b>	<b>914-1054</b>	Dark Red	0.054" / 1.40 mm
<b>14-1062</b>	<b>914-1062 - SLV</b>	<b>914-1062</b>	Green	0.062" / 1.57 mm
<b>14-1078</b>	<b>914-1078 - SLV</b>	<b>914-1078</b>	Pink	0.078" / 1.98 mm
<b>14-1093</b>	<b>914-1093 - SLV</b>	<b>914-1093</b>	Grey	0.093" / 2.36 mm
<b>14-1109</b>	<b>914-1109 - SLV</b>	<b>914-1109</b>	Beige	0.109" / 2.77 mm
<b>14-1125</b>	<b>914-1125 - SLV</b>	<b>914-1125</b>	Purple	0.125" / 3.17 mm
<b>14-1141</b>	<b>914-1141- SLV</b>	<b>914-1141</b>	Brown	0.141" / 3.58 mm
<b>14-1156</b>	<b>914-1156 - SLV</b>	<b>914-1156</b>	Dark Blue	0.156" / 3.96 mm
<b>14-1188</b>	<b>914-1188 - SLV</b>	<b>914-1188</b>	Red	0.188" / 4.77 mm

STERILE PACKAGING: 80 units/box (2 units per sterile pouch - 40 pouches per box)

NON-STERILE PACKAGING: 100 units/box

# BATRIK SAFE-SOAK™ SPAC

Our Safe-Soak™ SPAC floor suction disk quickly removes fluid waste from the operating room floor through the application of suction.



- Highly durable: it won't be damaged if stepped on
- Easily maneuverable by foot around operating room floor
- Latex free, disposable, and non-sterile
- The bottom opening permits fluid collection through a highly absorbent material that comes into contact with its surface
- The hose is securely affixed on its side to reduce risk of catching equipment
- Its highly absorbent base material allows for quick fluid wicking and removal when suction is applied

# BATRIK SAFE-SOAK™ ABSORBENT FLOOR MATS

Our Safe-Soak™ Absorbent Floor Mats are designed to absorb a variety of fluids from the operating room floor quickly and efficiently, ensuring the floor beneath remains slip-free, safe, and dry during surgical procedures.

- Highly absorbent
- Facilitates quick and easy clean up
- Latex free, single-use, disposable, and non-sterile
- Conforms to uneven surfaces
- Skid-resistant backing (blue absorbent mat only)




## PRODUCT SELECTION

ITEM #	DESCRIPTION	ABSORBENCY	SIZE	PACKAGING
<b>SS-20-1000</b>	Safe-Soak™ / White extra absorbent mat	4.5 L / mat	28" x 48" (71 cm x 122 cm)	20 / Box
<b>SS-20-2000</b>	Safe-Soak™ / Blue absorbent mat with skid-resistant backer	2.5 L / mat	28" x 40" (71 cm x 102 cm)	20 / Box
<b>SS-20-3000</b>	Safe-Soak™ / Yellow absorbent mat	3 L / mat	30" x 52" (76 cm x 132 cm)	15 / Box
<b>SS-21-1000</b>	Safe-Soak™ SPAC	N/A	8.625" (22 cm diameter)	12 / Box

# BATRIK MAGNETIC MAT

Our Magnetic Mat provides a sterile, neutral field for the safe placement and hands-free transfer of instruments. It is latex free, reusable, and autoclavable (up to 134° C).

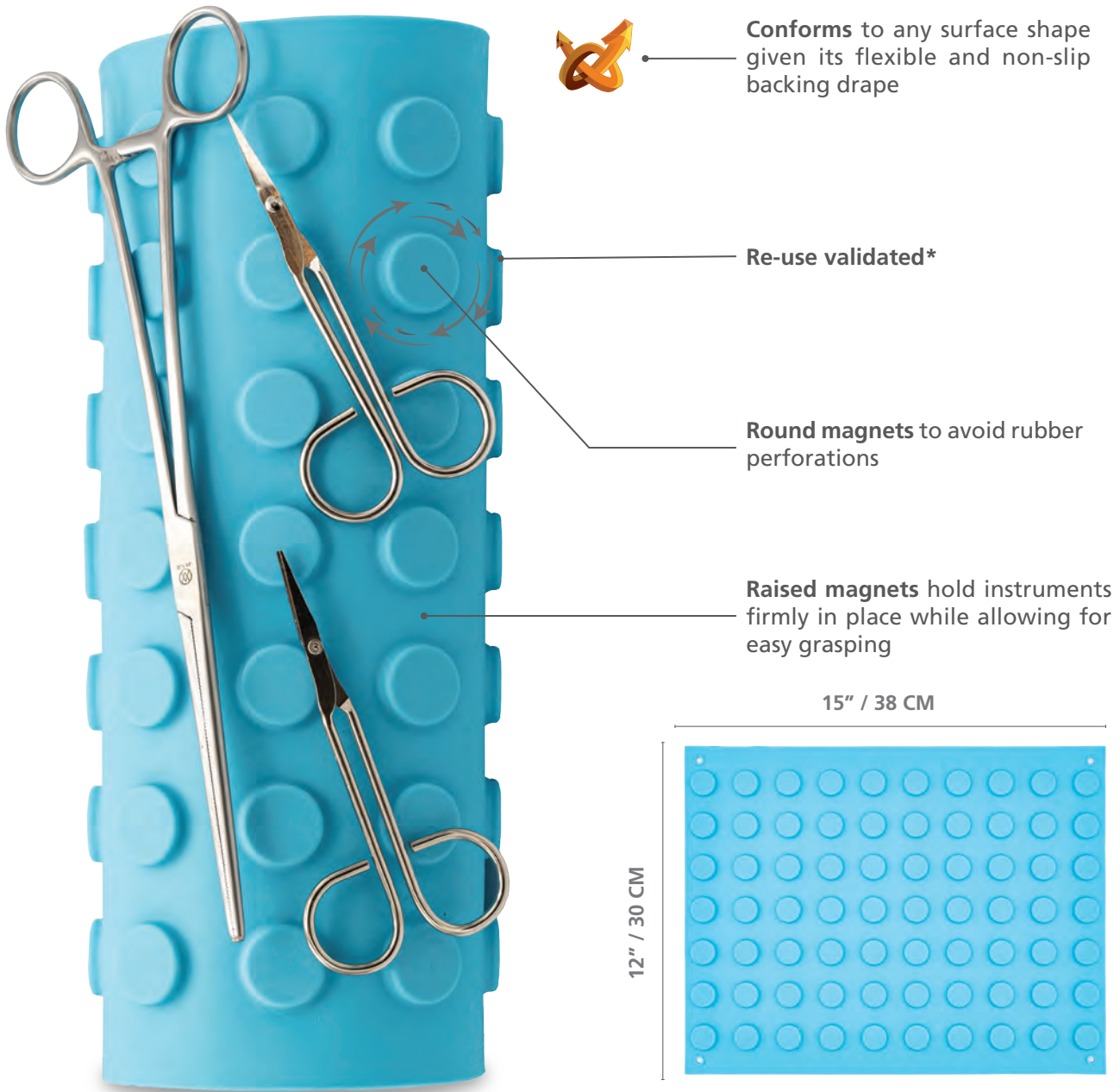


## PRODUCT SELECTION

ITEM #	DESCRIPTION	DIMENSIONS	PACKAGING
<b>MAG-02-1</b>	Magnetic Mat	12" x 15" / 30 cm x 38 cm	1 Unit
<b>MAG-02-5</b>	Magnetic Mat (Volume Packaging)	12" x 15" / 30 cm x 38 cm	5 Units

\*Independent laboratory tested to validate pre-vacuum steam sterilization efficacy for reprocessing of medical devices in a healthcare setting

# BATRIK TOURNY™ + TOURNIQUETS

Our Tourniquets are comfortable and easy to use on all patients, including the elderly, undernourished, children, and newborns.

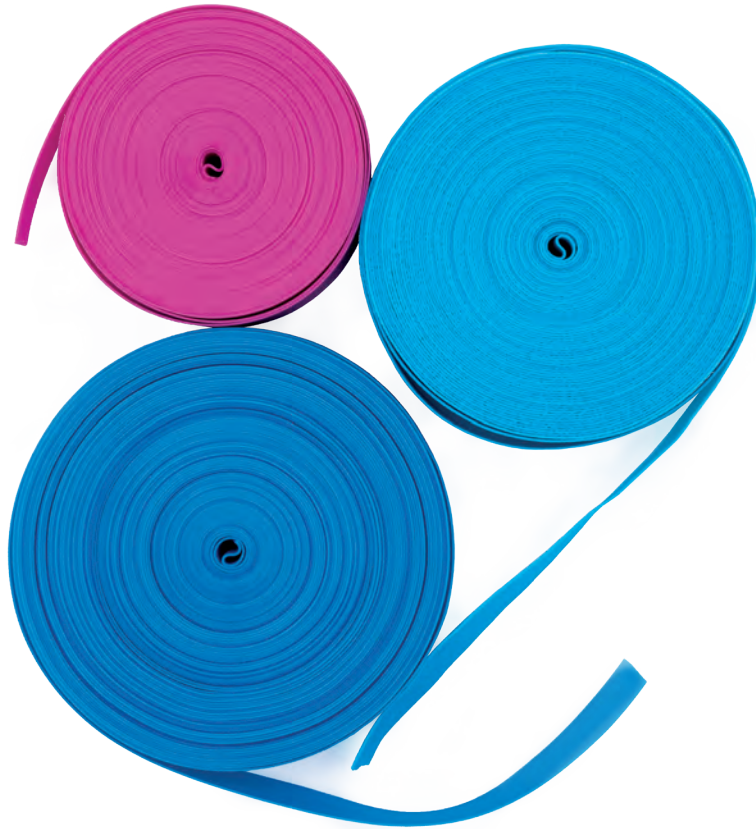

## AVAILABLE COLORS



- Available in different sizes including compact sizes, which are convenient for phlebotomy trays
- Innovative latex free material has a more comfortable texture for patients, puts less stress on skin, and safeguards hospital staff and patients from latex allergens

## SINGLE-USE BENEFITS

- Reduced risk of cross-contamination (patient to patient or patient to healthcare provider)
- Reduced risk of spread of nosocomial infections
- May improve your facility's compliance with OSHA's 1910-1030 (d)(1) "Blood-borne Pathogens Standard"

FOR YOUR CONVENIENCE:  
25 UNITS PER ROLL  
COME IN A DISPENSER BOX  
FOR EASY ACCESS



## PRODUCT SELECTION

ITEM #	DESCRIPTION	COLOR	DIMENSIONS	PACKAGING
<b>K-002</b>	TOURNY™+ Non-Sterile	Blue	1" x 18" / 2.5 cm x 45 cm	1 / Bag
<b>K-004</b>	TOURNY™+ Non-Sterile	Blue	1" x 18" / 2.5 cm x 45 cm	100 / Box
<b>K-009</b>	TOURNY™+ Non-Sterile	Blue	0.200" x 22" / 0.5 cm x 56 cm	20 / Box
<b>K-010</b>	TOURNY™+ Non-Sterile	White	0.200" x 22" / 0.5 cm x 56 cm	20 / Box
<b>K-400</b>	TOURNY™+ Non-Sterile	Blue	1" x 18" / 2.5 cm x 45 cm	16 rolls 400 / Box
<b>K-400PB</b>	TOURNY™+ Non-Sterile (Pebbled)	Light Blue	1" x 18" / 2.5 cm x 45 cm	16 rolls 400 / Box
<b>K-100</b>	TOURNY™+ Non-Sterile	Assorted (3 rolls of each) Blue, Orange, Green, and Pink	1" x 18" / 2.5 cm x 45 cm	12 rolls 300 / Box
<b>K-254</b>	TOURNY™+ Non-Sterile (Pebbled)	Light Blue	1" x 18" / 2.5 cm x 45 cm	250 / Box
<b>K-2-18RL</b>	TOURNY™+ Non-Sterile	Blue	1" x 18" / 2.5 cm x 45 cm	1 Roll
<b>K-2-18PBRL</b>	TOURNY™+ Non-Sterile (Pebbled)	Light Blue	1" x 18" / 2.5 cm x 45 cm	1 Roll
<b>K-3-18RL</b>	TOURNY™+ Non-Sterile	Green	1" x 18" / 2.5 cm x 45 cm	1 Roll
<b>K-6-18RL</b>	TOURNY™+ Non-Sterile	Pink	1" x 18" / 2.5 cm x 45 cm	1 Roll
<b>K-7-18RL</b>	TOURNY™+ Non-Sterile	Orange	1" x 18" / 2.5 cm x 45 cm	1 Roll
<b>K-2-29RL</b>	Long TOURNY™+ Non-Sterile	Blue	1" x 29.5" / 2.5 cm x 75 cm	1 Roll

# SHOOZCOVERS™ | SOFT LINE™


Batrik ShoozCovers™ Soft Line™ is manufactured in Canada from high quality 35 g/m<sup>2</sup> polypropylene fabric and 100% latex free elastics for performance and durability you can rely on. Our anti-skid tread reduces the risk of slipping and our smart design delivers easy on-and-off convenience and great fit.

## SHOOZ-REPELLENT™ TECHNOLOGY:

The Soft Line™ integrates Shooz-Repellent™ technology which creates a protective barrier that repels fluids, delivering basic incidental fluid contact repellent. Shoe covers are not waterproof or water-resistant.

## ANTI-SKID TREAD:

Our anti-skid tread reduces the risk of slipping with triple ultrasonic welds that increases tear strength and delivers easy on-and-off convenience.

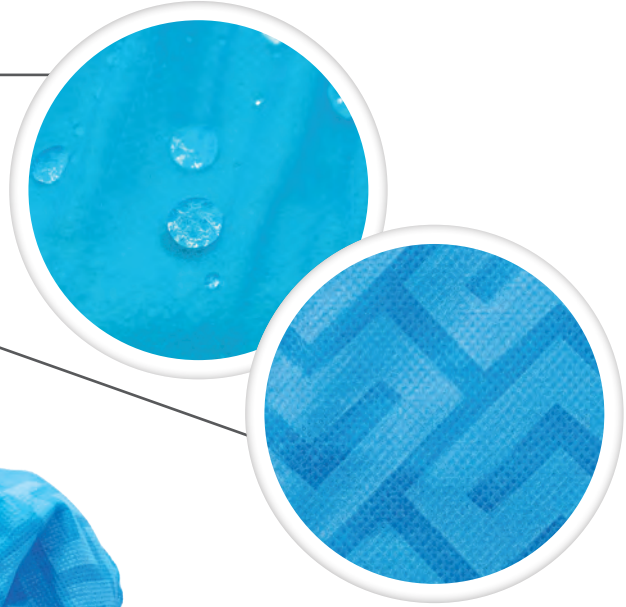

MADE IN CANADA

## OUR SOFT LINE™ IS PERFECT FOR:



MEDICAL ESTABLISHMENTS



PHARMACEUTICAL LABS



GMP SITES




ELECTRONIC MANUFACTURING



### WHY YOU'LL LOVE THE SOFT LINE™:

- Both shoe covers are manufactured in Canada with top-quality 35 g/m<sup>2</sup> polypropylene fabric and 100% latex-free USA-made elastics<sup>1</sup>.
- This line is particle count tested and complies with IEST-RP-CC003.4.
- Other colors available upon requestCC003.4.

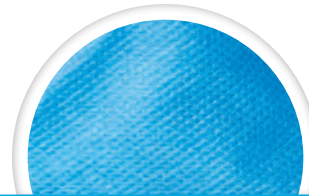

UN XL

AVAILABLE IN TWO VERSATILE SIZES

### BUBBLE GUM PINK | SOFT LINE™

### SUPPORT A CAUSE, WEAR PINK:

- A portion of our Bubble Gum Pink shoe cover sales will be donated to charity to help further breast cancer research. Join us in making a meaningful contribution by choosing Pink ShoozCovers™!



UN

AVAILABLE IN UNIVERSAL

### PRODUCT SELECTION:

ITEM #	DESCRIPTION	COLOR	SIZE	PACKAGING
<b>SC-1210-PB</b>	ShoozCovers™ - Soft Line™ 35 g/m <sup>2</sup> polypropylene non-woven fabric	Lolly Blue	Universal Up to shoe size 12 ~ 42 cm x 16.5 cm / 16.5" x 6.5"	Case of 1,000 (10 bags of 100 units)
<b>SC-1216-PB</b>	ShoozCovers™ - Soft Line™ 35 g/m <sup>2</sup> polypropylene non-woven fabric	Lolly Blue	Universal Up to shoe size 12 ~ 42 cm x 16.5 cm / 16.5" x 6.5"	Case of 600 (6 bags of 100 units)
<b>SC-1226-PB</b>	ShoozCovers™ - Soft Line™ 35 g/m <sup>2</sup> polypropylene non-woven fabric	Lolly Blue	XL Shoe size 12-15 ~ 47 cm x 19 cm / 18.5" x 7.5"	Case of 600 (6 bags of 100 units)
<b>SC-1220-PB</b>	ShoozCovers™ - Soft Line™ 35 g/m <sup>2</sup> polypropylene non-woven fabric	Lolly Blue	XL Shoe size 12-15 ~ 47 cm x 19 cm / 18.5" x 7.5"	Case of 1,000 (10 bags of 100 units)
<b>SC-1616-PB</b>	ShoozCovers™ - Soft Line™ 35 g/m <sup>2</sup> polypropylene non-woven fabric	Bubble Gum Pink	Universal Up to shoe size 12 ~ 42 cm x 16.5 cm / 16.5" x 6.5"	Case of 600 (6 bags of 100 units)

ALTERNATIVE PACKAGING: Add "-PDB" to end of code for a dispenser box instead of bags

<sup>1</sup>Many suppliers claim their products are latex free, however this often only relates to the fabric and not the elastics used. We source our ShoozCovers™ elastics from the USA to ensure that our components conform fully as a 100% latex free product.

# SHOOZCOVERS™ | TOUGH LINE™

## EXTRA-LARGE FOR EXTRA-TOUGH WORK:


Our Batrik ShoozCovers™ Tough Line™ comes in an extra-large size that fits over a men's size 15 steel-toed boot, delivering optimal full coverage and protection while you work.<sup>1</sup>

## ENHANCED DURABILITY:

These shoe covers are formulated with our **thickest fabric yet!** Made of ~55 g/m<sup>2</sup> polypropylene fabric, these highly durable covers provide added protection, increase tear strength and reduce linting.



SHOE SIZES 12-15 (45-49.5 EU)  
~ 47 CM X 19 CM / 18.5" X 7.5"

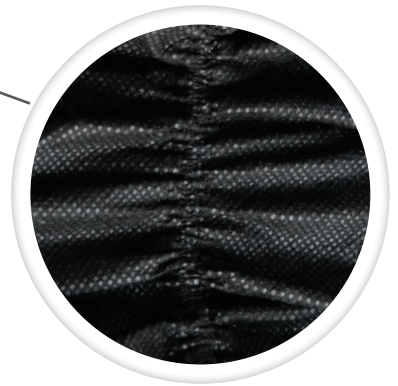

MADE IN CANADA



## OUR TOUGH LINE™ IS PERFECT FOR:



HOMES OR  
COMMERCIAL  
BUILDINGS



SHIELDING  
AGAINST MUD  
& GREASE



MEDICAL  
ESTABLISHMENTS



PHARMACEUTICAL  
LABS


<sup>1</sup> SHOOZCOVERS™ ARE MADE LONGER (42 & 47CM) FOR ACCURATE SHOE SIZE SELECTION.



## BLACK LIQUORICE | TOUGH LINE™

### WHY YOU'LL LOVE THE TOUGH LINE™:

- Triple ultrasonic welds that are long lasting to reduce waste
- High quality made in the USA latex-free elastics are triple reinforced to contract the shoe cover around your foot and deliver a strong and secure fit<sup>2</sup>
- Shooz-Repellent™ technology delivers basic incidental fluid contact repellency<sup>3</sup>
- Other colors available upon request

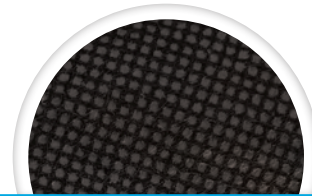

XL

AVAILABLE IN EXTRA LARGE

### ANY TIME, ANYWHERE:



RESIDENTIAL



COMMERCIAL



INDUSTRIAL

### PRODUCT SELECTION:

ITEM #	DESCRIPTION	COLOR	SIZE	PACKAGING
<b>SC-1922-PB</b>	ShoozCovers™ - Tough Line™ ~55 g/m <sup>2</sup> polypropylene non-woven fabric	Black Liquorice	XL Shoe size 12-15 ~ 47 cm x 19 cm / 18.5" x 7.5"	Case of 200 (2 bags of 100 units)

ALTERNATIVE PACKAGING: Add "-PDB" to end of code for a dispenser box instead of bags

<sup>2</sup> Many suppliers claim their products are latex free, however this often only relates to the fabric and not the elastics used. We source our ShoozCovers™ elastics from the USA to ensure that our components conform fully as a 100% latex free product. <sup>3</sup> Shooz-Repellent™ technology provides basic incidental fluid contact repellent. It is not waterproof nor water resistant.


# SHOOZCOVERS™ | COATED LINE™

## STAY DRY:

Our Batrik ShoozCovers™ Coated Line™ is made from a ~55 g/m<sup>2</sup> blend of polypropylene and polyethylene, providing high-quality water resistance.<sup>1</sup> Ideal for environments where maximum protection against fluids is required.

## TRIPLE WELDED TECHNOLOGY:

These shoe covers are reinforced with triple ultrasonic welds, a technique that joins materials together without the need for adhesives or stitching. This technique increase tear strength and reduce linting, for long lasting, resilient covers.

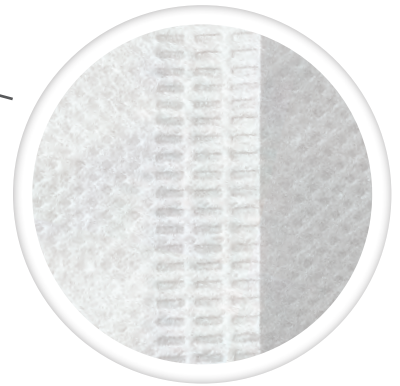

MADE IN CANADA

## OUR COATED LINE™ IS PERFECT FOR:



MEDICAL ESTABLISHMENTS



PHARMACEUTICAL LABS



GMP SITES



ELECTRONIC MANUFACTURING

<sup>1</sup>OUR SHOOZCOVERS™ COATED LINE™ IS RESISTANT AGAINST ALL NON PETROLEUM BASED LIQUIDS.

## MARSHMALLOW WHITE | COATED LINE™

### WHY YOU'LL LOVE THE COATED LINE™:



- **Superior protection using our unique blend of fabrics** which protects against not only water but all non-petroleum-based liquids, making them the ideal choice for environments where spills and leaks are a concern
- **Elastics around the ankle create a strong and secure fit** which is especially important when working in hazardous environments where a misplaced shoe cover could lead to contamination or injury
- **Smart design** with a large opening, providing a better fit and easy on-and-off convenience
- **Other colors available upon request**



UN XL

AVAILABLE IN TWO VERSATILE SIZES

### PRODUCT SELECTION:

ITEM #	DESCRIPTION	COLOR	SIZE	PACKAGING
<b>SC-2212-NPB</b>	ShoozCovers™ - Coated Line™ <i>Polypropylene + polyethylene non-woven fabric, waterproof</i>	Marshmallow White	Universal Up to shoe size 12 ~ 42 cm x 17.5 cm / 16.5" x 6.88"	Case of 300 (2 bags of 150 units)
<b>SC-2212-NPDB</b>	ShoozCovers™ - Coated Line™ <i>Polypropylene + polyethylene non-woven fabric, waterproof</i>	Marshmallow White	Universal Up to shoe size 12 ~ 42 cm x 17.5 cm / 16.5" x 6.88"	Case of 300 (2 boxes of 150 units)
<b>SC-2122-NPB</b>	ShoozCovers™ - Coated Line™ <i>Polypropylene + polyethylene non-woven fabric, waterproof</i>	Marshmallow White	XL Shoe size 12-15 ~ 47 cm x 19 cm / 18.5" x 6.88"	Case of 240 (2 bags of 120 units)
<b>SC-2122-NPDB</b>	ShoozCovers™ - Coated Line™ <i>Polypropylene + polyethylene non-woven fabric, waterproof</i>	Marshmallow White	XL Shoe size 12-15 ~ 47 cm x 19 cm / 18.5" x 6.88"	Case of 240 (2 boxes of 120 units)

ALTERNATIVE PACKAGING: Add "-PDB" to end of code for a dispenser box instead of bags

# BATRIK WTFOG™ by G.O.L.F.F.™ REUSABLE ANTI-FOG CLOTHS

## “WHAT THE FOG?!”

We’ve all experienced it.

Eyewear fogging up while wearing a mask is not only unpleasant and highly frustrating, but it can be dangerous, too.

## INTRODUCING THE FOG FIX:

Our WTFog™ reusable anti-fog cloths prevent fogging on all types of eyewear. Ideal for use when wearing a mask or face shield, WTFog™ takes you from “What The Fog?!” to confidently ready for your next patient in a snap.

## WHY YOU’LL LOVE IT:

- Eliminates fogging on all glass, polycarbonate, and plastic eyewear lenses
- Each reusable nano-infused cloth lasts up to 600 wipes or 6 months
- Microfiber cloths are infused with our proprietary heavy duty anti-fog treatment
- Quick and easy application - just wipe and wear
- Lint-free, scratch-free, no greasy smearing
- Long lasting protective coating provides effective protection against fog, humidity, sweat, and fast temperature changes.
- 100% Natural Titanium Dioxide (TiO<sub>2</sub>) Nano Formula developed in Germany meets strict REACH standards
- Non-Sterile



“To be able to complete a whole appointment without adjusting my eye protection makes maintaining COVID PPE protocols that much easier.”

**DR. MOLLY FRIEDMAN**

HOLLAND BLOORVIEW KIDS REHABILITATION HOSPITAL  
TORONTO, ON

## PERFECT WHEN:

- You wear glasses/goggles and a mask/face shield
- Your profession regularly exposes you to fogging and you need crystal clear vision without the risk or hassle of fogging up



### DIRECTIONS FOR USE:

Open resealable pouch and remove the WTFog™ reusable anti-fog cloth. Wipe cloth on clean surface (glasses/goggle lens) being careful not to remove the protective coating delivered by the cloth. There is no wait time required before wearing your protected item. Enjoy your fog-free experience!

### CLOTH STORAGE & CARE:

After use, keep cloth in the resealable pouch until next use, ensuring the seal is fully closed. If cloth becomes too dry, add up to 1 ml of water to pouch to reactivate cloth. For best results, store at ambient temperature/humidity. Dispose cloth after 600 uses or 6 months.

**External surface use only.**  
**Keep out of reach of children.**  
**Do not wash**

## PRODUCT SELECTION

ITEM #	DESCRIPTION	PACKAGING
<b>WTFOG-30</b>	30 WTFog™ in a dispenser box measuring: 3 7/16" x 5 2/16" x 6 2/16" (8.7 cm x 13 cm x 15.4 cm)	30 WTFog™ / dispenser box
<b>WTFOG-100</b>	100 WTFog™ - bulk packed	100 WTFog™ / box
<b>WTFOG-1200</b>	1200 WTFog™ - bulk packed	1200 WTFog™ / case (12 boxes of 100 bulk packed)

1 WTFog™ Cloth per Resealable Pouch  
 WTFog™ Cloth: 5.5" x 5.5" (13.9 cm x 13.9 cm)  
 WTFog™ Resealable Pouch: 3.5" x 5" (8.9 cm x 12.7 cm)

**IMPERIAL SURGICAL LTD.**  
**STACKABLE**  
**STEP STOOLS**

Made of lightweight cast aluminium and protected with an antimicrobial powder coating, Imperial Surgical Ltd.'s patented Stackable Step Stools have been designed with infection prevention as a prime focus.



**PATENTED**

- Weight of stool: 8.3 lbs (3.8 kg)
- Dimensions: 17.4" L x 13" W x 4.7" H (44.2 cm L x 33 cm W x 11.9 cm H)
- Antimicrobial coating prevents the spread of dangerous micro-organisms and suppresses the growth of mold, bacteria, yeast, and mildew
- Cart washable. Model OR-36364 (without mat) can be sterilized
- Anti-slip surface for enhanced stability
- Rubber-free corners avoids fluid collection
- No welds or large hollow openings to trap debris or fluids
- Transport cart available
- 300 lb (136 kg) weight capacity
- 7 year warranty covers all parts (step stool only)

**PRODUCT SELECTION**

ITEM #	DESCRIPTION	COLOR
<b>OR-36363</b>	Stackable Step Stool with mat	Black
<b>OR-36363-02</b>	Stackable Step Stool with mat	Dark Blue
<b>OR-36363-03</b>	Stackable Step Stool with mat	Green
<b>OR-36363-07</b>	Stackable Step Stool with mat	Orange
<b>OR-36364</b>	Stackable Step Stool without mat	N/A

**PACKAGING:** Sold by Unit

**OPTIONAL ACCESSORIES**

ITEM #	DESCRIPTION	PACKAGING
<b>OR-3693</b>	Transport cart	1 Unit
<b>809-1705</b>	Reusable Anti-Fatigue Mat with adhesive strips to secure it in place	5 / Box
<b>899-1710</b>	Single-Use Anti-Fatigue Mat with adhesive strips to secure it in place	10 / Box

**REPLACEMENT PARTS**

ITEM #	DESCRIPTION	PACKAGING
<b>OR-36363-1A</b>	Replacement mat for OR-36363	1 Unit
<b>OR-36363-1B</b>	Replacement Foot Pegs	1 Unit


**TIP: COLOR-CODE YOUR STEP STOOLS BY DEPARTMENT**



# BATRIK BLUE CLOUD-GEL™ ANTI-FATIGUE MATS

Our Blue Cloud-Gel™ Anti-Fatigue Mats have been scientifically engineered for use in the O.R. using patented gel technology and premium foam to provide the ultimate comfort and support for surgeons, nurses, and technicians particularly when movement is restricted for prolonged durations.

- Proprietary 3/4" thickness Blue Cloud-Gel™ foam technology will not "bottom-out" like conventional foam mats
- Conforms to and supports feet promoting efficient circulation, reducing physical discomfort and fatigue
- Promotes proper posture while diminishing back, knee, and foot discomfort
- Latex free, water resistant, and nonabsorbent
- Reusable, easy to clean, and Cart Washer compatible (contact us for more details)
- Excellent friction properties when damp or dry surgical booties are worn
- Completely sealed with an RF weld for excellent protection against moisture absorption
- US Standards: CALL 117 and NFPA 260 Compliant
- 1 Year Warranty

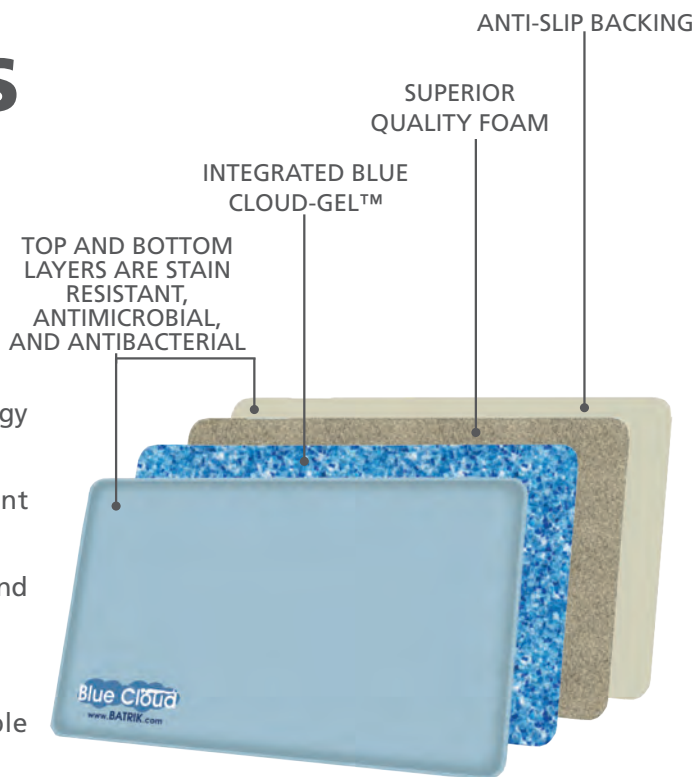

DID YOU KNOW? STUDIES HAVE DEMONSTRATED THAT ANTI-FATIGUE MATS IMPROVE PERFORMANCE, ENERGY, AND POST-OPERATIVE DISCOMFORT.

## PRODUCT SELECTION

ITEM #	DIMENSIONS	PACKAGING
<b>800-2401</b>	18" X 24" / 46 cm x 61 cm	1 / Box
<b>800-2405</b>	18" X 24" / 46 cm x 61 cm	5 / Box
<b>800-3201</b>	20" X 32" / 51 cm x 81 cm	1 / Box
<b>800-3203</b>	20" X 32" / 51 cm x 81 cm	3 / Box
<b>800-4801</b>	20" X 48" / 51 cm x 122 cm	1 / Box
<b>800-4802</b>	20" X 48" / 51 cm x 122 cm	2 / Box
<b>800-7201</b>	30" X 72" / 76 cm x 183 cm (Special Order)	1 / Box



# BATRIK BLUE-CLOUD NIMBUS™ ANTI-FATIGUE MATS

Our Blue-Cloud Nimbus™ Anti-Fatigue Mats are designed for general hospital use. Made with proprietary high-density “energy-return” polyurethane materials, they provide unmatched comfort, relief, and ergonomic support.

- Our Energy-Return properties will help reduce leg, foot, lower back discomfort, and work-related injuries
- Longer lasting support as 3/4” thickness proprietary foam blend will not bottom out like most other foam based anti-fatigue mats
- Antimicrobial properties blended into the mat provides added protection
- Highly water resistant surface will not absorb water or dirt
- Lightweight and easy to clean
- Textured design on surface and high-traction base ensures greater safety. (US Standard NFSI 101-C)
- Beveled edging reduces tripping hazards and anti-slip backing enhances safety
- Durable design is ideal for commercial applications, yet still ergonomically and aesthetically pleasing for home or office environments
- Flammability/Burn resistant per CFR 16 Parts 1630/1631
- Made with eco-friendly materials, is non-toxic, and latex free
- 5 Year Warranty (non-industrial usage), 4 Year Warranty (industrial usage)

GOOD TO KNOW: IF WITHIN 30 DAYS YOU ARE NOT 100% SATISFIED, CONTACT YOUR DEALER FOR A REFUND (EXCLUDES TRANSPORT, HANDLING, AND TAXES)



Available in  
Taupe or Black

## PRODUCT SELECTION

ITEM #	DIMENSIONS INCHES / CM	PACKAGING
<b>BLACK COLOR</b>		
<b>809-3201</b>	20" X 32" / 51 cm x 81 cm	1 / Box
<b>809-3204</b>	20" X 32" / 51 cm x 81 cm	4 / Box
<b>809-4801</b>	20" X 48" / 51 cm x 122 cm	1 / Box
<b>809-4804</b>	20" X 48" / 51 cm x 122 cm	4 / Box
<b>809-3601</b>	24" X 36" / 61 cm x 91 cm	1 / Box
<b>809-7201</b>	20" X 72" / 51 cm x 183 cm	1 / Box
<b>809-6001</b>	36" X 60" / 91 cm x 152 cm	1 / Box
<b>TAUPE COLOR</b>		
<b>810-3201</b>	20" X 32" / 51 cm x 81 cm	1 / Box
<b>810-3204</b>	20" X 32" / 51 cm x 81 cm	4 / Box
<b>810-4801</b>	20" X 48" / 51 cm x 122 cm	1 / Box
<b>810-4804</b>	20" X 48" / 51 cm x 122 cm	4 / Box
<b>810-3601</b>	24" X 36" / 61 cm x 91 cm	1 / Box
<b>810-7201</b>	20" X 72" / 51 cm x 183 cm	1 / Box
<b>810-6001</b>	36" X 60" / 91 cm x 152 cm	1 / Box



**NEW!**



## Introducing PandaPatch™: Innovative Adhesive Eye Protection for Safer Patient Care

PandaPatch™ offers an effective alternative to traditional tape for protecting patients' eyes during surgery. These adhesive eye protectors are designed to reduce the risk of corneal injuries and ophthalmic complications for patients under general anesthesia.

With PandaPatch™, you ensure superior patient outcome by preventing avoidable eye damage associated with the use of standard surgical tape.

### KEY FEATURES:

- **Ergonomic Design:** Contours comfortably around the eyes.
- **Post-Surgery Friendly:** Particularly effective for children, the elderly, and individuals with sensitive or delicate skin.
- **Latex-Free & Hypoallergenic Adhesive:** Made from medical-grade material, gentle on the skin and designed to minimize irritation.
- **Non-Adhesive Colored Tabs:** Allowing quick and hassle-free application and removal, facilitating effortless monitoring of pupil dilation during surgery.



**70-04ADU**  
Adult Sensitive Eye Protector

### HERE'S WHY IT MATTERS:

- During prolonged surgery under general anesthesia, **up to 59%<sup>1</sup>** of patients may not fully close their eyes, risking corneal abrasions.
- **PandaPatch™** gentle adhesive eye protectors ensure eyelids remain safely and comfortably closed, minimizing the risks typically associated with the use of standard surgical tape.



#### REFERENCES

1. Anson, Jonathan, M.D. "Perioperative Corneal Abrasions: Etiology, Prevention, and Management." Clinical Updates, Penn State University College of Medicine, 2021.




**PANDAPATCH™ RANGE:**

Three adhesive eye protectors to accommodate different age groups and skin sensitivities. Providing optimal protection for every patient.

- Adult Extra Sensitive Adhesive (Blue Tabs)
- Adult Adhesive (Red Tabs)
- Pediatric Adhesive (Red Tabs)



**70-04ADU**  
Adult Adhesive Eye Protector

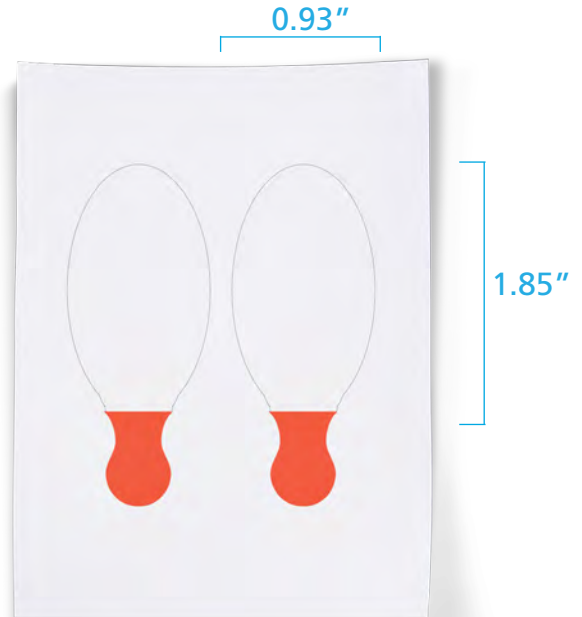

**70-04PED**  
Pediatric Adhesive Eye Protector

**PRODUCT SELECTION**

ITEM #	DESCRIPTION	DIMENSIONS	PACKAGING
<b>70-02ADU</b>	ADULT EXTRA SENSITIVE ADHESIVE EYE PROTECTORS (BLUE)	2.15" x 1.35"	50 / Box
<b>70-04ADU</b>	ADULT ADHESIVE EYE PROTECTORS (RED)	2.15" x 1.35"	50 / Box
<b>70-04PED</b>	PEDIATRIC ADHESIVE EYE PROTECTORS (RED)	1.85" x 0.93"	50 / Box

# BATRIK G.O.L.F.F.™ ANTI-FOG SOLUTION

Our G.O.L.F.F.™ Anti-Fog Solution solves the fogging problem often encountered during Laparoscopic surgical procedures.

- Prevents the occurrence of fog on lenses
- Ensures clearer vision and avoids unnecessary interruptions during surgery
- Set includes an easy-squeeze 6 mL bottle and an adhesive radiopaque sponge for secure placement
- OEM (available)
- Sterile
- Latex free



## PRODUCT SELECTION

ITEM #	DESCRIPTION	PACKAGING
<b>G.O.L.F.F.</b>	Anti-Fog Solution & Adhesive Sponge set	20 / Box


# BATRIK ENDOSCOPY BO2™ BITE BLOCKS WITH NASAL OXYGEN DELIVERY

Our BO2™ Bite Blocks significantly enhance patient safety by providing supplemental oxygenation during gastrointestinal endoscopy, preventing hypoxia in sedated patients, and ensuring superior patient comfort.

- Saves time and money by eliminating the need for separate nasal prongs and oxygen masks
- Designed to be free of sharp edges to ensure patient comfort and safety
- Available exclusively in Single-Use
- Optional retention strap is available to ensure a secure fit
- The Bite Blocks and optional retention straps are all latex free, lead free and DEHP free

**BO2**  
(FRONT SIDE)

**BO2**  
(BACK SIDE)

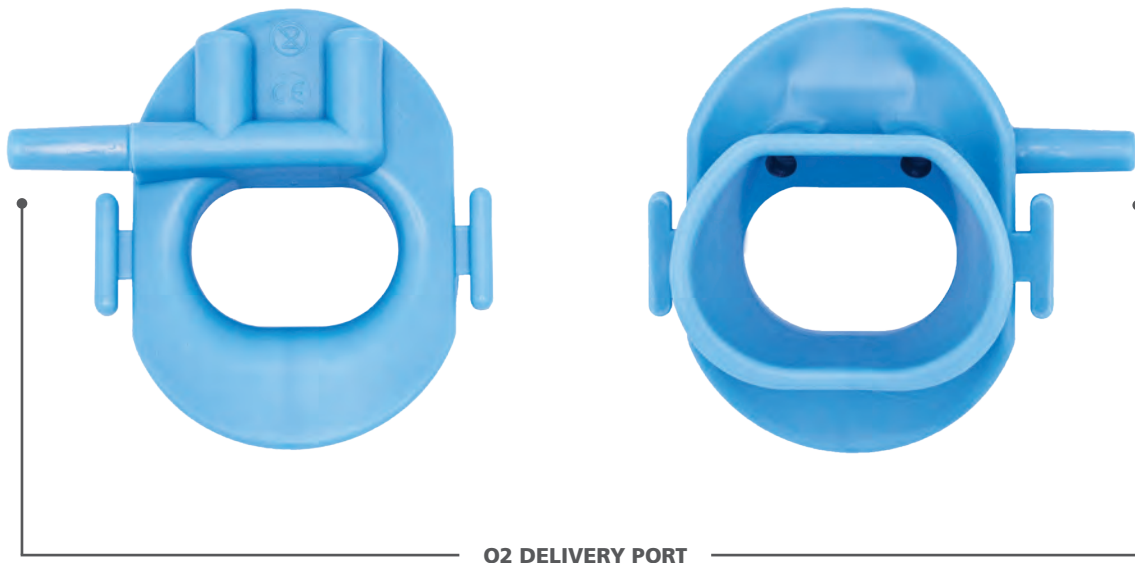

## PRODUCT SELECTION

ITEM #	TYPE	DESCRIPTION	PACKAGING
<b>BO2</b>	Single-Use	18mm / 54FR, includes strap	100 / Box
<b>BO2-N</b>	Single-Use	18mm / 54FR, without strap	100 / Box

# BATRIK ENDOSCOPY BITESTRONG™ BITE BLOCKS

## ADULT & PEDIATRIC BITE BLOCKS

Our Bitestrong™ Bite Blocks provide a safe passage through the patient's mouth during TEE, protecting the patient's teeth and endoscope.

Our molding process ensures a smooth finish with no sharp edges (flashing) that can often irritate a patient's mouth. All Bite Blocks are individually wrapped in sealed, easy to open PE bags and are latex free, lead free, and DEHP free. Additional straps are available, if required

### SINGLE USE

- Available without dental retention

### REUSABLE

- Designed for repeated usage and will withstand commonly practiced chemical and steam sterilization cycles (134°C / 275°F).
- Available with or without an adjustable retention strap

## PRODUCT SELECTION

ITEM #	TYPE	DESCRIPTION	PACKAGING
<b>B1</b>	Single-Use <a href="#">Pediatric</a>	16 mm / 48FR, includes strap	100 / Box
<b>B1-25</b>	Single-Use <a href="#">Pediatric</a>	16 mm / 48FR, includes strap	25 / Box
<b>B1-N</b>	Single-Use <a href="#">Pediatric</a>	16 mm / 48FR, without strap	100 / Box
<b>B3</b>	Single-Use	20 mm / 60FR, includes strap	100 / Box
<b>B3-N</b>	Single-Use	20 mm / 60FR, without strap	100 / Box
<b>B3-SA</b>	Single-Use	20 mm / 60FR, with strap attached to one side	100 / Box
<b>B3-NDR</b>	Single-Use	20 mm / 60FR, no dental retention rim	100 / Box
<b>B3-NDR-N</b>	Single-Use	Same as B3-N, no dental retention rim, without strap.	100 / Box
<b>B4</b>	Single-Use	19 mm / 60FR	100 / Box
<b>B4T</b>	Single-Use	19 mm / 60FR, includes strap	100 / Box
<b>B4T-N</b>	Single-Use	19 mm / 60FR, without strap	100 / Box
<b>B6</b>	Reusable	19 mm / 60FR	100 / Box
<b>B6T</b>	Reusable	19 mm / 60FR, includes strap	100 / Box
<b>B6T-N</b>	Reusable	19 mm / 60FR, without strap	100 / Box
<b>B6-25</b>	Reusable	19 mm / 60FR	25 / Box
<b>B6T-25</b>	Reusable	19 mm / 60FR, includes strap	25 / Box
<b>B6T-N-25</b>	Reusable	19 mm / 60FR, without strap	25 / Box
<b>BS</b>	Single-Use	Bite-Block straps 21.7" x 0.79" / 55 cm x 2 cm	100 / Box



**B5**



**B1**



**B3**



**B4**



**B6**

# BATRIK COUVRETTE™ VALVE COVERS

Our Couvrette™ Valve Covers eliminate the need for manual cleaning and reprocessing of reusable valves. Our single use Valve Covers effectively reduce cross contamination between patients and also reduce fluid leakage from the parts during procedures.

- Easily and securely fit to all valves
- Reduce the potential for errors
- Single-Use
- Non-Sterile
- Latex free



TIP: REUSABLE VALVE COVERS ARE NOTORIOUSLY DIFFICULT TO CLEAN.

SAVE TIME AND REDUCE RISK WITH COUVRETTE™.

COMING SOON!



## PRODUCT SELECTION

ITEM #	DESCRIPTION	COLOR	PACKAGING
<b>16-0002</b>	Biopsy Valve Cover - Olympus®	Blue	100 / Box
<b>16-0004</b>	Biopsy Valve Cover - Pentax®	Red	100 / Box
<b>16-0022</b>	Biopsy Valve Cover - Olympus® & Fujinon®	Light Blue	100 / Box <b>COMING SOON!</b>
<b>16-0044</b>	Biopsy Valve Cover - Pentax®	Pink	100 / Box <b>COMING SOON!</b>


NOTE: Contact us to verify compatibility with your scope

# BATRIK DAYBRUSH™ CHANNEL BRUSHES

## REUSABLE

Our Reusable DayBrush™ Channel Brushes can be used for cleaning endoscopes for an entire day. In fact, it's been proven to effectively clean up to 10 scopes before the need to discard.\*

- Soft tip ensures easy manipulation within the channel with no risk of damage
- The back-end is sealed and softened to prevent fluid from getting trapped inside the tube. The softening process also prevents damage if the cleaning brush is used in reverse.
- Can be steam sterilized.
- Bristle volume is 20mm in length – nearly double most brushes on the market. A large brush provides more surface contact during the cleaning process.
- Significant per use savings when compared to other reusable or disposable channel cleaning brushes.

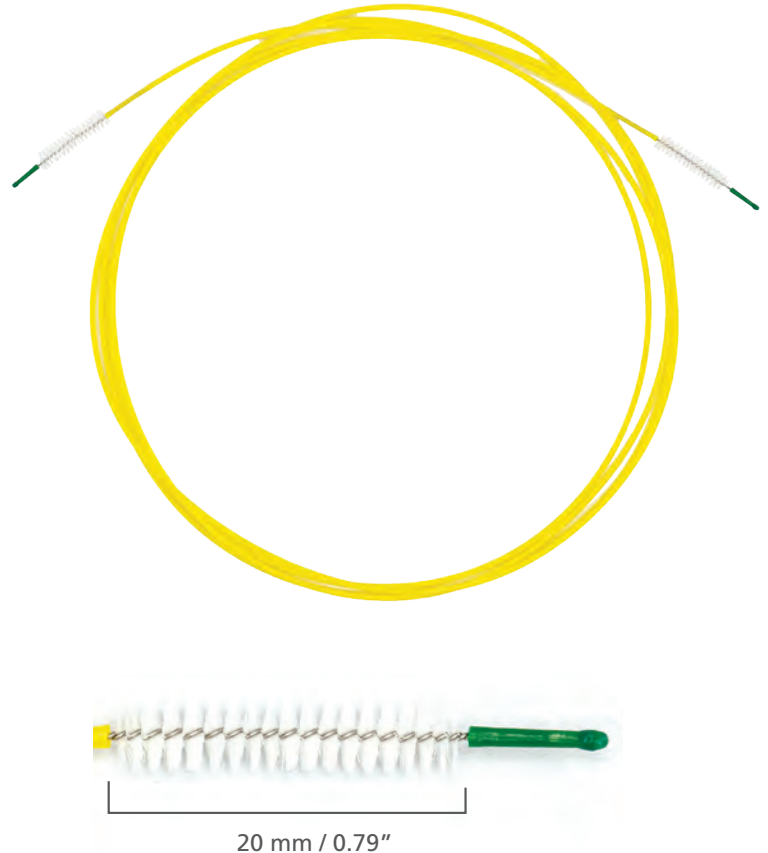

*(\*) Discard after cleaning 10 scopes or sooner if any brush degradation is noticed. After ten applications, the brush may have the appearance that it can continue to clean effectively, but we discourage prolonged use due to risk of cross contamination. Actual number of uses will vary based on handling and processing.*

## PRODUCT SELECTION

### 5 MM CYLINDRICAL BRUSHES – SHEATH Ø 1.75 MM

ITEM #	COLOUR	LENGTH	BRUSH LENGTH	BRUSH (Ø)	SHEATH (Ø)	BRUSH ENDS
<b>30-18050</b>	White	70" / 180 cm	0.79" / 20 mm	0.20" / 5 mm	0.07" / 1.75 mm	Single
<b>30-18050-T</b>	Yellow	70" / 180 cm	0.79" / 20 mm	0.20" / 5 mm	0.07" / 1.75 mm	Single
<b>30-18050-2T</b>	Yellow	70" / 180 cm	0.79" / 20 mm	0.20" / 5 mm	0.07" / 1.75 mm	Double
<b>30-23050</b>	White	90" / 230 cm	0.79" / 20 mm	0.20" / 5 mm	0.07" / 1.75 mm	Single
<b>30-23050-T</b>	Yellow	90" / 230 cm	0.79" / 20 mm	0.20" / 5 mm	0.07" / 1.75 mm	Single
<b>30-23050-2T</b>	Yellow	90" / 230 cm	0.79" / 20 mm	0.20" / 5 mm	0.07" / 1.75 mm	Double

PACKAGING : 50 / BOX

### OPTIONAL KIT ADD-ONS

ADD THE SUFFIX:

"F" to receive a kit with our Flat Sponge (Packaging: Box of 40). Applicable on certain models only, call for details.


"K" to receive a kit with the single use Valve / Control Cleaning Brush 5/11 Double Brushes (item # 27-57787). Applicable on all models.

"FK" to receive both the sponge and valve cleaning brushes.

# BATRIK ARC™ CHANNEL BRUSHES SINGLE-USE

BATRIK offers one of the largest selection of Endoscopy Channel Cleaning Brushes in the industry. Our brushes are available in a variety of lengths, materials, and styles to accommodate the widest range of endoscopic devices used in the healthcare industry today.

No matter what your requirements, you'll be able to find the brush that suits your specific needs. If you don't see the brush you require, simply call us. We have over 50 different brushes available and can custom fabricate to your precise specifications.

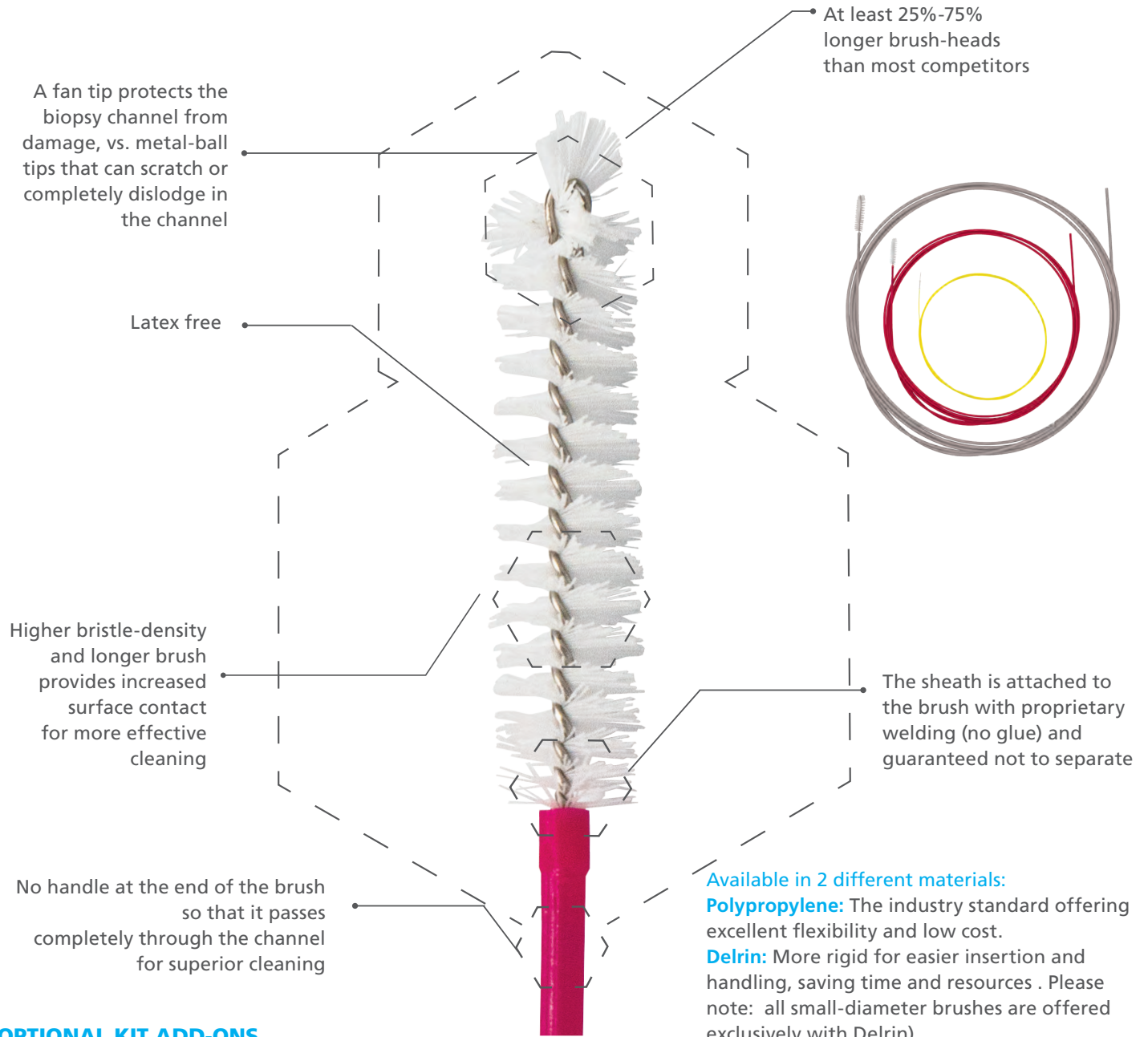

## OPTIONAL KIT ADD-ONS

ADD THE SUFFIX:

"F" to receive a kit with our Flat Sponge (Packaging: Box of 40). Applicable on certain models only, call for details.

"K" to receive a kit with the single use Valve / Control Cleaning Brush 5/11 Double Brushes (item # 27-57787). Applicable on all models.

"FK" to receive both the sponge and valve cleaning brushes.



## SERIES 311 & 333



### PRODUCT SELECTION

#### 5 MM CYLINDRICAL BRUSHES – SHEATH Ø 1.85 MM - PUNCHY PINK POLYPROPYLENE

ITEM #	LENGTH	BRUSH LENGTH	BRUSH (Ø)	SHEATH (Ø)	BRUSH ENDS
<b>311-10050</b>	39" / 100 cm	0.94" / 24 mm	0.20" / 5 mm	0.073" / 1.85 mm	Single
<b>311-18050</b>	70" / 180 cm	0.94" / 24 mm	0.20" / 5 mm	0.073" / 1.85 mm	Single
<b>311-23050</b>	90" / 230 cm	0.94" / 24 mm	0.20" / 5 mm	0.073" / 1.85 mm	Single
<b>311-23050-2</b>	90" / 230 cm	0.94" / 24 mm	0.20" / 5 mm	0.073" / 1.85 mm	Double

PACKAGING: 50 / BOX

#### 5/6 MM TAPERED BRUSHES – SHEATH Ø 1.85 MM - PUNCHY PINK POLYPROPYLENE

ITEM #	LENGTH	BRUSH LENGTH	BRUSH (Ø)	SHEATH (Ø)	BRUSH ENDS
<b>333-23056</b>	90" / 230 cm	0.94" / 24 mm	0.20", 0.24" / 5, 6 mm	0.073" / 1.85 mm	Single
<b>333-23056-2</b>	90" / 230 cm	0.94" / 24 mm	0.20", 0.24" / 5, 6 mm	0.073" / 1.85 mm	Double

PACKAGING: 50 / Box

## SERIES 300 & 322



### PRODUCT SELECTION

#### 5 MM CYLINDRICAL BRUSHES – SHEATH Ø 1.75 MM - METALLIC GREY DELRIN

ITEM #	LENGTH	BRUSH LENGTH	BRUSH (Ø)	SHEATH (Ø)	BRUSH ENDS
<b>300-18050</b>	70" / 180 cm	0.79" / 20 mm	0.20" / 5 mm	0.07" / 1.75 mm	Single
<b>300-23050</b>	90" / 230 cm	0.79" / 20 mm	0.20" / 5 mm	0.07" / 1.75 mm	Single
<b>300-23050-2</b>	90" / 230 cm	0.79" / 20 mm	0.20" / 5 mm	0.07" / 1.75 mm	Double


PACKAGING: 50 / Box

#### 5/6 MM TAPERED BRUSHES – SHEATH Ø 1.75 MM - METALLIC GREY DELRIN

ITEM #	LENGTH	BRUSH LENGTH	BRUSH (Ø)	SHEATH (Ø)	BRUSH ENDS
<b>322-23056</b>	90" / 230 cm	0.79" / 20 mm	0.20", 0.24" / 5, 6 mm	0.07" / 1.75 mm	Single
<b>322-23056-2</b>	90" / 230 cm	0.79" / 20 mm	0.20", 0.24" / 5, 6 mm	0.07" / 1.75 mm	Double

PACKAGING: 50 / Box

## SERIES 344

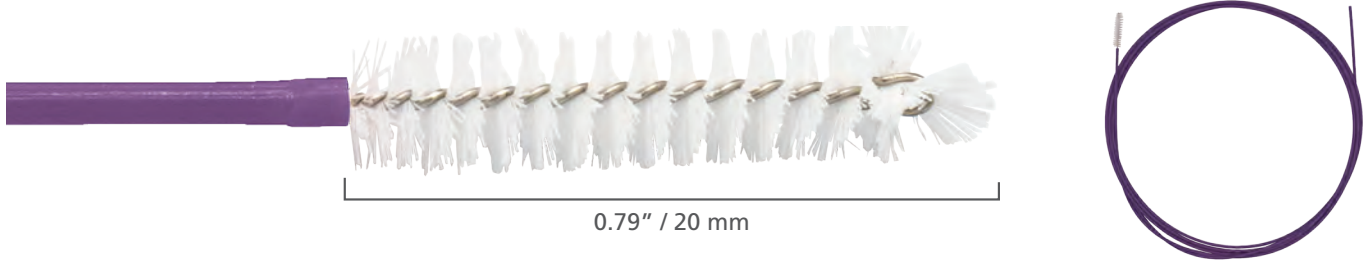


### PRODUCT SELECTION

#### 3 MM CYLINDRICAL BRUSHES – SHEATH Ø 1.55 MM - METALLIC PURPLE DELRIN

ITEM #	LENGTH	BRUSH LENGTH	BRUSH (Ø)	SHEATH (Ø)	BRUSH ENDS
<b>344-10030</b>	39" / 100 cm	0.79" / 20 mm	0.12" / 3 mm	0.06" / 1.55 mm	Single
<b>344-18030</b>	70" / 180 cm	0.79" / 20 mm	0.12" / 3 mm	0.06" / 1.55 mm	Single
<b>344-23030</b>	90" / 230 cm	0.79" / 20 mm	0.12" / 3 mm	0.06" / 1.55 mm	Single

PACKAGING: 50 / Box

## SERIES 35

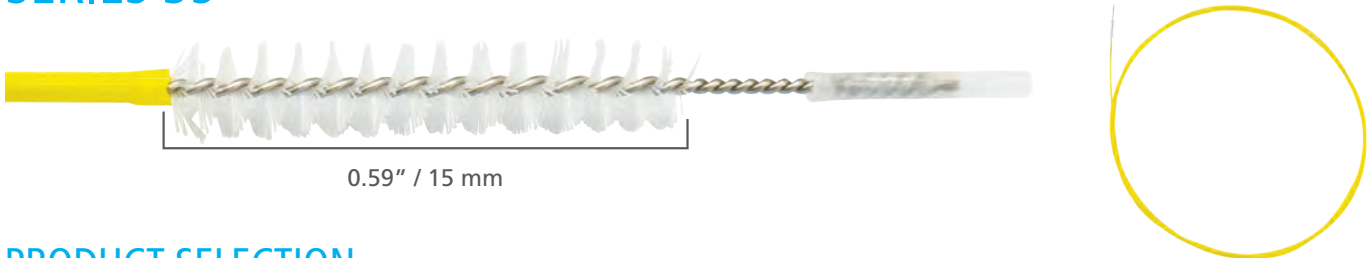

### PRODUCT SELECTION

#### 2.5 MM CYLINDRICAL BRUSHES – SHEATH Ø 1.00 MM - YELLOW DELRIN

ITEM #	LENGTH	BRUSH LENGTH	BRUSH (Ø)	SHEATH (Ø)	BRUSH ENDS
<b>35-12025</b>	47.2" / 120 cm	0.59" / 15 mm	0.098" / 2.5 mm	0.039" / 1.00 mm	Single
<b>35-23025</b>	90" / 230 cm	0.59" / 15 mm	0.098" / 2.5 mm	0.039" / 1.00 mm	Single
<b>35-23025-2</b>	90" / 230 cm	0.59" / 15 mm	0.098" / 2.5 mm	0.039" / 1.00 mm	Double

PACKAGING: 50 / Box

#### 2 MM CYLINDRICAL BRUSHES – SHEATH Ø 1.00 MM - YELLOW DELRIN

ITEM #	LENGTH	BRUSH LENGTH	BRUSH (Ø)	SHEATH (Ø)	BRUSH ENDS
<b>35-12020</b>	47.2" / 120 cm	0.59" / 15 mm	0.079" / 2.0 mm	0.039" / 1.00 mm	Single
<b>35-18020</b>	70.8" / 180 cm	0.59" / 15 mm	0.079" / 2.0 mm	0.039" / 1.00 mm	Single
<b>35-23020</b>	90" / 230 cm	0.59" / 15 mm	0.079" / 2.0 mm	0.039" / 1.00 mm	Single
<b>35-23020-2</b>	90" / 230 cm	0.59" / 15 mm	0.079" / 2.0 mm	0.039" / 1.00 mm	Double


PACKAGING: 50 / Box

# BATRIK VALVE / CONTROL BRUSHES


27-05197



27-57787

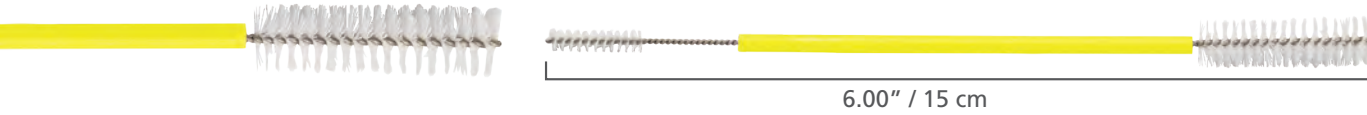

40-18511

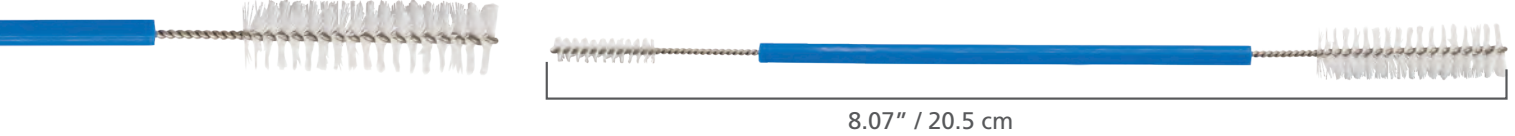

36-10010



27-15311




## PRODUCT SELECTION

ITEM #	DESCRIPTION	LENGTH	BRUSH 1 LENGTH	BRUSH 2 LENGTH	BRUSH 1 (Ø)	BRUSH 2 (Ø)
<b>27-05197</b>	Double Trumpet 5/7 White Color	5.00" / 13 cm	0.20" / 5 mm	1.00" / 25 mm	0.20" / 5 mm	0.28" / 7 mm
<b>PACKAGING:</b> Bag of 3 (single-use recommended). Volume Packaging: Bag of 50 use ITEM # 27-05197-50						
<b>27-57787</b>	Double. Universal 5/11 Yellow Color	6.00" / 15 cm	1.18" / 30 mm	0.79" / 20 mm	0.43" / 11 mm	0.20" / 5 mm
<b>PACKAGING:</b> Bag of 3 (single-use recommended). Volume Packaging: Bag of 50 use ITEM # 27-57787-50						
<b>40-18511</b>	Double Universal 5/11 Extra Long Blue Color	8.07" / 20.5 cm	1.57" / 40 mm	0.79" / 20 mm	0.43" / 11 mm	0.20" / 5 mm
<b>PACKAGING:</b> Bag of 3 (reusable and can be steam sterilized). Volume Packaging: Bag of 50 use ITEM # 40-18511-50						
<b>36-10010</b>	Single end large Red Color	4.3" / 11 cm	1.57" / 40 mm	N/A	0.39" / 10 mm	N/A
<b>PACKAGING:</b> Bag of 3 (single-use recommended). Volume Packaging: Bag of 50 use ITEM # 36-10010-50						
<b>27-15311</b>	Single end large Blue Color	4" / 10.3 cm	1.57" / 40 mm	N/A	0.39" / 10 mm	N/A
<b>PACKAGING:</b> Bag of 3 (single-use recommended). Volume Packaging: Bag of 50 use ITEM # 27-15311-50						

# BATRIK V-BRUSHES EXTERNAL CLEANING BRUSHES

Our External Cleaning Brushes for endoscopes / optical instruments provide excellent cleaning capabilities. With their unique U-shape bristles and innovative design, our brushes are the ultimate tool for cleaning the exterior of your endoscopic instruments. Free from sharp or metallic ends that can scratch or cause damage to valuable equipment, both single-use and reusable models are latex free and feature ergonomic handles.

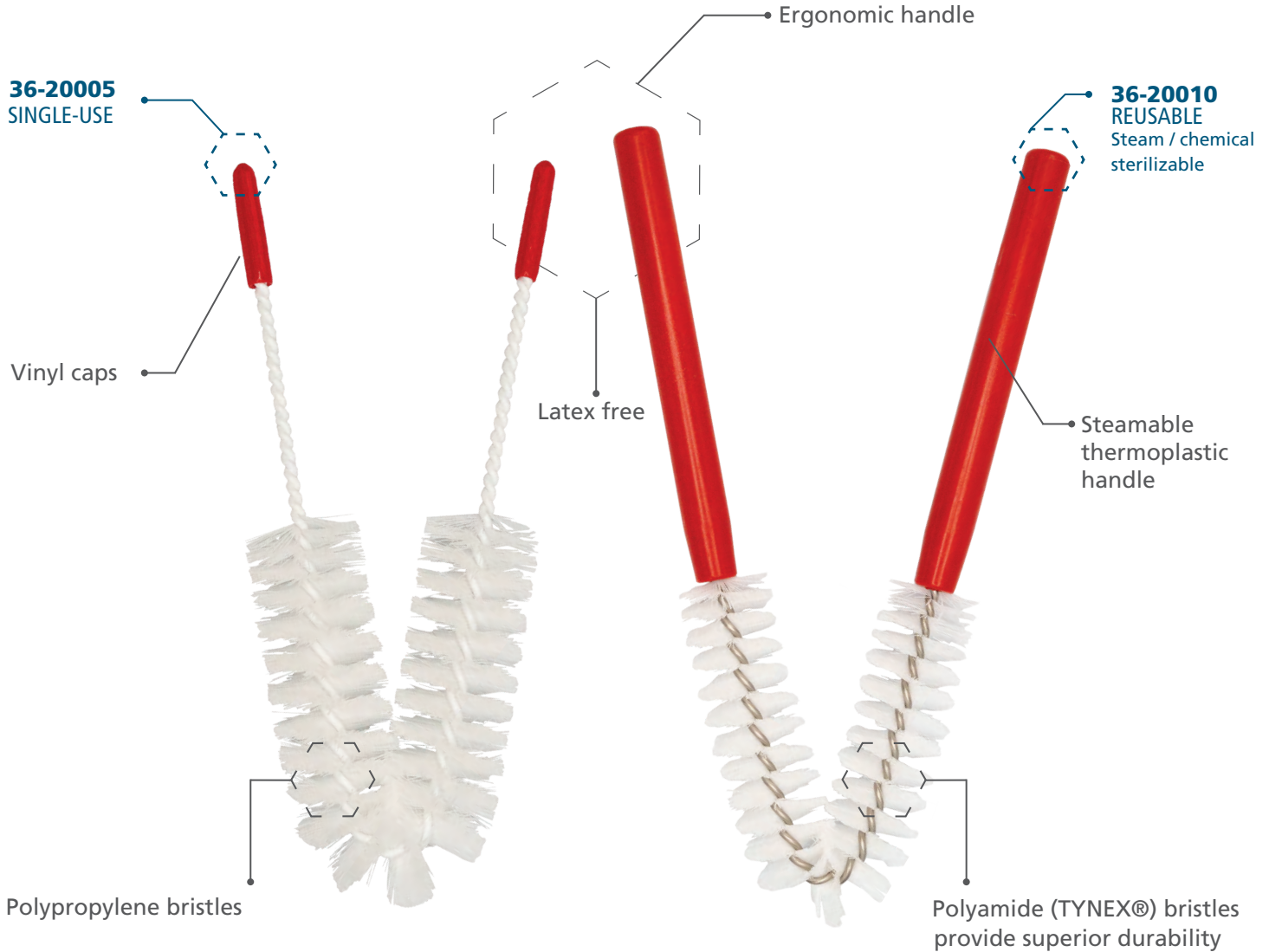

## PRODUCT SELECTION

ITEM #	DESCRIPTION	BRUSH LENGTH / (Ø)	HANDLE LENGTH / (Ø)	TOTAL LENGTH	PACKAGING
<b>36-20005</b>	Single-Use External Cleaning Brush	70 mm / 25 mm	80 mm / 3 mm	150 mm	50 / Box
<b>36-20010</b>	Reusable External Cleaning Brush	90 mm / 25 mm	120 mm / 12 mm	210 mm	3 / Bag

**NEW!**  
**BATRIK**

# TAG-IT™ | NOTES

## WHEN SIMPLY BROKEN DOESN'T CUT IT—TAG IT, TELL IT, OWN IT!

Tag-It™ Notes offer the freedom to **identify and document instructions** for instruments, trays, and carts before and after **sterilization**. Tag-It™ Notes are **tear-resistant, withstand extreme chemicals, high temperature** and ensures **reliable communication**.

### WHY YOU'LL LOVE THEM:

- Tag-It™ Notes gives you the **freedom** to clearly document important information **without the clutter**.
- Robust and easy-to-use **tie wraps** designed to **securely attach** to instruments (134°C & chemical tested).
- Includes a **fast-drying, heat-resistant permanent marker** for writing your important notes.
- Choose from other **customizable cards** for personalized needs or **pre-made cards** for fast, convenient use.
- When you need more than a "repair" or "sharpen" seal to identify an issue. Instead, **Tag it, Tell it, Own it.**

**WATERPROOF**  
The synthetic paper withstands spills and liquid exposure without any damage.

**ECO-FRIENDLY**  
Made from sustainable materials.

**HEAT RESISTANT**  
Can resist high temperatures of up to 273.2°F (134°C). Compatible with Gas Plasma Sterilization (STERAD®100 NX Sterilization).

**DURABLE & TEAR-PROOF**  
Handles rigorous cleaning, tough conditions, and resists tearing.

### PRODUCT SELECTION

ITEM #	COLOR	SIZE	PACKAGING	
TI-NOTE-R0	Yellow	NOTES	3" x 5" / 8 cm x 13cm	100 / PKG


All Tag-It™ are single-use. PACKAGING: 1 pack of 100 tags, 1 permanent marker, and 100 ties.

TAG-IT™

Tag-It™ Cards are an effortless way to add information with either pre-printed or blank formats. These **durable** tags streamline your workflow, **withstand tough conditions**, and feature **(optional) indicator dots** for quick verification. They help keep **instruments, trays, and carts** organized, ensuring efficiency before and after sterilization/decontamination.

KEY PRODUCT FEATURES:

- **Special Permanent Marker:** The ink will not run or appear blurry after exposure to chemicals or high temperature (max. 260°C/500°F) sterilization. *(Regular pens or markers can be used only if the tags are not going to be exposed to any sterilization/decontamination).*
- **Tie Wraps:** Assist in securing the tags to almost any instrument or container. *(Tie wraps have been tested for chemical and steam sterilization).*
- **Our Tag-It™:** Will not shrink, absorb water/chemicals, or lose its intended integrity after chemical or steam sterilization.

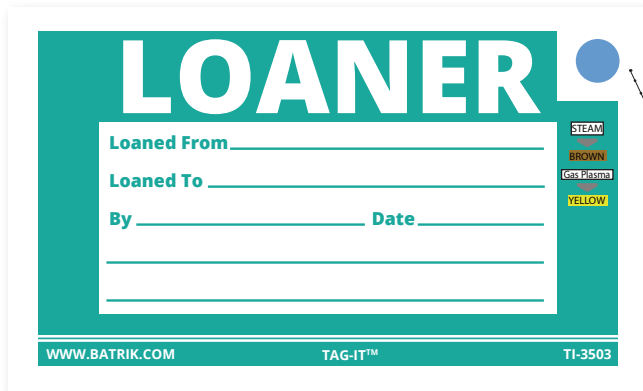

**(Optional)  
All Tags Available  
with Indicator Dots**  
Provides visual confirmation that  
the tag was processed.

CLICK BELOW TO  
WATCH OUR VIDEO



PRODUCT SELECTION:

ITEM #	COLOR	DESCRIPTION	SIZE	PACKAGING
TI-3502	Blue	VENDOR	3" x 5" / 8 cm x 13cm	100 / PKG
TI-3503	Green	LOANER	3" x 5" / 8 cm x 13cm	100 / PKG
TI-3504	Red	PRIORITY	3" x 5" / 8 cm x 13cm	100 / PKG
TI-3505	Yellow	RE-RUN	3" x 5" / 8 cm x 13cm	100 / PKG
TI-3506	Pink	MISSING	3" x 5" / 8 cm x 13cm	100 / PKG
TI-3507	Orange	BROKEN	3" x 5" / 8 cm x 13cm	100 / PKG
TI-3509	Black	BLANK	3" x 5" / 8 cm x 13cm	100 / PKG
TI-35B1	White	BLANK	3" x 5" / 8 cm x 13cm	100 / PKG
TI-35B4	Red	BROKEN	3" x 5" / 8 cm x 13cm	100 / PKG
TI-35C22	Light Blue	VENDOR CHECK IN + OUT	3" x 5" / 8 cm x 13cm	100 / PKG
TI-35L3	Green	LOANER CHECK IN + OUT	3" x 5" / 8 cm x 13cm	100 / PKG
TI-35T2	Blue	QUICK TURNOVER	3" x 5" / 8 cm x 13cm	100 / PKG
TI-PM09	Black	FINE PERMANENT MARKERS	1.5" Ø x 5.5" / 3.81cm Ø x 12.7 cm	36 / Bag

PACKAGING: 1 pack of 100 tags, 100 tie wraps, and 1 permanent marker.

OPTIONAL | PROCESSING INDICATOR DOTS: Add suffix "ST" to end of code for steam dots, "PS" for plasma dots.

Available in other languages or customize your own.

PRODUCT SELECTION: (customization also available)


**VENDOR / BLUE**

ITEM #	SIZE
TI-3502	3" x 5" (8 cm x 13cm)



**NOTES / YELLOW**

ITEM #	SIZE
TI-NOTE-RO	3" x 5" (8 cm x 13cm)

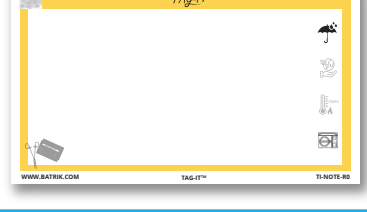

**LOANER / GREEN**

ITEM #	SIZE
TI-3503	3" x 5" (8 cm x 13cm)




**BLANK / WHITE**

ITEM #	SIZE
TI-35B1	3" x 5" (8 cm x 13cm)




**PRIORITY / RED**

ITEM #	SIZE
TI-3504	3" x 5" (8 cm x 13cm)

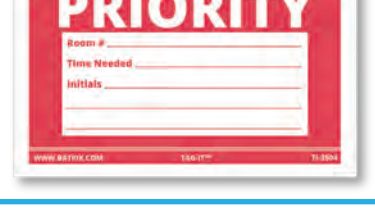


**BROKEN / RED**

ITEM #	SIZE
TI-35B4	3" x 5" (8 cm x 13cm)

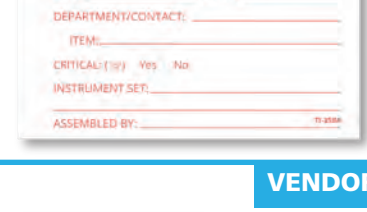

**RE-RUN / YELLOW**

ITEM #	SIZE
TI-3505	3" x 5" (8 cm x 13cm)

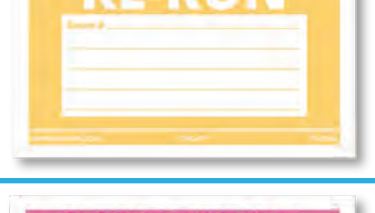


**VENDOR IN + OUT / LIGHT BLUE**

ITEM #	SIZE
TI-35C22	5" x 3" (13 cm x 8 cm)



**MISSING / PINK**

ITEM #	SIZE
TI-3506	3" x 5" (8 cm x 13cm)

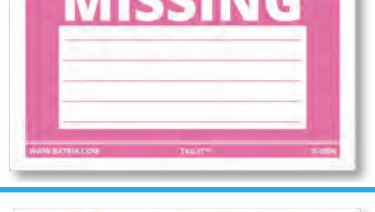


**LOANER IN + OUT / GREEN**

ITEM #	SIZE
TI-35L3	5" x 3" (13 cm x 8 cm)



**BROKEN / ORANGE**

ITEM #	SIZE
TI-3507	3" x 5" (8 cm x 13cm)

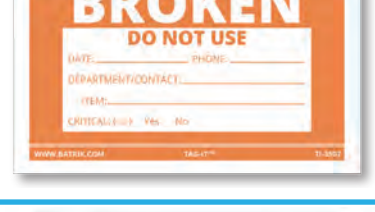

**BLANK / BLACK**

ITEM #	SIZE
TI-3509	3" x 5" (8 cm x 13cm)



**QUICK TURNOVER/ BLUE**

ITEM #	SIZE
TI-35T2	3" x 5" (8 cm x 13cm)



All information tags are single-use. PACKAGING: 1 pack of 100 tags, 100 tie wraps, and 1 permanent marker.

NEW!

BATRIK

# TAG-IT™ | SCOPE

Batrik Scope Tags are tear-resistant, double-sided tags that will prove indispensable throughout the reprocessing journey. Our tags make it easy to **record important information**, including processing dates, with a user-friendly feature that allows you to select the month, day, and date. Crafted from **waterproof synthetic paper**, these tags feature a **self-looping design**, making it effortless to attach them to the scope without the need for additional fastening materials. Elevate your reprocessing protocols with Batrik!

### KEY FEATURES

- The paper, tie wraps, and smudge-proof ink are designed to withstand high temperatures of up to 273.2°F (134°C) and are compatible with Gas plasma sterilization (STERRAD®100 NX Sterilization System).
- Tags include the month, day, and date selection on one side, making it convenient to mark the processing date accurately



### COLOR-CODED AND CUSTOMIZABLE

- Available in seven colors
- The tags facilitate color coding for easy identification

### SELF LOOPING

- The self-looping design allows for easy attachment to the scope
- Eliminating the need for additional fastening materials


### SINGLE-USE AND LATEX FREE

- Batrik Scope Tags are designed for single use and are latex-free
- Promoting hygiene and reducing the risk of cross-contamination

### PRODUCT SELECTION

ITEM #	COLOR	SIZE	PACKAGING
TI-ST05	Yellow	11.8" x 1.55" / 33.8 cm x 4 cm	100 / PKG
TI-ST02	Blue	11.8" x 1.55" / 33.8 cm x 4 cm	100 / PKG
TI-ST03	Green	11.8" x 1.55" / 33.8 cm x 4 cm	100 / PKG
TI-ST07	Orange	11.8" x 1.55" / 33.8 cm x 4 cm	100 / PKG
TI-ST06	Purple	11.8" x 1.55" / 33.8 cm x 4 cm	100 / PKG
TI-ST04	Red	11.8" x 1.55" / 33.8 cm x 4 cm	100 / PKG
TI-ST01	White	11.8" x 1.55" / 33.8 cm x 4 cm	100 / PKG

All information tags are single-use. PACKAGING: 1 pack of 100 tags and 1 permanent marker.

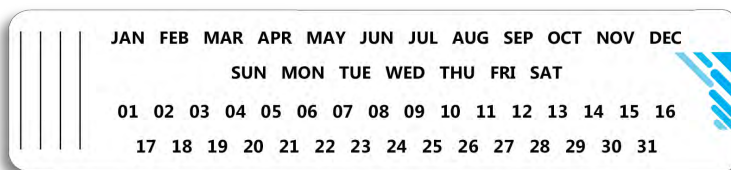




**NEW!**  
**BATRIK**

# TAG-IT™ | SCOPE REPROCESSING

Batrik Scope Reprocessing Tags are durable and can withstand handling during reprocessing **without tearing**, ensuring that information remains intact. These tags are vital in recording important information such as **date to reprocess on, scope model and serial number, sterilizer type, meeting sterilizer parameter conditions, alcohol flush, and hang date**. These tags have been meticulously designed to elevate **operational efficiency** and contribute to the delivery of the highest standard of patient care.



**TI-RST**

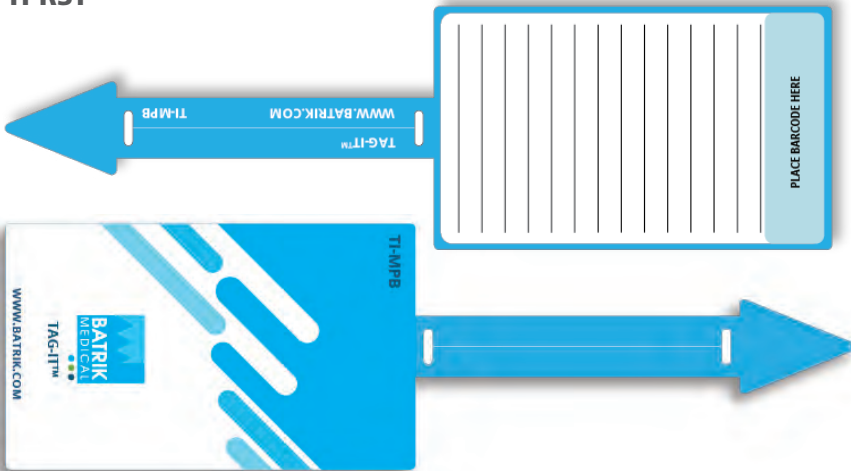

**TI-MPB**

\*Suitable for various information tagging purposes

### LOGGING IMPORTANT INFORMATION

- Easy for staff to indicate when the scope was processed
- Key scope, technician, and cycle-related information easily reported and identified.

### SMART DESIGN

- The self-looping design allows for easy attachment to the scope
- No requirement for extra fastening materials
- Secure fastening as the loop is tightened around the instrument
- Features an integrated barcode area



## PRODUCT SELECTION



ITEM #	DESCRIPTION	DIMENSIONS	PACKAGING
<b>TI-RST</b>	REPROCESSING SCOPE TAG	3" x 10.5" / 8 cm x 27 cm	100 / PKG
<b>TI-MPB</b>	MULTI-PURPOSE BLANK TAG	3" x 10.5" / 8 cm x 27 cm	100 / PKG

All information tags are single-use. **PACKAGING: 1 pack of 100 tags and 1 permanent marker.**


# BATRIK SILICONE PIN MAT

Our Silicone Pin Mat cushions, protects, and organizes delicate surgical instruments for sterilization. The elevated nubs provide greater air-flow ventilation and drying, significantly reducing residual moisture collection.

- Reusable
- Autoclavable
- Latex free

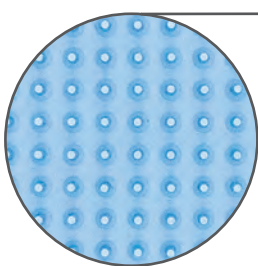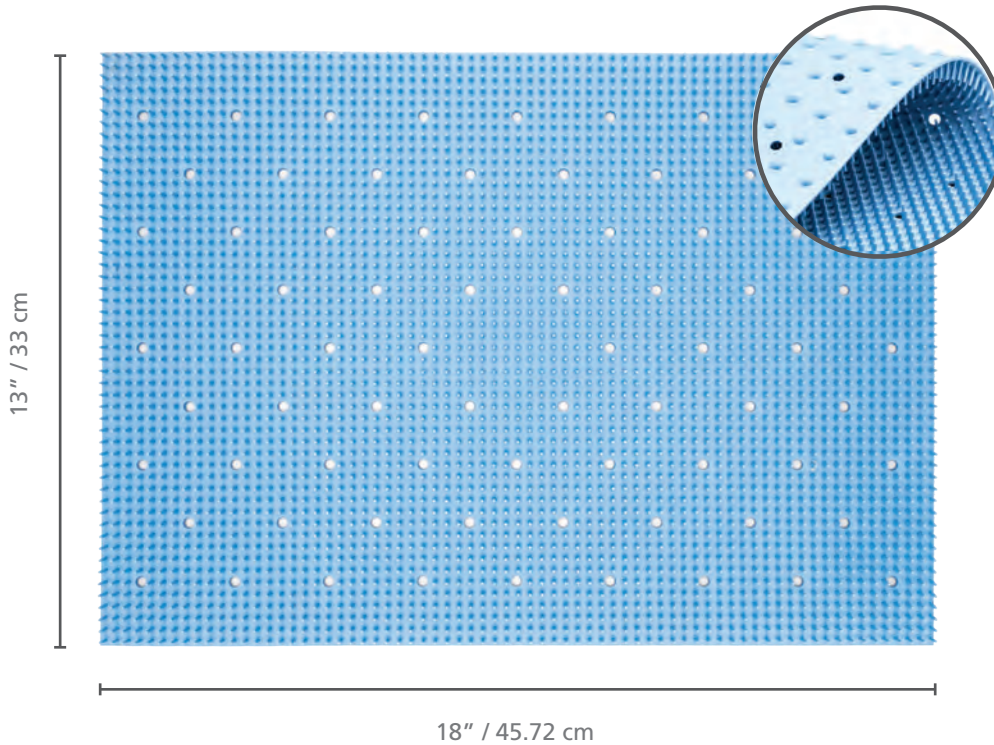

**FRONT SIDE**  
Pins keep instruments securely in place

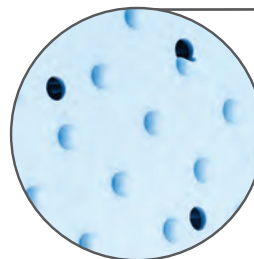

**BACK SIDE**  
Numerous holes and elevated nubs provide superior venting and reduce water collection

## PRODUCT SELECTION

ITEM #	DIMENSIONS	PACKAGING
<b>13-2001</b>	13" x 18" (33 cm x 45.72 cm) 11.97 mm high from base of mat to end of pin 15.53 mm high from base of elevated nub to end of pin	1 unit
<b>13-2005</b>	13" x 18" (33 cm x 45.72 cm) 11.97 mm high from base of mat to end of pin 15.53 mm high from base of elevated nub to end of pin	5 units


GOOD TO KNOW:  
CUSTOM SIZES  
ARE AVAILABLE  
UPON REQUEST

# BATRIK CORNER POCKETS™


With three models to choose from, our Corner Pockets™ portfolio is a truly comprehensive line of products designed to prevent tearing associated with sharp corners or feet on surgical trays.



18-CP02

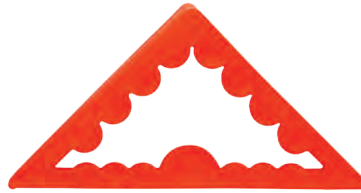

18-CP07

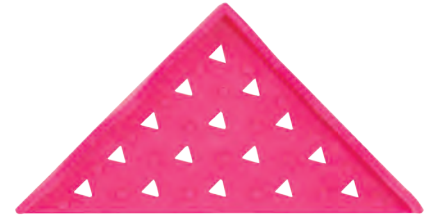


18-CP06

COLOUR	BLUE	ORANGE	PINK
USAGE	SINGLE-USE	REUSABLE <i>withstands repeated sterilization up to 134°C</i>	REUSABLE <i>withstands repeated sterilization up to 134°C</i>
DIMENSIONS	2.8" x 4" (7.3 cm x 10.3 cm)	2.8" x 4" (7.3 cm x 10.3 cm)	3.4" x 4.7" (8.5 cm x 12 cm)
LATEX-FREE	✓	✓	✓
NO MOISTURE RETENTION	✓	✓	✓
FABRICATION	Flexible, silicone-like plastic	Flexible, silicone-like plastic	Our most flexible plastic for easiest placement and most secure protection

# SOFT POCKETS™

Our Soft Pockets™ Foam Corner Protectors are made from high density foam to effectively cushion and protect the corners of surgical trays to prevent tearing. They are especially useful with heavier instrument sets.

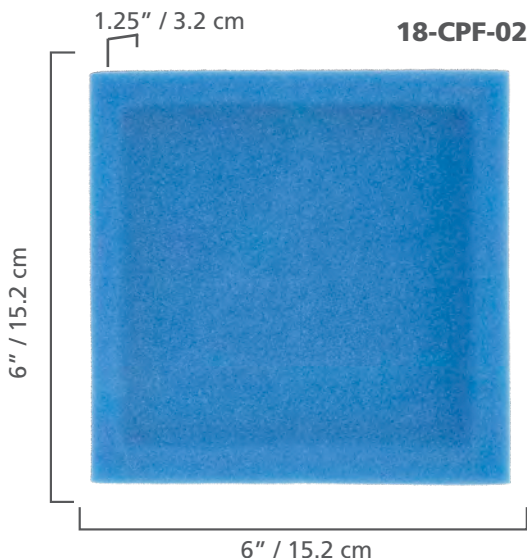

18-CPF-02

## PRODUCT SELECTION

ITEM #	DESCRIPTION	PACKAGING
18-CP02**	Corner Pockets™ <i>Single-Use, Blue</i>	240 / Bag
18-EVAL2**	Corner Pockets™ <i>Single-Use, Blue</i>	24 / Bag
18-CP07	Corner Pockets™ <i>Reusable, Orange</i>	80 / Bag
18-CORNER	Corner Pockets™ <i>Reusable, Orange</i>	320 / Box (4 Bags of 80 pcs)
18-EVAL7	Corner Pockets™ <i>Reusable, Orange</i>	24 / Bag
18-CP06	Corner Pockets™ <i>Reusable, Pink</i>	80 / Bag
18-CPF-02	Soft Pockets™ <i>Single-Use Foam Corner Protectors</i>	80 / Box (24 Boxes / Pallet)

\*\* Although the single use Corner Pockets™ appear similar to the reusable models they are blended differently and have only been validated as single-use. Re-using them may compromise sterility.

BATRIK

# LILA'S INSTRUMENT STAIN REMOVER

Lila's special blend of naturally sourced powders provides outstanding removal of stubborn stains.

Convenient-to-use formula, you can achieve spotless results every time!

It can be used on surgical/dental instruments, autoclaves, ultrasonic cleaners, stainless steel countertops, trays, bowls and more.

Specifically designed for use on stainless steel, chrome, nickel, brass, copper, and silver.



LSR-8501


## EASE TO USE AS 1,2,3...

1. Dampen a cloth or grab a Lila's Low Linting Dura-Towel and dab into the jar.
2. Gently rub the area that is stained.
3. Rinse thoroughly for a flawless shine.<sup>2</sup>

## A little Goes a long Way!

Even great to tackle stubborn stains on steel pots and pans in your kitchen with ease.

Lila's - If you have used Spectra-Scrub or Integra's instrument stain remover then you will love Lila's!

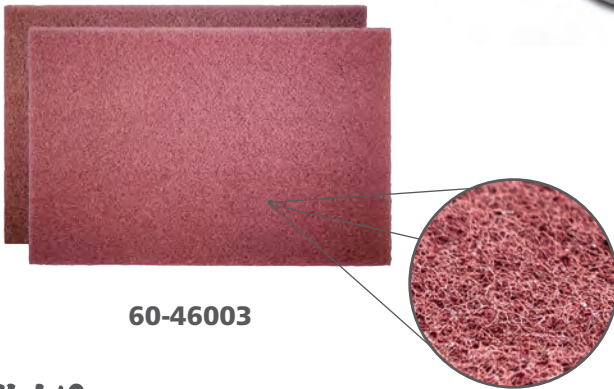

60-46003



17-DT1213C

## LILA'S SCRUB PAD

- Say goodbye to extra stubborn stains with Lila's Scrub Pad without leaving any metal residue.
- Simply add a small amount of Lila's Stain Remover with one of our Scrub Pad's for very large areas.
- The gentle yet effective scrub is the perfect choice for Ophthalmic Instruments.<sup>1</sup>

## LILA'S LOW-LINTING DURA-TOWEL

- Lila's Dura-Towel has long lasting strength and has a minimal lint composition.
- These single-use towels are clean room certified and are specifically designed for instrument wiping.
- Its ideal dimensions and unique attributes guarantee optimal performance, setting it apart from ordinary paper towels.

## PRODUCT SELECTION




ITEM #	DESCRIPTION	DIMENSIONS	PACKAGING
<b>LSR-8501</b>	LILA'S INSTRUMENT STAIN REMOVER 3 oz (85g) plastic jar with a convenient screw-on lid	Ø 2.8" / 70 mm	1 Jar / Bag
<b>LSR-8512</b>	LILA'S INSTRUMENT STAIN REMOVER 3 oz (85g) plastic jar with a convenient screw-on lid	Ø 2.8" / 70 mm	12 Jars / Box
<b>17-DT1213C</b>	LILA'S LOW LINTING DURA-TOWEL 1080 towels	12" x 13" / 31 cm x 33 cm Folded: 3" x 3.25" / 7.5 cm x 8.25 cm	90 / Bag 12 Bags / Box
<b>60-46003</b>	LILA'S SCRUB PAD	4" x 6" / 10 cm x 15 cm	3 / Bag
<b>LSR-BUNDLE</b>	LILA'S STAIN REMOVER + LILA'S DURA-TOWEL 3 oz (85g) plastic jar + 90 towels	Ø 2.8" / 70 mm 12" x 13" / 31 cm x 33 cm	1 Jar / Bag 90 / Bag

1 : Although Lila's is completely safe, it is not recommended for use on ophthalmic instruments. Use Lila's Scrub Pad on its own.  
2: The stain remover should not be mixed with bleach and is for external use only to remove instrument stains.

# BATRIK INSTRUMENT PROTECTORS

Our Instrument Protectors protect the tips and edges of instruments to prevent damage during handling and sterilization. Available in a wide array of color coded sizes for almost every application, including osteotome, retractor, double, and perforated.

- Steam and gas permeable
- Perforated protectors reduce vacuum seal and moisture build-up
- Validated for pre-vacuum sterilization\*
- Multiple color-coded sizes
- Latex Free
- Non-Sterile

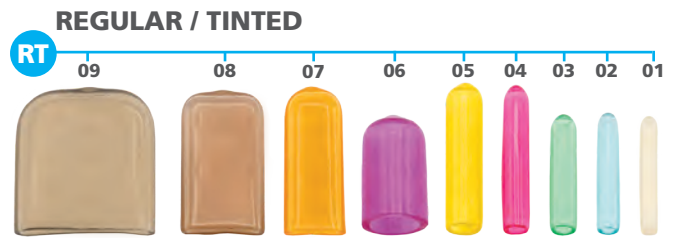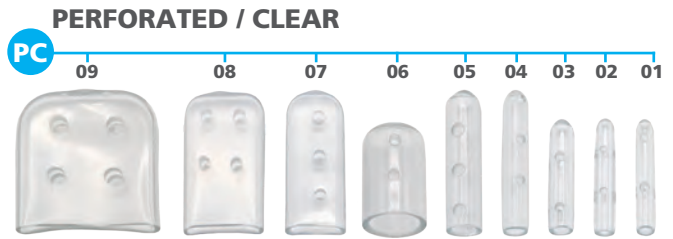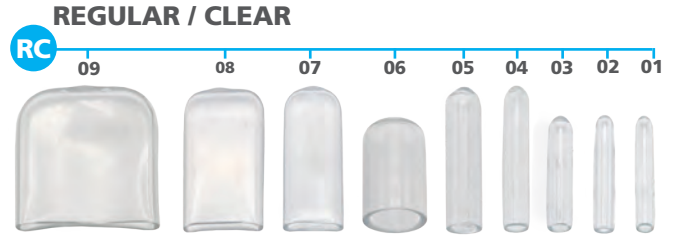

\*Independent laboratory tested to validate pre-vacuum steam sterilization efficacy in a healthcare setting


## PRODUCT SELECTION

	ITEM # REGULAR	ITEM # PERFORATED	COLOR	DIMENSIONS
<b>REGULAR &amp; PERFORATED / CLEAR</b>				
 <b>RC-07</b>  <b>PC-07</b>	<b>RC-01</b>	<b>PC-01</b>	Clear	0.062" x 0.75" (1.57 mm x 19 mm)
	<b>RC-02</b>	<b>PC-02</b>	Clear	0.075" x 0.75" (2 mm x 19 mm)
	<b>RC-03</b>	<b>PC-03</b>	Clear	0.109" x 0.75" (2.77 mm x 19 mm)
	<b>RC-04</b>	<b>PC-04</b>	Clear	0.125" x 1" (3.18 mm x 25.4 mm)
	<b>RC-05</b>	<b>PC-05</b>	Clear	0.187" x 1" (4.75 mm x 25.4 mm)
	<b>RC-06</b>	<b>PC-06</b>	Clear	0.406" x 0.75" (10.3 mm x 19 mm)
	<b>RC-07</b>	<b>PC-07</b>	Clear	0.062" x 0.375" x 1" (1.57 mm x 9.53 mm x 25.4 mm)
	<b>RC-08</b>	<b>PC-08</b>	Clear	0.078" x 0.625" x 1" (2 mm x 15.875 mm x 25.4 mm)
	<b>RC-09</b>	<b>PC-09</b>	Clear	0.125" x 1" x 1" (3.18 mm x 25.4 mm x 25.4 mm)
	<b>RC-10</b>	<b>PC-10</b>	Clear	Assorted
PACKAGING: Bag of 100				
<b>REGULAR &amp; PERFORATED / TINTED</b>				
 <b>RT-07</b>  <b>PT-07</b>	<b>RT-01</b>	<b>PT-01</b>	Tinted White	0.062" x 0.75" (1.57 mm x 19 mm)
	<b>RT-02</b>	<b>PT-02</b>	Tinted Blue	0.075" x 0.75" (2 mm x 19 mm)
	<b>RT-03</b>	<b>PT-03</b>	Tinted Green	0.109" x 0.75" (2.77 mm x 19 mm)
	<b>RT-04</b>	<b>PT-04</b>	Tinted Red	0.125" x 1" (3.18 mm x 25.4 mm)
	<b>RT-05</b>	<b>PT-05</b>	Tinted Yellow	0.187" x 1" (4.75 mm x 25.4 mm)
	<b>RT-06</b>	<b>PT-06</b>	Tinted Maroon	0.406" x 0.75" (10.3 mm x 19 mm)
	<b>RT-07</b>	<b>PT-07</b>	Tinted Orange	0.062" x 0.375" x 1" (1.57 mm x 9.53 mm x 25.4 mm)
	<b>RT-08</b>	<b>PT-08</b>	Tinted Brown	0.078" x 0.625" x 1" (2 mm x 15.875 mm x 25.4 mm)
	<b>RT-09</b>	<b>PT-09</b>	Tinted Black	0.125" x 1" x 1" (3.18 mm x 25.4 mm x 25.4 mm)
	<b>RT-10</b>	<b>PT-10</b>	Tinted Assorted	Assorted
PACKAGING: Bag of 100				
<b>REGULAR &amp; PERFORATED / SOLID</b>				
 <b>RS-07</b>  <b>PS-07</b>	<b>RS-01</b>	<b>PS-01</b>	Solid White	0.062" x 0.75" (1.57 mm x 19 mm)
	<b>RS-02</b>	<b>PS-02</b>	Solid Blue	0.075" x 0.75" (2 mm x 19 mm)
	<b>RS-03</b>	<b>PS-03</b>	Solid Green	0.109" x 0.75" (2.77 mm x 19 mm)
	<b>RS-04</b>	<b>PS-04</b>	Solid Red	0.125" x 1" (3.18 mm x 25.4 mm)
	<b>RS-05</b>	<b>PS-05</b>	Solid Yellow	0.187" x 1" (4.75 mm x 25.4 mm)
	<b>RS-06</b>	<b>PS-06</b>	Solid Maroon	0.406" x 0.75" (10.3 mm x 19 mm)
	<b>RS-07</b>	<b>PS-07</b>	Solid Orange	0.062" x 0.375" x 1" (1.57 mm x 9.53 mm x 25.4 mm)
	<b>RS-08</b>	<b>PS-08</b>	Solid Brown	0.078" x 0.625" x 1" (2 mm x 15.875 mm x 25.4 mm)
	<b>RS-09</b>	<b>PS-09</b>	Solid Black	0.125" x 1" x 1" (3.18 mm x 25.4 mm x 25.4 mm)
	<b>RS-10</b>	<b>PS-10</b>	Solid Assorted	Assorted
PACKAGING: Bag of 100				

ITEM #	DIMENSIONS
OP-11	0.25" x 0.75" (6.35 mm x 19.05 mm)
OP-12	0.375" x 1" (9.53 mm x 25.4 mm)
OP-13	0.5" x 1" (12.7 mm x 25.4 mm)
OP-14	0.75" x 1" (19.05 mm x 25.4 mm)
OP-15	1" x 1" (25.4 mm x 25.4 mm)
OP-16	1.25" x 1.25" (31.75 mm x 31.75 mm)
OP-17	1.5" x 1.5" (38.1 mm x 38.1 mm)
OP-18	1.75" x 1.25" (44.45 mm x 31.75 mm)
OP-20	Assorted

PACKAGING: Bag of 50

**OSTEOTOME PROTECTORS**

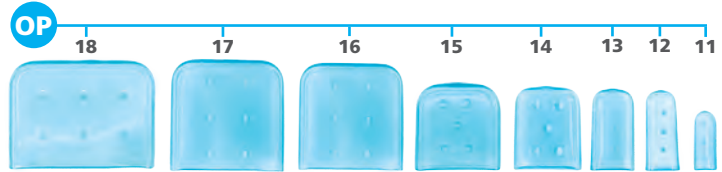

ITEM #	DIMENSIONS	PACKAGING
RGS-01	1.5" (38 mm)	55 / Bag
RGS-02	2" (51 mm)	45 / Bag
RGS-03	3" (76.2 mm)	40 / Bag


**SINGLE-TIP RETRACTOR PROTECTORS**



ITEM #	COLOR	SIZE
RGM-01	Pink	Small
RGM-02	Clear	Small
RGM-03	Green	Medium
RGM-04	Clear	Medium

PACKAGING: Bag of 50


**WEITLANER RETRACTOR PROTECTORS**

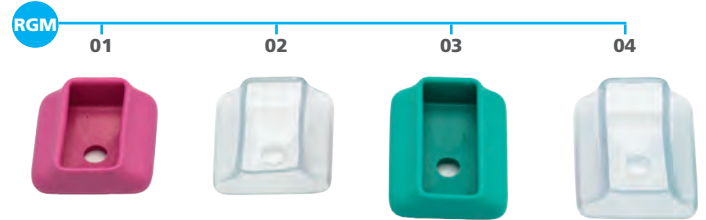

ITEM #	COLOR	DIMENSIONS
DG-01	Clear	0.0625" x 0.75" (1.59 mm x 19.05 mm)
DG-02	Blue	0.078" x 0.75" (1.98 mm x 19.05 mm)
DG-03	Green	0.109" x 0.75" (2.78 mm x 19.05 mm)
DG-04	Red	0.125" x 1" (3.18 mm x 25.4 mm)
DG-05	Yellow	0.187" x 1" (4.76 mm x 25.4 mm)
DG-07	Orange	0.062" x 0.375" x 1" (1.59 mm x 9.53 mm x 25.4 mm)
DG-10	Assorted	Assorted

PACKAGING: Bag of 50

**PERFORATED DOUBLE PROTECTORS**

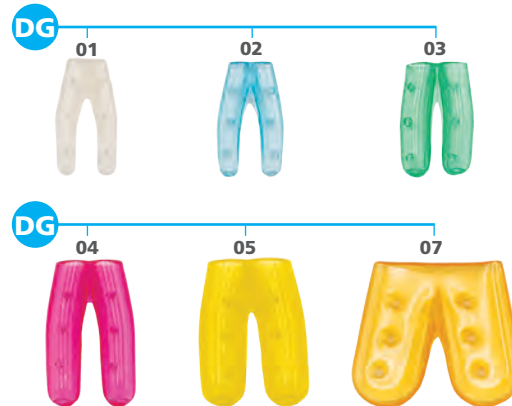

# BATRIK PrAmLi Power<sup>3</sup> LUBE™ Instrument Lubricant

Specially formulated to lubricate surgical instruments, hinges, box locks, and general equipment during autoclaving. PrAmLi Power<sup>3</sup> LUBE™ is uniquely formulated with an all natural corrosion inhibitor.

- Forms a protective coating which can extend the instrument's, useful life
- Bacteriostatic
- Neutral pH 6.5 - 7.5
- Prevents corrosion

**TIP:**  
FOR MILK BATH SOLUTION,  
DILUTE 1 PART SOLUTION  
IN 6 PARTS WATER




## PRODUCT SELECTION

ITEM #	DILUTION LUBRICANT / WATER	PACKAGING
SS-21-6001	1 to 4 ml / L (0.125 to 0.5 oz / Gallon)	1 x 3.8 L (1 Gallon)
SS-21-6004	1 to 4 ml / L (0.125 to 0.5 oz / Gallon)	4 x 3.8 L (4 Gallons)
SS-21-6025	1 to 4 ml / L (0.125 to 0.5 oz / Gallon)	2 x 9.5 L (2 x 2.5 Gallons)
SS-21-6032	No dilution required	1 x 946 ml (32 oz) Spray bottle
SS-21-6012	No dilution required	12 X 946 ml (32 oz) Spray Bottles



# BATRIK S.T.I.C.S.

Our Simple Tracking Identification Color System, or S.T.I.C.S., is an easy to use, latex free, tool that allows you to quickly identify your instruments by department and/or procedure. The tape adheres securely and every tape sheet comes with a protective cover to ensure it lasts as long as you need.

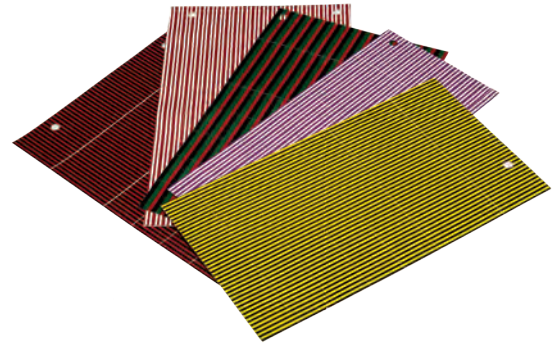

Each 8.5" x 11" (216 mm x 279 mm) sheet contains the following adhesive strips:

- 1" (25.4 mm) x 220
- 1.5" (38.1 mm) x 44
- 2" (50.8 mm) x 44
- Total tape per sheet 374" (9500 mm)

## PRODUCT SELECTION

SOLID COLORS	
ITEM #	COLOR
T-850-2-01	White
T-850-2-02	Royal Blue
T-850-2-03	Kelly Green
T-850-2-04	Red
T-850-2-05	Yellow
T-850-2-06	Purple
T-850-2-07	Orange
T-850-2-08	Brown
T-850-2-09	Black
T-850-2-10	Gray
T-850-2-11	Lavender
T-850-2-12	Beige
T-850-2-13	Light Blue
T-850-2-14	Lime Green
T-850-2-15	Pink
T-850-2-17	Plum
T-850-2-18	Mint Green

BLACK STRIPES	
ITEM #	COLOR
T-850-3-01	Black/White
T-850-3-03	Black/Kelly Green
T-850-3-04	Black/Red
T-850-3-05	Black/Yellow
T-850-3-07	Black/Orange
T-850-3-10	Black/Gray
T-850-3-11	Black/Lavender
T-850-3-12	Black/Beige
T-850-3-13	Black/Light Blue
T-850-3-14	Black/Lime Green
T-850-3-15	Black/Pink
T-850-3-17	Black/Plum

TRI-COLOR	
ITEM #	COLOR
T-850-5-51	Red/White/Blue
T-850-5-54	Red/Black/Green

SPECIALTY	
ITEM #	COLOR
T-850-7-04	Yellow/Red


WHITE STRIPES	
ITEM #	COLOR
T-850-4-02	White/Royal Blue
T-850-4-03	White/Kelly Green
T-850-4-04	White/Red
T-850-4-05	White/Yellow
T-850-4-07	White/Orange
T-850-4-08	White/Brown
T-850-4-11	White/Lavender
T-850-4-13	White/Light Blue
T-850-4-14	White/Lime Green
T-850-4-15	White/Pink
T-850-4-17	White/Plum

FOUR-COLOR	
ITEM #	COLOR
T-850-2-31	Red/Yellow/ Purple/Orange
T-850-2-32	Blue/White/ Black/Green

# BATRIK ID TAPE



### AVAILABLE WIDTHS

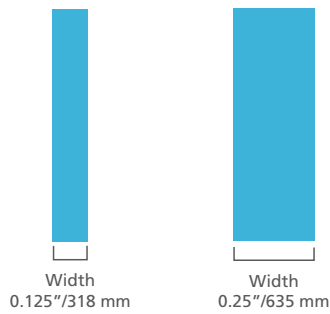

A simple solution for marking surgical instruments with the aid of color-coding, our ID Tape is durable and will not chip, flake, or peel. It can be used in autoclaves, withstanding repeated exposure to sterilization. Independent validation tests demonstrated the efficiency of sterilization of marked instruments, even underneath the Batrik ID Tape.



Quickly identifies instruments by department or procedure



Easy to apply



Withstands repeated exposure to sterilization



Has more per tape roll (300" / 7620 mm)



Available in 9 colors and 2 widths



Conveniently packaged in a clear box



Non-Sterile

## PRODUCT SELECTION

ITEM # 0.125" / 318 MM	ITEM# 0.25" / 635 MM	COLOR
T-125-01	T-250-01	White
T-125-02	T-250-02	Blue
T-125-03	T-250-03	Green
T-125-04	T-250-04	Red
T-125-05	T-250-05	Yellow
T-125-06	T-250-06	Purple
T-125-07	T-250-07	Orange
T-125-08	T-250-08	Brown
T-125-09	T-250-09	Black

# BATRIK TAPE-B-GONE™ TAPE REMOVER

Tape-B-Gone™ is a uniquely formulated and ready-to-use proprietary formula that easily removes adhesives from medical equipment and instruments. Unlike many other brands, Tape-B-Gone™ does not contain highly dangerous or caustic chemicals, is not flammable or combustible, and is safe for use on metal, glass, plastic, and rubber without any special equipment or ventilation.



## PRODUCT SELECTION

ITEM #	DESCRIPTION	PACKAGING
<b>T-08RMV</b>	8 oz (237 ml) bottle	1 Unit
<b>T-32RMV</b>	8 oz (237 ml) bottle	4 / Pack
<b>T-SAMP-RMV</b>	0.5 oz (15 ml) bottle	1 Unit


# BATRIK PADLOCKS SECURITY TAGS

Our Padlock Security Tags for instrument containers and general hospital applications are specifically designed for the healthcare market. Our material selection, fabrication process, and engineering ensure consistent quality and high temperature compliance with the various steam sterilization temperatures in worldwide hospital applications.


Available with optional reactive ink indicator dots that provide visual sterility indication. Light blue dots turn brown/dark brown when steam sterilized.

## AVAILABLE IN THREE FORMATS


### 700-SERIES

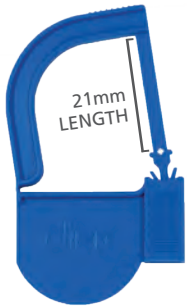

### 800-SERIES

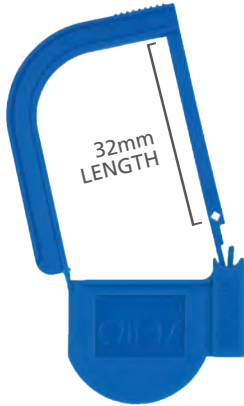

### 900-SERIES

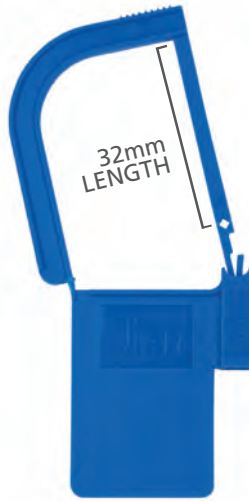

Made of polypropylene



Latex and DEHP Free



Will not melt or soften with any established steam sterilization protocol




Fully open design (45°) facilitates insertion into instrument containers



Easy twist breaking design



Large base for easier handling

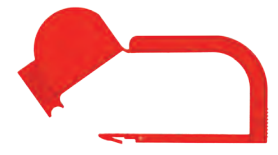

Available in different size



Custom Colors and OEM available upon request

## PRINTED & CUSTOM NUMBERED PADLOCKS

Available in 700, 800 & 900 SERIES

### PRODUCT SELECTION

COLOR	PRINTED WORD
<b>PRINTED PADLOCKS</b>	
WHITE	CLEAN
BLUE	STERILE
RED	REPAIR
YELLOW	SHARPEN
YELLOW	LOANER
BLACK	DISCARD
<b>CUSTOM NUMBERED</b>	
ALL COLORS	up to 8 numbers



**CUSTOM NUMBERED**  
Custom printed padlocks (up to 8 characters)



**PRINTED**  
6 different printed labels available

NOTE: Printed & Custom Numbered Padlocks are available in 700, 800, and 900 Series and can be ordered in 100, 200, and 500 unit quantities. Please contact us for ordering details.

INDICATOR DOTS



REGULAR



# 700-SERIES PADLOCKS

The 700 series padlocks are perfect for "Aesculap" containers and those with a short, large opening.

## PRODUCT SELECTION

### REGULAR PADLOCKS

COLOR	100 UNITS	200 UNITS	500 UNITS
WHITE	15-701	15-711	15-710
BLUE	15-702	15-712	15-720
GREEN	15-703	15-713	15-730
RED	15-704	15-714	15-740
YELLOW	15-705	15-715	15-750
PINK	15-706	15-716	15-760
ORANGE	15-707	15-717	15-770
BLACK	15-709	15-719	15-790
LIGHT BLUE	15-722	15-732	15-742

### WITH INDICATORS DOTS

COLOR	100 UNITS	200 UNITS	500 UNITS
WHITE	15-701-IND	15-711-IND	15-710-IND
BLUE	15-702-IND	15-712-IND	15-720-IND
GREEN	15-703-IND	15-713-IND	15-730-IND
RED	15-704-IND	15-714-IND	15-740-IND
YELLOW	15-705-IND	15-715-IND	15-750-IND
PINK	15-706-IND	15-716-IND	15-760-IND
ORANGE	15-707-IND	15-717-IND	15-770-IND
BLACK	15-709-IND	15-719-IND	15-790-IND
LIGHT BLUE	15-722-IND	15-732-IND	15-742-IND

Double locking hooks can be inserted in a large guide hole



The grooves on the padlock facilitate handling



INDICATOR DOTS

REGULAR



# 800-SERIES PADLOCKS

The 800 series padlocks have a longer arm (32 mm) and have a single, thinner hook than the 700 series. They are perfect for almost any type of container, even those with long and thin channels.

## PRODUCT SELECTION

### REGULAR PADLOCKS

COLOR	100 UNITS	200 UNITS	500 UNITS
WHITE	15-801	15-811	15-810
BLUE	15-802	15-812	15-820
GREEN	15-803	15-813	15-830
RED	15-804	15-814	15-840
YELLOW	15-805	15-815	15-850
PINK	15-806	15-816	15-860
ORANGE	15-807	15-817	15-870
BLACK	15-809	15-819	15-890
LIGHT BLUE	15-822	15-832	15-842

### WITH INDICATORS DOTS

COLOR	100 UNITS	200 UNITS	500 UNITS
WHITE	15-801-IND	15-811-IND	15-810-IND
BLUE	15-802-IND	15-812-IND	15-820-IND
GREEN	15-803-IND	15-813-IND	15-830-IND
RED	15-804-IND	15-814-IND	15-840-IND
YELLOW	15-805-IND	15-815-IND	15-850-IND
PINK	15-806-IND	15-816-IND	15-860-IND
ORANGE	15-807-IND	15-817-IND	15-870-IND
BLACK	15-809-IND	15-819-IND	15-890-IND
LIGHT BLUE	15-822-IND	15-832-IND	15-842-IND

Single locking hook can be inserted in a large guide hole




The locking system allows for easy closure



**INDICATOR  
DOTS**

**REGULAR**

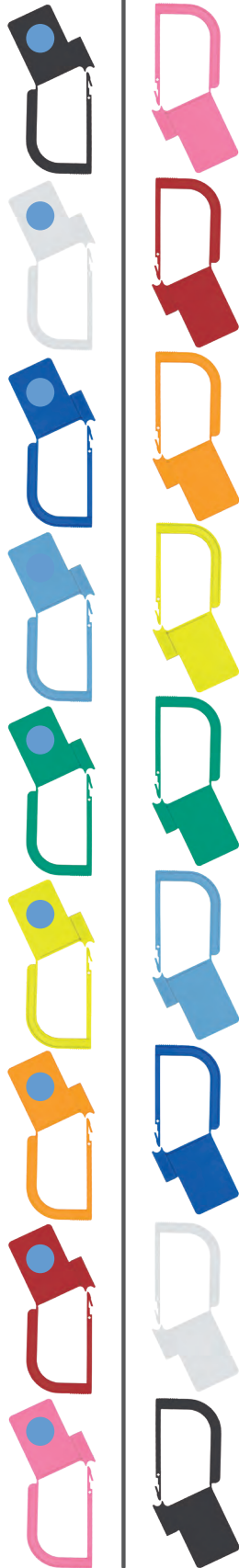

# 900-SERIES PADLOCKS

The 900 series padlocks are identical to the 800 series, but with a larger rectangular base which allows for the placement of a 20 mm x 30 mm tracking label.

## PRODUCT SELECTION

### REGULAR PADLOCKS

COLOR	100 UNITS	200 UNITS	500 UNITS
WHITE	15-901	15-911	15-910
BLUE	15-902	15-912	15-920
GREEN	15-903	15-913	15-930
RED	15-904	15-914	15-940
YELLOW	15-905	15-915	15-950
PINK	15-906	15-916	15-960
ORANGE	15-907	15-917	15-970
BLACK	15-909	15-919	15-990
LIGHT BLUE	15-922	15-932	15-942

### WITH INDICATORS DOTS

COLOR	100 UNITS	200 UNITS	500 UNITS
WHITE	15-901-IND	15-911-IND	15-910-IND
BLUE	15-902-IND	15-912-IND	15-920-IND
GREEN	15-903-IND	15-913-IND	15-930-IND
RED	15-904-IND	15-914-IND	15-940-IND
YELLOW	15-905-IND	15-915-IND	15-950-IND
PINK	15-906-IND	15-916-IND	15-960-IND
ORANGE	15-907-IND	15-917-IND	15-970-IND
BLACK	15-909-IND	15-919-IND	15-990-IND
LIGHT BLUE	15-922-IND	15-932-IND	15-942-IND

**Single locking hook**  
for use with smaller  
diameters



Ideal for adding labels/  
indicators indicators or  
markings to the large  
rectangular base



# SAVEFACE™ DRAPED SHIELDS

**REDUCE PERSONAL INJURY**  
**SUPERIOR INSTRUMENT DECONTAMINATION PROTECTION**

## SAVEFACE DRAPED - FACE SHIELD

In the Sterile Processing Department (SPD)/Central Sterile Supply Department (CSSD), **personal safety** is paramount.

Our Draped Face Shield provides essential **protection** during high-risk instrument decontamination.

## EFFECTIVE PROTECTION IS NON-NEGOTIABLE

SAVEFACE Draped dynamically reduces harmful exposure to caustic chemicals as well as the risk of cross-contamination. Reduce costly injury to hospital staff.

## PPE CATEGORY 2 CERTIFIED

SAVEFACE DRAPED is tested to protect. Wear it with confidence.



PATENT PENDING



INCLUDES CLOSING SYSTEM TO  
 REDUCE RISK OF FLUID CONTACT  
 FROM BELOW THE SHIELD

## SAVEFACE: KEY FEATURES

### ANSI/AAMI PB70 LEVEL 3 COMPLIANT DRAPE:

- Provides a high level of fluid penetration protection.
- Tuck into clothing or pinch together.

### EXTRA-LONG 9"/23CM WRAP-AROUND SHIELD:

- Ensures full coverage, reducing the risk of exposure from different angles.

### COMFORTABLE:

- The soft, self-shaping foam headpiece molds comfortably to your forehead, allowing for extended wear without discomfort.

### OPTICAL CLARITY WITH ANTI-FOG AND ANTI-GLARE



## DON'T MISS THE OPPORTUNITY TO EXPLORE OUR BROADER RANGE OF PROTECTIVE SHIELDS.

Visit our website [www.surgmed.com](http://www.surgmed.com) or consult our detailed literature to find the shield that best fits your needs.

### VERY IMPORTANT

#### IS YOUR PRESENT FACE SHIELD CATEGORY 2 CERTIFIED?

Batrik face shields are independently certified to comply with the PPE Category 2 standard. Many facilities are now requesting compliance to this standard.



SF-240



SF-8240



### PRODUCT SELECTION

ITEM #	DESCRIPTION	DIMENSIONS	TOTAL LENGTH	PACKAGING
SF-240	BATRIK SAVEFACE™ no drape	Shield Length: 9" / 23 cm Shield Width: 13" / 33 cm	9" (23 cm)	240 / Case (8 Boxes of 30 Pcs)
SF-8240	BATRIK SAVEFACE™ with ANSI/AAMI PB70 Level 3 compliant protection drape	Shield Length: 9" / 23 cm Shield Width: 13" / 33 cm Drape Length: 6" / 15 cm	15" (38 cm)	240 / Case (8 Boxes of 30 Pcs)

# BATRIK SCRUB PAD

Our Scrub Pad is ideal for pre-cleaning your surgical instruments without leaving any metal residue that often occurs when using steel wool or metal sponges. The abrasiveness is perfectly balanced to provide you with efficient cleaning without damaging your instruments. Its fine fibers allow the user to polish, clean, and debur all types of material surfaces.

- Abrasive Material: Aluminum Oxide
- Grit: 400
- Grade: Very Fine
- Style: Flexible
- Pad material: Non-Woven
- Color: Maroon
- Type: Hand Pad



## PRODUCT SELECTION

ITEM #	DESCRIPTION	DIMENSIONS	PACKAGING
<b>60-46050</b>	Scrub Pad	4" x 6" / 10 cm x 15 cm	150 / Box (50 bags of 3 units each)
<b>60-46003</b>	Scrub Pad	4" x 6" / 10 cm x 15 cm	3 / Bag

Your First ™

# ENZYMATIC & NON-ENZYMATIC CLEANERS

Our Enzymatic and Non-Enzymatic Cleaners are specially formulated to remove organic material from medical equipment and instrumentation.

Our **Enzymatic Cleaner** formats are formulated with PrAmLi Power3 Enzymes™, our proprietary unique blend of enzymes, and a surfactant to efficiently remove all proteins, blood, fatty tissue, and more. They also remove stains from surgical scrubs, linens, and medical clothing.

Our **Non-Enzymatic Cleaner** formats are formulated with a free-rinsing, neutral pH detergent blend for effective cleaning.



Both formats are powerful solutions that break down all bioburden for effective surface and instrument pre-cleaning\*

## PRODUCT SELECTION

ENZYMATIC CLEANER	NON-ENZYMATIC CLEANER	DILUTION DETERGENT / WATER	PACKAGING
<b>LOW FOAM</b>			
<b>SS-21-1001</b>	<b>SS-721-1001</b>	4 to 8 ml / L (0.5 to 1 oz / Gal)	1 x 3.8 L (1 Gal)
<b>SS-21-1004</b>	<b>SS-721-1004</b>	4 to 8 ml / L (0.5 to 1 oz / Gal)	4 x 3.8 L (4 Gal)
<b>MILD FOAM</b>			
<b>SS-21-2001</b>	<b>SS-721-2001</b>	4 to 8 ml / L (0.5 to 1 oz / Gal)	1 x 3.8 L (1 Gal)
<b>SS-21-2004</b>	<b>SS-721-2004</b>	4 to 8 ml / L (0.5 to 1 oz / Gal)	4 x 3.8 L (4 Gal)
<b>SUPER CONCENTRATE</b>			
<b>SS-21-3001</b>	<b>SS-721-3001</b>	1 to 2 ml / L (0.125 to 0.25 oz / Gal)	1 x 3.8 L (1 Gal)
<b>SS-21-3004</b>	<b>SS-721-3004</b>	1 to 2 ml / L (0.125 to 0.25 oz / Gal)	4 x 3.8 L (4 Gal)

### GOOD TO KNOW:

OUR LOW FOAM & SUPER CONCENTRATE DETERGENTS (ENZYMATIC & NON-ENZYMATIC) MAY BE USED IN WASHERS / ULTRASONIC CLEANERS.

WE RECOMMEND CONFIRMING SUGGESTED DILUTION RATIO AND WATER TEMPERATURE WITH THE MANUFACTURER OF YOUR WASHER / ULTRASONIC CLEANER AS MACHINES / SETTINGS VARY.


NOTE: Our Low Foam & Super Concentrate Detergent formulas are available in a 2.5 Gallon format based on a minimum order quantity. Contact us for more information.

\*Caution: We recommend consulting your hospital pre-cleaning/cleaning protocols for various equipment/instrument cleaning detergent preferences.

Your First ™

# ENVIRO-SACHET™

## ENZYMATIC GEL & NON-ENZYMATIC CLEANER

Your First ™ Enviro-Sachet™, our unique cleaning solutions now available in Enzymatic & Non-Enzymatic formats.



- Our Enzymatic Gel format is formulated with the same triple action **PrAmLi Power<sup>3</sup> Enzymes™** that you know and trust.
- Our Non-Enzymatic Cleaner format is formulated with a free-rinsing, neutral pH detergent blend.

Both formats are powerful solutions that break down all bioburden for effective surface and instrument pre-cleaning.

### PRODUCT SELECTION

ENZYMATIC GEL	NON-ENZYMATIC CLEANER	DESCRIPTION	PACKAGING
<b>17-166</b>	<b>717-166</b>	1 Dry Tubular Sponge & Enviro-Sachet™ set	25 / Box 4 Boxes / Case
<b>17-266</b>	<b>717-266</b>	1 Dry Flat Sponge & Enviro-Sachet™ set	50 / Box 4 Boxes / Case
<b>17-177</b>	<b>717-177</b>	1 Dry Tubular Sponge, Enviro-Sachet™ & Stainless Steel Bowl kit	1 Kit
<b>17-277</b>	<b>717-277</b>	1 Dry Flat Sponge, Enviro-Sachet™ & Stainless Steel Bowl kit	1 Kit
<b>17-EZSAMP</b>	<b>717-EZSAMP</b>	SAMPLE KIT INCLUDES: 3 Dry Flat Sponges with Enviro-Sachet™ sets 1 Dry Tubular Sponge with Enviro-Sachet™ set 1 Reusable Stainless Steel Bowl	1 Kit
<b>H-EZ500</b>		Reusable Stainless Steel Bowl 4.3" x 5.3" x 3" / 11 cm x 13.5 cm x 7.5 cm	1 Unit

Your First ™

# FOAMING SPRAY

## READY-TO-USE ENZYME FOAMING SPRAY

Our Foaming Spray is ideal for instrument pre-cleaning as the thick foam spray enzymatically breaks down all bioburden to allow for effective cleaning and sterilization.

Formulated with **PrAmLi Power<sup>3</sup> Enzymes™**, a unique blend of three powerful enzymes along with a surfactant to provide the ultimate in cleaning efficacy



### INSTRUCTIONS

1. Shake well
2. Spray immediately after procedure
3. Keep wet until the decontamination process is started

### PRODUCT SELECTION

ITEM #	PACKAGING
<b>SS-21-4001</b>	1 x 946 ml (32 oz) Spray Bottle
<b>SS-21-4012</b>	12 x 946 ml (32 oz) Spray Bottles

Important: In addition to our enzymatic formulation, our spray bottle contains as much as 25% more solution when compared to the competition.

Your First ™

# FIRST WIPES

Our cleanroom wipes are designed to do the first wipe-down of an endoscope after a procedure. They are made from hydroentangled non-woven material composed of 45% polyester and 55% cellulose. These wipes are manufactured and packaged in a cleanroom environment.



- Can be used for a cleanroom environment of ISO Class 5 through 8
- Production process guarantees extremely low particle generation and chemical extractables
- Compatible with virtually any solvent or solution for cleaning and wiping (such as **PrAmLi Power<sup>3</sup> Enzymes™**)
- Double-bag packaging
- Available in two versatile sizes

## PRODUCT SELECTION

ITEM #	DESCRIPTION	DIMENSIONS	PACKAGING
<b>17-A0909</b>	Non-woven Cleanroom Wipes	9" x 9" (23 cm x 23 cm)	300 / Bag 14 Bags / Case
<b>17-A1212</b>	Non-woven Cleanroom Wipes	12" x 12" (31 cm x 31 cm)	150 / Bag 10 Bags / Case

Your First ™

# DRY SPONGES

Our sponges are specifically fabricated for use with all medical instrumentation. Unlike most competitor sponges, our sponges are white to better evidence the effectiveness of the cleaning process. Dark colored sponges easily hide contaminants, but with our white sponges you can visually evaluate the pre-clean efficacy.



OUR WHITE SPONGES MAY APPEAR SLIGHTLY YELLOWED. THIS IS ONLY DUE TO UV EXPOSURE, IS NOT A REFLECTION OF ITS AGE.

## PRODUCT SELECTION

ITEM #	DESCRIPTION	DIMENSIONS	PACKAGING
<b>17-100</b>	DRY TUBULAR SPONGE <a href="#">White Dispenser Box</a>	2.0" ø x 4.25" (2.7 cm ø x 10.8 cm)	25 / Box (4 Boxes / Case)
<b>17-101</b>	DRY TUBULAR SPONGE <a href="#">Bulk</a>	2.0" ø x 4.25" (2.7 cm ø x 10.8 cm)	560 / Case
<b>17-200</b>	DRY FLAT SPONGE <a href="#">White Dispenser Box</a>	4" x 5" x 0.5" (10.2 cm x 12.7 cm x 1.3 cm)	50 / Box (4 Boxes / Case)
<b>17-201</b>	DRY FLAT SPONGE <a href="#">Bulk</a>	4" x 5" x 0.5" (10.2 cm x 12.7 cm x 1.3 cm)	1,000 / Case
ITEM #	DESCRIPTION	DIMENSIONS	PACKAGING
<b>17-600</b>	DRY TUBULAR SPONGE <a href="#">Purple Dispenser Box</a>	2.0" ø x 4.25" (2.7 cm ø x 10.8 cm)	25 / Box (4 Boxes / Case)
<b>17-620</b>	DRY FLAT SPONGE <a href="#">Purple Dispenser Box</a>	4" x 5" x 0.5" (10.2 cm x 12.7 cm x 1.3 cm)	50 / Box (4 Boxes / Case)

Your First ™

**ENZYMATIC BEDSIDE SPONGE KITS**

# RightSwipe™

Introducing RightSwipe™, our newly launched Enzymatic Bedside Sponge Kits now feature sponges treated with an enhanced formula for a more effective pre-clean or surface clean solution. Our improved three-enzyme blend ensures thorough and safe cleaning and a simplified application process making RightSwipe™ the perfect match for your endoscope care needs. The tubular sponge design is compatible with diverse endoscope models, providing versatility alongside the superior benefits of our enzymatic solution.



FLAT SPONGE



TUBULAR SPONGE



FREESTANDING POUCH



ENVIRO-SACHET™  
ENZYMATIC GEL KIT



DISPENSER BOX

## PRODUCT SELECTION

ITEM #	DESCRIPTION	PACKAGING
17-111	PRE-TREATED TUBULAR SPONGE Add 100-500 ml of water to activate detergent solution. Comes in plastic bag. <sup>1</sup>	25 / Box (4 Boxes / Case)
17-211	PRE-TREATED FLAT SPONGE Add 100-500 ml of water to activate detergent solution. Comes in plastic bag. <sup>1</sup>	50 / Box (4 Boxes / Case)
17-122	PRE-TREATED TUBULAR SPONGE Comes pre-treated in 100 ml enzymatic solution in freestanding pouch. Simply add 500 ml of water to the pouch.	15 / Box (4 Boxes / Case)
17-222	FLAT SPONGE Comes pre-treated in 100 ml enzymatic solution in freestanding pouch. Simply add 500 ml of water to the pouch.	25 / Box (4 Boxes / Case)
17-133	READY TO USE TUBULAR SPONGE Comes pre-treated in 500 ml of an enzymatic solution in freestanding pouch.	24 / Case
17-233	READY TO USE FLAT SPONGE Comes pre-treated in 500 ml of an enzymatic solution in freestanding pouch.	24 / Case
17-144	READY TO USE DUAL PURPOSE TUBULAR SPONGE Comes pre-treated in 100 ml of an enzymatic solution in freestanding pouch.	12 / Box (4 Boxes / Case)
17-244	READY TO USE DUAL PURPOSE FLAT SPONGE Comes pre-treated in 100 ml of an enzymatic solution in freestanding pouch.	12 / Box (4 Boxes / Case)
17-155	READY TO USE TUBULAR SPONGE Comes pre-treated in 220 ml of an enzymatic solution in freestanding pouch.	24 / Case
17-255	READY TO USE FLAT SPONGE Comes pre-treated in 220 ml of an enzymatic solution in freestanding pouch.	24 / Case
17-277	ENZYMATIC GEL Comes with 1 Dry Flat Sponge, Enviro-Sachet™ & Stainless Steel Bowl kit	1 Kit


Our stainless steel bowl can be sold separately and used as a basin for easy and convenient sponge activation<sup>1</sup>



Your First ™

## NON-ENZYMATIC DETERGENT BEDSIDE SPONGE KITS

# RightSwipe™

Presenting our enhanced Non-Enzymatic Bedside Sponge Kits. Each kit includes sponges pre-treated with a free-rinsing, neutral pH detergent blend from RightSwipe™, ensuring an effective pre-cleaning and surface cleaning solution for medical instruments. This kit stands as Your First ™ step in the medical instrument cleaning process, showcasing the advanced benefits for superior instrument care. Our enhanced bedside kits will provide you with a streamlined workflow.



FREESTANDING POUCH



FLAT SPONGE



TUBULAR SPONGE



ENVIRO-SACHET™  
NON-ENZYMATIC CLEANER



DISPENSER BOX

## PRODUCT SELECTION

ITEM #	DESCRIPTION	PACKAGING
717-122	PRE-TREATED TUBULAR SPONGE Add 100-500 ml of water to activate non-enzymatic detergent. Comes in plastic bag.1	15 / Box (4 Boxes / Case)
717-222	PRE-TREATED FLAT SPONGE Add 100-500 ml of water to activate non-enzymatic detergent. Comes in plastic bag.1	25 / Box (4 Boxes / Case)
717-133	READY TO USE TUBULAR SPONGE Comes pre-treated in 500 ml of non-enzymatic solution in freestanding pouch.	24 / Case
717-233	READY TO USE FLAT SPONGE Comes pre-treated in 500 ml of non-enzymatic solution in freestanding pouch.	24 / Case
717-144	READY TO USE DUAL PURPOSE TUBULAR SPONGE Comes pre-treated in 100 ml of non-enzymatic solution in freestanding pouch.	12 / Box (4 Boxes / Case)
717-244	READY TO USE DUAL PURPOSE FLAT SPONGE Comes pre-treated in 100 ml of non-enzymatic solution in freestanding pouch.	12 / Box (4 Boxes / Case)
717-155	READY TO USE TUBULAR SPONGE Comes pre-treated in 220 ml of non-enzymatic solution in freestanding pouch.	24 / Case
717-255	READY TO USE FLAT SPONGE Comes pre-treated in 220 ml of non-enzymatic solution in freestanding pouch.	24 / Case
717-277	NON-ENZYMATIC CLEANER Comes with 1 Dry Flat Sponge, Enviro-Sachet™ & Stainless Steel Bowl kit.	1 Kit

Our stainless steel bowl can be sold separately and used as a basin for easy and convenient sponge activation<sup>1</sup>

Your First ™

# TRANSPORT BAGS FOR ENDOSCOPE AND MEDICAL INSTRUMENTS

Thoughtfully designed to eliminate touch-points during endoscope transport from point-of-care to the reprocessing area, and back, focused on safety of patients and medical staff.

- Made of recyclable LDPE material known to be quite flexible, highly resilient, and tough
- Latex and DEHP-free
- Excellent resistance to dilute and concentrated acids, alcohols, bases and esters\*
- Good resistance to aldehydes, ketones, and vegetable-based oils\*
- Can withstand high temperatures of up to 90 °C (194°F)
- Tamper-evident seal for added safety, also aids in keeping the instrument moist to avoid drying of bio-debris
- Clearly marked bags with color-coded visuals facilitates immediate identification of clean vs. used instruments
- Users can write important information directly on the bag in spaces marked with visual icons
- Available in two convenient sizes

GOOD TO KNOW: USING **Your First ™** SEALABLE TRANSPORT BAGS REDUCES USER ERROR, PREVENTS CROSS-CONTAMINATION, AND PROMOTES PATIENT AND MEDICAL STAFF SAFETY.

## PRODUCT SELECTION


ITEM #	DESCRIPTION	DIMENSIONS	PACKAGING
<b>17-HAP24</b>	Clean Instruments Bag	24" x 22" / 61.0 cm x 56.0 cm	200 / Case
<b>17-SAD24</b>	Soiled Instruments Bag	24" x 22" / 61.0 cm x 56.0 cm	200 / Case
<b>17-HAP30</b>	Clean Instruments Bag	30.5" x 26" / 77.0 cm x 66.0 cm	200 / Case
<b>17-SAD30</b>	Soiled Instruments Bag	30.5" x 26" / 77.0 cm x 66.0 cm	200 / Case

\*Reference on file

### CLEAN INSTRUMENTS BAG



### SOILED INSTRUMENTS BAG

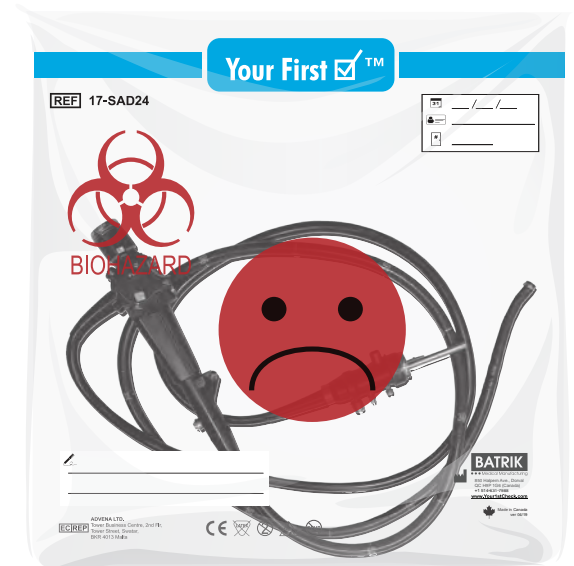

MADE IN CANADA

# BATRIK ARC™ WIRE STYLE BRUSHES

## SINGLE-USE AND REUSABLE

You've got options with our ARC™ Wire Brushes. Our single-use and reusable brushes are available in both flexible and rigid varieties and come in a multitude of sizes. ARC™ Wire Brushes are guaranteed to last longer than twisted wire and preserve the integrity of your valuable instruments.

### BRUSH SELECTION TIPS:

#### INSTRUMENT SPECIFICATIONS

##### What is the maximum length of the instrument(s) to be cleaned?

Add at least 2" (5cm) to that length, so you have an adequate length to work with.

##### What is the inner diameter of the instrument?

Add at most 0.125" (3mm) to determine the diameter of the bristle to ensure effective cleaning. The bristle diameter has to slightly exceed the inner tube diameter. Bristles that are too large will bend too much and not clean effectively.

#### BRUSH SELECTION

Simply find the section in our catalogue that best corresponds to the overall length and diameter of your instrument. As a rule, if your present brush is 8" (20.32cm) in length and we only offer 10" (25.4cm) length with your required diameter, that is entirely acceptable.

In many cases, our ARC™ Wire Brushes will have more bristles on the brush. As an example, many 5mm (0.19") diameter brushes are only 1" (2.5cm) in length. At BATRIK, we strived to provide the greatest amount of bristles to ensure the most effective cleaning.

#### CUSTOMIZABLE WIRE BRUSHES


A unique aspect of single-use ARC™ Wire Brushes is that they can be customized to any length:

- Up to 2000mm for Wire Brushes on plastic semi-rigid sheaths, rigid brushes from 5 to 20mm in diameter.
- No limit in length for Wire Brushes on flexible plastic sheaths, flexible brushes from 2-5mm in diameter

## RIGID



## FLEXIBLE

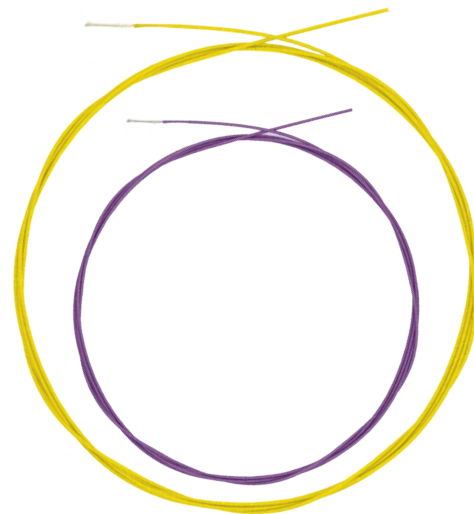

**12" / 30 CM SINGLE-END WIRE BRUSH - DELRIN PLASTIC**

	ITEM # SINGLE-USE	ITEM # REUSABLE	COLOR	TOTAL LENGTH	BRUSH (Ø)	BRUSH LENGTH	SHEATH (Ø)
<b>FLEXIBLE</b> BRUSH & SHEATH	<b>41-31502-S</b>	<b>41-31502</b>	White	12.401" / 31.5 cm	0.078" / 2 mm	0.59" / 15 mm	0.39" / 1 mm
	<b>41-32503-S</b>	<b>41-32503</b>	Yellow	12.795" / 32.5 cm	0.118" / 3 mm	0.98" / 25 mm	0.39" / 1 mm
	<b>411-32504-S</b>	<b>411-32504</b>	Violet	12.795" / 32.5 cm	0.157" / 4 mm	0.98" / 25 mm	0.61" / 1.55 mm
	<b>411-32505-S</b>	<b>411-32505</b>	Grey	12.795" / 32.5 cm	0.19" / 5 mm	0.98" / 25 mm	0.68" / 1.75 mm
<b>RIGID</b> BRUSH & SHEATH	<b>41-35005-S</b>	<b>41-35005</b>	White	13.779" / 35 cm	0.19" / 5 mm	1.968" / 50 mm	0.118" / 3 mm
	<b>41-35006-S</b>	<b>41-35006</b>	Green	13.779" / 35 cm	0.23" / 6 mm	1.968" / 50 mm	0.118" / 3 mm
	<b>41-35008-S</b>	<b>41-35008</b>	Blue	13.779" / 35 cm	0.31" / 8 mm	1.968" / 50 mm	0.157" / 4 mm
	<b>41-35010-S</b>	<b>41-35010</b>	Red	13.779" / 35 cm	0.39" / 10 mm	1.968" / 50 mm	0.157" / 4 mm
	<b>41-35012-S</b>	<b>41-35012</b>	Yellow	13.779" / 35 cm	0.49" / 12.5 mm	1.968" / 50 mm	0.157" / 4 mm
	<b>41-35015-S</b>	<b>41-35015</b>	Green	13.779" / 35 cm	0.59" / 15 mm	1.968" / 50 mm	0.228" / 5.8 mm
	<b>41-35020-S</b>	<b>41-35020</b>	Blue	13.779" / 35 cm	0.78" / 20 mm	1.968" / 50 mm	0.228" / 5.8 mm

PACKAGING Single-Use: Bag of 20 / Reusable: Bag of 3

**18" / 45 CM SINGLE-END WIRE BRUSH - DELRIN PLASTIC**

	ITEM # SINGLE-USE	ITEM # REUSABLE	COLOR	TOTAL LENGTH	BRUSH (Ø)	BRUSH LENGTH	SHEATH (Ø)
<b>FLEXIBLE</b> BRUSH & SHEATH	<b>41-46502-S</b>	<b>41-46502</b>	White	18.307" / 46.5 cm	0.078" / 2 mm	0.59" / 15 mm	0.39" / 1 mm
	<b>41-47503-S</b>	<b>41-47503</b>	Yellow	18.7" / 47.5 cm	0.118" / 3 mm	0.98" / 25 mm	0.39" / 1 mm
	<b>41-47505-S</b>	<b>N/A</b>	White	18.7" / 47.5 cm	0.19" / 5 mm	0.98" / 25 mm	0.68" / 1.75 mm
	<b>411-47504-S</b>	<b>411-47504</b>	Violet	18.7" / 47.5 cm	0.157" / 4 mm	0.98" / 25 mm	0.61" / 1.55 mm
	<b>411-47505-S</b>	<b>411-47505</b>	Grey	18.7" / 47.5 cm	0.19" / 5 mm	0.98" / 25 mm	0.68" / 1.75 mm
<b>RIGID</b> BRUSH & SHEATH	<b>41-50005-S</b>	<b>41-50005</b>	White	19.685" / 50 cm	0.19" / 5 mm	1.968" / 50 mm	0.118" / 3 mm
	<b>41-50006-S</b>	<b>41-50006</b>	Green	19.685" / 50 cm	0.23" / 6 mm	1.968" / 50 mm	0.118" / 3 mm
	<b>41-50008-S</b>	<b>41-50008</b>	Blue	19.685" / 50 cm	0.31" / 8 mm	1.968" / 50 mm	0.157" / 4 mm
	<b>41-50010-S</b>	<b>41-50010</b>	Red	19.685" / 50 cm	0.39" / 10 mm	1.968" / 50 mm	0.157" / 4 mm
	<b>41-50012-S</b>	<b>41-50012</b>	Yellow	19.685" / 50 cm	0.49" / 12.5 mm	1.968" / 50 mm	0.157" / 4 mm
	<b>41-50015-S</b>	<b>41-50015</b>	Green	19.685" / 50 cm	0.59" / 15 mm	1.968" / 50 mm	0.228" / 5.8 mm
	<b>41-50020-S</b>	<b>41-50020</b>	Blue	19.685" / 50 cm	0.78" / 20 mm	1.968" / 50 mm	0.228" / 5.8 mm

PACKAGING: Single -Use: Bag of 20 / Reusable: Bag of 3

**24" / 60 CM SINGLE-END WIRE BRUSH - DELRIN PLASTIC**

	ITEM # SINGLE-USE	ITEM # REUSABLE	COLOR	TOTAL LENGTH	BRUSH (Ø)	BRUSH LENGTH	SHEATH (Ø)
<b>FLEXIBLE</b> BRUSH & SHEATH	<b>41-61502-S</b>	<b>41-61502</b>	White	24.212" / 61.5 cm	0.078" / 2 mm	0.59" / 15 mm	0.39" / 1mm
	<b>41-62503-S</b>	<b>41-62503</b>	Yellow	24.6" / 62.5 cm	0.118" / 3 mm	0.98" / 25 mm	0.39" / 1mm
	<b>41-62505-S</b>	<b>N/A</b>	White	24.6" / 62.5 cm	0.19" / 5 mm	0.98" / 25 mm	0.68" / 1.75mm
	<b>411-62504-S</b>	<b>411-62504</b>	Violet	24.6" / 62.5 cm	0.157" / 4 mm	0.98" / 25 mm	0.61" / 1.55mm
	<b>411-62505-S</b>	<b>411-62505</b>	Grey	24.6" / 62.5 cm	0.19" / 5 mm	0.98" / 25 mm	0.68" / 1.75mm
<b>RIGID</b> BRUSH & SHEATH	<b>41-65005-S</b>	<b>41-65005</b>	White	25.59" / 65 cm	0.19" / 5 mm	1.968" / 50 mm	0.118" / 3mm
	<b>41-65006-S</b>	<b>41-65006</b>	Green	25.59" / 65 cm	0.23" / 6 mm	1.968" / 50 mm	0.118" / 3mm
	<b>41-65008-S</b>	<b>41-65008</b>	Blue	25.59" / 65 cm	0.31" / 8 mm	1.968" / 50 mm	0.157" / 4 mm
	<b>41-65010-S</b>	<b>41-65010</b>	Red	25.59" / 65 cm	0.39" / 10 mm	1.968" / 50 mm	0.157" / 4 mm
	<b>41-65012-S</b>	<b>41-65012</b>	Yellow	25.59" / 65 cm	0.49" / 12.5 mm	1.968" / 50 mm	0.157" / 4 mm
	<b>41-65015-S</b>	<b>41-65015</b>	Green	25.59" / 65 cm	0.59" / 15 mm	1.968" / 50 mm	0.228" / 5.8 mm
	<b>41-65020-S</b>	<b>41-65020</b>	Blue	25.59" / 65 cm	0.78" / 20 mm	1.968" / 50 mm	0.228" / 5.8 mm



PACKAGING: Single-Use: Bag of 20 / Reusable: Bag of 3

**30" / 75 CM SINGLE-END WIRE BRUSH - DELRIN PLASTIC**

	ITEM # SINGLE-USE	ITEM # REUSABLE	COLOR	TOTAL LENGTH	BRUSH (Ø)	BRUSH LENGTH	SHEATH (Ø)
<b>FLEXIBLE</b> BRUSH & SHEATH	<b>41-76502-S</b>	<b>41-76502</b>	White	30.118" / 76.5 cm	0.078" / 2 mm	0.59" / 15 mm	0.39" / 1 mm
	<b>41-77503-S</b>	<b>41-77503</b>	Yellow	30.511" / 77.5 cm	0.118" / 3 mm	0.98" / 25 mm	0.39" / 1 mm
	<b>411-77504-S</b>	<b>411-77504</b>	Violet	30.511" / 77.5 cm	0.157" / 4 mm	0.98" / 25 mm	0.61" / 1.55 mm
	<b>411-77505-S</b>	<b>411-77505</b>	Grey	30.511" / 77.5 cm	0.19" / 5 mm	0.98" / 25 mm	0.68" / 1.75 mm
	<b>41-80005-S</b>	<b>41-80005</b>	White	31.496" / 80 cm	0.19" / 5 mm	1.968" / 50 mm	0.118" / 3 mm
<b>RIGID</b> BRUSH & SHEATH	<b>41-80006-S</b>	<b>41-80006</b>	Green	31.496" / 80 cm	0.23" / 6 mm	1.968" / 50 mm	0.118" / 3 mm
	<b>41-80008-S</b>	<b>41-80008</b>	Blue	31.496" / 80 cm	0.31" / 8 mm	1.968" / 50 mm	0.157" / 4 mm
	<b>41-80010-S</b>	<b>41-80010</b>	Red	31.496" / 80 cm	0.39" / 10 mm	1.968" / 50 mm	0.157" / 4 mm
	<b>41-80012-S</b>	<b>41-80012</b>	Yellow	31.496" / 80 cm	0.49" / 12.5 mm	1.968" / 50 mm	0.157" / 4 mm
	<b>41-80015-S</b>	<b>41-80015</b>	Green	31.496" / 80 cm	0.59" / 15 mm	1.968" / 50 mm	0.228" / 5.8 mm
	<b>41-80020-S</b>	<b>41-80020</b>	Blue	31.496" / 80 cm	0.78" / 20 mm	1.968" / 50 mm	0.228" / 5.8 mm


PACKAGING: Single-Use: Bag of 20 / Reusable: Bag of 3

BRUSHES

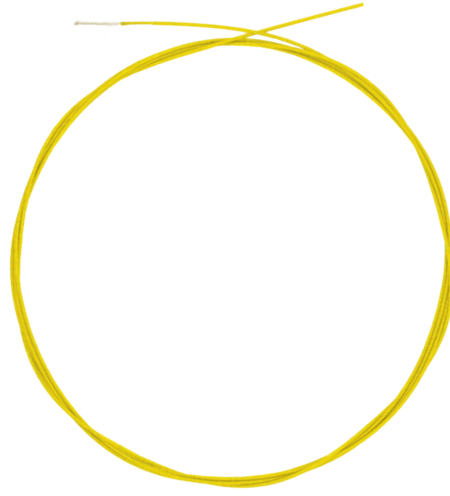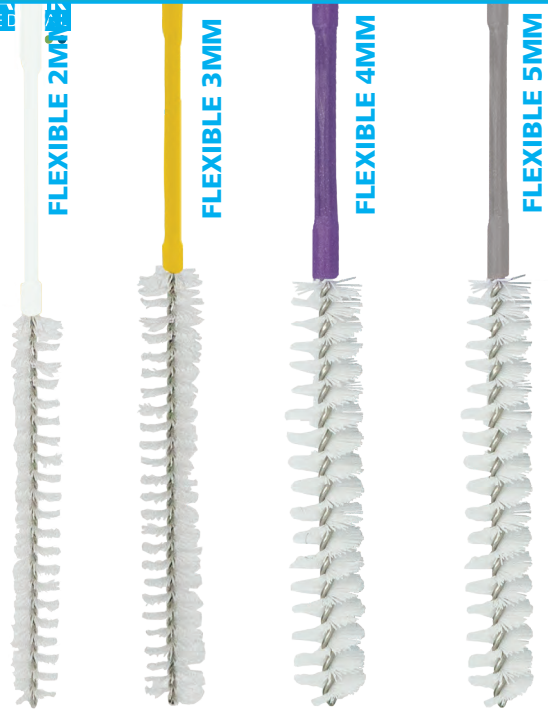

FLEXIBLE

Flexible brush diameter (Ø) ranging from 2 mm to 5 mm

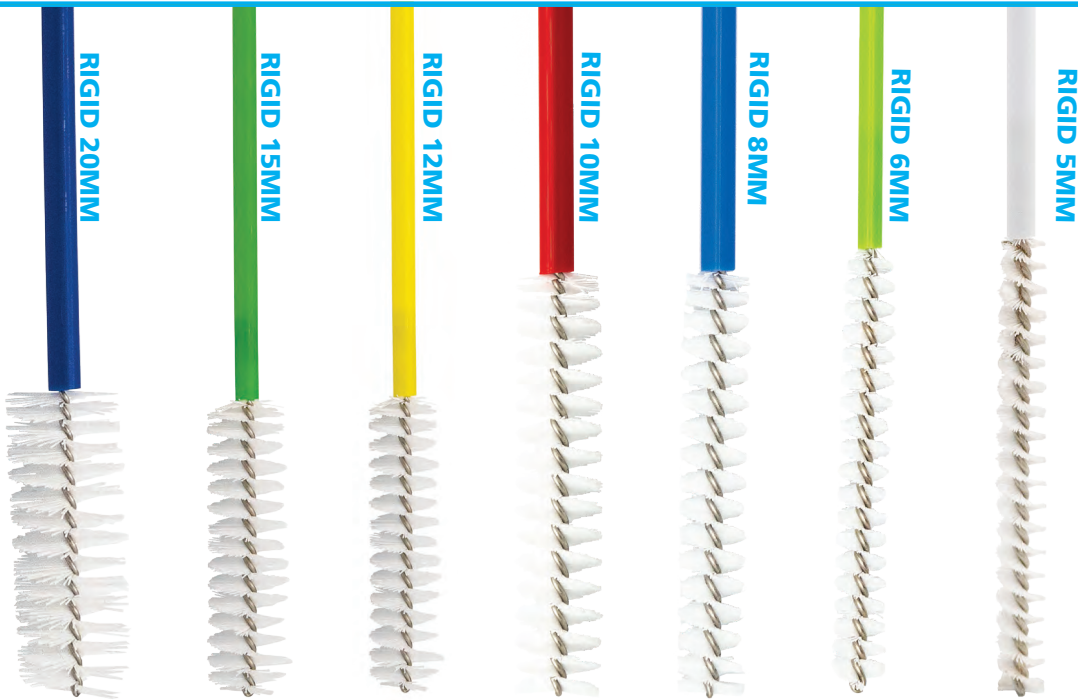

RIGID

Rigid brush diameter (Ø) ranging from 5 mm to 20 mm

BRUSH SELECTION SYSTEM (Actual Size)

2mm 0.078"	3 mm 0.118"	4 mm 0.157"	5 mm 0.196"	6 mm 0.236"	8 mm 0.314"	10 mm 0.394"	12 mm 0.472"	15 mm 0.591"	20 mm 0.787"

# BATRIK REUSABLE DOUBLE ARC™ WIRE BRUSHES

## ASYMMETRIC

### 14" / 35 CM ASYMMETRIC DOUBLE-END WIRE BRUSHES

ITEM #	COLOR	TOTAL LENGTH	BRUSH 1 (Ø)	BRUSH 2 (Ø)	BRUSH LENGTH	SHEATH (Ø)
40-35056	White	13.779" / 35 cm	0.19" / 5 mm	0.23" / 6 mm	1.968" / 50 mm	0.118" / 3 mm
40-35112	Yellow	13.779" / 35 cm	0.39" / 10 mm	0.49" / 12.5 mm	1.968" / 50 mm	0.157" / 4 mm
40-35152	Green	13.779" / 35 cm	0.59" / 15 mm	0.78" / 20 mm	1.968" / 50 mm	0.228" / 5.8 mm
40-35810	Blue	13.779" / 35 cm	0.31" / 8 mm	0.39" / 10 mm	1.968" / 50 mm	0.157" / 4 mm

PACKAGING: Bag of 3

### 20" / 50 CM ASYMMETRIC DOUBLE-END WIRE BRUSHES

ITEM #	COLOR	TOTAL LENGTH	BRUSH 1 (Ø)	BRUSH 2 (Ø)	BRUSH LENGTH	SHEATH (Ø)
40-50056	White	19.685" / 50 cm	0.19" / 5 mm	0.23" / 6 mm	1.968" / 50 mm	0.118" / 3 mm
40-50112	Yellow	19.685" / 50 cm	0.39" / 10 mm	0.49" / 12.5 mm	1.968" / 50 mm	0.157" / 4 mm
40-50152	Green	19.685" / 50 cm	0.59" / 15 mm	0.78" / 20 mm	1.968" / 50 mm	0.228" / 5.8 mm
40-50810	Blue	19.685" / 50 cm	0.31" / 8 mm	0.39" / 10 mm	1.968" / 50 mm	0.157" / 4 mm

PACKAGING: Bag of 3

### 26" / 65 CM ASYMMETRIC DOUBLE-END WIRE BRUSHES

ITEM #	COLOR	TOTAL LENGTH	BRUSH 1 (Ø)	BRUSH 2 (Ø)	BRUSH LENGTH	SHEATH (Ø)
40-65056	White	25.59" / 65 cm	0.19" / 5 mm	0.23" / 6 mm	1.968" / 50 mm	0.118" / 3 mm
40-65112	Yellow	25.59" / 65 cm	0.39" / 10 mm	0.49" / 12.5 mm	1.968" / 50 mm	0.157" / 4 mm
40-65152	Green	25.59" / 65 cm	0.59" / 15 mm	0.78" / 20 mm	1.968" / 50 mm	0.228" / 5.8 mm
40-65810	Blue	25.59" / 65 cm	0.31" / 8 mm	0.39" / 10 mm	1.968" / 50 mm	0.157" / 4 mm

PACKAGING: Bag of 3

### 31" / 80 CM ASYMMETRIC DOUBLE-END WIRE BRUSHES

ITEM #	COLOR	TOTAL LENGTH	BRUSH 1 (Ø)	BRUSH 2 (Ø)	BRUSH LENGTH	SHEATH (Ø)
40-80056	White	31.496" / 80 cm	0.19" / 5 mm	0.23" / 6 mm	1.968" / 50 mm	0.118" / 3 mm
40-80112	Yellow	31.496" / 80 cm	0.39" / 10 mm	0.49" / 12.5 mm	1.968" / 50 mm	0.157" / 4 mm
40-80152	Green	31.496" / 80 cm	0.59" / 15 mm	0.78" / 20 mm	1.968" / 50 mm	0.228" / 5.8 mm
40-80810	Blue	31.496" / 80 cm	0.31" / 8 mm	0.39" / 10 mm	1.968" / 50 mm	0.157" / 4 mm

PACKAGING: Bag of 3

# BATRIK REUSABLE DOUBLE ARC™ WIRE BRUSHES SYMMETRIC

## 14" / 35 CM SYMMETRIC DOUBLE-END WIRE BRUSHES

ITEM #	COLOR	TOTAL LENGTH	BRUSH 1 (Ø)	BRUSH LENGTH	SHEATH (Ø)
<b>40-35005</b>	White	13.779" / 35 cm	0.19" / 5 mm	1.968" / 50 mm	0.118" / 3 mm
<b>40-35006</b>	Green	13.779" / 35 cm	0.23" / 6 mm	1.968" / 50 mm	0.118" / 3 mm
<b>40-35008</b>	Blue	13.779" / 35 cm	0.31" / 8 mm	1.968" / 50 mm	0.157" / 4 mm
<b>40-35010</b>	Red	13.779" / 35 cm	0.39" / 10 mm	1.968" / 50 mm	0.157" / 4 mm
<b>40-35012</b>	Yellow	13.779" / 35 cm	0.49" / 12.5 mm	1.968" / 50 mm	0.157" / 4 mm
<b>40-35015</b>	Green	13.779" / 35 cm	0.59" / 15 mm	1.968" / 50 mm	0.228" / 5.8 mm
<b>40-35020</b>	Blue	13.779" / 35 cm	0.78" / 20 mm	1.968" / 50 mm	0.228" / 5.8 mm

PACKAGING: Bag of 3

## 20" / 50 CM SYMMETRIC DOUBLE-END WIRE BRUSHES

ITEM #	COLOR	TOTAL LENGTH	BRUSH 1 (Ø)	BRUSH LENGTH	SHEATH (Ø)
<b>40-50005</b>	White	19.685" / 50 cm	0.19" / 5 mm	1.968" / 50 mm	0.118" / 3 mm
<b>40-50006</b>	Green	19.685" / 50 cm	0.23" / 6 mm	1.968" / 50 mm	0.118" / 3 mm
<b>40-50008</b>	Blue	19.685" / 50 cm	0.31" / 8 mm	1.968" / 50 mm	0.157" / 4 mm
<b>40-50010</b>	Red	19.685" / 50 cm	0.39" / 10 mm	1.968" / 50 mm	0.157" / 4 mm
<b>40-50012</b>	Yellow	19.685" / 50 cm	0.49" / 12.5 mm	1.968" / 50 mm	0.157" / 4 mm
<b>40-50015</b>	Green	19.685" / 50 cm	0.59" / 15 mm	1.968" / 50 mm	0.228" / 5.8 mm
<b>40-50020</b>	Blue	19.685" / 50 cm	0.78" / 20 mm	1.968" / 50 mm	0.228" / 5.8 mm

PACKAGING: Bag of 3

## 26" / 65 CM SYMMETRIC DOUBLE-END WIRE BRUSHES

ITEM #	COLOR	TOTAL LENGTH	BRUSH 1 (Ø)	BRUSH LENGTH	SHEATH (Ø)
<b>40-65005</b>	White	25.59" / 65 cm	0.19" / 5 mm	1.968" / 50 mm	0.118" / 3 mm
<b>40-65006</b>	Green	25.59" / 65 cm	0.23" / 6 mm	1.968" / 50 mm	0.118" / 3 mm
<b>40-65008</b>	Blue	25.59" / 65 cm	0.31" / 8 mm	1.968" / 50 mm	0.157" / 4 mm
<b>40-65010</b>	Red	25.59" / 65 cm	0.39" / 10 mm	1.968" / 50 mm	0.157" / 4 mm
<b>40-65012</b>	Yellow	25.59" / 65 cm	0.49" / 12.5 mm	1.968" / 50 mm	0.157" / 4 mm
<b>40-65015</b>	Green	25.59" / 65 cm	0.59" / 15 mm	1.968" / 50 mm	0.228" / 5.8 mm
<b>40-65020</b>	Blue	25.59" / 65 cm	0.78" / 20 mm	1.968" / 50 mm	0.228" / 5.8 mm

PACKAGING: Bag of 3



### 31" / 80 CM SYMMETRIC DOUBLE-END WIRE BRUSHES

ITEM #	COLOR	TOTAL LENGTH	BRUSH 1 (Ø)	BRUSH LENGTH	SHEATH (Ø)
40-80005	White	31.496" / 80 cm	0.19" / 5 mm	1.968" / 50 mm	0.118" / 3 mm
40-80006	Green	31.496" / 80 cm	0.23" / 6 mm	1.968" / 50 mm	0.118" / 3 mm
40-80008	Blue	31.496" / 80 cm	0.31" / 8 mm	1.968" / 50 mm	0.157" / 4 mm
40-80010	Red	31.496" / 80 cm	0.39" / 10 mm	1.968" / 50 mm	0.157" / 4 mm
40-80012	Yellow	31.496" / 80 cm	0.49" / 12.5 mm	1.968" / 50 mm	0.157" / 4 mm
40-80015	Green	31.496" / 80 cm	0.59" / 15 mm	1.968" / 50 mm	0.228" / 5.8 mm
40-80020	Blue	31.496" / 80 cm	0.78" / 20 mm	1.968" / 50 mm	0.228" / 5.8 mm

PACKAGING: Bag of 3

## ASYMMETRIC BRUSHES

5 MM

6 MM



8 MM

10 MM



10 MM

12.5 MM



15 MM

20 MM




BATRIK

# REUSABLE IN-BRUSH™ ANTIMICROBIAL TWISTED WIRE BRUSHES

Our In-Brush™ Antimicrobial Reusable Twisted Wire Brushes are designed for use with Frazier/Baron suction tubes. The twisted wires have a fine/medium diameter yet are rigid enough to ensure excellent cleaning capabilities.

Our wide assortment of Nylon Twisted Wire Brushes are ideal for large bore, uniquely shaped, and endoscopic instruments.

GOOD TO KNOW: OUR TWISTED WIRE BRUSHES FEATURE A PROTECTIVE BLUE TIP FOR ADDED SAFETY.

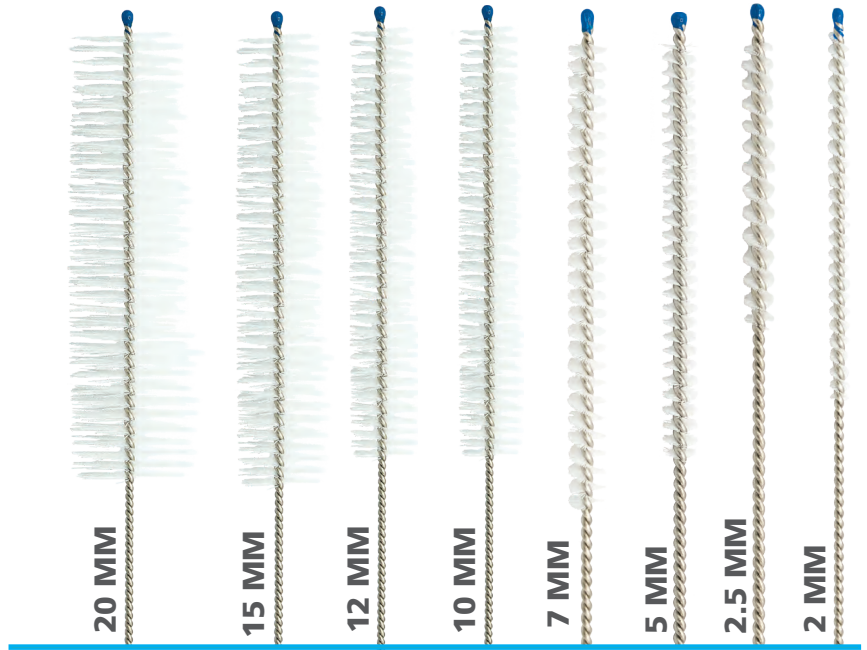

## 10" / 250 MM FRAZIER /BARON SUCTION TUBE BRUSHES

ITEM #	SIZE	BRUSH LENGTH	BRUSH (Ø)
<b>21-10039</b>	3 FR	0.5" / 12.7 mm	0.039" / 1 mm
<b>21-10067</b>	5 FR	0.5" / 12.7 mm	0.065" / 1.67 mm
<b>20-10059</b>	5 FR	1" / 25.4 mm	0.059" / 1.50 mm
<b>20-10079</b>	6 FR	1" / 25.4 mm	0.079" / 2 mm
<b>20-10090</b>	7 FR	1" / 25.4 mm	0.090" / 2.3 mm
<b>20-10106</b>	8 FR	1" / 25.4 mm	0.106" / 2.7 mm
<b>20-10118</b>	9 FR	1" / 25.4 mm	0.118" / 3 mm
<b>20-10129</b>	10 FR	1.5" / 38 mm	0.129" / 3.3 mm
<b>20-10146</b>	11 FR	1.5" / 38 mm	0.0146" / 3.7 mm
<b>20-10157</b>	12 FR	1.5" / 38 mm	0.157" / 4 mm

PACKAGING: Bag of 3.  
Also available in bags of 50, add the suffix "-50".

## 12" / 305 MM NYLON TWISTED WIRE BRUSHES

ITEM #	BRUSH LENGTH	BRUSH (Ø)
<b>22-12078</b>	2" / 50.8 mm	0.078" / 2 mm
<b>22-12118</b>	2" / 50.8 mm	0.118" / 3 mm
<b>22-12157</b>	3" / 76 mm	0.157" / 4 mm
<b>22-12197</b>	3" / 76 mm	0.197" / 5 mm
<b>22-12236</b>	3" / 76 mm	0.236" / 6 mm
<b>22-12393</b>	3" / 76 mm	0.394" / 10 mm

PACKAGING: Bag of 3.  
Also available in bags of 50, add the suffix "- 50".

### 16" / 406 MM NYLON TWISTED WIRE BRUSHES

ITEM #	BRUSH LENGTH	BRUSH (Ø)
<b>23-16197</b>	2.5" / 63.5 mm	0.197" / 5 mm
<b>23-16236</b>	2.5" / 63.5 mm	0.236" / 6 mm
<b>23-16276</b>	3" / 76 mm	0.276" / 7 mm
<b>23-16394</b>	3" / 76 mm	0.394" / 10 mm
<b>23-16472</b>	3" / 76 mm	0.472" / 12 mm
<b>23-16591</b>	3" / 76 mm	0.591" / 15 mm
<b>23-16787</b>	3" / 76 mm	0.787" / 20 mm

PACKAGING: Bag of 3.  
 Also available in bags of 50, add the suffix "-50".

### 18" / 457 MM NYLON TWISTED WIRE BRUSHES

ITEM #	BRUSH LENGTH	BRUSH (Ø)
<b>19-18010</b>	0.47" / 12 mm	0.039" / 1 mm
<b>19-18015</b>	0.5" / 12 mm	0.06" / 1.6 mm

PACKAGING: Bag of 3.  
 Also available in bags of 50, add the suffix "-50".

### 24" / 610 MM NYLON TWISTED WIRE BRUSHES

ITEM #	BRUSH LENGTH	BRUSH (Ø)
<b>24-24098</b>	2.5" / 63.5 mm	0.098" / 2.5 mm
<b>24-24197</b>	2.5" / 63.5 mm	0.197" / 5 mm
<b>24-24276</b>	3" / 76 mm	0.276" / 7 mm
<b>24-24394</b>	3" / 76 mm	0.394" / 10 mm
<b>24-24472</b>	3" / 76 mm	0.472" / 12 mm

PACKAGING: Bag of 3.  
 Also available in bags of 50, add the suffix "-50".



### BRONCHO BRUSHES

ITEM #	DESCRIPTION	BRUSH LENGTH	BRUSH (Ø) INCHES / MM
<b>26-06175</b>	Small	6" / 152 mm	1.75" / 44.5 mm
<b>26-08325</b>	Large	8" / 203 mm	3.25" / 82.6 mm

PACKAGING: Bag of 3.  
 Also available in bags of 50, add the suffix "-50".

### TRACHEOSTOMY BRUSHES

ITEM #	DESCRIPTION	TOTAL LENGTH	BRUSH LENGTH	BRUSH (Ø)
<b>25-08500</b>	Large 6-12.3 3 / 4"	8" / 203 mm	3.75" / 95 mm	0.5" / 12.7 mm
<b>25-37550</b>	Small 000-5.1 7 / 8"	4" / 93 mm	1.88" / 48 mm	0.5" / 12.7 mm

PACKAGING: Bag of 3. Also available in bags of 50, add the suffix "-50".

# BATRIK STAINLESS STEEL ARC™ MICRO BRUSHES

Our Stainless Steel ARC™ Micro Brushes are some of the smallest brushes in the industry. They are designed to clean channels with an inner diameter ranging from 0.5 mm to 1.40 mm.

- Flexible stainless steel sheath will not bend like the regular metal twisted wire brushes
- Polyamide bristles prevents or slows bacterial growth between sterilizations
- Reusable\*
- Available in two lengths

TIP: USING THE PULL-THROUGH TECHNIQUE IS HIGHLY RECOMMENDED WITH THESE BRUSHES. SIMPLY INSERT THE BLUNT, NON-BRUSH END FIRST.



## PRODUCT SELECTION

ITEM #	TOTAL LENGTH	BRUSH SIZE (Ø)	BRUSH LENGTH	SHEATH (Ø)	PACKAGING
<b>35-03015</b>	12" / 30 cm	0.06" / 1.5 mm	0.6" / 15 mm	0.02" / 0.40 mm	6 / Tube
<b>35-04515</b>	18" / 45 cm	0.06" / 1.5 mm	0.6" / 15 mm	0.02" / 0.40 mm	6 / Tube

\*Although all materials used in the fabrication of these brushes can be re-processed repeatedly, due to the delicate nature of these brushes re-use is conditional on their handling and application.

# BATRIK INSTRUMENT CLEANING FLAT BRUSHES

Our Instrument Cleaning Flat Brushes are specifically designed for hospital use in the pre-clean and preparation of instruments.

- Long ergonomic handle permits ease-of-use
- Extra long length of brush head (60 [±5] mm in Single-Use and 75 [±5] mm in Reusable)
- Transparent brush fibers allow for quick verification of its cleanliness
- The plastic handle will not release any toxic components
- Available in 6 models (4 reusable/autoclavable and 2 single-use)



## FLEXIBLE

This brush is designed for large instruments that require a more gentle cleaning application such as optics, plastics, or hard to reach spot.



## EXTRA FLEXIBLE

The very long, flexible bristles eliminate any stress during sanitation, and makes this the ideal brush for fragile instruments as well as both pliable and rigid endoscopes. The length of the bristles assist in cleaning areas that would otherwise be difficult to reach, such as the accessory channel or tip deflection of pliable endoscopes.



## RIGID

The rigidity of this brush, along with the length of the bristles and handle, improves and simplifies instrument cleaning.



## EXTRA RIGID

Instrument safe polyamide fibers offer steel-like rigidity and allows deep cleaning of soiled instruments, particularly in grooved or uneven surfaces.



50-40919



50-40921



50-40920



50-40918



50-40922



50-40917

## PRODUCT SELECTION

### DELTRIN PLASTIC


ITEM #	BRUSH TYPE	DESCRIPTION	COLOR	BRUSH LENGTH	BRUSH SIZE	PACKAGING
<b>50-40917</b>	Extra Rigid	Reusable Autoclavable Flat Brush	Blue	220 mm	75 mm x 20 mm x 10 mm	Bag of 3
<b>50-40918</b>	Rigid	Reusable Autoclavable Flat Brush	Green	220 mm	75 mm x 20 mm x 15 mm	Bag of 3
<b>50-40919</b>	Flexible	Reusable Autoclavable Flat Brush	Yellow	220 mm	75 mm x 20 mm x 15 mm	Bag of 3
<b>50-40920</b>	Extra Flexible	Reusable Autoclavable Flat Brush	Red	220 mm	75 mm x 20 mm x 25 mm	Bag of 3
<b>50-40921</b>	Flexible	Single-Use Flat Brush	Orange	220 mm	60 mm x 20 mm x 15 mm	Bag of 20
<b>50-40922</b>	Rigid	Single-Use Flat Brush	Light Blue	220 mm	60 mm x 20 mm x 15 mm	Bag of 20

Reusable models are also available in bags of 50. Add the suffix "- 50".

# BATRIK REUSABLE ORTHOPEDIC RASPS CLEANING BRUSH

Our Reusable Orthopedic Rasps Cleaning Brush is used to clean orthopedic rasps and other similar instruments used in orthopedic surgery.

- Long lasting and can withstand many sterilization cycles
- Ergonomic steamable thermoplastic handle
- Antibacterial polyamide bristles
- Reusable
- Latex free

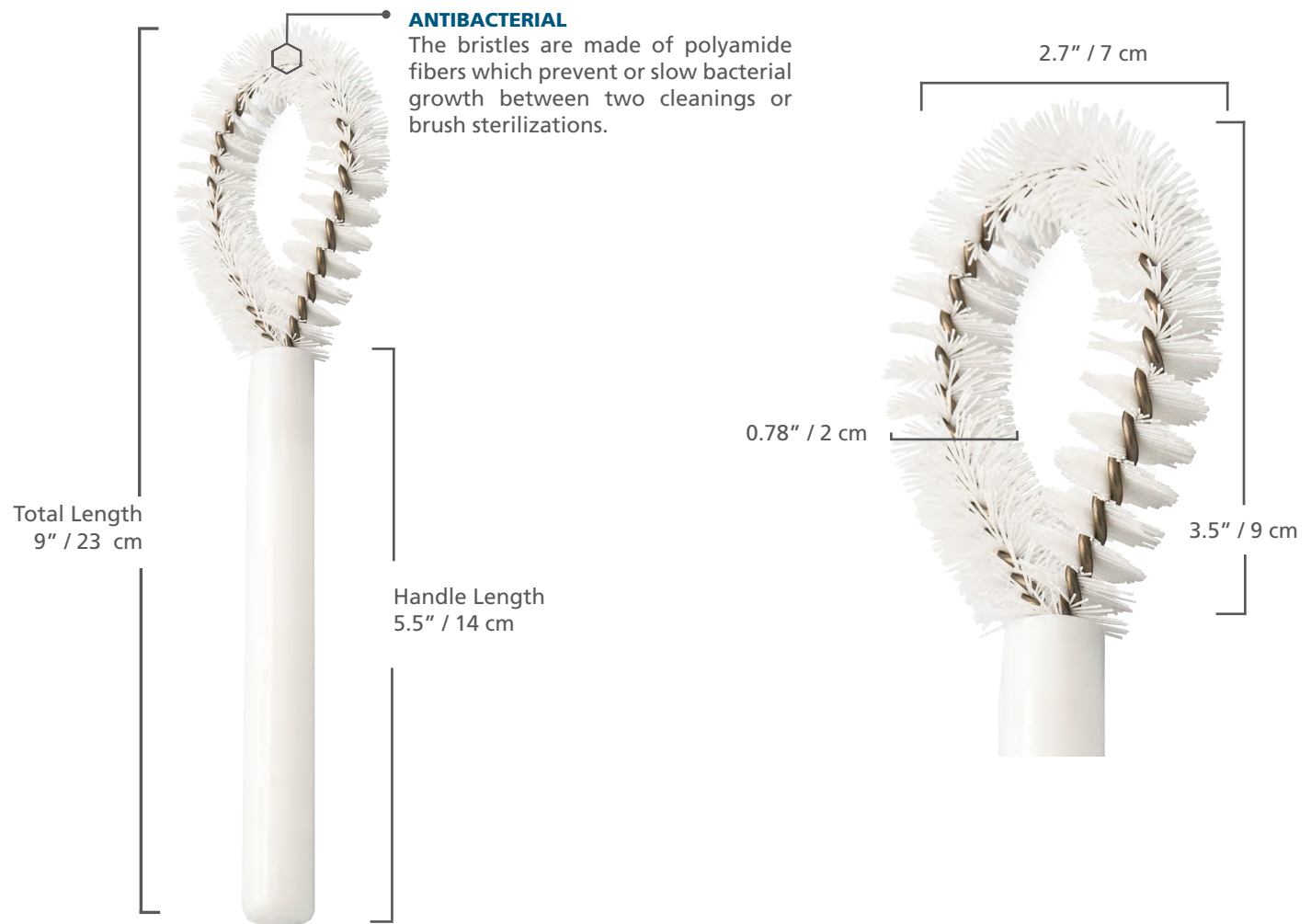

## PRODUCT DESCRIPTION

ITEM #	TOTAL LENGTH	HANDLE LENGTH	BRUSH SIZE	PACKAGING
50-RASP	9" / 23 cm	5.5" / 14 cm	0.78" x 2.7" x 3.5" (2 cm x 7 cm x 9 cm)	3 / Bag


BATRIK

# REUSABLE PRECISION BRUSH

## INSTRUMENT CLEANING

Our Reusable Precision Instrument Cleaning Brush is the ideal brush for the quick removal of difficult debris. It cleans externally as well as deep into instrument jaws. Although the firm and durable bristles feel like stainless steel, there is absolutely no risk of instrument damage or personal injury.

- Bristles have an extra firm density and superior durability
- Polyamide (TYNEX®) bristles
- Steamable thermoplastic ergonomic handle
- Steam / chemical sterilizable
- Latex free

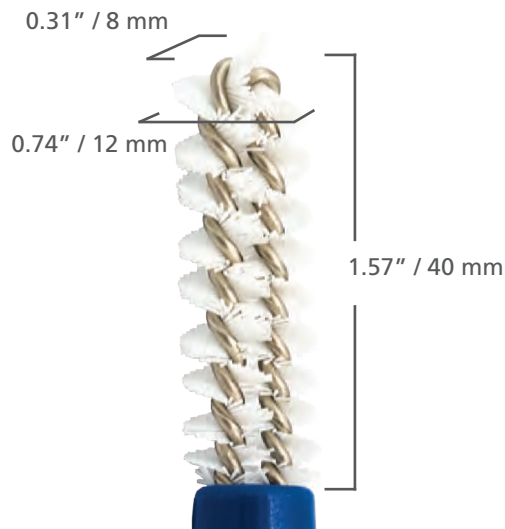

GOOD TO KNOW: ELIMINATES THE NEED FOR THE STAINLESS STEEL TOOTHBRUSH-STYLE BRUSH, WHICH CAN DAMAGE VALUABLE INSTRUMENTS!

## PRODUCT SELECTION

### BLUE DELRIN PLASTIC

ITEM #	DESCRIPTION	TOTAL LENGTH
50-40800	Reusable Precision External Cleaning Brush	6.29" / 160 mm

PACKAGING: Bag of 3




# BATRIK REUSABLE INSTRUMENT CLEANING BRUSHES

Designed specifically for hospital applications, our Reusable Instrument Cleaning Brushes feature long and flexible bristles to maximize the effectiveness of large-surface cleaning.

- High quality polyamide clear bristles
- Easy to identify debris
- Autoclavable
- Latex free

**REUSABLE JAR/BOTTLE BRUSH**

**STERILIZER DRAIN BRUSH**

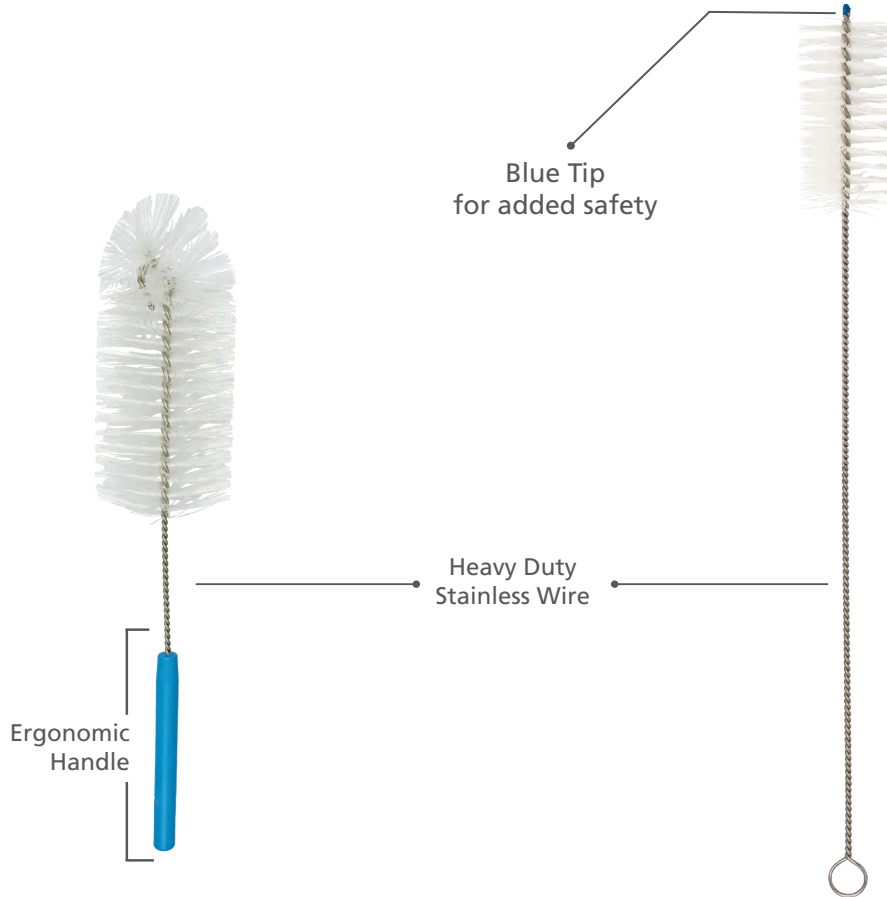

## PRODUCT SELECTION

### BLUE DELRIN PLASTIC

ITEM #	DESCRIPTION	HANDLE LENGTH	BRUSH LENGTH	BRUSH (Ø)	TOTAL LENGTH
<b>50-10295</b>	Reusable Jar / Bottle Brush	120 mm ± 2 mm	120 mm ± 5 mm	7 mm ± 2 mm	13.19" / 355 mm
<b>50-10224</b>	Sterilizer Drain Cleaning Brush	N/A	4" / 10.1 cm	2" / 5.08 cm	24" / 609 mm

PACKAGING: Bag of 3



BATRIK

# ACETABULAR REAMER BRUSH

Our double-ended cleaning brush is manufactured exclusively for cleaning acetabular bone reamers. The "Pig Tail" rounded style circular head brush is designed to effectively clean inside the reamer, while the rigid straight brush on the reverse side is used for cleaning individual cutting holes. Special fusion of the brush heads to the hand piece along with its high quality antimicrobial medical grade bristles ensures a durable, reusable, and highly effective cleaning tool.



STRAIGHT BRUSH



1.18" / 30 mm length

ROUND BRUSH



2" / 50 mm length

## PRODUCT SELECTION

### GREEN DELRIN PLASTIC

ITEM #	DESCRIPTION	HANDLE LENGTH	STRAIGHT BRUSH	ROUND BRUSH HEAD
<b>50-ACETAB</b>	Acetabular Reamer Cleaning Brush	4.72" / 120 mm ø 0.47" / 12 mm	1.18" / 30 mm ø 0.3125" / 8 mm	2" / 50 mm ø 1.75" / 45 mm

PACKAGING: Bag of 3

NOTE: Specifications subject to small variances

# BATRIK ALL PURPOSE REUSABLE BRUSHES

## Plastic Handle Brushes

We offer a wide selection of brushes for instrument cleaning applications. These brushes are ideal for large bore, uniquely shaped, and endoscopic instruments.

### In-Brush™ Plastic Handle

Available in regular toothbrush style (nylon or stainless) and double-ended brushes.



28-72550



28-69394



29-72550



29-69394



29-70050



50-40911

## PRODUCT SELECTION

ITEM #	DESCRIPTION	LENGTH TOTAL	BRUSH LENGTH	BRUSH HEIGHT
<b>28-72550</b>	Plastic Single End Nylon Reusable Handle Brush	7.25" / 184 mm	1.5" / 38 mm	0.5" / 12.7 mm
<b>29-72550</b>	Plastic Single End Stainless Steel Reusable Handle Brush	7.25" / 184 mm	1.5" / 38 mm	0.5" / 12.7 mm
<b>29-70050</b>	Plastic Single End Brass Reusable Handle Brush	7.25" / 184 mm	1.5" / 38 mm	0.5" / 12.7 mm
<b>28-69394</b>	Plastic Double End Nylon Reusable Handle Brush	6.9" / 175 mm	1.38" / 35 mm 1" / 25 mm	0.39" / 10 mm 0.2" / 5 mm
<b>29-69394</b>	Plastic Double End Stainless Steel Reusable Handle Brush	6.9" / 175 mm	1.38" / 35 mm 1" / 25 mm	0.39" / 10 mm 0.2" / 5 mm
<b>50-40911</b>	Hand Scrub Nylon Brush Reusable autoclavable	4.09" / 104 mm	3.74" / 95 mm	1.1" / 28 mm

**PACKAGING:** Bag of 3 - For a bag of 50, add "-50" following the item code.

BATRIK

# AUTOCLAVE CLEANING BRUSH KITS



## ROUND BRUSH

- Features 3 holes for adjustable pole placement
- Adjustable brush placement (parallel or perpendicular) and longer brush length provides easy access to hard-to-reach corners
- Made with FDA approved polyester bristles

## HAND BRUSH

- Ergonomic handle for better hand comfort
- Longer, FDA approved polyester bristles
- Easy to clean (no cavity)

## PRODUCT SELECTION

### REGULAR AUTOCLAVE CLEANING BRUSH KIT

ITEM # (KIT)	DESCRIPTION	ITEM # (UNIT)	DESCRIPTION	DIMENSIONS
<b>50-CLAVEREG</b>	Autoclave REGULAR Brush Cleaning Kit	<b>50-35900</b>	1 Hand Brush	Length: 3.5" x 8" / 8.89 cm x 20.32 cm
		<b>50-55800</b>	1 Rounded Corner Brush (Head Brush)	Head size: 5.5 x 10" / 14 x 25.4 cm
		<b>50-EXTPOLE</b>	1 Anodized Aluminum Extension Pole	2-4 ft / 60-122 cm

### SUPERIOR AUTOCLAVE CLEANING BRUSH KIT

ITEM # (KIT)	DESCRIPTION	ITEM # (UNIT)	DESCRIPTION	DIMENSIONS
<b>50-CLAVESUP</b>	Autoclave SUPERIOR Brush Cleaning Kit	<b>50-35900</b>	1 Hand Brush	Length: 3.5" x 8" / 11.5 cm x 20.32cm
		<b>50-55800</b>	1 Rounded Corner Brush (Head Brush)	Head size: 5.5 x 10" / 14 x 25.4 cm
		<b>50-EXTPOLE</b>	1 Anodized Aluminum Extension Pole	2-4 ft / 60-122 cm
		<b>SF-100</b>	2 SAVEFACE™ Face Shields	9" / 23 cm

## REPLACEMENT COMPONENTS



**HAND BRUSH**  
**50-35900**

Packaging: 1



**BRUSH HEAD**  
**50-55800**

Packaging: 1



**EXTENSION POLE**  
**50-EXTPOLE**

Packaging: 1



**SAVEFACE™**  
**FACE SHIELDS SF-100**

Packaging: 100 / Box



# SURGMED

GROUP 

Proudly Canadian, **The Surgmed™ Group of Companies** began operations in 1992 as a single distribution company known as **Advanced Surgi-Pharm Inc.** The company quickly became known as an industry leader and specialty distributor of value-added surgical products in the Canadian hospital market. Five years later, **Surgmed expanded upon its success with the launch of Batrik Medical Manufacturing Inc.** in 1997 and the subsequent acquisition of **Imperial Surgical Ltd.** one year later in 1998.




**Together, our three unique companies serve the Canadian and global healthcare markets with high quality and competitively priced products, all backed by world-class customer care.**



Founded in 1992, Advanced Surgi-Pharm Inc. is the distribution arm of Surgmed™ in Canada, representing leading US, international and domestic manufacturers for a wide range of specialty products with unique benefits for applications in the Operating Room, Women's Health, Labor & Delivery, ICU, and Sterile Processing areas. Equipped with a direct sales team, ASP provides effective nationwide coverage as the Canadian distributor of Batrik, Vygon, Rocket Medical, Cooper Surgical, SunMed, and other industry leading products including a wide-ranging portfolio of top-rated PPE products.

Established in 1997 and sold in 85 countries globally, Batrik Medical Manufacturing Inc. specializes in the manufacturing of endoscopy accessories, enzymatic solutions, infection control, central sterilization, and surgical devices. Batrik's in-house controlled environment production facility adheres to strictly enforced quality and regulatory control systems. They are registered as a device manufacturing establishment with the United States FDA and are in full compliance with CE and Health Canada Regulations. All Batrik products are rigorously tested, liability insured, and meet/exceed high international quality standards.

Established in 1935, Imperial Surgical Ltd. is one of Canada's oldest hospital capital equipment manufacturers, well-known for their stainless steel healthcare equipment for Sterile Processing, Operating Room, Examination and Treatment, Patient Handling, and Mortuary/ Pathology areas. Their growing portfolio of products also has applications in Infection Control & Prevention, an important area of focus for their future.

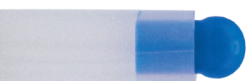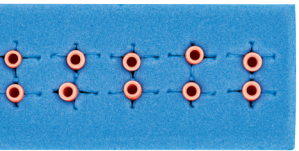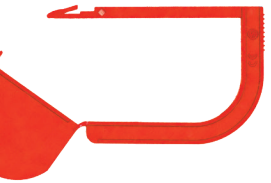

SURGMED GROUP

2165 46th Ave., Lachine, Quebec H8T 2P1  
Phone: +1 514-631-7988 | Toll Free: 800-661-5432 | Fax.: +1 514-631-9083  
Email: [info@surgmed.com](mailto:info@surgmed.com) | Website: [www.batrik.com](http://www.batrik.com)