



COMPLETE PRODUCT CATALOGUE

Established in 1997. Proudly Canadian. 🇨🇦
Sold in over 85 countries worldwide.



Proudly Canadian and established in 1997, we specialize in the manufacturing of enzymatic solutions offered in different formats as well as the fabrication of silicone and plastic products. Our products are sold in over 80 countries worldwide and used primarily in surgery, endoscopy, and central sterilization areas.

Our in-house controlled environment production facility in Montreal, Canada adheres to strictly enforced quality and regulatory control systems. We are registered as a device manufacturing establishment with the United States FDA and we are in full compliance with CE and Health Canada Regulations. We continually invest in upgrading our plant, machinery, and processes to ensure high quality products at high volumes. Most products are available off-the-shelf.

At Batrik, we are committed to innovation and providing our customers with the finest quality products backed by world-class customer care.

CONTENTS

SURGICAL PRODUCTS & GENERAL ACCESSORIES

04	SUTURE CLAMP TAGS	
05	VASCULAR SILICONE TIES	
06	CLAMP JAW COVERS	
07	G.O.L.F.F.™ ANTI-FOG SOLUTION	
08	SKINSTM Surgical MULTI-TIP Markers	
09	PANDAPATCH™ EYE SAFETY TAPE	NEW!
11	MAGNETIC MAT	
12	PIN GUARDS	
13	CAUTERY TIP CLEANER	
14	SAFE-SOAK™ SPAC & MATS	
15	BLUE CLOUD-GEL™ ANTI-FATIGUE MATS	
16	BLUE-CLOUD NIMBUS™ ANTI-FATIGUE MATS	
17	STACKABLE STEP STOOLS	
18	TOURNY™+ TOURNIQUETS	
20	SHOOZCOVERS™ SOFT LINE™	
22	SHOOZCOVERS™ GRIPPEEZ™	NEW!
24	SHOOZCOVERS™ COATED LINE™	
26	SHOOZCOVERS™ TOUGH LINE™	

ENDOSCOPY PRODUCTS

28	BITE BLOCKS Endoscopy Bo2™ / Bitestrong™	
29	ENDOSCOPY BITESTRONG™ BITE BLOCKS	
30	COUVRETTE™ VALVE COVERS	
31	DAYBRUSH™ CHANNEL BRUSHES	
32	ARC™ CHANNEL BRUSHES	
35	VALVE / CONTROL BRUSHES	
36	V-BRUSHES	

INFECTION CONTROL & STERILE PROCESSING

37	TAG-IT™	NEW!
42	SILICONE PIN MAT	
43	CORNER POCKETS™ & SOFT POCKETS™	
44	LILA'S INSTRUMENT STAIN REMOVER	
45	INSTRUMENT PROTECTORS	
48	PrAmLi Power ³ LUBE™ INSTRUMENT LUBRICANT	
49	S.T.I.C.S.	
50	ID TAPE	
51	TAPE-B-GONE™ TAPE REMOVER	

52	PADLOCKS SECURITY TAGS	
56	SAVEFACE™ DRAPEP SHIELDS	
58	SCRUB PAD	
59	Your First <input checked="" type="checkbox"/> ™ ENZYMATIC & NON-ENZYMATIC CLEANERS	
60	Your First <input checked="" type="checkbox"/> ™ ENVIRO-SACHET™	
61	Your First <input checked="" type="checkbox"/> ™ FOAMING SPRAY	
62	Your First <input checked="" type="checkbox"/> ™ F1RST WIPES	
63	Your First <input checked="" type="checkbox"/> ™ DRY SPONGES	
64	Your First <input checked="" type="checkbox"/> ™ RIGHTSWIPE™ ENZYMATIC SPONGES BEDSIDE SPONGE KITS	NEW!
65	Your First <input checked="" type="checkbox"/> ™ RIGHTSWIPE™ NON-ENZYMATIC DETERGENT BEDSIDE SPONGE KITS	NEW!
66	Your First <input checked="" type="checkbox"/> ™ TRANSPORT BAGS FOR ENDOSCOPE AND MEDICAL INSTRUMENTS	

BRUSHES


67	ARC™ WIRE STYLE BRUSHES	
71	DOUBLE ARC™ WIRE BRUSHES Asymmetrics / Symmetric	
74	IN-BRUSH™ ANTIMICROBIAL TWISTED WIRE BRUSHES	
76	STAINLESS STEEL ARC™ MICRO BRUSHES	
77	INSTRUMENT CLEANING FLAT BRUSHES	
78	PRECISION BRUSH	
79	TRIFECTA™ BRUSH	NEW!
80	ALL PURPOSE REUSABLE BRUSHES	
81	ACETABULAR REAMER BRUSH	
82	ORTHOPEDIC RASPS CLEANING BRUSH	
83	BOTTLE & DRAIN CLEANING BRUSHES	
84	AUTOCLAVE CLEANING BRUSH KITS	

BATRIK SUTURE CLAMP TAGS


Non-toxic, radiopaque, and latex free, our Suture Clamp Tags are designed to protect and safely hold sutures without any damage.

Additionally, its foam base, affixed to a sterile area with the clamps, allows for easy insertion.

- Packed in a convenient foam base with easy peel-off adhesive backing, even while wearing gloves
- Available in sterile and non-sterile
- Available in 2, 3, or 5 pairs
- Pre-counted

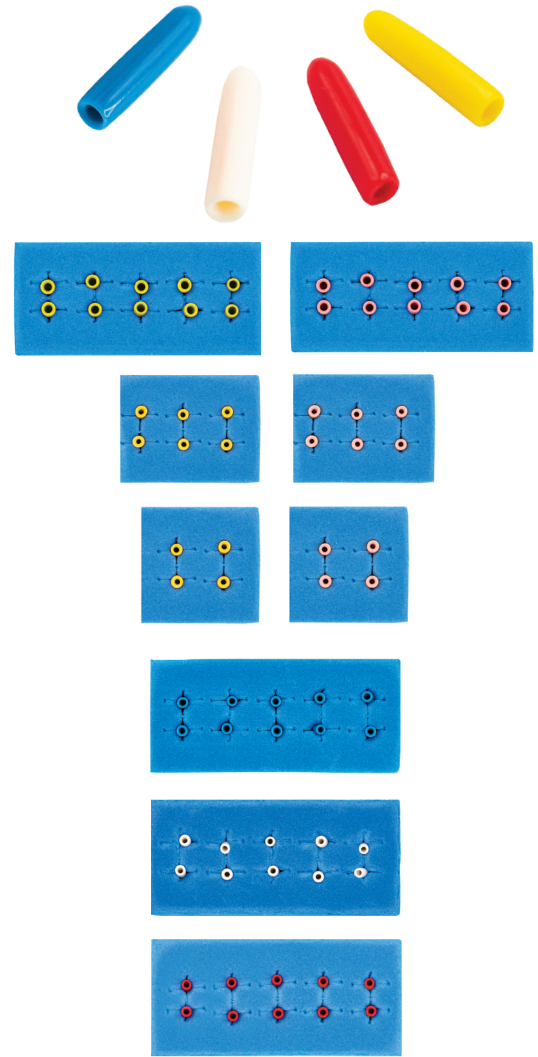

OUR "EASY PEEL CORNER" IS A FEATURE UNIQUE TO BATRIK

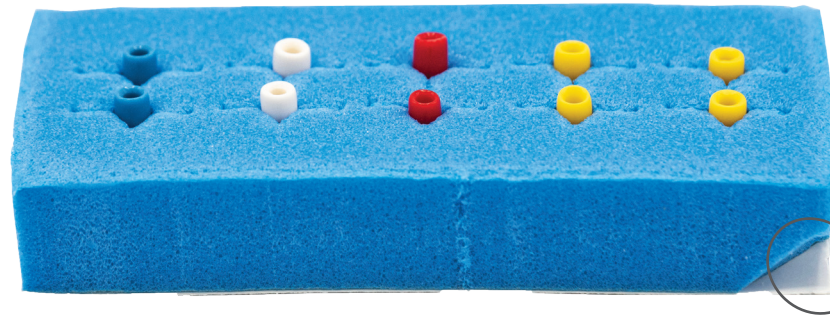

PRODUCT SELECTION


ITEM # STERILE	ITEM# NON-STERILE	COLOR	DIMENSIONS	PACKAGING
11-1002	911-1002	Yellow	Standard \varnothing (0.06" x 0.5" / 0.15 cm x 1.27 cm)	2 pairs / foam base
11-1003	911-1003	Yellow	Standard \varnothing (0.06" x 0.5" L / 0.15 cm x 1.27 cm)	3 pairs / foam base
11-1005	911-1005	Yellow	Standard \varnothing (0.06" x 0.5" L / 0.15 cm x 1.27 cm)	5 pairs / foam base
11-1025	911-1025	Assorted	Standard \varnothing (0.06" x 0.5" L / 0.15 cm x 1.27 cm)	5 pairs / foam base
11-1045	911-1045	White	Standard \varnothing (0.06" x 0.5" L / 0.15 cm x 1.27 cm)	5 pairs / foam base
11-1065	911-1065	Blue	Standard \varnothing (0.06" x 0.5" L / 0.15 cm x 1.27 cm)	5 pairs / foam base
11-1085	911-1085	Red	Standard \varnothing (0.06" x 0.5" L / 0.15 cm x 1.27 cm)	5 pairs / foam base
11-3002	911-3002	Pink	Large \varnothing (0.09" x 0.75" L / 0.23 cm x 1.91 cm)	2 pairs / foam base
11-3003	911-3003	Pink	Large \varnothing (0.09" x 0.75" L / 0.23 cm x 1.91 cm)	3 pairs / foam base
11-3005	911-3005	Pink	Large \varnothing (0.09" x 0.75" L / 0.23 cm x 1.91 cm)	5 pairs / foam base

STERILE PACKAGING: 1 Foam Base / Pouch – 5 Pouches / Box

NON-STERILE PACKAGING: Sold per foam base

BATRIK VASCULAR SILICONE TIES

AVAILABLE IN 4 SIZES
AND 4 COLORS

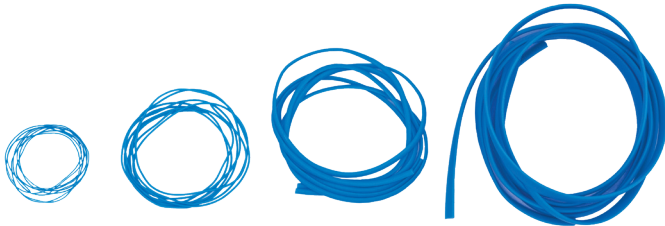


PRODUCT SELECTION

ITEM #	DIMENSIONS	COLOR
10-1000	0.075" X 16" (0.2 cm x 40 cm)	Mini-Red
10-1001	0.075" X 16" (0.2 cm x 40 cm)	Mini-Blue
10-1002	0.075" X 16" (0.2 cm x 40 cm)	Mini-Yellow
10-1003	0.075" X 16" (0.2 cm x 40 cm)	Mini-White
10-1010	0.120" x 16" (0.3 cm x 40 cm)	Maxi-Red
10-1011	0.120" x 16" (0.3 cm x 40 cm)	Maxi-Blue
10-1012	0.120" x 16" (0.3 cm x 40 cm)	Maxi-Yellow
10-1013	0.120" x 16" (0.3 cm x 40 cm)	Maxi-White
10-1101	0.021" x 16" (0.5 cm x 40 cm)	Extra-Mini-Blue
10-1111	0.200" x 22" (0.5 cm x 56 cm)	Extra-Maxi-Blue
10-1113	0.200" x 22" (0.5 cm x 56 cm)	Extra-Maxi-White
10-1800	0.075" x 18" / 0.2 cm x 46 cm	Mini-Red
10-1801	0.075" x 18" / 0.2 cm x 46 cm	Mini-Blue
10-1802	0.075" x 18" / 0.2 cm x 46 cm	Mini-Yellow
10-1803	0.075" x 18" / 0.2 cm x 46 cm	Mini-White
10-1810	0.12" x 18" / 0.3 cm x 46 cm	Maxi-Red
10-1811	0.12" x 18" / 0.3 cm x 46 cm	Maxi-Blue
10-1812	0.12" x 18" / 0.3 cm x 46 cm	Maxi-Yellow
10-1813	0.12" x 18" / 0.3 cm x 46 cm	Maxi-White

PACKAGING: 2 ties per sterile pouch – 10 pouches per box

Our Vascular Silicone Ties permit safe and smooth retraction, occlusion and identification of vessels during surgical procedures.

Unlike clamps and retractors which can pinch the vessels flat, or even crack them, the BATRIK elliptical shape allows for the gentle handling of arteries, veins, nerves, tendons, spermatic cords, ureters, or whenever delicate retraction / occlusion is required.

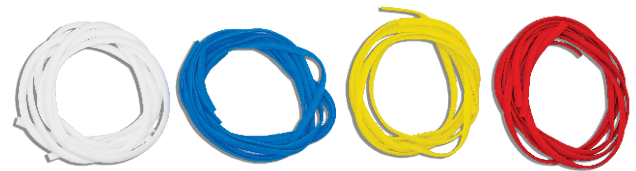

- Manufactured of medical grade silicone
- Do not twist on themselves when unpacked
- Easy to separate, eliminating the need for
- Sterile, radiopaque, and latex free cardboard inserts

PRODUCT SELECTION

ITEM #	DIMENSIONS	COLOR
10-1831	0.04" x 18" / 0.1 cm x 46 cm	Extra-Mini-Blue
10-1841	0.20" x 18" / 0.5 cm x 46 cm	Extra-Maxi-Blue
10-1843	0.20" x 18" / 0.5 cm x 46 cm	Extra-Maxi-White
10-3010	0.075" x 32" (0.2 cm x 81 cm)	Maxi-Red
10-3202	0.075" x 32" (0.2 cm x 81 cm)	Mini-Yellow
10-3200	0.075" x 32" (0.2 cm x 81 cm)	Mini-Red
10-3201	0.075" x 32" (0.2 cm x 81 cm)	Mini-Blue
10-3202	0.075" x 32" (0.2 cm x 81 cm)	Mini-Yellow
10-3203	0.075" x 32" (0.2 cm x 81 cm)	Mini-White
10-3210	0.120" x 32" (0.3 cm x 81 cm)	Maxi-Red
10-3211	0.120" x 32" (0.3 cm x 81 cm)	Maxi-Blue
10-3212	0.120" x 32" (0.3 cm x 81 cm)	Maxi-Yellow
10-3213	0.120" x 32" (0.3 cm x 81 cm)	Maxi-White

PACKAGING: 2 ties per sterile pouch – 10 pouches per box

BATRIK CLAMP JAW COVERS

Manufactured from medical grade silicone, our Clamp Jaw Covers reduce trauma to vessels requiring occlusion, particularly when a hemostatic clamp is utilized. The thread-like indents allow for gentle, constant and non-slip tension.



GOOD TO KNOW: OUR CLAMP JAW COVERS CONVERT ALMOST ANY INSTRUMENT TO ATRAUMATIC



- Gamma sterilized
- Available in 3 colors and 3 widths to fit most clamp sizes
- Radiopaque and latex free

PRODUCT SELECTION

ITEM # STERILE	ITEM # NON-STERILE	DIMENSIONS	COLOR
12-1000	912-1000	0.1" X 5" / 0.25 cm X 12.7 cm	Red
12-1001	912-1001	0.2" X 5" / 0.51 cm X 12.7 cm	White
12-1002	912-1002	0.3" X 5" / 0.76 cm X 12.7 cm	Blue

STERILE PACKAGING: 2 clamps per sterile pouch - 10 pouches per box

NON-STERILE PACKAGING: 20 clamps per bag (also available 100 units / bag)

BATRIK G.O.L.F.F.™ ANTI-FOG SOLUTION

Our G.O.L.F.F.™ Anti-Fog Solution solves the fogging problem often encountered during Laparoscopic surgical procedures.

- Prevents the occurrence of fog on lenses
- Ensures clearer vision and avoids unnecessary interruptions during surgery
- Set includes an easy-squeeze 6 mL bottle and an adhesive radiopaque sponge for secure placement
- OEM (available)
- Sterile
- Latex free



PRODUCT SELECTION

ITEM #	DESCRIPTION	PACKAGING	CASE DIMENSIONS
G.O.L.F.F.	Anti-Fog Solution & Adhesive Sponge set	20 / Box	20" x 16" x 18"

BATRIK SKINS™ SURGICAL MULTI-TIP SKIN MARKERS

Our SKINS™ Surgical Multi-Tip Skin Markers are designed with a beveled tip to enable the creation of fine or bold tracing lines by holding the pen at different angles.

- Available with or without 6 labels and a flexible 6" calibrated (in/mm) ruler
- Available in non-sterile for pre-operative use, or sterile for use in the operating room
- Surgical prep resistant
- Made of lasting gentian violet ink that is non-toxic and non-tattooing
- Latex free





PRODUCT SELECTION

ITEM #	TYPE	DESCRIPTION	PACKAGING
S-3000	Sterile	Multi-Tip fine marker with ruler and label set	50 / Box
S-3001	Sterile	Multi-Tip fine marker	50 / Box
S-30005	Sterile	Multi-Tip fine marker with ruler and label set (Volume Packaging)	500 / Box
S-30015	Sterile	Multi-Tip fine marker (Volume Packaging)	500 / Box
S-93000	Non-Sterile	Multi-Tip fine marker with ruler and label set	1000 / Box
S-93001	Non-Sterile	Multi-Tip fine marker	1000 / Box
S-94000	Non-Sterile	Ruler	100 / Bag
S-94001	Non-Sterile	Ruler (Bulk)	20,000 / Box
S-95000	Non-Sterile	Set of 3 labels	100 / Bag

SAVE by ordering our SKINS™ MARKERS in volume packaging. Cost effective and reduces waste.

NEW!

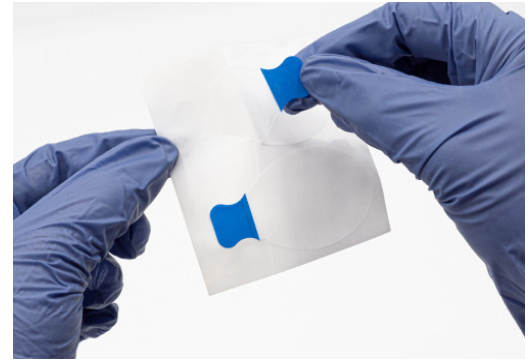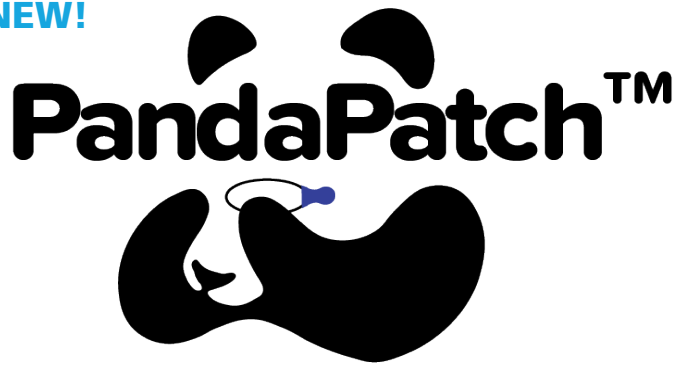

Introducing PandaPatch™: Innovative Adhesive Eye Protection for Safer Patient Care

PandaPatch™ offers an effective alternative to traditional tape for protecting patients' eyes during surgery. These adhesive eye protectors are designed to reduce the risk of corneal injuries and ophthalmic complications for patients under general anesthesia.

With PandaPatch™, you ensure superior patient outcome by preventing avoidable eye damage associated with the use of standard surgical tape.

KEY FEATURES:

- **Ergonomic Design:** Contours comfortably around the eyes.
- **Post-Surgery Friendly:** Particularly effective for children, the elderly, and individuals with sensitive or delicate skin.
- **Latex-Free & Hypoallergenic Adhesive:** Made from medical-grade material, gentle on the skin and designed to minimize irritation.
- **Non-Adhesive Colored Tabs:** Allowing quick and hassle free application and removal, facilitating effortless monitoring of pupil dilation during surgery.



70-04ADU
Adult Sensitive Eye Protector

HERE'S WHY IT MATTERS:

- During prolonged surgery under general anesthesia, up to 59%¹ of patients may not fully close their eyes, risking corneal abrasions.
- PandaPatch™ gentle adhesive eye protectors ensure eyelids remain safely and comfortably closed, minimizing the risks typically associated with the use of standard surgical tape.



REFERENCES


1. Anson, Jonathan, M.D. "Perioperative Corneal Abrasions: Etiology, Prevention, and Management." Clinical Updates, Penn State University College of Medicine, 2021.




PANDAPATCH™ RANGE:

Three adhesive eye protectors to accommodate different age groups and skin sensitivities. Providing optimal protection for every patient.

- Adult Extra Sensitive Adhesive (Blue Tabs)
- Adult Adhesive (Red Tabs)
- Pediatric Adhesive (Red Tabs)

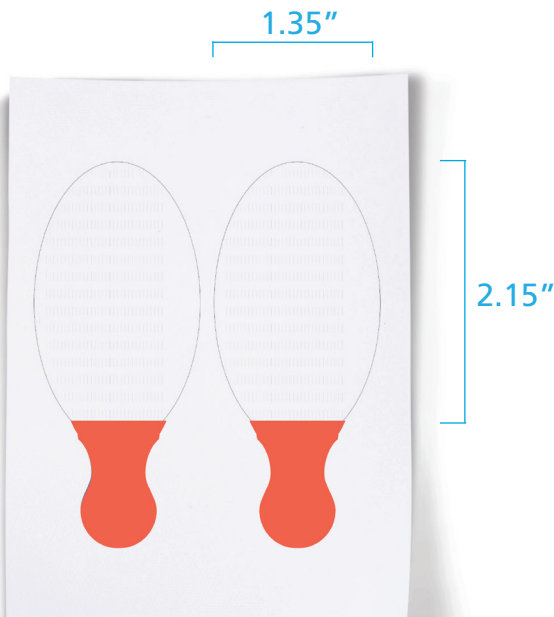

70-04ADU
Adult Adhesive Eye Protector

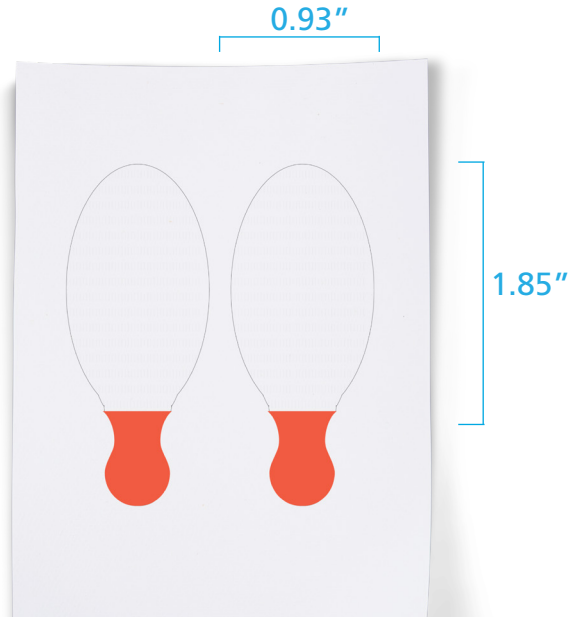


70-04PED
Pediatric Adhesive Eye Protector

PRODUCT SELECTION

ITEM #	DESCRIPTION	DIMENSIONS	PACKAGING
70-02ADU	ADULT EXTRA SENSITIVE ADHESIVE EYE PROTECTORS (BLUE)	2.15" x 1.35"	50 / Box
70-04ADU	ADULT ADHESIVE EYE PROTECTORS (RED)	2.15" x 1.35"	50 / Box
70-04PED	PEDIATRIC ADHESIVE EYE PROTECTORS (RED)	1.85" x 0.93"	50 / Box

BATRIK MAGNETIC MAT

Our Magnetic Mat provides a sterile, neutral field for the safe placement and hands-free transfer of instruments. It is latex free, reusable, and autoclavable (up to 134° C).

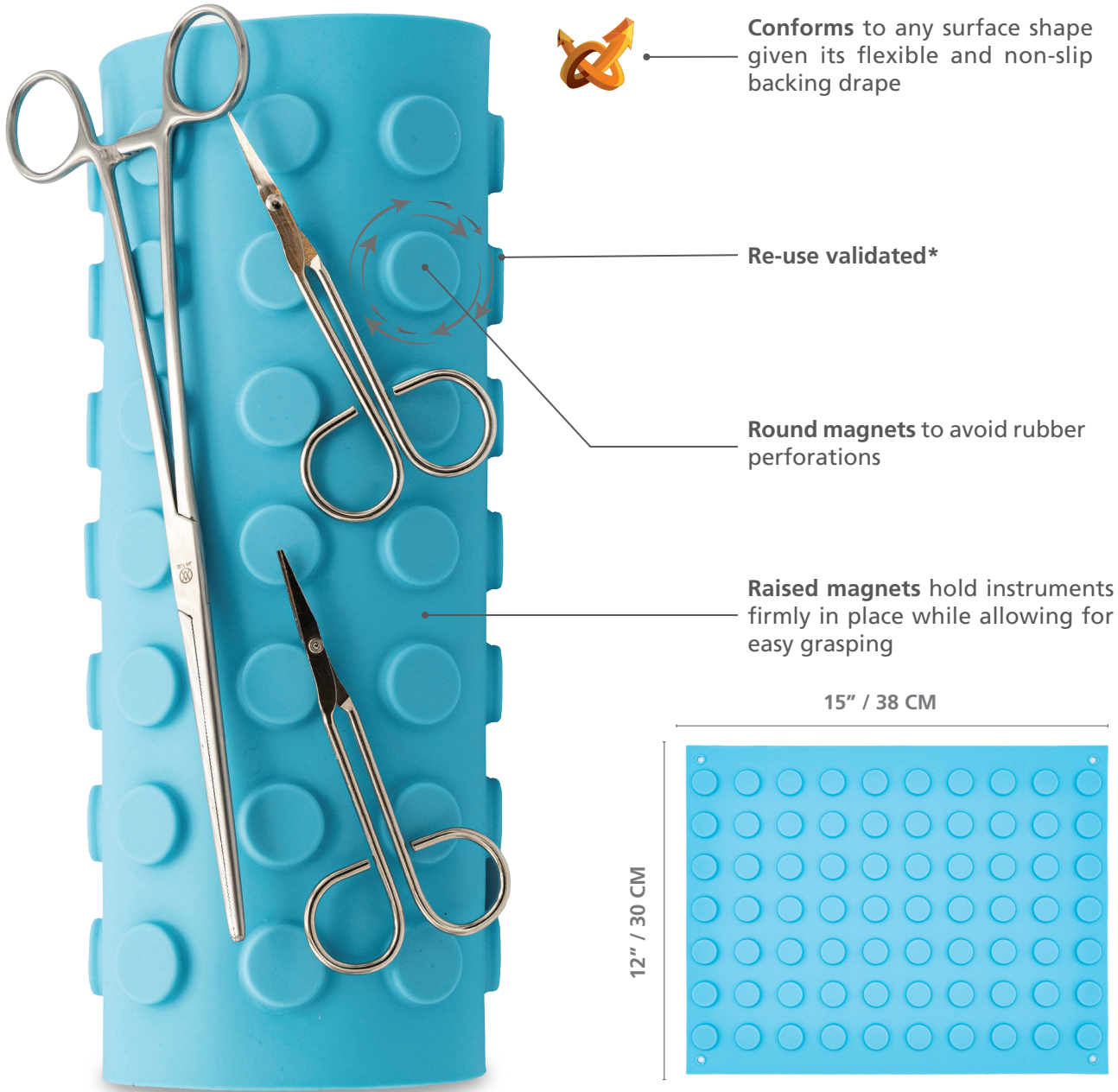

PRODUCT SELECTION


ITEM #	DESCRIPTION	DIMENSIONS	PACKAGING
MAG-02-1	Magnetic Mat	12" x 15" / 30 cm x 38 cm	1 Unit
MAG-02-5	Magnetic Mat (Volume Packaging)	12" x 15" / 30 cm x 38 cm	5 Units

*Independent laboratory tested to validate pre-vacuum steam sterilization efficacy for reprocessing of medical devices in a healthcare setting

BATRIK PIN GUARDS



- Safely cover the sharp ends of your K-Wire and Steinmann Pins
- Color coded by wire size which eliminates the risk of selecting of the wrong size
- A removable sleeve is provided to avoid injury and facilitate application (sterile only)
- Non-sterile model is available with or without plastic sleeves
- Packaging designed to reduce waste
- Latex free

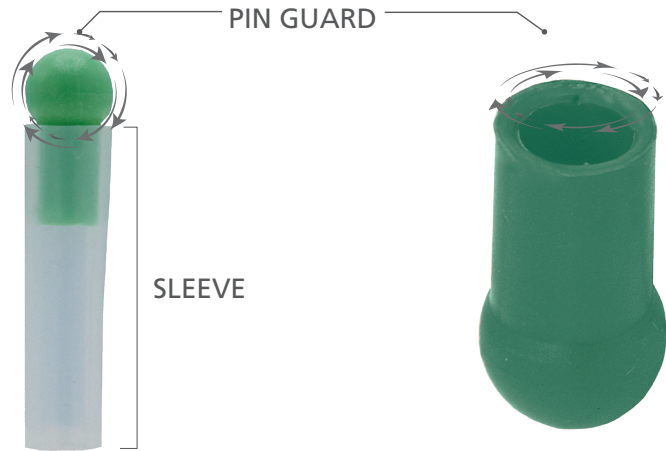

PRODUCT SELECTION

ITEM # STERILE PIN GUARDS + SLEEVE	ITEM # NON-STERILE PIN GUARDS + SLEEVE	ITEM # NON-STERILE PIN GUARDS (NO SLEEVE)	COLOR	DIMENSIONS (WIRE SIZE)
14-1028	914-1028 - SLV	914-1028	Yellow	0.028" / 0.71 mm
14-1035	914-1035 - SLV	914-1035	Light Blue	0.035" / 0.89 mm
14-1045	914-1045 - SLV	914-1045	White	0.045" / 1.14 mm
14-1054	914-1054 - SLV	914-1054	Dark Red	0.054" / 1.40 mm
14-1062	914-1062 - SLV	914-1062	Green	0.062" / 1.57 mm
14-1078	914-1078 - SLV	914-1078	Pink	0.078" / 1.98 mm
14-1093	914-1093 - SLV	914-1093	Grey	0.093" / 2.36 mm
14-1109	914-1109 - SLV	914-1109	Beige	0.109" / 2.77 mm
14-1125	914-1125 - SLV	914-1125	Purple	0.125" / 3.17 mm
14-1141	914-1141- SLV	914-1141	Brown	0.141" / 3.58 mm
14-1156	914-1156 - SLV	914-1156	Dark Blue	0.156" / 3.96 mm
14-1188	914-1188 - SLV	914-1188	Red	0.188" / 4.77 mm

STERILE PACKAGING: 80 units/box (2 units per sterile pouch - 40 pouches per box)

NON-STERILE PACKAGING: 100 units/box

BATRIK CAUTERY TIP CLEANER

Our Cautery Tip Cleaner assists in removing material from cautery tips for both mono-polar and bi-polar electrodes.

- Foam pad with integrated abrasive material
- Adheres to disposable and non-disposable drapes
- Radiopaque
- Available exclusively in Non-Sterile
- Latex free



PRODUCT SELECTION

ITEM #	DESCRIPTION	DIMENSIONS	PACKAGING
60-10000	Cautery Tip Cleaner Non-Sterile	1.8" x 1.8" / 4.5 cm x 4.5 cm	2,200 /Box

BATRIK SAFE-SOAK™ SPAC

Our Safe-Soak™ SPAC floor suction disk quickly removes fluid waste from the operating room floor through the application of suction.



- Highly durable: it won't be damaged if stepped on
- Easily maneuverable by foot around operating room floor
- Latex free, disposable, and non-sterile
- The bottom opening permits fluid collection through a highly absorbent material that comes into contact with its surface
- The hose is securely affixed on its side to reduce risk of catching equipment
- Its highly absorbent base material allows for quick fluid wicking and removal when suction is applied

BATRIK SAFE-SOAK™ ABSORBENT FLOOR MATS

Our Safe-Soak™ Absorbent Floor Mats are designed to absorb a variety of fluids from the operating room floor quickly and efficiently, ensuring the floor beneath remains slip-free, safe, and dry during surgical procedures.

- Highly absorbent
- Facilitates quick and easy clean up
- Latex free, single-use, disposable, and non-sterile
- Conforms to uneven surfaces
- Skid-resistant backing (blue absorbent mat only)



PRODUCT SELECTION

ITEM #	DESCRIPTION	ABSORBENCY	SIZE	PACKAGING
SS-20-1000	Safe-Soak™ / White extra absorbent mat	4.5 L / mat	28" x 48" (71 cm x 122 cm)	20 / Box
SS-20-2000	Safe-Soak™ / Blue absorbent mat with skid-resistant backer	2.5 L / mat	28" x 40" (71 cm x 102 cm)	20 / Box
SS-20-3000	Safe-Soak™ / Yellow absorbent mat	3 L / mat	30" x 52" (76 cm x 132 cm)	15 / Box
SS-21-1000	Safe-Soak™ SPAC	N/A	8.625" (22 cm diameter)	12 / Box

BATRIK BLUE CLOUD-GEL™ ANTI-FATIGUE MATS


Our Blue Cloud-Gel™ Anti-Fatigue Mats have been scientifically engineered for use in the O.R. using patented gel technology and premium foam to provide the ultimate comfort and support for surgeons, nurses, and technicians particularly when movement is restricted for prolonged durations.

- Proprietary 3/4" thickness Blue Cloud-Gel™ foam technology will not "bottom-out" like conventional foam mats
- Conforms to and supports feet promoting efficient circulation, reducing physical discomfort and fatigue
- Promotes proper posture while diminishing back, knee, and foot discomfort
- Latex free, water resistant, and nonabsorbent
- Reusable, easy to clean, and Cart Washer compatible (contact us for more details)
- Excellent friction properties when damp or dry surgical booties are worn
- Completely sealed with an RF weld for excellent protection against moisture absorption
- US Standards: CALL 117 and NFPA 260 Compliant
- 1 Year Warranty

DID YOU KNOW? STUDIES HAVE DEMONSTRATED THAT ANTI-FATIGUE MATS IMPROVE PERFORMANCE, ENERGY, AND POST-OPERATIVE DISCOMFORT.

PRODUCT SELECTION

ITEM #	DIMENSIONS	PACKAGING
800-2401	18" X 24" / 46 cm x 61 cm	1 / Box
800-2405	18" X 24" / 46 cm x 61 cm	5 / Box
800-3201	20" X 32" / 51 cm x 81 cm	1 / Box
800-3203	20" X 32" / 51 cm x 81 cm	3 / Box
800-4801	20" X 48" / 51 cm x 122 cm	1 / Box
800-4802	20" X 48" / 51 cm x 122 cm	2 / Box
800-7201	30" X 72" / 76 cm x 183 cm (Special Order)	1 / Box

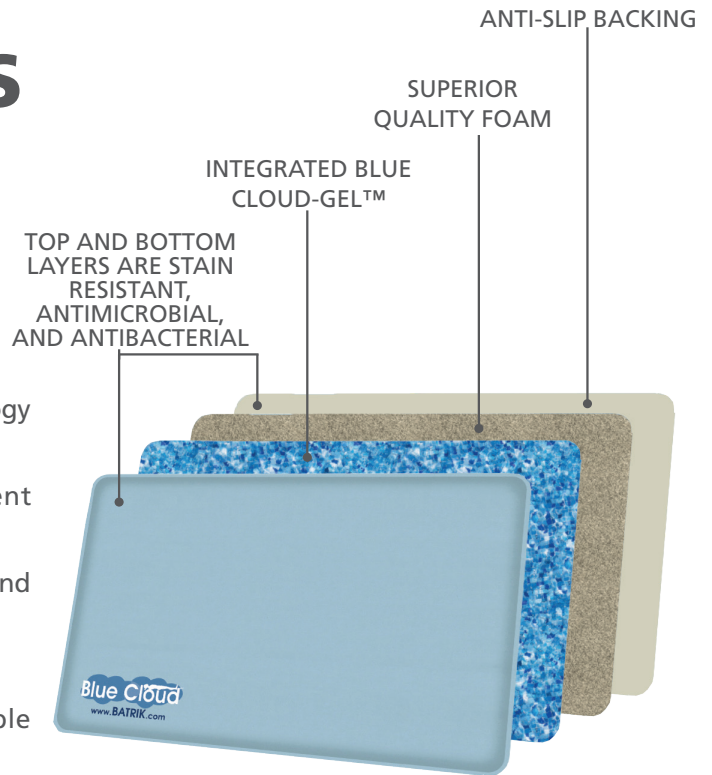


BATRIK BLUE-CLOUD NIMBUS™ ANTI-FATIGUE MATS

Our Blue-Cloud Nimbus™ Anti-Fatigue Mats are designed for general hospital use. Made with proprietary high-density “energy-return” polyurethane materials, they provide unmatched comfort, relief, and ergonomic support.

- Our Energy-Return properties will help reduce leg, foot, lower back discomfort, and work-related injuries
- Longer lasting support as 3/4” thickness proprietary foam blend will not bottom out like most other foam based anti-fatigue mats
- Antimicrobial properties blended into the mat provides added protection
- Highly water resistant surface will not absorb water or dirt
- Lightweight and easy to clean
- Textured design on surface and high-traction base ensures greater safety. (US Standard NFSI 101-C)
- Beveled edging reduces tripping hazards and anti-slip backing enhances safety
- Durable design is ideal for commercial applications, yet still ergonomically and aesthetically pleasing for home or office environments
- Flammability/Burn resistant per CFR 16 Parts 1630/1631
- Made with eco-friendly materials, is non-toxic, and latex free
- 5 Year Warranty (non-industrial usage), 4 Year Warranty (industrial usage)

GOOD TO KNOW: IF WITHIN 30 DAYS YOU ARE NOT 100% SATISFIED, CONTACT YOUR DEALER FOR A REFUND (EXCLUDES TRANSPORT, HANDLING, AND TAXES)

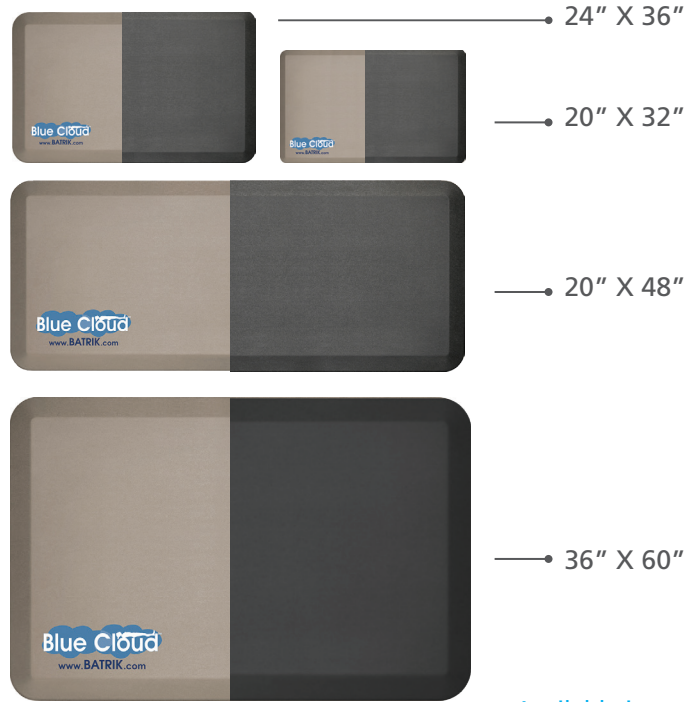

Available in
Taupe or Black

PRODUCT SELECTION

ITEM #	DIMENSIONS INCHES / CM	PACKAGING
BLACK COLOR		
809-3201	20" X 32" / 51 cm x 81 cm	1 / Box
809-3204	20" X 32" / 51 cm x 81 cm	4 / Box
809-4801	20" X 48" / 51 cm x 122 cm	1 / Box
809-4804	20" X 48" / 51 cm x 122 cm	4 / Box
809-3601	24" X 36" / 61 cm x 91 cm	1 / Box
809-7201	20" X 72" / 51 cm x 183 cm	1 / Box
809-6001	36" X 60" / 91 cm x 152 cm	1 / Box
TAUPE COLOR		
810-3201	20" X 32" / 51 cm x 81 cm	1 / Box
810-3204	20" X 32" / 51 cm x 81 cm	4 / Box
810-4801	20" X 48" / 51 cm x 122 cm	1 / Box
810-4804	20" X 48" / 51 cm x 122 cm	4 / Box
810-3601	24" X 36" / 61 cm x 91 cm	1 / Box
810-7201	20" X 72" / 51 cm x 183 cm	1 / Box
810-6001	36" X 60" / 91 cm x 152 cm	1 / Box

BATRIK

STACKABLE STEP STOOLS

Made of lightweight cast aluminium and protected with an antimicrobial powder coating, our patented Stackable Step Stools have been designed with infection prevention as a prime focus.



PRODUCT SELECTION

ITEM #	DESCRIPTION	COLOR	CASE DIMENSIONS
OR-36363	Stackable Step Stool with mat	Black	19" x 14" x 6"
OR-36363-02	Stackable Step Stool with mat	Dark Blue	19" x 14" x 6"
OR-36363-03	Stackable Step Stool with mat	Green	19" x 14" x 6"
OR-36363-07	Stackable Step Stool with mat	Orange	19" x 14" x 6"
OR-36364	Stackable Step Stool without mat	N/A	19" x 14" x 6"

PACKAGING: Sold by Unit

REPLACEMENT PARTS

ITEM #	DESCRIPTION	PACKAGING	CASE DIMENSIONS
OR-36363-1A	Replacement mat for OR-36363	1 Unit	15" x 18" x 5" (Ziplock)
OR-36363-1B	Replacement Foot Pegs	1 Unit	Ziplock Based on Quantity ordered

OPTIONAL ACCESSORIES

ITEM #	DESCRIPTION	PACKAGING	CASE DIMENSIONS
OR-3693	Transport cart	1 Unit	34" x 17" x 7"
809-1705	Reusable Anti-Fatigue Mat with adhesive strips to secure it in place	5 / Box	12" x 12" x 4"
899-1710	Single-Use Anti-Fatigue Mat with adhesive strips to secure it in place	10 / Box	12" x 12" x 8"

PRODUCT FEATURES

- Weight: 8.3 lbs (3.8 kg)
- Dimensions: 17.4" L x 13" W x 4.7" H (44.2 cm L x 33 cm W x 11.9 cm H)
- Transport Cart Available
- Weight Capacity: 300 lb (136 kg)
- Warranty: 7 years (step stool only)

ANTIMICROBIAL COATING:

Prevents the spread of hazardous microorganisms and suppresses the growth of mold, bacteria, yeast, and mildew

ANTI-SLIP SURFACE:

Provides enhanced stability.

RUBBER-FREE CORNERS:


Avoids fluid collection.

WASHABLE CART:

Model OR-36364 (without mat) can be sterilized.

HYGIENIC DESIGN:

No welds or large hollow openings to trap debris or fluids

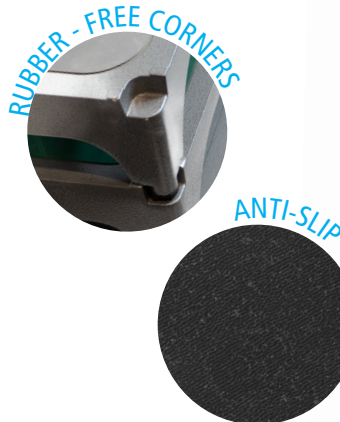


OR-3693

Ver. 2024-11-29

BATRIK TOURNY™ + TOURNIQUETS

Our Tourniquets are comfortable and easy to use on all patients, including the elderly, undernourished, children, and newborns.

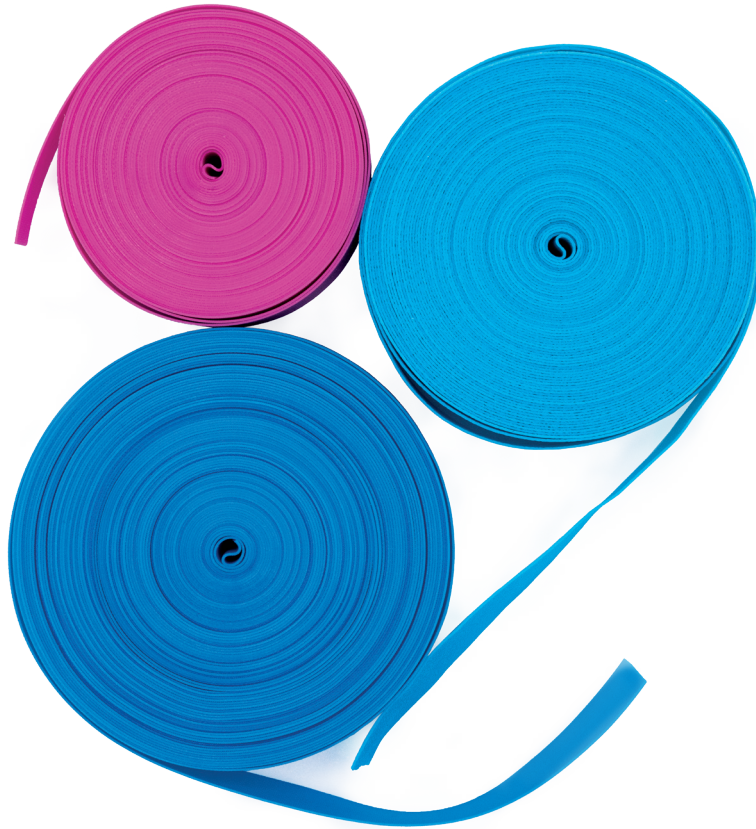

AVAILABLE COLORS



- Available in different sizes including compact sizes, which are convenient for phlebotomy trays
- Innovative latex free material has a more comfortable texture for patients, puts less stress on skin, and safeguards hospital staff and patients from latex allergens

SINGLE-USE BENEFITS

- Reduced risk of cross-contamination (patient to patient or patient to healthcare provider)
- Reduced risk of spread of nosocomial infections
- May improve your facility's compliance with OSHA's 1910-1030 (d)(1) "Blood-borne Pathogens Standard"

FOR YOUR CONVENIENCE:
25 UNITS PER ROLL
COME IN A DISPENSER BOX
FOR EASY ACCESS



PRODUCT SELECTION

ITEM #	DESCRIPTION	COLOR	DIMENSIONS	PACKAGING
K-002	TOURNY™+ Non-Sterile	Blue	1" x 18" / 2.5 cm x 45 cm	1 / Bag
K-004	TOURNY™+ Non-Sterile	Blue	1" x 18" / 2.5 cm x 45 cm	100 / Box
K-009	TOURNY™+ Non-Sterile	Blue	0.200" x 22" / 0.5 cm x 56 cm	20 / Box
K-010	TOURNY™+ Non-Sterile	White	0.200" x 22" / 0.5 cm x 56 cm	20 / Box
K-400	TOURNY™+ Non-Sterile	Blue	1" x 18" / 2.5 cm x 45 cm	16 rolls 400 / Box
K-400PB	TOURNY™+ Non-Sterile (Pebbled)	Light Blue	1" x 18" / 2.5 cm x 45 cm	16 rolls 400 / Box
K-100	TOURNY™+ Non-Sterile	Assorted (3 rolls of each) Blue, Orange, Green, and Pink	1" x 18" / 2.5 cm x 45 cm	12 rolls 300 / Box
K-254	TOURNY™+ Non-Sterile (Pebbled)	Light Blue	1" x 18" / 2.5 cm x 45 cm	250 / Box
K-2-18RL	TOURNY™+ Non-Sterile	Blue	1" x 18" / 2.5 cm x 45 cm	1 Roll
K-2-18PBRL	TOURNY™+ Non-Sterile (Pebbled)	Light Blue	1" x 18" / 2.5 cm x 45 cm	1 Roll
K-3-18RL	TOURNY™+ Non-Sterile	Green	1" x 18" / 2.5 cm x 45 cm	1 Roll
K-6-18RL	TOURNY™+ Non-Sterile	Pink	1" x 18" / 2.5 cm x 45 cm	1 Roll
K-7-18RL	TOURNY™+ Non-Sterile	Orange	1" x 18" / 2.5 cm x 45 cm	1 Roll
K-2-29RL	Long TOURNY™+ Non-Sterile	Blue	1" x 29.5" / 2.5 cm x 75 cm	1 Roll

SHOOZCOVERS™ | SOFT LINE™


Batrik ShoozCovers™ Soft Line™ is manufactured in Canada from high quality 35 g/m² polypropylene fabric and 100% latex free elastics for performance and durability you can rely on. Our anti-skid tread reduces the risk of slipping and our smart design delivers easy on-and-off convenience and great fit.

SHOOZ-REPELLENT™ TECHNOLOGY:

The Soft Line™ integrates Shooz-Repellent™ technology which creates a protective barrier that repels fluids, delivering basic incidental fluid contact repellent. Shoe covers are not waterproof or water-resistant.

ANTI-SKID TREAD:

Our anti-skid tread reduces the risk of slipping with triple ultrasonic welds that increases tear strength and delivers easy on-and-off convenience.

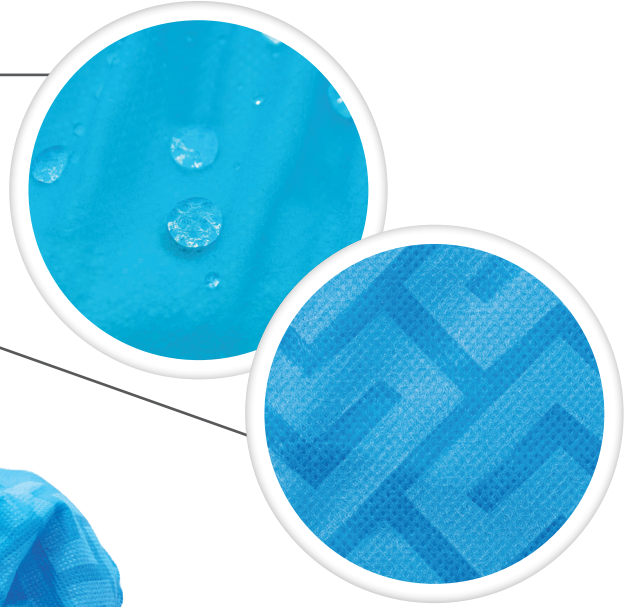

MADE IN CANADA

OUR SOFT LINE™ IS PERFECT FOR:



MEDICAL ESTABLISHMENTS



PHARMACEUTICAL LABS



GMP SITES




ELECTRONIC MANUFACTURING

BUBBLE GUM PINK | **SOFT LINE™**



WHY YOU'LL LOVE THE SOFT LINE™:

- Both shoe covers are manufactured in Canada with top-quality 35 g/m² polypropylene fabric and 100% latex-free USA-made elastics¹.
- This line is particle count tested and complies with IEST-RP-CC003.4.
- Other colors available upon request CC003.4.

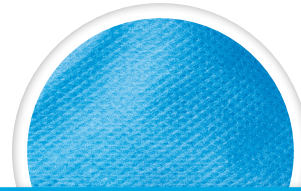

UN XL

AVAILABLE IN TWO VERSATILE SIZES

BUBBLE GUM PINK | **SOFT LINE™**



SUPPORT A CAUSE, WEAR PINK:

- A portion of our Bubble Gum Pink shoe cover sales will be donated to charity to help further breast cancer research. Join us in making a meaningful contribution by choosing Pink ShoozCovers™!



UN

AVAILABLE IN UNIVERSAL

PRODUCT SELECTION:

ITEM #	DESCRIPTION	COLOR	SIZE	PACKAGING	CASE DIMENSIONS
SC-1210-PB	ShoozCovers™ - Soft Line™ <i>polypropylene non-woven fabric</i>	Lolly Blue	Universal Up to shoe size 12 ~ 42 cm x 16.5 cm / 16.5" x 6.5"	Case of 1,000 (10 bags of 100 units)	20" x 16" x 16"
SC-1216-PB	ShoozCovers™ - Soft Line™ <i>polypropylene non-woven fabric</i>	Lolly Blue	Universal Shoe size 12 ~ 42 cm x 16.5 cm / 16.5" x 6.5"	Case of 600 (6 bags Dispenser Boxes of 100 units)	16.25" x 8.875" x 18.75"
SC-1216-PB	ShoozCovers™ - Soft Line™ <i>polypropylene non-woven fabric</i>	Lolly Blue	Universal Shoe size 12 ~ 42 cm x 16.5 cm / 16.5" x 6.5"	Case of 600 (6 bags of 100 units)	19.5" x 19.5" x 14.25"
SC-1226-PDB	ShoozCovers™ - Soft Line™ <i>polypropylene non-woven fabric</i>	Lolly Blue	XL Shoe size 12-15 ~ 47 cm x 19 cm / 18.5" x 7.5"	Case of 600 (6 bags Dispenser Boxes of 100 units)	17.25" x 9.375" x 21.75"
SC-1220-PB	ShoozCovers™ - Soft Line™ <i>polypropylene non-woven fabric</i>	Lolly Blue	XL Shoe size 12-15 ~ 47 cm x 19 cm / 18.5" x 7.5"	Case of 1,000 (10 bags of 100 units)	22" x 16" x 16"
SC-1616-PDB	ShoozCovers™ - Soft Line™ <i>polypropylene non-woven fabric</i>	Bubble Gum Pink	Universal Shoe size 12-15 ~ 42 cm x 16.5 cm / 16.5" x 6.5"	Case of 600 (6 bags Dispenser Boxes of 100 units)	16.25" x 8.875" x 18.75"
SC-1616-PB	ShoozCovers™ - Soft Line™ <i>polypropylene non-woven fabric</i>	Bubble Gum Pink	Universal Shoe size 12-15 ~ 42 cm x 16.6 cm / 16.5" x 6.5"	Case of 600 (6 bags of 100 units)	19.5" x 19.5" x 14.25"

ALTERNATIVE PACKAGING: Add "-PDB" to end of code for a dispenser box instead of

¹Many suppliers claim their products are latex free, however this often only relates to the fabric and not the elastics used. We source our ShoozCovers™ elastics from the USA to ensure that our components conform fully as a 100% latex free product.

GRIPPEEZ SHOOZ™


*RELIABLE, DURABLE, AND TRULY ANTI-SLIP SHOE COVERS
DESIGNED FOR THE MOST DEMANDING ENVIRONMENTS*

ShoozCovers™ Grippeez provide exceptional protection and traction for professionals working in demanding environments. Featuring waterproof protection, our shoe covers ensure safety and hygiene without compromising on quality. They are the ideal solution for tough jobs.

MADE IN CANADA 

WHY CHOOSE GRIPPEEZ

- **Grip You Can Count On:** High-performance rubberized coating ensures superior traction on wet, oily, and uneven surfaces.
- Shooz-Repellent™ technology delivers waterproof protection and basic incidental fluid contact repellency.
- High quality material made in the **USA latex-free elastics.**

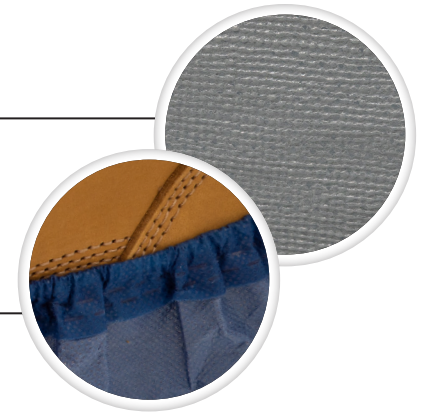

PERFECT FOR:



AEROSPACE



FOOD
PROCESSING
INDUSTRY



OIL, GAS AND
MINING SITES



CONSTRUCTION &
RENOVATION SITES



GRIPPEEZ SHOOZ™ ARE A PERFECT FIT FOR EVERYONE:

- **Fast, Reliable Fit:** Stretch material ensures a secure, easy-on fit. Coating provides protection against slippery surfaces.
- Comes in an **extra-large size** that fits over a men's size 15 steel-toed boot.¹
- **Enhanced Safety:** Anti-slip rubberized coating provides reliable grip on all surfaces, helping to prevent sliding in critical areas.
- Available in two colors: Gecko Navy and Rhino Grey.

RHINO GREY



GECKO NAVY



PRODUCT SELECTION:

ITEM #	DESCRIPTION	COLOR	SIZE	PACKAGING	CASE DIMENSIONS
SC-29312-NPDB	Grippez ShoozCovers™ Polypropylene + polyethylene non-woven fabric, waterproof	Rhino Grey	Universal Up to shoe size 12 ~ 42 cm x 19 cm / 16.5" x 6.88"	Case of 300 (2 Dispenser boxes of 150 Units)	16.62"x 14.875"x 8.75"
SC-29312-NPB	Grippez ShoozCovers™ Polypropylene + polyethylene non-woven fabric, waterproof	Rhino Grey	Universal Shoe size 12-15 ~ 42 cm x 19 cm / 16.5" x 6.88"	Case of 300 (2 bags of 150 Units)	16"x 13.375"x 10"
SC-29322-NPDB	Grippez ShoozCovers™ Polypropylene + polyethylene non-woven fabric, waterproof	Rhino Grey	XL Up to shoe size 12 ~ 47 cm x 19 cm / 16.5" x 6.88"	Case of 240 (2 Dispenser boxes of 120 Units)	16.62"x 14.875"x 8.75"
SC-29322-NPB	Grippez ShoozCovers™ Polypropylene + polyethylene non-woven fabric, waterproof	Rhino Grey	XL Up to shoe size 12 ~ 47 cm x 19 cm / 16.5" x 6.88"	Case of 240 (2 bags of 120 units)	16"x 13.375"x 10"
SC-22212-NPDB	Grippez ShoozCovers™ Polypropylene + polyethylene non-woven fabric, waterproof	Gecko Navy	Universal Shoe size 12-15 ~ 42 cm x 19 cm / 16.5" x 6.88"	Case of 300 (2 Dispenser boxes of 150 units)	16.62"x 14.875"x 8.75"
SC-22212-NPB	Grippez ShoozCovers™ Polypropylene + polyethylene non-woven fabric, waterproof	Gecko Navy	Universal Up to shoe size 12 ~ 42 cm x 19 cm / 16.5" x 6.88"	Case of 300 (2 bags of 150 units)	16"x 13.375"x 10"
SC-22212-NPDB	Grippez ShoozCovers™ Polypropylene + polyethylene non-woven fabric, waterproof	Gecko Navy	XL Up to shoe size 12 ~ 47 cm x 19 cm / 16.5" x 6.88"	Case of 240 (2 Dispenser boxes of 120 units)	16.62"x 14.875"x 8.75"
SC-22212-NPB	Grippez ShoozCovers™ Polypropylene + polyethylene non-woven fabric, waterproof	Gecko Navy	XL Up to shoe size 12 ~ 47 cm x 19 cm / 16.5" x 6.88"	Case of 240 (2 bags of 120 units)	16"x 13.375" x 10"

ALTERNATIVE PACKAGING: Add "-PDB" to end of code for a dispenser box instead of bag

SHOOZCOVERS™ | COATED LINE™


STAY DRY:

Our Batrik ShoozCovers™ Coated Line™ is made from a ~55 g/m² blend of polypropylene and polyethylene, providing high-quality water resistance.¹ Ideal for environments where maximum protection against fluids is required.



TRIPLE WELDED TECHNOLOGY:

These shoe covers are reinforced with triple ultrasonic welds, a technique that joins materials together without the need for adhesives or stitching. This technique increase tear strength and reduce linting, for long lasting, resilient covers.

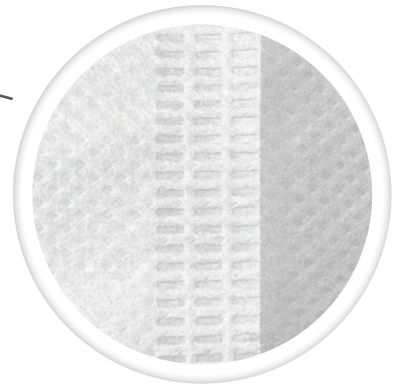

MADE IN CANADA

OUR COATED LINE™ IS PERFECT FOR:



MEDICAL ESTABLISHMENTS



PHARMACEUTICAL LABS



GMP SITES



ELECTRONIC MANUFACTURING

¹OUR SHOOZCOVERS™ COATED LINE™ IS RESISTANT AGAINST ALL NON PETROLEUM BASED LIQUIDS.

MARSHMALLOW WHITE | COATED LINE™

WHY YOU'LL LOVE THE COATED LINE™:



- **Superior protection using our unique blend of fabrics** which protects against not only water but all non-petroleum-based liquids, making them the ideal choice for environments where spills and leaks are a concern
- **Elastics around the ankle create a strong and secure fit** which is especially important when working in hazardous environments where a misplaced shoe cover could lead to contamination or injury
- **Smart design** with a large opening, providing a better fit and easy on-and-off convenience
- **Other colors available upon request**



UN XL

AVAILABLE IN TWO VERSATILE SIZES

PRODUCT SELECTION:

ITEM #	DESCRIPTION	COLOR	SIZE	PACKAGING	CASE DIMENSIONS
SC-2112-NPB	ShoozCovers™ - Coated Line™ <i>Polypropylene + polyethylene non-woven fabric, waterproof</i>	Marshmallow White	Universal Up to shoe size 12 ~ 42 cm x 17.5 cm / 16.5" x 6.88"	Case of 300 (2 bags of 150 units)	16" x 13.375" x 10"
SC-2112-NPDB	ShoozCovers™ - Coated Line™ <i>Polypropylene + polyethylene non-woven fabric, waterproof</i>	Marshmallow White	Universal Up to shoe size 12 ~ 42 cm x 17.5 cm / 16.5" x 6.88"	Case of 300 (2 boxes of 150 units)	16.62" x 14.875" x 8.75"
SC-2122-NPB	ShoozCovers™ - Coated Line™ <i>Polypropylene + polyethylene non-woven fabric, waterproof</i>	Marshmallow White	XL Shoe size 12-15 ~ 47 cm x 19 cm / 18.5" x 6.88"	Case of 240 (2 bags of 120 units)	16" x 13.375" x 10"
SC-2122-NPDB	ShoozCovers™ - Coated Line™ <i>Polypropylene + polyethylene non-woven fabric, waterproof</i>	Marshmallow White	XL Shoe size 12-15 ~ 47 cm x 19 cm / 18.5" x 6.88"	Case of 240 (2 boxes of 120 units)	16.62" x 14.875" x 8.75"

ALTERNATIVE PACKAGING: Add "-PDB" to end of code for a dispenser box instead of bags

SHOOZCOVERS™ | TOUGH LINE™

EXTRA-LARGE FOR EXTRA-TOUGH WORK:


Our Batrik ShoozCovers™ Tough Line™ comes in an extra-large size that fits over a men's size 15 steel-toed boot, delivering optimal full coverage and protection while you work.¹

ENHANCED DURABILITY:

These shoe covers are formulated with our **thickest fabric yet!** Made of ~55 g/m² polypropylene fabric, these highly durable covers provide added protection, increase tear strength and reduce linting.



SHOE SIZES 12-15 (45-49.5 EU)
~ 47 CM X 19 CM / 18.5" X 7.5"

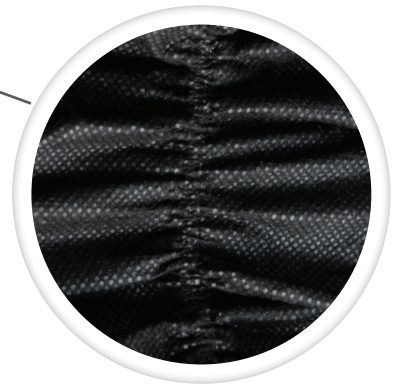

MADE IN CANADA



OUR TOUGH LINE™ IS PERFECT FOR:



HOMES OR
COMMERCIAL
BUILDINGS



SHIELDING
AGAINST MUD
& GREASE



MEDICAL
ESTABLISHMENTS



PHARMACEUTICAL
LABS

¹ SHOOZCOVERS™ ARE MADE LONGER (42 & 47CM) FOR ACCURATE SHOE SIZE SELECTION.

BLACK LIQUORICE | TOUGH LINE™

WHY YOU'LL LOVE THE TOUGH LINE™:

- Triple ultrasonic welds that are long lasting to reduce waste
- High quality made in the USA latex-free elastics are triple reinforced to contract the shoe cover around your foot and deliver a strong and secure fit²
- Shooz-Repellent™ technology delivers basic incidental fluid contact repellency³
- Other colors available upon request



XL

AVAILABLE IN EXTRA LARGE

ANY TIME, ANYWHERE:



RESIDENTIAL



COMMERCIAL



INDUSTRIAL

PRODUCT SELECTION:

ITEM #	DESCRIPTION	COLOR	SIZE	PACKAGING	CASE DIMENSIONS
SC-1922-PB	ShoozCovers™ - Tough Line™ Polypropylene + polyethylene non-woven fabric	Black Liquorice	XL Shoe size 12-15 ~ 47 cm x 19 cm / 18.5" x 7.5"	Case of 200 (2 bags of 100 units)	16" x 13.375" x 10"

ALTERNATIVE PACKAGING: Add "-PDB" to end of code for a dispenser box instead of bags

² Many suppliers claim their products are latex free, however this often only relates to the fabric and not the elastics used. We source our ShoozCovers™ elastics from the USA to ensure that our components conform fully as a 100% latex free product. ³ Shooz-Repellent™ technology provides basic incidental fluid contact repellent. It is not waterproof nor water resistant.


BATRIK ENDOSCOPY BO2™ BITE BLOCKS WITH NASAL OXYGEN DELIVERY

Our BO2™ Bite Blocks significantly enhance patient safety by providing supplemental oxygenation during gastrointestinal endoscopy, preventing hypoxia in sedated patients, and ensuring superior patient comfort.

- Saves time and money by eliminating the need for separate nasal prongs and oxygen masks
- Designed to be free of sharp edges to ensure patient comfort and safety
- Available exclusively in Single-Use
- Optional retention strap is available to ensure a secure fit
- The Bite Blocks and optional retention straps are all latex free, lead free and DEHP free

BO2
(FRONT SIDE)

BO2
(BACK SIDE)

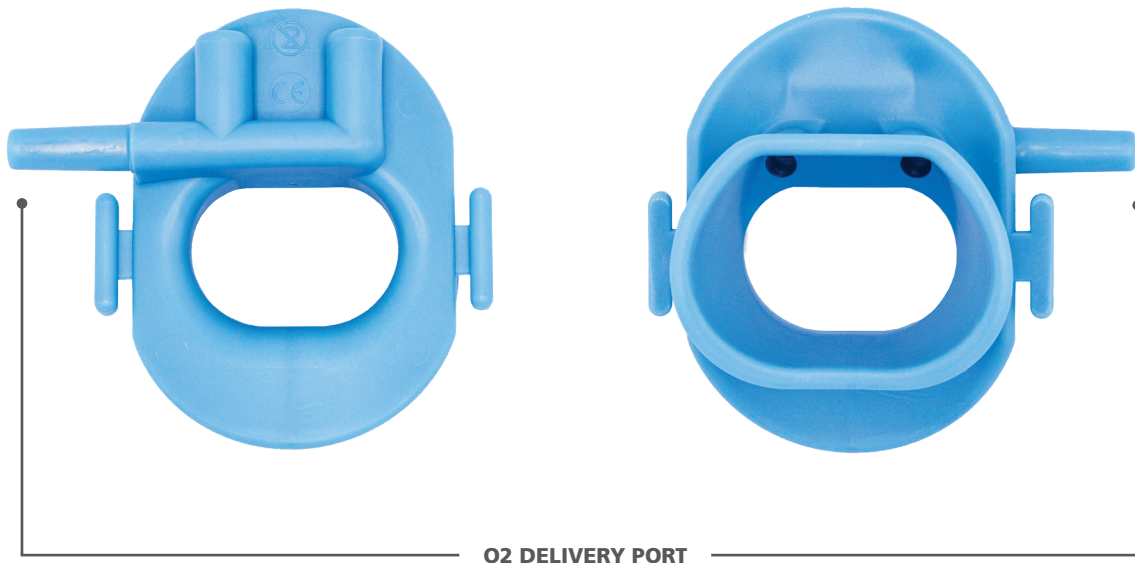

PRODUCT SELECTION

ITEM #	TYPE	DESCRIPTION	PACKAGING	CASE DIMENSION
BO2	Single-Use	18mm / 54FR, includes strap	100 / Box	10" x 10" x 10"
BO2-N	Single-Use	18mm / 54FR, without strap	100 / Box	10" x 10" x 10"

BATRIK ENDOSCOPY BITESTRONG™ BITE BLOCKS

ADULT & PEDIATRIC BITE BLOCKS

Our Bitestrong™ Bite Blocks provide a safe passage through the patient's mouth during TEE, protecting the patient's teeth and endoscope.

Our molding process ensures a smooth finish with no sharp edges (flashing) that can often irritate a patient's mouth. All Bite Blocks are individually wrapped in sealed, easy to open PE bags and are latex free, lead free, and DEHP free. Additional straps are available, if required

SINGLE USE

- Available without dental retention

REUSABLE

- Designed for repeated usage and will withstand commonly practiced chemical and steam sterilization cycles (134°C / 275°F).
- Available with or without an adjustable retention strap

PRODUCT SELECTION

ITEM #	TYPE	DESCRIPTION	PACKAGING	CASE DIMENSIONS
B1	Single-Use Pediatric	16 mm / 48FR, includes strap	100 / Box	10" x 10" x 10"
B1-25	Single-Use Pediatric	16 mm / 48FR, includes strap	25 / Box	10" x 7" x 7"
B1-N	Single-Use Pediatric	16 mm / 48FR, without strap	100 / Box	10" x 10" x 10"
B3	Single-Use	20 mm / 60FR, includes strap	100 / Box	10" x 10" x 10"
B3-N	Single-Use	20 mm / 60FR, without strap	100 / Box	10" x 10" x 10"
B3-SA	Single-Use	20 mm / 60FR, with strap <u>attached</u> to one side	100 / Box	10" x 10" x 10"
B3-NDR	Single-Use	20 mm / 60FR, no dental retention rim	100 / Box	10" x 10" x 10"
B3-NDR-N	Single-Use	Same as B3-N, no dental retention rim, without strap.	100 / Box	10" x 10" x 10"
B4	Single-Use	19 mm / 60FR	100 / Box	10" x 10" x 10"
B4T	Single-Use	19 mm / 60FR, includes strap	100 / Box	10" x 10" x 10"
B4T-N	Single-Use	19 mm / 60FR, without strap	100 / Box	10" x 10" x 10"
B6	Reusable	19 mm / 60FR	100 / Box	10" x 10" x 10"
B6T	Reusable	19 mm / 60FR, includes strap	100 / Box	10" x 10" x 10"
B6T-N	Reusable	19 mm / 60FR, without strap	100 / Box	10" x 10" x 10"
B6-25	Reusable	19 mm / 60FR	25 / Box	10" x 7" x 7"
B6T-25	Reusable	19 mm / 60FR, includes strap	25 / Box	10" x 7" x 7"
B6T-N-25	Reusable	19 mm / 60FR, without strap	25 / Box	10" x 7" x 7"
BS	Single-Use	Bite-Block straps 21.7" x 0.79" / 55 cm x 2 cm	100 / Box	10" x 13" x 1.5"



BS



B1



B3



B4



B6

BATRIK COUVRETTE™ VALVE COVERS

Our Couvrette™ Valve Covers eliminate the need for manual cleaning and reprocessing of reusable valves. Our single use Valve Covers effectively reduce cross contamination between patients and also reduce fluid leakage from the parts during procedures.

- Easily and securely fit to all valves
- Reduce the potential for errors
- Single-Use
- Non-Sterile
- Latex free



TIP: REUSABLE VALVE COVERS ARE NOTORIOUSLY DIFFICULT TO CLEAN.

SAVE TIME AND REDUCE RISK WITH COUVRETTE™.

COMING SOON!



PRODUCT SELECTION

ITEM #	DESCRIPTION	COLOR	PACKAGING
16-0002	Biopsy Valve Cover - Olympus®	Blue	100 / Box
16-0004	Biopsy Valve Cover - Pentax®	Red	100 / Box
16-0022	Biopsy Valve Cover - Olympus® & Fujinon®	Light Blue	100 / Box COMING SOON!
16-0044	Biopsy Valve Cover - Pentax®	Pink	100 / Box COMING SOON!


NOTE: Contact us to verify compatibility with your scope

BATRIK DAYBRUSH™ CHANNEL BRUSHES

REUSABLE

Our Reusable DayBrush™ Channel Brushes can be used for cleaning endoscopes for an entire day. In fact, it's been proven to effectively clean up to 10 scopes before the need to discard.*

- Soft tip ensures easy manipulation within the channel with no risk of damage
- The back-end is sealed and softened to prevent fluid from getting trapped inside the tube. The softening process also prevents damage if the cleaning brush is used in reverse.
- Can be steam sterilized.
- Bristle volume is 20mm in length – nearly double most brushes on the market. A large brush provides more surface contact during the cleaning process.
- Significant per use savings when compared to other reusable or disposable channel cleaning brushes.

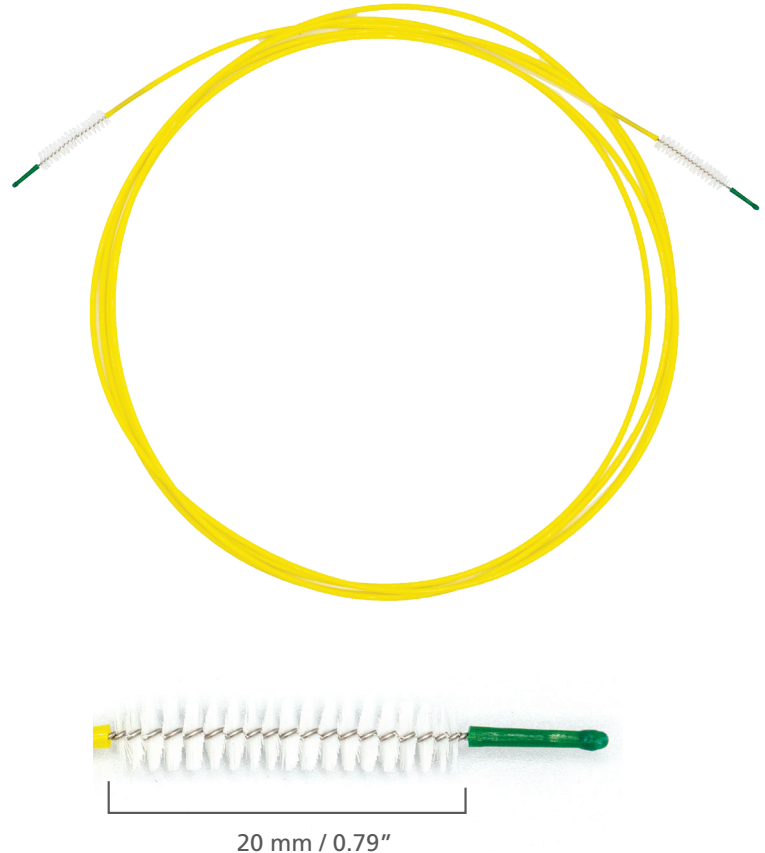

() Discard after cleaning 10 scopes or sooner if any brush degradation is noticed. After ten applications, the brush may have the appearance that it can continue to clean effectively, but we discourage prolonged use due to risk of cross contamination. Actual number of uses will vary based on handling and processing.*

PRODUCT SELECTION

5 MM CYLINDRICAL BRUSHES – SHEATH Ø 1.75 MM

ITEM #	COLOUR	LENGTH	BRUSH LENGTH	BRUSH (Ø)	SHEATH (Ø)	BRUSH ENDS
30-18050	White	70" / 180 cm	0.79" / 20 mm	0.20" / 5 mm	0.07" / 1.75 mm	Single
30-18050-T	Yellow	70" / 180 cm	0.79" / 20 mm	0.20" / 5 mm	0.07" / 1.75 mm	Single
30-18050-2T	Yellow	70" / 180 cm	0.79" / 20 mm	0.20" / 5 mm	0.07" / 1.75 mm	Double
30-23050	White	90" / 230 cm	0.79" / 20 mm	0.20" / 5 mm	0.07" / 1.75 mm	Single
30-23050-T	Yellow	90" / 230 cm	0.79" / 20 mm	0.20" / 5 mm	0.07" / 1.75 mm	Single
30-23050-2T	Yellow	90" / 230 cm	0.79" / 20 mm	0.20" / 5 mm	0.07" / 1.75 mm	Double

PACKAGING : 50 / BOX

OPTIONAL KIT ADD-ONS

ADD THE SUFFIX:

"F" to receive a kit with our Flat Sponge (Packaging: Box of 40). Applicable on certain models only, call for details.


"K" to receive a kit with the single use Valve / Control Cleaning Brush 5/11 Double Brushes (item # 27-57787). Applicable on all models.

"FK" to receive both the sponge and valve cleaning brushes.

BATRIK ARC™ CHANNEL BRUSHES SINGLE-USE

BATRIK offers one of the largest selection of Endoscopy Channel Cleaning Brushes in the industry. Our brushes are available in a variety of lengths, materials, and styles to accommodate the widest range of endoscopic devices used in the healthcare industry today.

No matter what your requirements, you'll be able to find the brush that suits your specific needs. If you don't see the brush you require, simply call us. We have over 50 different brushes available and can custom fabricate to your precise specifications.

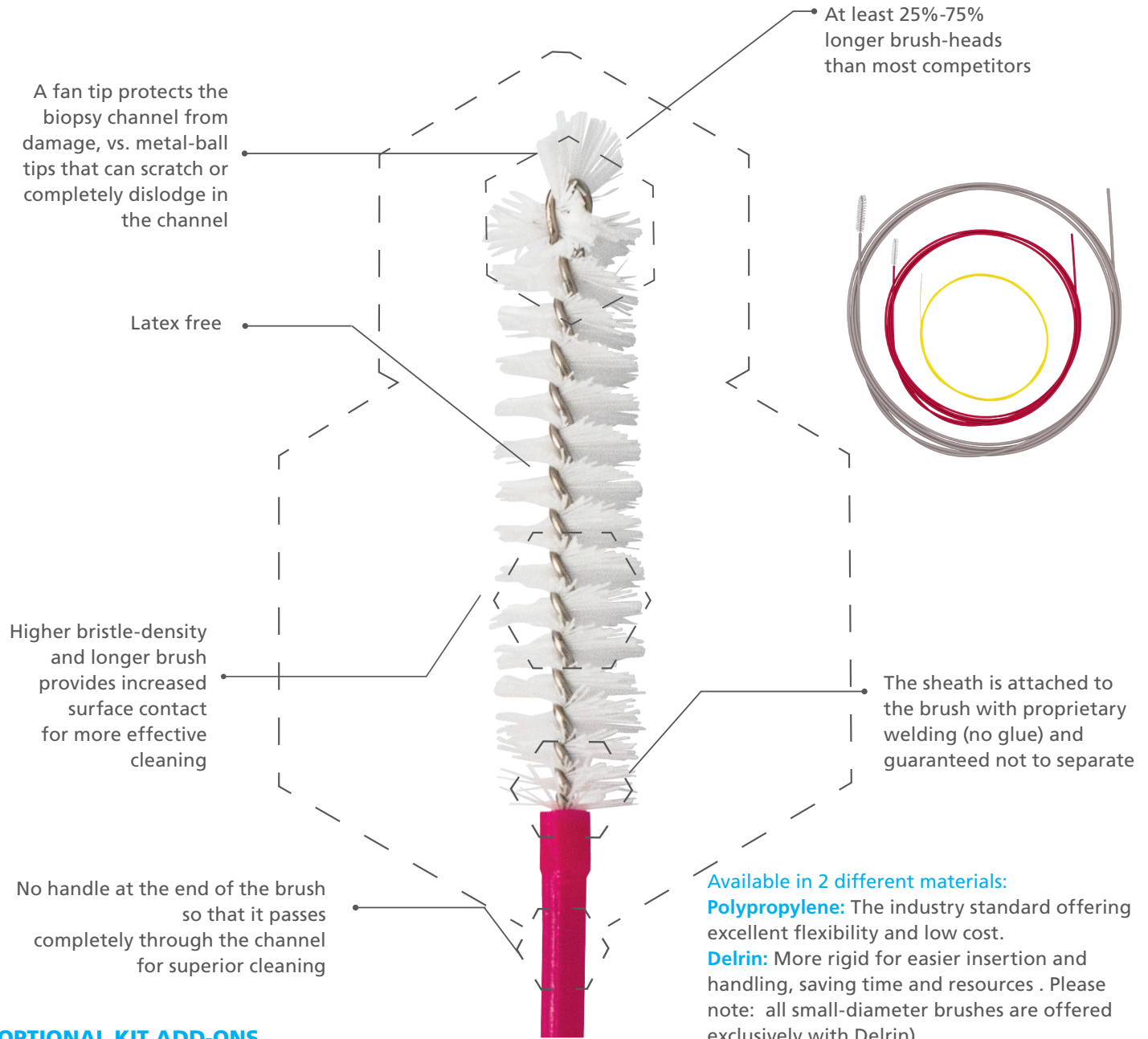

OPTIONAL KIT ADD-ONS


ADD THE SUFFIX:

"F" to receive a kit with our Flat Sponge (Packaging: Box of 40). Applicable on certain models only, call for details.

"K" to receive a kit with the single use Valve / Control Cleaning Brush 5/11 Double Brushes (item # 27-57787). Applicable on all models.

"FK" to receive both the sponge and valve cleaning brushes.

SERIES 311 & 333

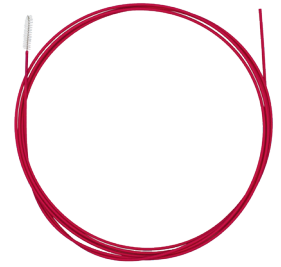

PRODUCT SELECTION

5 MM CYLINDRICAL BRUSHES – SHEATH Ø 1.85 MM - PUNCHY PINK POLYPROPYLENE

ITEM #	LENGTH	BRUSH LENGTH	BRUSH (Ø)	SHEATH (Ø)	BRUSH ENDS
311-10050	39" / 100 cm	0.94" / 24 mm	0.20" / 5 mm	0.073" / 1.85 mm	Single
311-18050	70" / 180 cm	0.94" / 24 mm	0.20" / 5 mm	0.073" / 1.85 mm	Single
311-23050	90" / 230 cm	0.94" / 24 mm	0.20" / 5 mm	0.073" / 1.85 mm	Single
311-23050-2	90" / 230 cm	0.94" / 24 mm	0.20" / 5 mm	0.073" / 1.85 mm	Double


PACKAGING: 50 / BOX

5/6 MM TAPERED BRUSHES – SHEATH Ø 1.85 MM - PUNCHY PINK POLYPROPYLENE

ITEM #	LENGTH	BRUSH LENGTH	BRUSH (Ø)	SHEATH (Ø)	BRUSH ENDS
333-23056	90" / 230 cm	0.94" / 24 mm	0.20", 0.24" / 5, 6 mm	0.073" / 1.85 mm	Single
333-23056-2	90" / 230 cm	0.94" / 24 mm	0.20", 0.24" / 5, 6 mm	0.073" / 1.85 mm	Double

PACKAGING: 50 / Box

SERIES 300 & 322

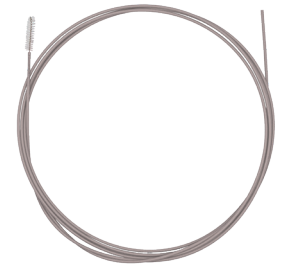

PRODUCT SELECTION

5 MM CYLINDRICAL BRUSHES – SHEATH Ø 1.75 MM - METALLIC GREY DELRIN

ITEM #	LENGTH	BRUSH LENGTH	BRUSH (Ø)	SHEATH (Ø)	BRUSH ENDS
300-18050	70" / 180 cm	0.79" / 20 mm	0.20" / 5 mm	0.07" / 1.75 mm	Single
300-23050	90" / 230 cm	0.79" / 20 mm	0.20" / 5 mm	0.07" / 1.75 mm	Single
300-23050-2	90" / 230 cm	0.79" / 20 mm	0.20" / 5 mm	0.07" / 1.75 mm	Double


PACKAGING: 50 / Box

5/6 MM TAPERED BRUSHES – SHEATH Ø 1.75 MM - METALLIC GREY DELRIN

ITEM #	LENGTH	BRUSH LENGTH	BRUSH (Ø)	SHEATH (Ø)	BRUSH ENDS
322-23056	90" / 230 cm	0.79" / 20 mm	0.20", 0.24" / 5, 6 mm	0.07" / 1.75 mm	Single
322-23056-2	90" / 230 cm	0.79" / 20 mm	0.20", 0.24" / 5, 6 mm	0.07" / 1.75 mm	Double

PACKAGING: 50 / Box

SERIES 344

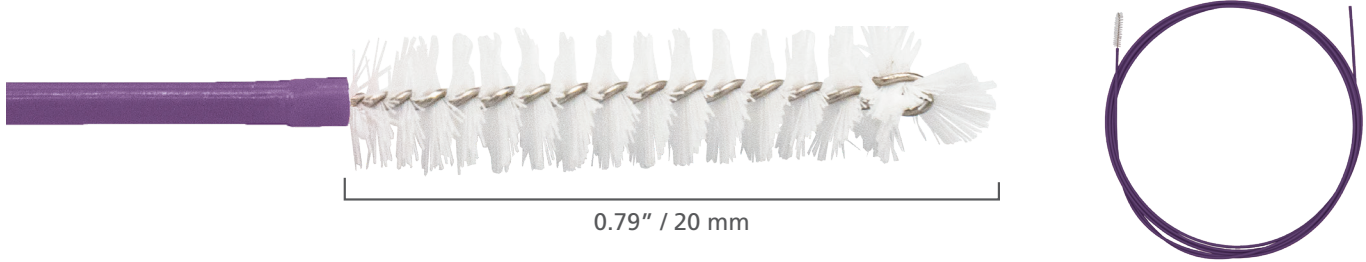


PRODUCT SELECTION

3 MM CYLINDRICAL BRUSHES – SHEATH Ø 1.55 MM - METALLIC PURPLE DELRIN

ITEM #	LENGTH	BRUSH LENGTH	BRUSH (Ø)	SHEATH (Ø)	BRUSH ENDS
344-10030	39" / 100 cm	0.79" / 20 mm	0.12" / 3 mm	0.06" / 1.55 mm	Single
344-18030	70" / 180 cm	0.79" / 20 mm	0.12" / 3 mm	0.06" / 1.55 mm	Single
344-23030	90" / 230 cm	0.79" / 20 mm	0.12" / 3 mm	0.06" / 1.55 mm	Single

PACKAGING: 50 / Box

SERIES 35

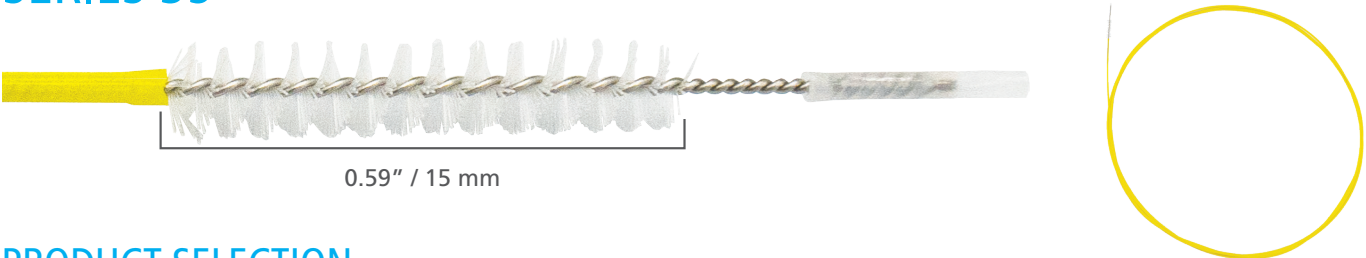

PRODUCT SELECTION

2.5 MM CYLINDRICAL BRUSHES – SHEATH Ø 1.00 MM - YELLOW DELRIN

ITEM #	LENGTH	BRUSH LENGTH	BRUSH (Ø)	SHEATH (Ø)	BRUSH ENDS
35-12025	47.2" / 120 cm	0.59" / 15 mm	0.098" / 2.5 mm	0.039" / 1.00 mm	Single
35-23025	90" / 230 cm	0.59" / 15 mm	0.098" / 2.5 mm	0.039" / 1.00 mm	Single
35-23025-2	90" / 230 cm	0.59" / 15 mm	0.098" / 2.5 mm	0.039" / 1.00 mm	Double

PACKAGING: 50 / Box

2 MM CYLINDRICAL BRUSHES – SHEATH Ø 1.00 MM - YELLOW DELRIN

ITEM #	LENGTH	BRUSH LENGTH	BRUSH (Ø)	SHEATH (Ø)	BRUSH ENDS
35-12020	47.2" / 120 cm	0.59" / 15 mm	0.079" / 2.0 mm	0.039" / 1.00 mm	Single
35-18020	70.8" / 180 cm	0.59" / 15 mm	0.079" / 2.0 mm	0.039" / 1.00 mm	Single
35-23020	90" / 230 cm	0.59" / 15 mm	0.079" / 2.0 mm	0.039" / 1.00 mm	Single
35-23020-2	90" / 230 cm	0.59" / 15 mm	0.079" / 2.0 mm	0.039" / 1.00 mm	Double


PACKAGING: 50 / Box

BATRIK VALVE / CONTROL BRUSHES


27-05197



27-57787

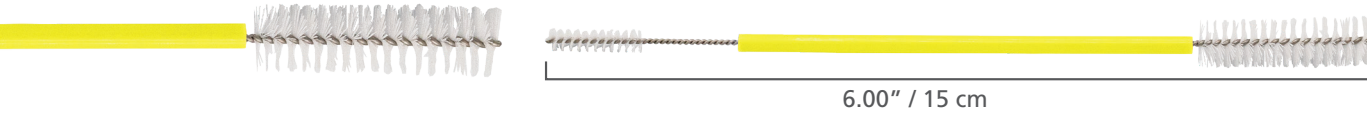

40-18511

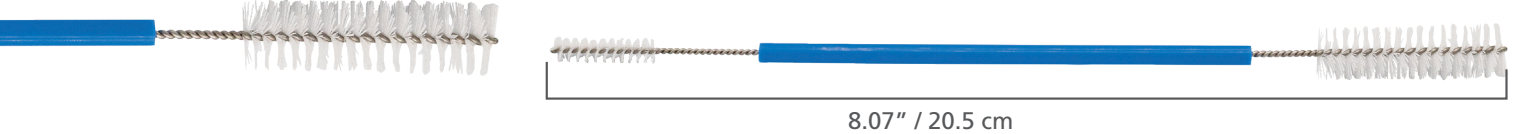

36-10010



27-15311



PRODUCT SELECTION


ITEM #	DESCRIPTION	LENGTH	BRUSH 1 LENGTH	BRUSH 2 LENGTH	BRUSH 1 (Ø)	BRUSH 2 (Ø)
27-05197	Double Trumpet 5/7 White Color	5.00" / 13 cm	0.20" / 5 mm	1.00" / 25 mm	0.20" / 5 mm	0.28" / 7 mm
PACKAGING: Bag of 3 (single-use recommended). Volume Packaging: Bag of 50 use ITEM # 27-05197-50						
27-57787	Double. Universal 5/11 Yellow Color	6.00" / 15 cm	1.18" / 30 mm	0.79" / 20 mm	0.43" / 11 mm	0.20" / 5 mm
PACKAGING: Bag of 3 (single-use recommended). Volume Packaging: Bag of 50 use ITEM # 27-57787-50						
40-18511	Double Universal 5/11 Extra Long Blue Color	8.07" / 20.5 cm	1.57" / 40 mm	0.79" / 20 mm	0.43" / 11 mm	0.20" / 5 mm
PACKAGING: Bag of 3 (reusable and can be steam sterilized). Volume Packaging: Bag of 50 use ITEM # 40-18511-50						
36-10010	Single end large Red Color	4.3" / 11 cm	1.57" / 40 mm	N/A	0.39" / 10 mm	N/A
PACKAGING: Bag of 3 (single-use recommended). Volume Packaging: Bag of 50 use ITEM # 36-10010-50						
27-15311	Single end large Blue Color	4" / 10.3 cm	1.57" / 40 mm	N/A	0.39" / 10 mm	N/A
PACKAGING: Bag of 3 (single-use recommended). Volume Packaging: Bag of 50 use ITEM # 27-15311-50						

BATRIK

V-BRUSHES

EXTERNAL CLEANING BRUSHES

Our External Cleaning Brushes for endoscopes / optical instruments provide excellent cleaning capabilities. With their unique U-shape bristles and innovative design, our brushes are the ultimate tool for cleaning the exterior of your endoscopic instruments. Free from sharp or metallic ends that can scratch or cause damage to valuable equipment, both single-use and reusable models are latex free and feature ergonomic handles.

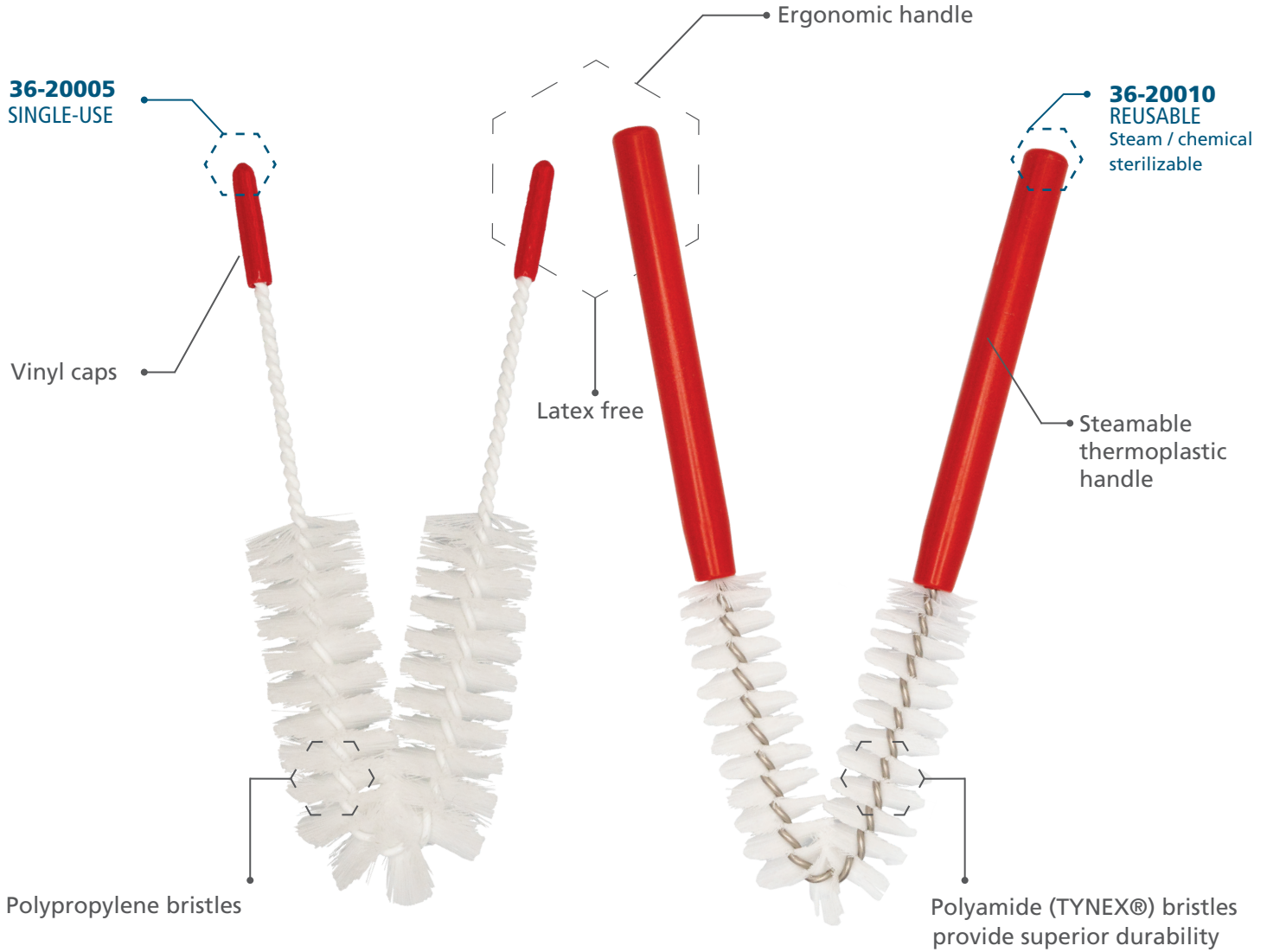

PRODUCT SELECTION

ITEM #	DESCRIPTION	BRUSH LENGTH / (Ø)	HANDLE LENGTH / (Ø)	TOTAL LENGTH	PACKAGING
36-20005	Single-Use External Cleaning Brush	70 mm / 25 mm	80 mm / 3 mm	150 mm	50 / Box
36-20010	Reusable External Cleaning Brush	90 mm / 25 mm	120 mm / 12 mm	210 mm	3 / Bag

NEW!
BATRIK


TAG-IT™ | NOTES

WHEN SIMPLY BROKEN DOESN'T CUT IT—TAG IT, TELL IT, OWN IT!

Tag-It™ Notes offer the freedom to **identify and document instructions for instruments, trays, and carts** before and after **sterilization**. Tag-It™ Notes are **tear-resistant, withstand extreme chemicals, high temperature** and ensures **reliable communication**.

WHY YOU'LL LOVE THEM:

- Tag-It™ Notes gives you the **freedom** to clearly document important information **without the clutter**.
- Robust and easy-to-use **tie wraps** designed to **securely attach** to instruments (134°C & chemical tested).
- Includes a **fast-drying, heat-resistant permanent marker** for writing your important notes.
- Choose from other **customizable cards** for personalized needs or **pre-made cards** for fast, convenient use.
- When you need more than a "repair" or "sharpen" seal to identify an issue. Instead, **Tag it, Tell it, Own it.**

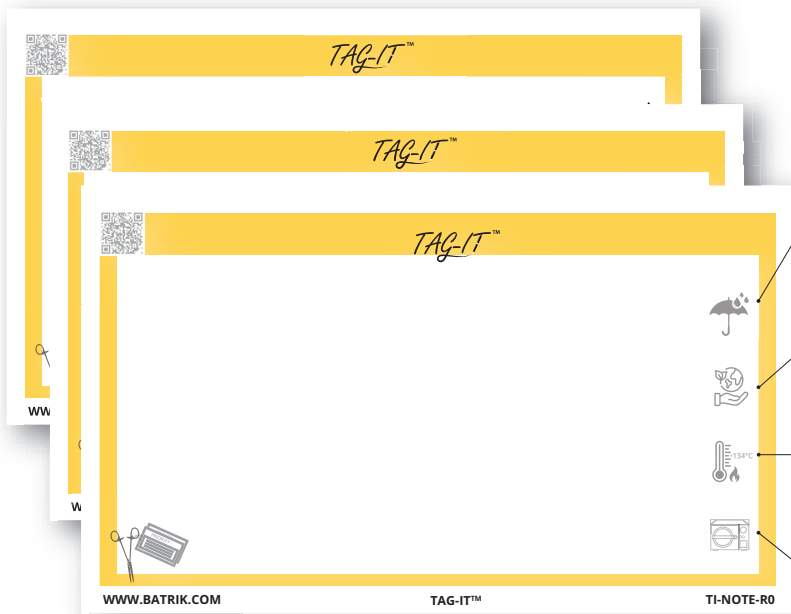


WATERPROOF

The synthetic paper withstands spills and liquid exposure without any damage.




ECO-FRIENDLY

Made from sustainable materials.

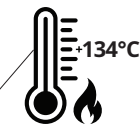

HEAT RESISTANT

Can resist high temperatures of up to 273.2°F (134°C). Compatible with Gas Plasma Sterilization (STERAD®100 NX Sterilization).

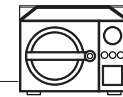

DURABLE & TEAR-PROOF

Handles rigorous cleaning, tough conditions, and resists tearing.

PRODUCT SELECTION

ITEM #	COLOR		SIZE	PACKAGING
TI-NOTE-R0	Yellow	NOTES	3" x 5" / 8 cm x 13cm	100 / PKG


All Tag-It™ are single-use. PACKAGING: 1 pack of 100 tags, 1 permanent marker, and 100 ties.

TAG-IT™

Tag-It™ Cards are an effortless way to add information with either pre-printed or blank formats. These **durable** tags streamline your workflow, **withstand tough conditions**, and feature **(optional) indicator dots** for quick verification. They help keep **instruments, trays, and carts** organized, ensuring efficiency before and after sterilization/decontamination.

KEY PRODUCT FEATURES:

- **Special Permanent Marker:** The ink will not run or appear blurry after exposure to chemicals or high temperature (max. 260°C/500°F) sterilization. *(Regular pens or markers can be used only if the tags are not going to be exposed to any sterilization/decontamination).*
- **Tie Wraps:** Assist in securing the tags to almost any instrument or container. *(Tie wraps have been tested for chemical and steam sterilization).*
- **Our Tag-It™:** Will not shrink, absorb water/chemicals, or lose its intended integrity after chemical or steam sterilization.

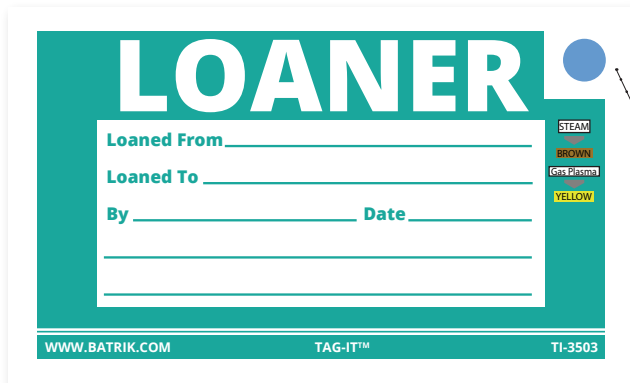

**(Optional)
All Tags Available
with Indicator Dots**
Provides visual confirmation that the tag was processed.

**CLICK BELOW TO
WATCH OUR VIDEO**



PRODUCT SELECTION:

ITEM #	COLOR	DESCRIPTION	SIZE	PACKAGING
TI-3502	Blue	VENDOR	3" x 5" / 8 cm x 13cm	100 / PKG
TI-3503	Green	LOANER	3" x 5" / 8 cm x 13cm	100 / PKG
TI-3504	Red	PRIORITY	3" x 5" / 8 cm x 13cm	100 / PKG
TI-3505	Yellow	RE-RUN	3" x 5" / 8 cm x 13cm	100 / PKG
TI-3506	Pink	MISSING	3" x 5" / 8 cm x 13cm	100 / PKG
TI-3507	Orange	BROKEN	3" x 5" / 8 cm x 13cm	100 / PKG
TI-3509	Black	BLANK	3" x 5" / 8 cm x 13cm	100 / PKG
TI-35B1	White	BLANK	3" x 5" / 8 cm x 13cm	100 / PKG
TI-35B4	Red	BROKEN	3" x 5" / 8 cm x 13cm	100 / PKG
TI-35C22	Light Blue	VENDOR CHECK IN + OUT	3" x 5" / 8 cm x 13cm	100 / PKG
TI-35L3	Green	LOANER CHECK IN + OUT	3" x 5" / 8 cm x 13cm	100 / PKG
TI-35T2	Blue	QUICK TURNOVER	3" x 5" / 8 cm x 13cm	100 / PKG
TI-PM09	Black	FINE PERMANENT MARKERS	1.5" Ø x 5.5" / 3.81cm Ø x 12.7 cm	36 / Bag

PACKAGING: 1 pack of 100 tags, 100 tie wraps, and 1 permanent marker.

OPTIONAL | PROCESSING INDICATOR DOTS: Add suffix "ST" to end of code for steam dots, "PS" for plasma dots.

Available in other languages or customize your own.

PRODUCT SELECTION: (customization also available)


VENDOR / BLUE

ITEM # **SIZE**
TI-3502 3" x 5"
 (8 cm x 13cm)




NOTES / YELLOW

ITEM # **SIZE**
TI-NOTE-R0 3" x 5"
 (8 cm x 13cm)

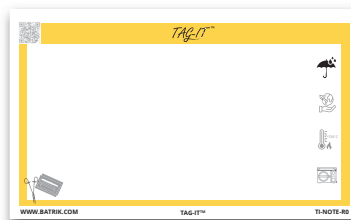

LOANER / GREEN

ITEM # **SIZE**
TI-3503 3" x 5"
 (8 cm x 13cm)

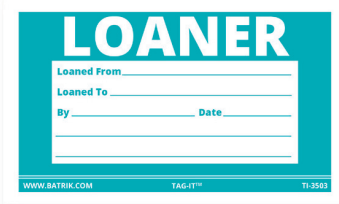


BLANK / WHITE

ITEM # **SIZE**
TI-35B1 3" x 5"
 (8 cm x 13cm)



PRIORITY / RED

ITEM # **SIZE**
TI-3504 3" x 5"
 (8 cm x 13cm)

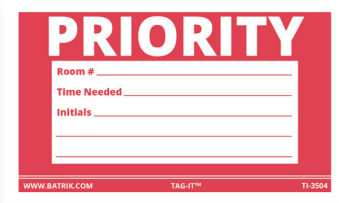

BROKEN / RED

ITEM # **SIZE**
TI-35B4 3" x 5"
 (8 cm x 13cm)



RE-RUN / YELLOW

ITEM # **SIZE**
TI-3505 3" x 5"
 (8 cm x 13cm)



VENDOR IN + OUT / LIGHT BLUE

ITEM #
TI-35C22
SIZE
 5" x 3"
 (13 cm x 8 cm)



MISSING / PINK

ITEM # **SIZE**
TI-3506 3" x 5"
 (8 cm x 13cm)



LOANER IN + OUT / GREEN

ITEM #
TI-35L3
SIZE
 5" x 3"
 (13 cm x 8 cm)




BROKEN/ ORANGE

ITEM # **SIZE**
TI-3507 3" x 5"
 (8 cm x 13cm)



BLANK / BLACK

ITEM # **SIZE**
TI-3509 3" x 5"
 (8 cm x 13cm)

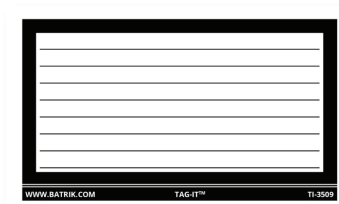

QUICK TURNOVER/ BLUE

ITEM # **SIZE**
TI-35T2 3" x 5"
 (8 cm x 13cm)



All information tags are single-use. PACKAGING: 1 pack of 100 tags, 100 tie wraps, and 1 permanent marker.

NEW!

BATRIK

TAG-IT™ | SCOPE

Batrik Scope Tags are tear-resistant, double-sided tags that will prove indispensable throughout the reprocessing journey. Our tags make it easy to **record important information**, including processing dates, with a user-friendly feature that allows you to select the month, day, and date. Crafted from **waterproof synthetic paper**, these tags feature a **self-looping design**, making it effortless to attach them to the scope without the need for additional fastening materials. Elevate your reprocessing protocols with Batrik!

KEY FEATURES

- The paper, tie wraps, and smudge-proof ink are designed to withstand high temperatures of up to 273.2°F (134°C) and are compatible with Gas plasma sterilization (STERRAD®100 NX Sterilization System).
- Tags include the month, day, and date selection on one side, making it convenient to mark the processing date accurately



COLOR-CODED AND CUSTOMIZABLE

- Available in seven colors
- The tags facilitate color coding for easy identification

SELF LOOPING

- The self-looping design allows for easy attachment to the scope
- Eliminating the need for additional fastening materials


SINGLE-USE AND LATEX FREE

- Batrik Scope Tags are designed for single use and are latex-free
- Promoting hygiene and reducing the risk of cross-contamination

PRODUCT SELECTION

ITEM #	COLOR	SIZE	PACKAGING
TI-ST05	Yellow	11.8" x 1.55" / 33.8 cm x 4 cm	100 / PKG
TI-ST02	Blue	11.8" x 1.55" / 33.8 cm x 4 cm	100 / PKG
TI-ST03	Green	11.8" x 1.55" / 33.8 cm x 4 cm	100 / PKG
TI-ST07	Orange	11.8" x 1.55" / 33.8 cm x 4 cm	100 / PKG
TI-ST06	Purple	11.8" x 1.55" / 33.8 cm x 4 cm	100 / PKG
TI-ST04	Red	11.8" x 1.55" / 33.8 cm x 4 cm	100 / PKG
TI-ST01	White	11.8" x 1.55" / 33.8 cm x 4 cm	100 / PKG


All information tags are single-use. PACKAGING: 1 pack of 100 tags and 1 permanent marker.

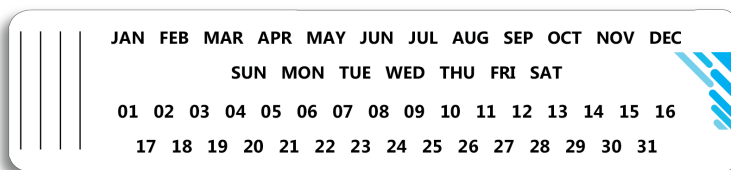

NEW!
BATRIK

TAG-IT™ | SCOPE REPROCESSING

Batrik Scope Reprocessing Tags are durable and can withstand handling during reprocessing **without tearing**, ensuring that information remains intact. These tags are vital in recording important information such as **date to reprocess on, scope model and serial number, sterilizer type, meeting sterilizer parameter conditions, alcohol flush, and hang date**. These tags have been meticulously designed to elevate **operational efficiency** and contribute to the delivery of the highest standard of patient care.

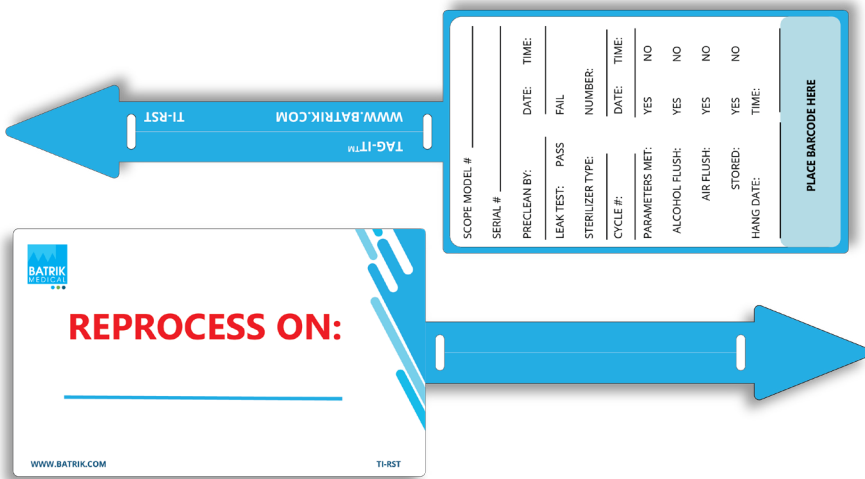


LOGGING IMPORTANT INFORMATION

- Easy for staff to indicate when the scope was processed
- Key scope, technician, and cycle-related information easily reported and identified.

SMART DESIGN

- The self-looping design allows for easy attachment to the scope
- No requirement for extra fastening materials
- Secure fastening as the loop is tightened around the instrument
- Features an integrated barcode area

TI-RST

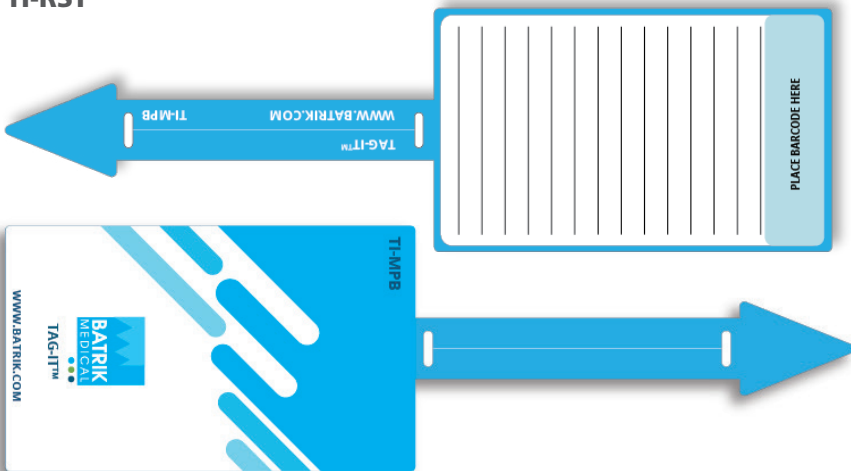

TI-MPB



*Suitable for various information tagging purposes

PRODUCT SELECTION



ITEM #	DESCRIPTION	DIMENSIONS	PACKAGING
TI-RST	REPROCESSING SCOPE TAG	3" x 10.5" / 8 cm x 27 cm	100 / PKG
TI-MPB	MULTI-PURPOSE BLANK TAG	3" x 10.5" / 8 cm x 27 cm	100 / PKG

All information tags are single-use. **PACKAGING: 1 pack of 100 tags and 1 permanent marker.**


BATRIK SILICONE PIN MAT

Our Silicone Pin Mat cushions, protects, and organizes delicate surgical instruments for sterilization. The elevated nubs provide greater air-flow ventilation and drying, significantly reducing residual moisture collection.

- Reusable
- Autoclavable
- Latex free

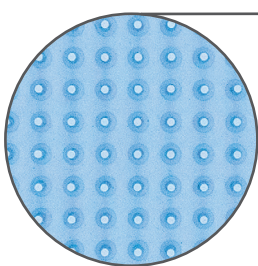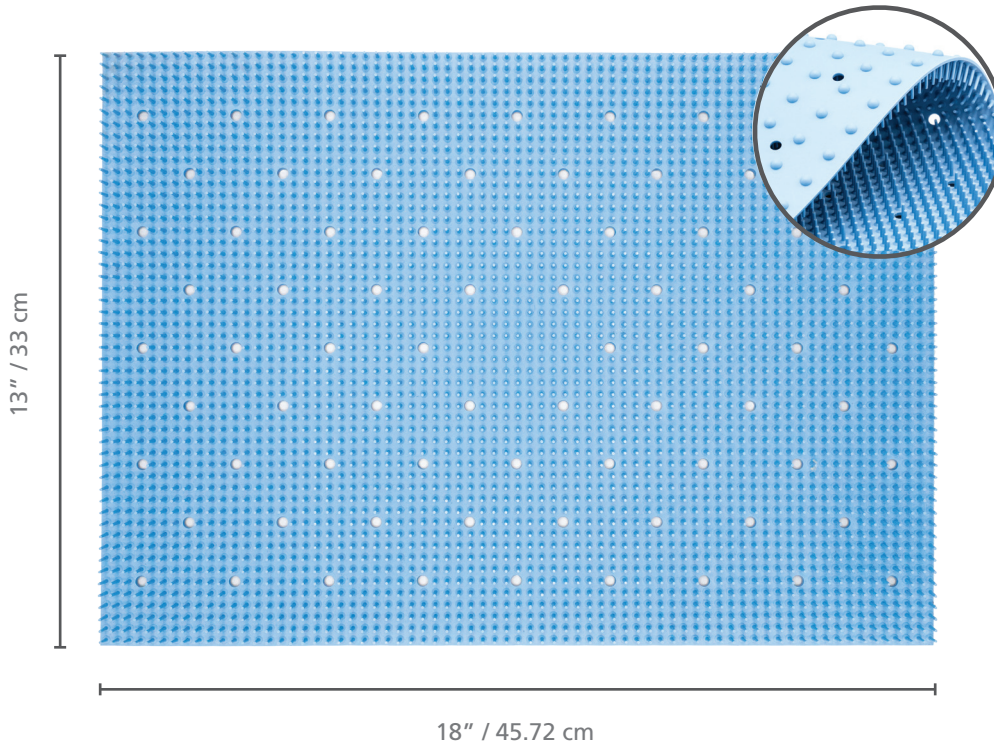

FRONT SIDE
Pins keep instruments securely in place

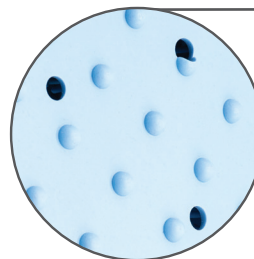

BACK SIDE
Numerous holes and elevated nubs provide superior venting and reduce water collection


PRODUCT SELECTION

ITEM #	DIMENSIONS	PACKAGING
13-2001	13" x 18" (33 cm x 45.72 cm) 11.97 mm high from base of mat to end of pin 15.53 mm high from base of elevated nub to end of pin	1 unit
13-2005	13" x 18" (33 cm x 45.72 cm) 11.97 mm high from base of mat to end of pin 15.53 mm high from base of elevated nub to end of pin	5 units

GOOD TO KNOW:
CUSTOM SIZES
ARE AVAILABLE
UPON REQUEST

BATRIK CORNER POCKETS™

With three models to choose from, our Corner Pockets™ portfolio is a truly comprehensive line of products designed to prevent tearing associated with sharp corners or feet on surgical trays.

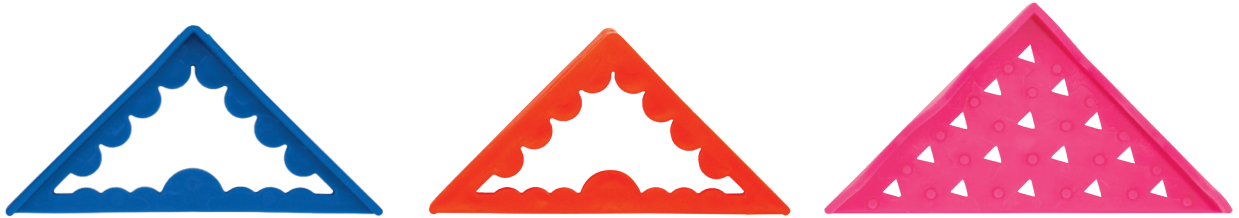

18-CP02


18-CP07

18-CP06

COLOUR	BLUE	ORANGE	PINK
USAGE	SINGLE-USE	REUSABLE <i>withstands repeated sterilization up to 134°C</i>	REUSABLE <i>withstands repeated sterilization up to 134°C</i>
DIMENSIONS	2.8" x 4" (7.3 cm x 10.3 cm)	2.8" x 4" (7.3 cm x 10.3 cm)	3.4" x 4.7" (8.5 cm x 12 cm)
LATEX-FREE	✓	✓	✓
NO MOISTURE RETENTION	✓	✓	✓
FABRICATION	Flexible, silicone-like plastic	Flexible, silicone-like plastic	Our most flexible plastic for easiest placement and most secure protection

SOFT POCKETS™

Our Soft Pockets™ Foam Corner Protectors are made from high density foam to effectively cushion and protect the corners of surgical trays to prevent tearing. They are especially useful with heavier instrument sets.

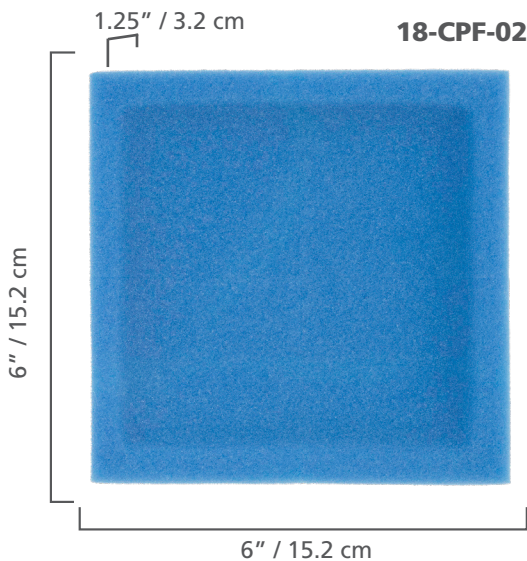

PRODUCT SELECTION

ITEM #	DESCRIPTION	PACKAGING	CASE DIMENSIONS
18-CP02**	Corner Pockets™ Single-Use, Blue	240 / Bag	80 / Box (24 Boxes / Pallet)
18-EVAL2**	Corner Pockets™ Single-Use, Blue	24 / Bag	-
18-CP07	Corner Pockets™ Reusable, Orange	80 / Bag	-
18-CORNER	Corner Pockets™ Reusable, Orange	320 / Box (4 Bags of 80 pcs)	-
18-EVAL7	Corner Pockets™ Reusable, Orange	24 / Bag	-
18-CP06	Corner Pockets™ Reusable, Pink	80 / Bag	-

** Although the single use Corner Pockets™ appear similar to the reusable models they are blended differently and have only been validated as single-use. Re-using them may compromise sterility.

NEW!

BATRIK

LILA'S INSTRUMENT STAIN REMOVER

Lila's special blend of naturally sourced powders provides outstanding removal of stubborn stains.

Convenient-to-use formula, you can achieve spotless results every time!

It can be used on surgical/dental instruments, autoclaves, ultrasonic cleaners, stainless steel countertops, trays, bowls and more.

Specifically designed for use on stainless steel, chrome, nickel, brass, copper, and silver.




EASE TO USE AS 1,2,3...

1. Dampen a cloth or grab a Lila's Low Linting Dura-Towel and dab into the jar.
2. Gently rub the area that is stained.¹
3. Rinse thoroughly for a flawless shine.²

A Little Goes a long Way!

Even great to tackle stubborn stains on steel pots and pans in your kitchen with ease.

Lila's - If you have used Spectra-Scrub or Integra's Instrument stain remover then you will love Lila's!

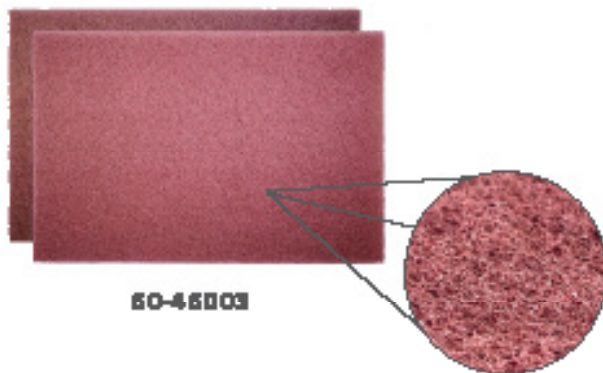

60-46003

LSR-BSD1



17-DT1213C

LILA'S SCRUB PAD

- Say goodbye to extra stubborn stains with Lila's Scrub Pad without leaving any metal residue.
- Simply add a small amount of Lila's Stain Remover with one of our Scrub Pad's for very large areas.
- The gentle yet effective scrub is the perfect choice for Ophthalmic Instruments.¹

LILA'S LOW-LINTING DURA-TOWEL

- Lila's Dura-Towel has long lasting strength and has a minimal lint composition.
- These single-use towels are clean room certified and are specifically designed for instrument wiping.
- Its ideal dimensions and unique attributes guarantee optimal performance, setting it apart from ordinary paper towels.

PRODUCT SELECTION




ITEM #	DESCRIPTION	DIMENSIONS	PACKAGING
LSR-BSD1	LILA'S INSTRUMENT STAIN REMOVER 3 oz (85g) plastic jar with a convenient screw-on lid	Ø 2.8" / 70 mm	1 Jar / Bag
LSR-B512	LILA'S INSTRUMENT STAIN REMOVER 3 oz (85g) plastic jar with a convenient screw-on lid	Ø 2.8" / 70 mm	12 Jars / Box
17-DT1213C	LILA'S LOW LINTING DURA-TOWEL 100% towels	12" x 13" / 31 cm x 33 cm Folded: 3" x 3.25" / 7.5 cm x 8.25 cm	80 / Bag 12 Bags / Box
60-46003	LILA'S SCRUB PAD	4" x 6" / 10 cm x 15 cm	3 / Bag

1: Although Lila's is completely safe, it is not recommended for use on ophthalmic instruments. Use Lila's Scrub Pad on its own.
2: The stain remover should not be mixed with bleach and is for external use only to remove instrument stains.

BATRIK INSTRUMENT PROTECTORS

Our Instrument Protectors protect the tips and edges of instruments to prevent damage during handling and sterilization. Available in a wide array of color coded sizes for almost every application, including osteotome, retractor, double, and perforated.

- Steam and gas permeable
- Perforated protectors reduce vacuum seal and moisture build-up
- Validated for pre-vacuum sterilization*
- Multiple color-coded sizes
- Latex Free
- Non-Sterile

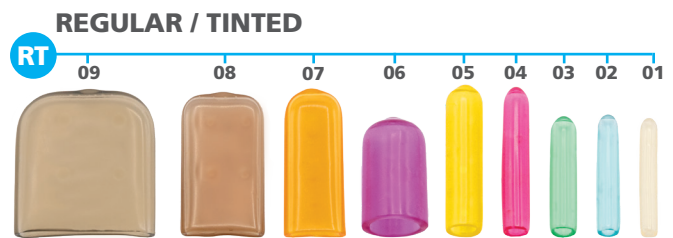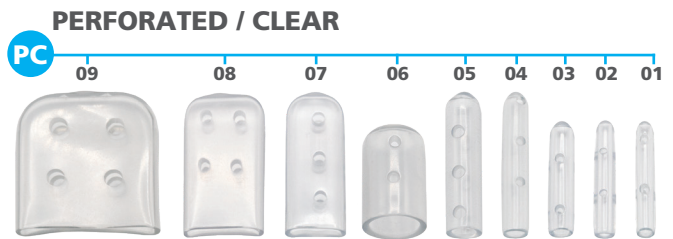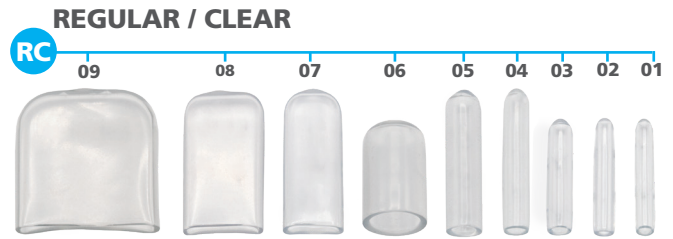

*Independent laboratory tested to validate pre-vacuum steam sterilization efficacy in a healthcare setting

PRODUCT SELECTION

REGULAR & PERFORATED / CLEAR



RC-07



PC-07

ITEM # REGULAR	ITEM # PERFORATED	COLOR	DIMENSIONS
RC-01	PC-01	Clear	0.062" x 0.75" (1.57 mm x 19 mm)
RC-02	PC-02	Clear	0.075" x 0.75" (2 mm x 19 mm)
RC-03	PC-03	Clear	0.109" x 0.75" (2.77 mm x 19 mm)
RC-04	PC-04	Clear	0.125" x 1" (3.18 mm x 25.4 mm)
RC-05	PC-05	Clear	0.187" x 1" (4.75 mm x 25.4 mm)
RC-06	PC-06	Clear	0.406" x 0.75" (10.3 mm x 19 mm)
RC-07	PC-07	Clear	0.062" x 0.375" x 1" (1.57 mm x 9.53 mm x 25.4 mm)
RC-08	PC-08	Clear	0.078" x 0.625" x 1" (2 mm x 15.875 mm x 25.4 mm)
RC-09	PC-09	Clear	0.125" x 1" x 1" (3.18 mm x 25.4 mm x 25.4 mm)
RC-10	PC-10	Clear	Assorted

PACKAGING: Bag of 100

REGULAR & PERFORATED / TINTED



RT-07



PT-07

RT-01	PT-01	Tinted White	0.062" x 0.75" (1.57 mm x 19 mm)
RT-02	PT-02	Tinted Blue	0.075" x 0.75" (2 mm x 19 mm)
RT-03	PT-03	Tinted Green	0.109" x 0.75" (2.77 mm x 19 mm)
RT-04	PT-04	Tinted Red	0.125" x 1" (3.18 mm x 25.4 mm)
RT-05	PT-05	Tinted Yellow	0.187" x 1" (4.75 mm x 25.4 mm)
RT-06	PT-06	Tinted Maroon	0.406" x 0.75" (10.3 mm x 19 mm)
RT-07	PT-07	Tinted Orange	0.062" x 0.375" x 1" (1.57 mm x 9.53 mm x 25.4 mm)
RT-08	PT-08	Tinted Brown	0.078" x 0.625" x 1" (2 mm x 15.875 mm x 25.4 mm)
RT-09	PT-09	Tinted Black	0.125" x 1" x 1" (3.18 mm x 25.4 mm x 25.4 mm)
RT-10	PT-10	Tinted Assorted	Assorted

PACKAGING: Bag of 100

REGULAR & PERFORATED / SOLID



RS-07



PS-07


RS-01	PS-01	Solid White	0.062" x 0.75" (1.57 mm x 19 mm)
RS-02	PS-02	Solid Blue	0.075" x 0.75" (2 mm x 19 mm)
RS-03	PS-03	Solid Green	0.109" x 0.75" (2.77 mm x 19 mm)
RS-04	PS-04	Solid Red	0.125" x 1" (3.18 mm x 25.4 mm)
RS-05	PS-05	Solid Yellow	0.187" x 1" (4.75 mm x 25.4 mm)
RS-06	PS-06	Solid Maroon	0.406" x 0.75" (10.3 mm x 19 mm)
RS-07	PS-07	Solid Orange	0.062" x 0.375" x 1" (1.57 mm x 9.53 mm x 25.4 mm)
RS-08	PS-08	Solid Brown	0.078" x 0.625" x 1" (2 mm x 15.875 mm x 25.4 mm)
RS-09	PS-09	Solid Black	0.125" x 1" x 1" (3.18 mm x 25.4 mm x 25.4 mm)
RS-10	PS-10	Solid Assorted	Assorted

PACKAGING: Bag of 100

ITEM #	DIMENSIONS
OP-11	0.25" x 0.75" (6.35 mm x 19.05 mm)
OP-12	0.375" x 1" (9.53 mm x 25.4 mm)
OP-13	0.5" x 1" (12.7 mm x 25.4 mm)
OP-14	0.75" x 1" (19.05 mm x 25.4 mm)
OP-15	1" x 1" (25.4 mm x 25.4 mm)
OP-16	1.25" x 1.25" (31.75 mm x 31.75 mm)
OP-17	1.5" x 1.5" (38.1 mm x 38.1 mm)
OP-18	1.75" x 1.25" (44.45 mm x 31.75 mm)
OP-20	Assorted

PACKAGING: Bag of 50

OSTEOTOME PROTECTORS

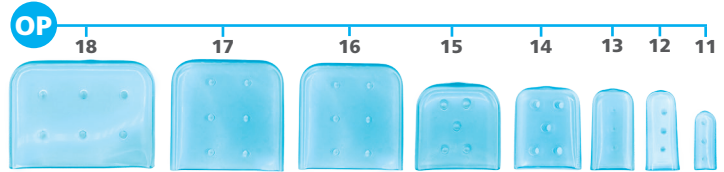

ITEM #	DIMENSIONS	PACKAGING
RGS-01	1.5" (38 mm)	55 / Bag
RGS-02	2" (51 mm)	45 / Bag
RGS-03	3" (76.2 mm)	40 / Bag


SINGLE-TIP RETRACTOR PROTECTORS



ITEM #	COLOR	SIZE
RGM-01	Pink	Small
RGM-02	Clear	Small
RGM-03	Green	Medium
RGM-04	Clear	Medium

PACKAGING: Bag of 50


WEITLANER RETRACTOR PROTECTORS

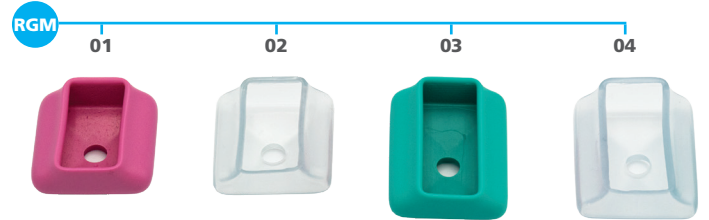

ITEM #	COLOR	DIMENSIONS
DG-01	Clear	0.0625" x 0.75" (1.59 mm x 19.05 mm)
DG-02	Blue	0.078" x 0.75" (1.98 mm x 19.05 mm)
DG-03	Green	0.109" x 0.75" (2.78 mm x 19.05 mm)
DG-04	Red	0.125" x 1" (3.18 mm x 25.4 mm)
DG-05	Yellow	0.187" x 1" (4.76 mm x 25.4 mm)
DG-07	Orange	0.062" x 0.375" x 1" (1.59 mm x 9.53 mm x 25.4 mm)
DG-10	Assorted	Assorted

PACKAGING: Bag of 50

PERFORATED DOUBLE PROTECTORS

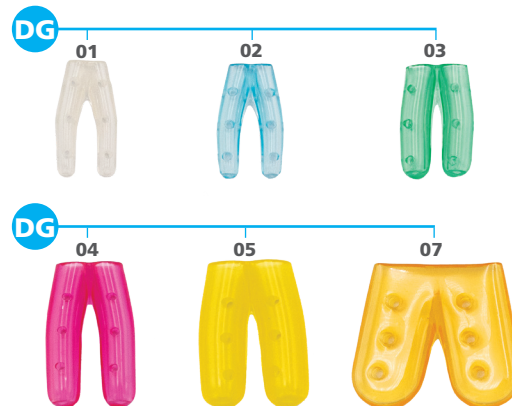

BATRIK PrAmLi Power³ LUBE™ Instrument Lubricant

Specially formulated to lubricate surgical instruments, hinges, box locks, and general equipment during autoclaving. PrAmLi Power³ LUBE™ is uniquely formulated with an all natural corrosion inhibitor.

- Forms a protective coating which can extend the instrument's, useful life
- Bacteriostatic
- Neutral pH 6.5 - 7.5
- Prevents corrosion

TIP:
FOR MILK BATH SOLUTION,
DILUTE 1 PART SOLUTION
IN 6 PARTS WATER

Formulated with
PrAmLi Power³ Enzymes™

ITEM #	DIMENSIONS
OP-11	0.25" x 0.75" (6.35 mm x 19.05 mm)
OP-12	0.375" x 1" (9.53 mm x 25.4 mm)
OP-13	0.5" x 1" (12.7 mm x 25.4 mm)
OP-14	0.75" x 1" (19.05 mm x 25.4 mm)
OP-15	1" x 1" (25.4 mm x 25.4 mm)
OP-16	1.25" x 1.25" (31.75 mm x 31.75 mm)
OP-17	1.5" x 1.5" (38.1 mm x 38.1 mm)
OP-18	1.75" x 1.25" (44.45 mm x 31.75 mm)
OP-20	Assorted

PACKAGING: Bag of 50




PRODUCT SELECTION

ITEM #	DILUTION LUBRICANT / WATER	PACKAGING	CASE DIMENSIONS
SS-21-6001	1 to 4 ml / L (0.125 to 0.5 oz / Gallon)	1 x 3.8 L (1 Gallon)	6x6x13"
SS-21-6004	1 to 4 ml / L (0.125 to 0.5 oz / Gallon)	4 x 3.8 L (4 Gallons)	13x13x14"
SS-21-6025	1 to 4 ml / L (0.125 to 0.5 oz / Gallon)	2 x 9.5 L (2 x 2.5 Gallons)	14.5x9.25x15"
SS-21-6032	No dilution required	1 x 946 ml (32 oz) Spray bottle	10x4"
SS-21-6012	No dilution required	12 X 946 ml (32 oz) Spray Bottles	18x12x11"

BATRIK S.T.I.C.S.

Our Simple Tracking Identification Color System, or S.T.I.C.S., is an easy to use, latex free, tool that allows you to quickly identify your instruments by department and/or procedure. The tape adheres securely and every tape sheet comes with a protective cover to ensure it lasts as long as you need.

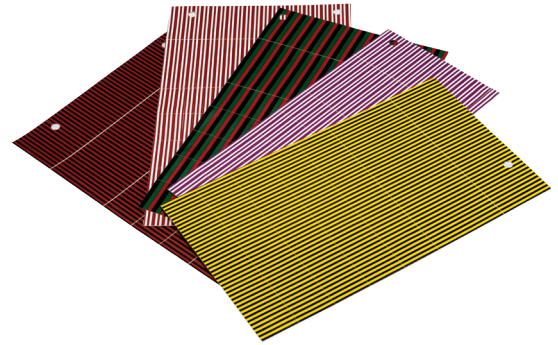

Each 8.5" x 11" (216 mm x 279 mm) sheet contains the following adhesive strips:

- 1" (25.4 mm) x 220
- 1.5" (38.1 mm) x 44
- 2" (50.8 mm) x 44
- Total tape per sheet 374" (9500 mm)

PRODUCT SELECTION

SOLID COLORS	
ITEM #	COLOR
T-850-2-01	White
T-850-2-02	Royal Blue
T-850-2-03	Kelly Green
T-850-2-04	Red
T-850-2-05	Yellow
T-850-2-06	Purple
T-850-2-07	Orange
T-850-2-08	Brown
T-850-2-09	Black
T-850-2-10	Gray
T-850-2-11	Lavender
T-850-2-12	Beige
T-850-2-13	Light Blue
T-850-2-14	Lime Green
T-850-2-15	Pink
T-850-2-17	Plum
T-850-2-18	Mint Green

BLACK STRIPES	
ITEM #	COLOR
T-850-3-01	Black/White
T-850-3-03	Black/Kelly Green
T-850-3-04	Black/Red
T-850-3-05	Black/Yellow
T-850-3-07	Black/Orange
T-850-3-10	Black/Gray
T-850-3-11	Black/Lavender
T-850-3-12	Black/Beige
T-850-3-13	Black/Light Blue
T-850-3-14	Black/Lime Green
T-850-3-15	Black/Pink
T-850-3-17	Black/Plum

TRI-COLOR	
ITEM #	COLOR
T-850-5-51	Red/White/Blue
T-850-5-54	Red/Black/Green

SPECIALTY	
ITEM #	COLOR
T-850-7-04	Yellow/Red

WHITE STRIPES	
ITEM #	COLOR
T-850-4-02	White/Royal Blue
T-850-4-03	White/Kelly Green
T-850-4-04	White/Red
T-850-4-05	White/Yellow
T-850-4-07	White/Orange
T-850-4-08	White/Brown
T-850-4-11	White/Lavender
T-850-4-13	White/Light Blue
T-850-4-14	White/Lime Green
T-850-4-15	White/Pink
T-850-4-17	White/Plum


FOUR-COLOR	
ITEM #	COLOR
T-850-2-31	Red/Yellow/ Purple/Orange
T-850-2-32	Blue/White/ Black/Green

BATRIK ID TAPE



A simple solution for marking surgical instruments with the aid of color-coding, our ID Tape is durable and will not chip, flake, or peel. It can be used in autoclaves, withstanding repeated exposure to sterilization. Independent validation tests demonstrated the efficiency of sterilization of marked instruments, even underneath the Batrik ID Tape.

AVAILABLE WIDTHS

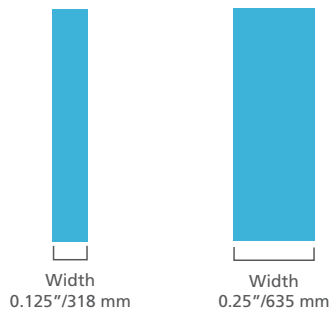

Quickly identifies instruments by department or procedure



Easy to apply



Withstands repeated exposure to sterilization



Has more per tape roll (300" / 7620 mm)



Available in 7 colors and 2 widths



Conveniently packaged in a clear box



Non-Sterile

PRODUCT SELECTION

ITEM #	ITEM#	COLOR
0.125" / 318 MM	0.25" / 635 MM	
T-125-01	T-250-01	White
T-125-02	T-250-02	Blue
T-125-03	T-250-03	Green
T-125-04	T-250-04	Red
T-125-05	T-250-05	Yellow
T-125-07	T-250-07	Orange
T-125-09	T-250-09	Black

BATRIK TAPE-B-GONE™ TAPE REMOVER

Tape-B-Gone™ is a uniquely formulated and ready-to-use proprietary formula that easily removes adhesives from medical equipment and instruments. Unlike many other brands, Tape-B-Gone™ does not contain highly dangerous or caustic chemicals, is not flammable or combustible, and is safe for use on metal, glass, plastic, and rubber without any special equipment or ventilation.



PRODUCT SELECTION

ITEM #	DESCRIPTION	PACKAGING
T-08RMV	8 oz (237 ml) bottle	1 Unit
T-32RMV	8 oz (237 ml) bottle	4 / Pack
T-SAMP-RMV	0.5 oz (15 ml) bottle	1 Unit


BATRIK PADLOCKS SECURITY TAGS

Our Padlock Security Tags for instrument containers and general hospital applications are specifically designed for the healthcare market. Our material selection, fabrication process, and engineering ensure consistent quality and high temperature compliance with the various steam sterilization temperatures in worldwide hospital applications.


Available with optional reactive ink indicator dots that provide visual sterility indication. Light blue dots turn brown/dark brown when steam sterilized.

AVAILABLE IN THREE FORMATS



700-SERIES

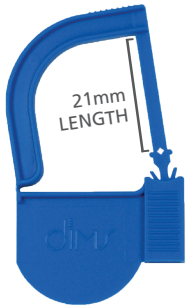


800-SERIES

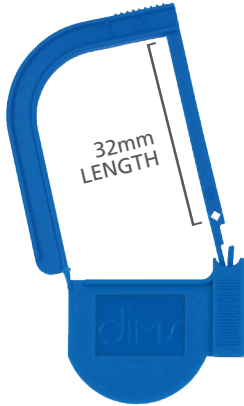


900-SERIES

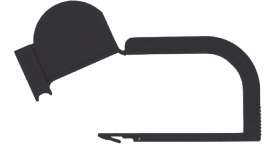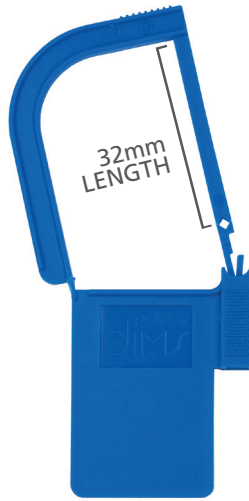

Made of polypropylene

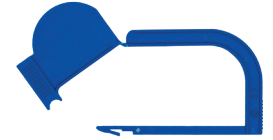

Latex and DEHP Free

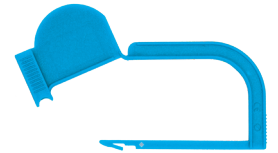

Will not melt or soften with any established steam sterilization protocol




Fully open design (45°) facilitates insertion into instrument containers



Easy twist breaking design



Large base for easier handling

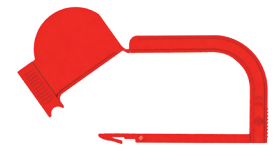

Available in different size



Custom Colors and OEM available upon request

PRINTED & CUSTOM NUMBERED PADLOCKS

Available in 700, 800 & 900 SERIES

PRODUCT SELECTION

COLOR	PRINTED WORD
PRINTED PADLOCKS	
WHITE	CLEAN
BLUE	STERILE
RED	REPAIR
YELLOW	SHARPEN
YELLOW	LOANER
BLACK	DISCARD
CUSTOM NUMBERED	
ALL COLORS	up to 8 numbers



CUSTOM NUMBERED
Custom printed padlocks (up to 8 characters)



PRINTED
6 different printed labels available

NOTE: Printed & Custom Numbered Padlocks are available in 700, 800, and 900 Series and can be ordered in 100, 200, and 500 unit quantities. Please contact us for ordering details.

INDICATOR DOTS



REGULAR



700-SERIES PADLOCKS

The 700 series padlocks are perfect for "Aesculap" containers and those with a short, large opening.

PRODUCT SELECTION

REGULAR PADLOCKS

COLOR	100 UNITS	200 UNITS	500 UNITS
WHITE	15-701	15-711	15-710
BLUE	15-702	15-712	15-720
GREEN	15-703	15-713	15-730
RED	15-704	15-714	15-740
YELLOW	15-705	15-715	15-750
PINK	15-706	15-716	15-760
ORANGE	15-707	15-717	15-770
BLACK	15-709	15-719	15-790
LIGHT BLUE	15-722	15-732	15-742

WITH INDICATORS DOTS

COLOR	100 UNITS	200 UNITS	500 UNITS
WHITE	15-701-IND	15-711-IND	15-710-IND
BLUE	15-702-IND	15-712-IND	15-720-IND
GREEN	15-703-IND	15-713-IND	15-730-IND
RED	15-704-IND	15-714-IND	15-740-IND
YELLOW	15-705-IND	15-715-IND	15-750-IND
PINK	15-706-IND	15-716-IND	15-760-IND
ORANGE	15-707-IND	15-717-IND	15-770-IND
BLACK	15-709-IND	15-719-IND	15-790-IND
LIGHT BLUE	15-722-IND	15-732-IND	15-742-IND

Double locking hooks can be inserted in a large guide hole



The grooves on the padlock facilitate handling



INDICATOR DOTS

REGULAR



800-SERIES PADLOCKS

The 800 series padlocks have a longer arm (32 mm) and have a single, thinner hook than the 700 series. They are perfect for almost any type of container, even those with long and thin channels.

PRODUCT SELECTION

REGULAR PADLOCKS

COLOR	100 UNITS	200 UNITS	500 UNITS
WHITE	15-801	15-811	15-810
BLUE	15-802	15-812	15-820
GREEN	15-803	15-813	15-830
RED	15-804	15-814	15-840
YELLOW	15-805	15-815	15-850
PINK	15-806	15-816	15-860
ORANGE	15-807	15-817	15-870
BLACK	15-809	15-819	15-890
LIGHT BLUE	15-822	15-832	15-842

WITH INDICATORS DOTS

COLOR	100 UNITS	200 UNITS	500 UNITS
WHITE	15-801-IND	15-811-IND	15-810-IND
BLUE	15-802-IND	15-812-IND	15-820-IND
GREEN	15-803-IND	15-813-IND	15-830-IND
RED	15-804-IND	15-814-IND	15-840-IND
YELLOW	15-805-IND	15-815-IND	15-850-IND
PINK	15-806-IND	15-816-IND	15-860-IND
ORANGE	15-807-IND	15-817-IND	15-870-IND
BLACK	15-809-IND	15-819-IND	15-890-IND
LIGHT BLUE	15-822-IND	15-832-IND	15-842-IND

Single locking hook can be inserted in a large guide hole




The locking system allows for easy closure



INDICATOR DOTS

REGULAR

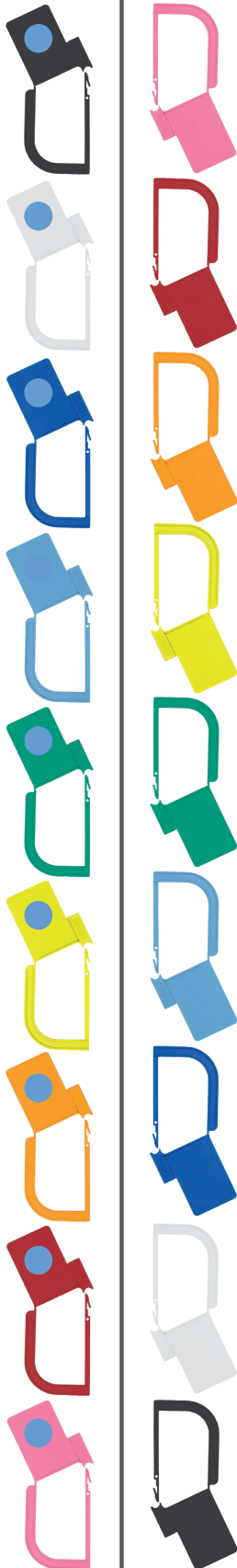

900-SERIES PADLOCKS

The 900 series padlocks are identical to the 800 series, but with a larger rectangular base which allows for the placement of a 20 mm x 30 mm tracking label.

PRODUCT SELECTION

REGULAR PADLOCKS

COLOR	100 UNITS	200 UNITS	500 UNITS
WHITE	15-901	15-911	15-910
BLUE	15-902	15-912	15-920
GREEN	15-903	15-913	15-930
RED	15-904	15-914	15-940
YELLOW	15-905	15-915	15-950
PINK	15-906	15-916	15-960
ORANGE	15-907	15-917	15-970
BLACK	15-909	15-919	15-990
LIGHT BLUE	15-922	15-932	15-942

WITH INDICATORS DOTS

COLOR	100 UNITS	200 UNITS	500 UNITS
WHITE	15-901-IND	15-911-IND	15-910-IND
BLUE	15-902-IND	15-912-IND	15-920-IND
GREEN	15-903-IND	15-913-IND	15-930-IND
RED	15-904-IND	15-914-IND	15-940-IND
YELLOW	15-905-IND	15-915-IND	15-950-IND
PINK	15-906-IND	15-916-IND	15-960-IND
ORANGE	15-907-IND	15-917-IND	15-970-IND
BLACK	15-909-IND	15-919-IND	15-990-IND
LIGHT BLUE	15-922-IND	15-932-IND	15-942-IND

Single locking hook
for use with smaller diameters



Ideal for adding labels/ indicators indicators or markings to the large rectangular base



SAVEFACE™ DRAPED SHIELDS

REDUCE PERSONAL INJURY
SUPERIOR INSTRUMENT DECONTAMINATION PROTECTION

SAVEFACE DRAPED-FACE SHIELD

In the Sterile Processing Department (SPD)/Central Sterile Supply Department (CSSD), personal safety is paramount.

Our Draped Face Shield provides essential protection during high-risk instrument decontamination.

EFFECTIVE PROTECTION IS NON-NEGOTIABLE

SAVEFACE Draped dynamically reduces harmful exposure to caustic chemicals as well as the risk of cross-contamination. Effective face protection is non-negotiable.

PPE CATEGORY 2 CERTIFIED

SAVEFACE DRAPED is tested to protect. Wear it with confidence.



PATENT PENDING



INCLUDES CLOSING SYSTEM TO REDUCE RISK OF FLUID CONTACT FROM BELOW THE SHIELD

SAVEFACE: KEY FEATURES

ANSI/AAMI PB70 LEVEL 3 COMPLIANT DRAPE:

- Provides a high level of fluid penetration protection.
- Tuck into clothing or pinch together.

EXTRA-LONG 9"/23CM WRAP-AROUND SHIELD:

- Ensures full coverage, reducing the risk of exposure from different angles.

COMFORTABLE:

- The soft, self-shaping foam headpiece molds comfortably to your forehead, allowing for extended wear without discomfort.

OPTICAL CLARITY WITH ANTI-FOG AND ANTI-GLARE

DON'T MISS THE OPPORTUNITY TO EXPLORE OUR BROADER RANGE OF PROTECTIVE SHIELDS.

Visit our website www.surgmed.com or consult our detailed literature to find the shield that best fits your needs.

VERY IMPORTANT

IS YOUR PRESENT FACE SHIELD CATEGORY 2 CERTIFIED?

Batrik face shields are independently certified to comply with the PPE Category 2 standard. Many facilities are now requesting compliance to this standard.

SHIELD LENGTH
9" / 23 CM

DRAPE LENGTH
6" / 15 CM



SF-240



SF-8240




PRODUCT SELECTION

ITEM #	DESCRIPTION	DIMENSIONS	TOTAL LENGTH	PACKAGING	CASE DIMENSIONS
SF-240	BATRIK SAVEFACE™ no drape	Shield Length: 9" / 23cm Shield width: 13' / 33cm	9" (23cm)	240 / Case (8 Boxes of 30 Pcs)	26.5" x 21" x 15"
SF-8240	BATRIK SAVEFACE™ with ANSI/AAMI PB70 Level 3 compliant protection drape	Shield Length: 9" / 23cm Shield width: 13' / 33cm Shield width: 6" / 15cm	58 cm / 23"	240 / Case (8 Boxes of 30 Pcs)	26.5" x 21" x 15"

BATRIK SCRUB PAD

Our Scrub Pad is ideal for pre-cleaning your surgical instruments without leaving any metal residue that often occurs when using steel wool or metal sponges. The abrasiveness is perfectly balanced to provide you with efficient cleaning without damaging your instruments. Its fine fibers allow the user to polish, clean, and debur all types of material surfaces.

- Abrasive Material: Aluminum Oxide
- Grit: 400
- Grade: Very Fine
- Style: Flexible
- Pad material: Non-Woven
- Color: Maroon
- Type: Hand Pad

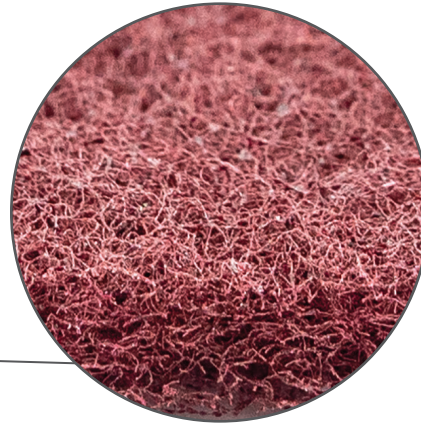

PRODUCT SELECTION

ITEM #	DESCRIPTION	DIMENSIONS	PACKAGING
60-46050	Scrub Pad	4" x 6" / 10 cm x 15 cm	150 / Box (50 bags of 3 units each)
60-46003	Scrub Pad	4" x 6" / 10 cm x 15 cm	3 / Bag

Your First ™

ENZYMATIC & NON-ENZYMATIC CLEANERS

Our Enzymatic and Non-Enzymatic Cleaners are specially formulated to remove organic material from medical equipment and instrumentation.

Our **Enzymatic Cleaner** formats are formulated with PrAmLi Power3 Enzymes™, our proprietary unique blend of enzymes, and a surfactant to efficiently remove all proteins, blood, fatty tissue, and more. They also remove stains from surgical scrubs, linens, and medical clothing.

Our **Non-Enzymatic Cleaner** formats are formulated with a free-rinsing, neutral pH detergent blend for effective cleaning.



Both formats are powerful solutions that break down all bioburden for effective surface and instrument pre-cleaning*

PRODUCT SELECTION

ENZYMATIC CLEANER	NON-ENZYMATIC CLEANER	DILUTION DETERGENT / WATER	PACKAGING	CASE DIMENSIONS
LOW FOAM				
SS-21-1001	SS-721-1001	4 to 8 ml / L (0.5 to 1 oz / Gal)	1 x 3.8 L (1 Gal)	6" x 6" x 13"
SS-21-1004	SS-721-1004	4 to 8 ml / L (0.5 to 1 oz / Gal)	4 x 3.8 L (4 Gal)	13" x 13" x 14"
MILD FOAM				
SS-21-2001	SS-721-2001	4 to 8 ml / L (0.5 to 1 oz / Gal)	1 x 3.8 L (1 Gal)	6" x 6" x 13"
SS-21-2004	SS-721-2004	4 to 8 ml / L (0.5 to 1 oz / Gal)	4 x 3.8 L (4 Gal)	13" x 13" x 14"
SUPER CONCENTRATE				
SS-21-3001	SS-721-3001	1 to 2 ml / L (0.125 to 0.25 oz / Gal)	1 x 3.8 L (1 Gal)	6" x 6" x 13"
SS-21-3004	SS-721-3004	1 to 2 ml / L (0.125 to 0.25 oz / Gal)	4 x 3.8 L (4 Gal)	13" x 13" x 14"

GOOD TO KNOW:

OUR LOW FOAM & SUPER CONCENTRATE DETERGENTS (ENZYMATIC & NON-ENZYMATIC) MAY BE USED IN WASHERS / ULTRASONIC CLEANERS.

WE RECOMMEND CONFIRMING SUGGESTED DILUTION RATIO AND WATER TEMPERATURE WITH THE MANUFACTURER OF YOUR WASHER / ULTRASONIC CLEANER AS MACHINES / SETTINGS VARY.

NOTE: Our Low Foam & Super Concentrate Detergent formulas are available in a 2.5 Gallon format based on a minimum order quantity. Contact us for more information.

*Caution: We recommend consulting your hospital pre-cleaning/cleaning protocols for various equipment/instrument cleaning detergent preferences.

Your First ™

ENVIRO-SACHET™

ENZYMATIC GEL & NON-ENZYMATIC CLEANER

Your First ™ Enviro-Sachet™, our unique cleaning solutions now available in Enzymatic & Non-Enzymatic formats.



- Our Enzymatic Gel format is formulated with the same triple action **PrAmLi Power³ Enzymes™** that you know and trust.
- Our Non-Enzymatic Cleaner format is formulated with a free-rinsing, neutral pH detergent blend.

Both formats are powerful solutions that break down all bioburden for effective surface and instrument pre-cleaning.

PRODUCT SELECTION

ENZYMATIC GEL	NON-ENZYMATIC CLEANER	DESCRIPTION	PACKAGING	CASE DIMENSIONS
17-166	717-166	1 Dry Tubular Sponge & Enviro-Sachet™ set	25 / Box 4 Boxes / Case	21" x 14" x 13.5"
17-266	717-266	1 Dry Flat Sponge & Enviro-Sachet™ set	50 / Box 4 Boxes / Case	21" x 14" x 13.5"
17-177	717-177	1 Dry Tubular Sponge, Enviro-Sachet™ & Stainless Steel Bowl kit	1 Kit	7" x 6" x 3"
17-277	717-277	1 Dry Flat Sponge, Enviro-Sachet™ & Stainless Steel Bowl kit	1 Kit	7" x 6" x 3"
17-EZSAMP	717-EZSAMP	SAMPLE KIT INCLUDES: 3 Dry Flat Sponges with Enviro-Sachet™ sets 1 Dry Tubular Sponge with Enviro-Sachet™ set 1 Reusable Stainless Steel Bowl	1 Kit	7" x 6" x 3"
H-EZ500		Reusable Stainless Steel Bowl 4.3" x 5.3" x 3" / 11 cm x 13.5 cm x 7.5 cm	1 Unit	7" x 6" x 3"

Your First ™

FOAMING SPRAY

READY-TO-USE ENZYME FOAMING SPRAY

Our Foaming Spray is ideal for instrument pre-cleaning as the thick foam spray enzymatically breaks down all bioburden to allow for effective cleaning and sterilization.

Formulated with **PrAmLi Power³ Enzymes™**, a unique blend of three powerful enzymes along with a surfactant to provide the ultimate in cleaning efficacy



INSTRUCTIONS



PRODUCT SELECTION

ITEM #	PACKAGING	CASE DIMENSIONS
SS-21-4001	1 x 946 ml (32 oz) Spray Bottle	10" x 4"
SS-21-4012	12 x 946 ml (32 oz) Spray Bottles	18" x 12" x 11"

Important: In addition to our enzymatic formulation, our spray bottle contains as much as 25% more solution when compared to the competition.

Your First ™

FIRST WIPES

Our cleanroom wipes are designed to do the first wipe-down of an endoscope after a procedure. They are made from hydroentangled non-woven material composed of 45% polyester and 55% cellulose. These wipes are manufactured and packaged in a cleanroom environment.



- Can be used for a cleanroom environment of ISO Class 5 through 8
- Production process guarantees extremely low particle generation and chemical extractables
- Compatible with virtually any solvent or solution for cleaning and wiping (such as **PrAmLi Power³ Enzymes™**)
- Double-bag packaging
- Available in two versatile sizes

PRODUCT SELECTION

ITEM #	DESCRIPTION	DIMENSIONS	PACKAGING	CASE DIMENSIONS
17-A0909	Non-woven Cleanroom Wipes	9" x 9" (23 cm x 23 cm)	300 / Bag 14 Bags / Case	18.5" x 10" x 25.25"
17-A1212	Non-woven Cleanroom Wipes	12" x 12" (31 cm x 31 cm)	150 / Bag 10 Bags / Case	13" x 13" x 19"

Your First ™

DRY SPONGES

Our sponges are specifically fabricated for use with all medical instrumentation. Unlike most competitor sponges, our sponges are white to better evidence the effectiveness of the cleaning process. Dark colored sponges easily hide contaminants, but with our white sponges you can visually evaluate the pre-clean efficacy.



Available in two formats, Flat and Tubular:

Will not flake or scratch expensive equipment

Eco-friendly

100% medical grade high-density foam

Compatible with enzymatic and disinfectant solutions

TIP: OUR TUBULAR SPONGE CAN ALSO BE USED AS A BUMPER TO PROTECT YOUR SCOPES

High absorbency retains more solution than conventional sponges or paper wipes

OUR WHITE SPONGES MAY APPEAR SLIGHTLY YELLOWED. THIS IS ONLY DUE TO UV EXPOSURE, IS NOT A REFLECTION OF ITS AGE.

PRODUCT SELECTION

ITEM #	DESCRIPTION	DIMENSIONS	PACKAGING
17-100	DRY TUBULAR SPONGE White Dispenser Box	2.0" ø x 4.25" (2.7 cm ø x 10.8 cm)	25 / Box (4 Boxes / Case)
17-101	DRY TUBULAR SPONGE Bulk	2.0" ø x 4.25" (2.7 cm ø x 10.8 cm)	560 / Case
17-200	DRY FLAT SPONGE White Dispenser Box	4" x 5" x 0.5" (10.2 cm x 12.7 cm x 1.3 cm)	50 / Box (4 Boxes / Case)
17-201	DRY FLAT SPONGE Bulk	4" x 5" x 0.5" (10.2 cm x 12.7 cm x 1.3 cm)	1,000 / Case

Your First ™

ENZYMATIC BEDSIDE SPONGE KITS

RightSwipe™

Introducing RightSwipe™, our newly launched Enzymatic Bedside Sponge Kits now feature sponges treated with an enhanced formula for a more effective pre-clean or surface clean solution. Our improved three-enzyme blend ensures thorough and safe cleaning and a simplified application process making RightSwipe™ the perfect match for your endoscope care needs. The tubular sponge design is compatible with diverse endoscope models, providing versatility alongside the superior benefits of our enzymatic solution.



FLAT SPONGE



TUBULAR SPONGE



FREESTANDING POUCH



ENVIRO-SACHET™
ENZYMATIC GEL KIT

PRODUCT SELECTION

ITEM #	DESCRIPTION	PACKAGING	CASE DIMENSIONS
17-111	PRE-TREATED TUBULAR SPONGE Add 100-500 ml of water to activate detergent solution. Comes in plastic bag. ¹	25 / Box (4 Boxes / Case)	21" x 14" x 13.5"
17-211	PRE-TREATED FLAT SPONGE Add 100-500 ml of water to activate detergent solution. Comes in plastic bag. ¹	50 / Box (4 Boxes / Case)	21" x 14" x 13.5"
17-122	PRE-TREATED TUBULAR SPONGE Comes in a freestanding pouch to use as a basin. Simply add 500 ml of water to the pouch.	15 / Box (4 Boxes / Case)	21" x 14" x 13.5"
17-222	PRE-TREATED FLAT SPONGE Comes in a freestanding pouch to use as a basin. Simply add 500 ml of water to the pouch.	25 / Box (4 Boxes / Case)	21" x 14" x 13.5"
17-133	READY TO USE TUBULAR SPONGE Comes pre-treated in 500 ml of an enzymatic solution in freestanding pouch.	24 / Case	14" x 12" x 9"
17-233	READY TO USE FLAT SPONGE Comes pre-treated in 500 ml of an enzymatic solution in freestanding pouch.	24 / Case	14" x 12" x 9"
17-144	READY TO USE DUAL PURPOSE TUBULAR SPONGE Comes pre-treated in 100 ml of an enzymatic solution in freestanding pouch.	12 / Box (4 Boxes / Case)	21" x 14" x 13.5"
17-244	READY TO USE DUAL PURPOSE FLAT SPONGE Comes pre-treated in 100 ml of an enzymatic solution in freestanding pouch.	12 / Box (4 Boxes / Case)	21" x 14" x 13.5"
17-155	READY TO USE TUBULAR SPONGE Comes pre-treated in 220 ml of an enzymatic solution in freestanding pouch.	24 / Case	14" x 12" x 9"
17-255	READY TO USE FLAT SPONGE Comes pre-treated in 220 ml of an enzymatic solution in freestanding pouch.	24 / Case	14" x 12" x 9"
17-277	ENZYMATIC GEL Comes with 1 Dry Flat Sponge, Tubular Sponge, Enviro-Sachet™ & Stainless Steel Bowl kit	1 Kit	1 Kit

Our stainless steel bowl can be sold separately and used as a basin for easy and convenient sponge activation.¹

Your First ™

NON-ENZYMATIC DETERGENT BEDSIDE SPONGE KITS

RightSwipe™

Presenting our enhanced Non-Enzymatic Bedside Sponge Kits. Each kit includes sponges pre-treated with a free-rinsing, neutral pH detergent blend from RightSwipe™, ensuring an effective pre-cleaning and surface cleaning solution for medical instruments. This kit stands as Your First ™ step in the medical instrument cleaning process, showcasing the advanced benefits for superior instrument care. Our enhanced bedside kits will provide you with a streamlined workflow.



FREESTANDING POUCH



FLAT SPONGE



TUBULAR SPONGE



ENVIRO-SACHET™
NON-ENZYMATIC CLEANER



DISPENSER BOX

PRODUCT SELECTION

ITEM #	DESCRIPTION	PACKAGING	CASE DIMENSIONS
717-122	PRE-TREATED TUBULAR SPONGE Add 100-500 ml of water to activate non-enzymatic detergent. Comes in plastic bag.1	15 / Box (4 Boxes / Case)	21" x 14" x 13.5"
717-222	PRE-TREATED FLAT SPONGE Add 100-500 ml of water to activate non-enzymatic detergent. Comes in plastic bag.1	25 / Box (4 Boxes / Case)	21" x 14" x 13.5"
717-133	READY TO USE TUBULAR SPONGE Comes pre-treated in 500 ml of non-enzymatic solution in freestanding pouch.	24 / Case	14" x 12" x 9"
717-233	READY TO USE FLAT SPONGE Comes pre-treated in 500 ml of non-enzymatic solution in freestanding pouch.	24 / Case	14" x 12" x 9"
717-144	READY TO USE DUAL PURPOSE TUBULAR SPONGE Comes pre-treated in 100 ml of non-enzymatic solution in freestanding pouch.	12 / Box (4 Boxes / Case)	21" x 14" x 13.5"
717-244	READY TO USE DUAL PURPOSE FLAT SPONGE Comes pre-treated in 100 ml of non-enzymatic solution in freestanding pouch.	12 / Box (4 Boxes / Case)	21" x 14" x 13.5"
717-155	READY TO USE TUBULAR SPONGE Comes pre-treated in 220 ml of non-enzymatic solution in freestanding pouch.	24 / Case	14" x 12" x 9"
717-255	READY TO USE FLAT SPONGE Comes pre-treated in 220 ml of non-enzymatic solution in freestanding pouch.	24 / Case	14" x 12" x 9"
717-277	NON-ENZYMATIC CLEANER Comes with 1 Dry Flat Sponge, Enviro-Sachet™ & Stainless Steel Bowl kit.	1 Kit	1 Kit

Our stainless steel bowl can be sold separately and used as a basin for easy and convenient sponge activation.1

Your First ™

TRANSPORT BAGS FOR ENDOSCOPE AND MEDICAL INSTRUMENTS

Thoughtfully designed to eliminate touch-points during endoscope transport from point-of-care to the reprocessing area, and back, focused on safety of patients and medical staff.


- Made of recyclable LDPE material known to be quite flexible, highly resilient, and tough
- Latex and DEHP-free
- Excellent resistance to dilute and concentrated acids, alcohols, bases and esters*
- Good resistance to aldehydes, ketones, and vegetable-based oils*
- Can withstand high temperatures of up to 90 °C (194°F)
- Tamper-evident seal for added safety, also aids in keeping the instrument moist to avoid drying of bio-debris
- Clearly marked bags with color-coded visuals facilitates immediate identification of clean vs. used instruments
- Users can write important information directly on the bag in spaces marked with visual icons
- Available in two convenient sizes

GOOD TO KNOW: USING **Your First ™** SEALABLE TRANSPORT BAGS REDUCES USER ERROR, PREVENTS CROSS-CONTAMINATION, AND PROMOTES PATIENT AND MEDICAL STAFF SAFETY.

CLEAN INSTRUMENTS BAG



SOILED INSTRUMENTS BAG

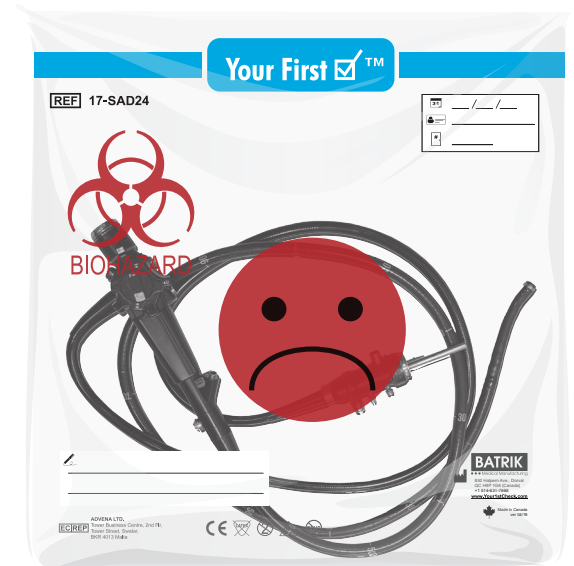

MADE IN CANADA

PRODUCT SELECTION

ITEM #	DESCRIPTION	DIMENSIONS	PACKAGING	CASE DIMENSIONS
17-HAP24	Clean Instruments Bag	24" x 22" / 61.0 cm x 56.0 cm	200 / Case	18" x 14" x 5"
17-SAD24	Soiled Instruments Bag	24" x 22" / 61.0 cm x 56.0 cm	200 / Case	18" x 14" x 5"
17-HAP30	Clean Instruments Bag	30.5" x 26" / 77.0 cm x 66.0 cm	200 / Case	21" x 19" x 7"
17-SAD30	Soiled Instruments Bag	30.5" x 26" / 77.0 cm x 66.0 cm	200 / Case	21" x 19" x 7"

*Reference on file

BATRIK ARC™ WIRE STYLE BRUSHES

SINGLE-USE AND REUSABLE

You've got options with our ARC™ Wire Brushes. Our single-use and reusable brushes are available in both flexible and rigid varieties and come in a multitude of sizes. ARC™ Wire Brushes are guaranteed to last longer than twisted wire and preserve the integrity of your valuable instruments.

BRUSH SELECTION TIPS:

INSTRUMENT SPECIFICATIONS

What is the maximum length of the instrument(s) to be cleaned?

Add at least 2" (5cm) to that length, so you have an adequate length to work with.

What is the inner diameter of the instrument?

Add at most 0.125" (3mm) to determine the diameter of the bristle to ensure effective cleaning. The bristle diameter has to slightly exceed the inner tube diameter. Bristles that are too large will bend too much and not clean effectively.

BRUSH SELECTION

Simply find the section in our catalogue that best corresponds to the overall length and diameter of your instrument. As a rule, if your present brush is 8" (20.32cm) in length and we only offer 10" (25.4cm) length with your required diameter, that is entirely acceptable.

In many cases, our ARC™ Wire Brushes will have more bristles on the brush. As an example, many 5mm (0.19") diameter brushes are only 1" (2.5cm) in length. At Batrik, we strived to provide the greatest amount of bristles to ensure the most effective cleaning.

CUSTOMIZABLE WIRE BRUSHES


A unique aspect of single-use ARC™ Wire Brushes is that they can be customized to any length:

- Up to 2000mm for Wire Brushes on plastic semi-rigid sheaths, rigid brushes from 5 to 20mm in diameter.
- No limit in length for Wire Brushes on flexible plastic sheaths, flexible brushes from 2-5mm in diameter

RIGID



FLEXIBLE

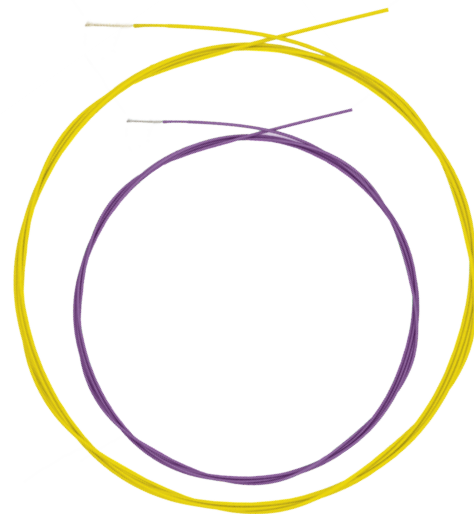

12" / 30 CM SINGLE-END WIRE BRUSH - DELRIN PLASTIC

	ITEM # SINGLE-USE	ITEM # REUSABLE	COLOR	TOTAL LENGTH	BRUSH (Ø)	BRUSH LENGTH	SHEATH (Ø)
FLEXIBLE BRUSH & SHEATH	41-31502-S	41-31502	White	12.401" / 31.5 cm	0.078" / 2 mm	0.59" / 15 mm	0.39" / 1 mm
	41-32503-S	41-32503	Yellow	12.795" / 32.5 cm	0.118" / 3 mm	0.98" / 25 mm	0.39" / 1 mm
	411-32504-S	411-32504	Violet	12.795" / 32.5 cm	0.157" / 4 mm	0.98" / 25 mm	0.61" / 1.55 mm
	411-32505-S	411-32505	Grey	12.795" / 32.5 cm	0.19" / 5 mm	0.98" / 25 mm	0.68" / 1.75 mm
RIGID BRUSH & SHEATH	41-35005-S	41-35005	White	13.779" / 35 cm	0.19" / 5 mm	1.968" / 50 mm	0.118" / 3 mm
	41-35006-S	41-35006	Green	13.779" / 35 cm	0.23" / 6 mm	1.968" / 50 mm	0.118" / 3 mm
	41-35008-S	41-35008	Blue	13.779" / 35 cm	0.31" / 8 mm	1.968" / 50 mm	0.157" / 4 mm
	41-35010-S	41-35010	Red	13.779" / 35 cm	0.39" / 10 mm	1.968" / 50 mm	0.157" / 4 mm
	41-35012-S	41-35012	Yellow	13.779" / 35 cm	0.49" / 12.5 mm	1.968" / 50 mm	0.157" / 4 mm
	41-35015-S	41-35015	Green	13.779" / 35 cm	0.59" / 15 mm	1.968" / 50 mm	0.228" / 5.8 mm
	41-35020-S	41-35020	Blue	13.779" / 35 cm	0.78" / 20 mm	1.968" / 50 mm	0.228" / 5.8 mm

PACKAGING Single-Use: Bag of 20 / Reusable: Bag of 3

18" / 45 CM SINGLE-END WIRE BRUSH - DELRIN PLASTIC

	ITEM # SINGLE-USE	ITEM # REUSABLE	COLOR	TOTAL LENGTH	BRUSH (Ø)	BRUSH LENGTH	SHEATH (Ø)
FLEXIBLE BRUSH & SHEATH	41-46502-S	41-46502	White	18.307" / 46.5 cm	0.078" / 2 mm	0.59" / 15 mm	0.39" / 1 mm
	41-47503-S	41-47503	Yellow	18.7" / 47.5 cm	0.118" / 3 mm	0.98" / 25 mm	0.39" / 1 mm
	41-47505-S	N/A	White	18.7" / 47.5 cm	0.19" / 5 mm	0.98" / 25 mm	0.68" / 1.75 mm
	411-47504-S	411-47504	Violet	18.7" / 47.5 cm	0.157" / 4 mm	0.98" / 25 mm	0.61" / 1.55 mm
	411-47505-S	411-47505	Grey	18.7" / 47.5 cm	0.19" / 5 mm	0.98" / 25 mm	0.68" / 1.75 mm
RIGID BRUSH & SHEATH	41-50005-S	41-50005	White	19.685" / 50 cm	0.19" / 5 mm	1.968" / 50 mm	0.118" / 3 mm
	41-50006-S	41-50006	Green	19.685" / 50 cm	0.23" / 6 mm	1.968" / 50 mm	0.118" / 3 mm
	41-50008-S	41-50008	Blue	19.685" / 50 cm	0.31" / 8 mm	1.968" / 50 mm	0.157" / 4 mm
	41-50010-S	41-50010	Red	19.685" / 50 cm	0.39" / 10 mm	1.968" / 50 mm	0.157" / 4 mm
	41-50012-S	41-50012	Yellow	19.685" / 50 cm	0.49" / 12.5 mm	1.968" / 50 mm	0.157" / 4 mm
	41-50015-S	41-50015	Green	19.685" / 50 cm	0.59" / 15 mm	1.968" / 50 mm	0.228" / 5.8 mm
	41-50020-S	41-50020	Blue	19.685" / 50 cm	0.78" / 20 mm	1.968" / 50 mm	0.228" / 5.8 mm

PACKAGING: Single-Use: Bag of 20 / Reusable: Bag of 3

24" / 60 CM SINGLE-END WIRE BRUSH - DELRIN PLASTIC

	ITEM # SINGLE-USE	ITEM # REUSABLE	COLOR	TOTAL LENGTH	BRUSH (Ø)	BRUSH LENGTH	SHEATH (Ø)
FLEXIBLE BRUSH & SHEATH	41-61502-S	41-61502	White	24.212" / 61.5 cm	0.078" / 2 mm	0.59" / 15 mm	0.39" / 1mm
	41-62503-S	41-62503	Yellow	24.6" / 62.5 cm	0.118" / 3 mm	0.98" / 25 mm	0.39" / 1mm
	41-62505-S	N/A	White	24.6" / 62.5 cm	0.19" / 5 mm	0.98" / 25 mm	0.68" / 1.75mm
	411-62504-S	411-62504	Violet	24.6" / 62.5 cm	0.157" / 4 mm	0.98" / 25 mm	0.61" / 1.55mm
	411-62505-S	411-62505	Grey	24.6" / 62.5 cm	0.19" / 5 mm	0.98" / 25 mm	0.68" / 1.75mm
RIGID BRUSH & SHEATH	41-65005-S	41-65005	White	25.59" / 65 cm	0.19" / 5 mm	1.968" / 50 mm	0.118" / 3mm
	41-65006-S	41-65006	Green	25.59" / 65 cm	0.23" / 6 mm	1.968" / 50 mm	0.118" / 3mm
	41-65008-S	41-65008	Blue	25.59" / 65 cm	0.31" / 8 mm	1.968" / 50 mm	0.157" / 4 mm
	41-65010-S	41-65010	Red	25.59" / 65 cm	0.39" / 10 mm	1.968" / 50 mm	0.157" / 4 mm
	41-65012-S	41-65012	Yellow	25.59" / 65 cm	0.49" / 12.5 mm	1.968" / 50 mm	0.157" / 4 mm
	41-65015-S	41-65015	Green	25.59" / 65 cm	0.59" / 15 mm	1.968" / 50 mm	0.228" / 5.8 mm
	41-65020-S	41-65020	Blue	25.59" / 65 cm	0.78" / 20 mm	1.968" / 50 mm	0.228" / 5.8 mm



PACKAGING: Single-Use: Bag of 20 / Reusable: Bag of 3

30" / 75 CM SINGLE-END WIRE BRUSH - DELRIN PLASTIC

	ITEM # SINGLE-USE	ITEM # REUSABLE	COLOR	TOTAL LENGTH	BRUSH (Ø)	BRUSH LENGTH	SHEATH (Ø)
FLEXIBLE BRUSH & SHEATH	41-76502-S	41-76502	White	30.118" / 76.5 cm	0.078" / 2 mm	0.59" / 15 mm	0.39" / 1 mm
	41-77503-S	41-77503	Yellow	30.511" / 77.5 cm	0.118" / 3 mm	0.98" / 25 mm	0.39" / 1 mm
	411-77504-S	411-77504	Violet	30.511" / 77.5 cm	0.157" / 4 mm	0.98" / 25 mm	0.61" / 1.55 mm
	411-77505-S	411-77505	Grey	30.511" / 77.5 cm	0.19" / 5 mm	0.98" / 25 mm	0.68" / 1.75 mm
	41-80005-S	41-80005	White	31.496" / 80 cm	0.19" / 5 mm	1.968" / 50 mm	0.118" / 3 mm
RIGID BRUSH & SHEATH	41-80006-S	41-80006	Green	31.496" / 80 cm	0.23" / 6 mm	1.968" / 50 mm	0.118" / 3 mm
	41-80008-S	41-80008	Blue	31.496" / 80 cm	0.31" / 8 mm	1.968" / 50 mm	0.157" / 4 mm
	41-80010-S	41-80010	Red	31.496" / 80 cm	0.39" / 10 mm	1.968" / 50 mm	0.157" / 4 mm
	41-80012-S	41-80012	Yellow	31.496" / 80 cm	0.49" / 12.5 mm	1.968" / 50 mm	0.157" / 4 mm
	41-80015-S	41-80015	Green	31.496" / 80 cm	0.59" / 15 mm	1.968" / 50 mm	0.228" / 5.8 mm
	41-80020-S	41-80020	Blue	31.496" / 80 cm	0.78" / 20 mm	1.968" / 50 mm	0.228" / 5.8 mm


PACKAGING: Single-Use: Bag of 20 / Reusable: Bag of 3

BRUSHES

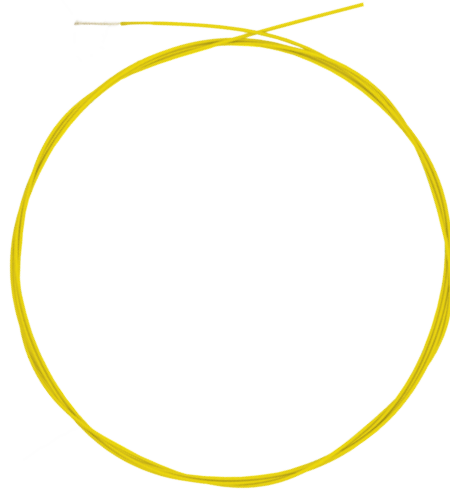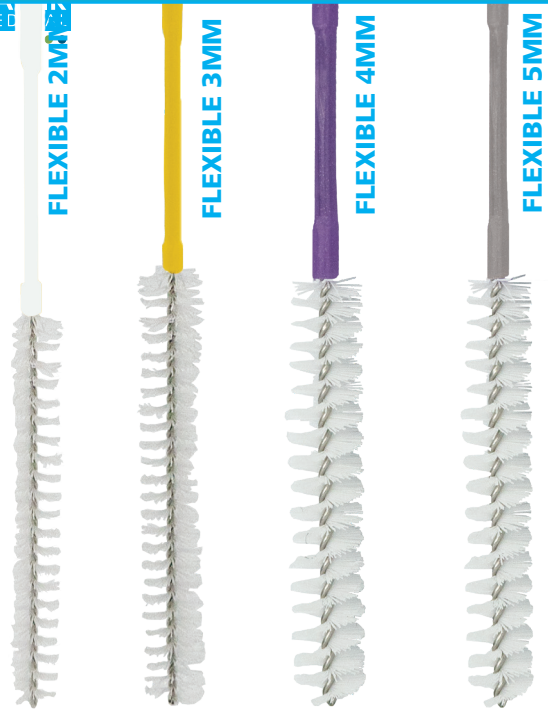

FLEXIBLE

Flexible brush diameter (Ø) ranging from 2 mm to 5 mm

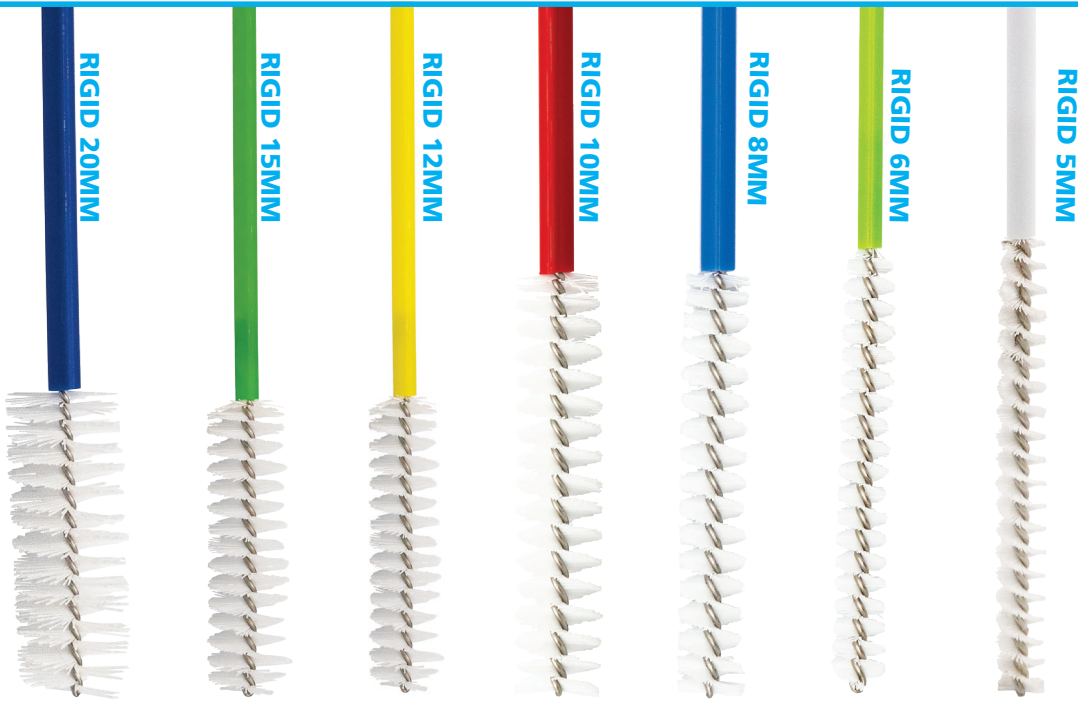

RIGID

Rigid brush diameter (Ø) ranging from 5 mm to 20 mm

BRUSH SELECTION SYSTEM (Actual Size)

2mm 0.078"	3 mm 0.118"	4 mm 0.157"	5 mm 0.196"	6 mm 0.236"	8 mm 0.314"	10 mm 0.394"	12 mm 0.472"	15 mm 0.591"	20 mm 0.787"

BATRIK REUSABLE DOUBLE ARC™ WIRE BRUSHES

ASYMMETRIC

14" / 35 CM ASYMMETRIC DOUBLE-END WIRE BRUSHES

ITEM #	COLOR	TOTAL LENGTH	BRUSH 1 (Ø)	BRUSH 2 (Ø)	BRUSH LENGTH	SHEATH (Ø)
40-35056	White	13.779" / 35 cm	0.19" / 5 mm	0.23" / 6 mm	1.968" / 50 mm	0.118" / 3 mm
40-35112	Yellow	13.779" / 35 cm	0.39" / 10 mm	0.49" / 12.5 mm	1.968" / 50 mm	0.157" / 4 mm
40-35152	Green	13.779" / 35 cm	0.59" / 15 mm	0.78" / 20 mm	1.968" / 50 mm	0.228" / 5.8 mm
40-35810	Blue	13.779" / 35 cm	0.31" / 8 mm	0.39" / 10 mm	1.968" / 50 mm	0.157" / 4 mm

PACKAGING: Bag of 3

20" / 50 CM ASYMMETRIC DOUBLE-END WIRE BRUSHES

ITEM #	COLOR	TOTAL LENGTH	BRUSH 1 (Ø)	BRUSH 2 (Ø)	BRUSH LENGTH	SHEATH (Ø)
40-50056	White	19.685" / 50 cm	0.19" / 5 mm	0.23" / 6 mm	1.968" / 50 mm	0.118" / 3 mm
40-50112	Yellow	19.685" / 50 cm	0.39" / 10 mm	0.49" / 12.5 mm	1.968" / 50 mm	0.157" / 4 mm
40-50152	Green	19.685" / 50 cm	0.59" / 15 mm	0.78" / 20 mm	1.968" / 50 mm	0.228" / 5.8 mm
40-50810	Blue	19.685" / 50 cm	0.31" / 8 mm	0.39" / 10 mm	1.968" / 50 mm	0.157" / 4 mm

PACKAGING: Bag of 3

26" / 65 CM ASYMMETRIC DOUBLE-END WIRE BRUSHES

ITEM #	COLOR	TOTAL LENGTH	BRUSH 1 (Ø)	BRUSH 2 (Ø)	BRUSH LENGTH	SHEATH (Ø)
40-65056	White	25.59" / 65 cm	0.19" / 5 mm	0.23" / 6 mm	1.968" / 50 mm	0.118" / 3 mm
40-65112	Yellow	25.59" / 65 cm	0.39" / 10 mm	0.49" / 12.5 mm	1.968" / 50 mm	0.157" / 4 mm
40-65152	Green	25.59" / 65 cm	0.59" / 15 mm	0.78" / 20 mm	1.968" / 50 mm	0.228" / 5.8 mm
40-65810	Blue	25.59" / 65 cm	0.31" / 8 mm	0.39" / 10 mm	1.968" / 50 mm	0.157" / 4 mm

PACKAGING: Bag of 3

31" / 80 CM ASYMMETRIC DOUBLE-END WIRE BRUSHES

ITEM #	COLOR	TOTAL LENGTH	BRUSH 1 (Ø)	BRUSH 2 (Ø)	BRUSH LENGTH	SHEATH (Ø)
40-80056	White	31.496" / 80 cm	0.19" / 5 mm	0.23" / 6 mm	1.968" / 50 mm	0.118" / 3 mm
40-80112	Yellow	31.496" / 80 cm	0.39" / 10 mm	0.49" / 12.5 mm	1.968" / 50 mm	0.157" / 4 mm
40-80152	Green	31.496" / 80 cm	0.59" / 15 mm	0.78" / 20 mm	1.968" / 50 mm	0.228" / 5.8 mm
40-80810	Blue	31.496" / 80 cm	0.31" / 8 mm	0.39" / 10 mm	1.968" / 50 mm	0.157" / 4 mm

PACKAGING: Bag of 3

BATRIK REUSABLE DOUBLE ARC™ WIRE BRUSHES SYMMETRIC

14" / 35 CM SYMMETRIC DOUBLE-END WIRE BRUSHES

ITEM #	COLOR	TOTAL LENGTH	BRUSH 1 (Ø)	BRUSH LENGTH	SHEATH (Ø)
40-35005	White	13.779" / 35 cm	0.19" / 5 mm	1.968" / 50 mm	0.118" / 3 mm
40-35006	Green	13.779" / 35 cm	0.23" / 6 mm	1.968" / 50 mm	0.118" / 3 mm
40-35008	Blue	13.779" / 35 cm	0.31" / 8 mm	1.968" / 50 mm	0.157" / 4 mm
40-35010	Red	13.779" / 35 cm	0.39" / 10 mm	1.968" / 50 mm	0.157" / 4 mm
40-35012	Yellow	13.779" / 35 cm	0.49" / 12.5 mm	1.968" / 50 mm	0.157" / 4 mm
40-35015	Green	13.779" / 35 cm	0.59" / 15 mm	1.968" / 50 mm	0.228" / 5.8 mm
40-35020	Blue	13.779" / 35 cm	0.78" / 20 mm	1.968" / 50 mm	0.228" / 5.8 mm

PACKAGING: Bag of 3

20" / 50 CM SYMMETRIC DOUBLE-END WIRE BRUSHES

ITEM #	COLOR	TOTAL LENGTH	BRUSH 1 (Ø)	BRUSH LENGTH	SHEATH (Ø)
40-50005	White	19.685" / 50 cm	0.19" / 5 mm	1.968" / 50 mm	0.118" / 3 mm
40-50006	Green	19.685" / 50 cm	0.23" / 6 mm	1.968" / 50 mm	0.118" / 3 mm
40-50008	Blue	19.685" / 50 cm	0.31" / 8 mm	1.968" / 50 mm	0.157" / 4 mm
40-50010	Red	19.685" / 50 cm	0.39" / 10 mm	1.968" / 50 mm	0.157" / 4 mm
40-50012	Yellow	19.685" / 50 cm	0.49" / 12.5 mm	1.968" / 50 mm	0.157" / 4 mm
40-50015	Green	19.685" / 50 cm	0.59" / 15 mm	1.968" / 50 mm	0.228" / 5.8 mm
40-50020	Blue	19.685" / 50 cm	0.78" / 20 mm	1.968" / 50 mm	0.228" / 5.8 mm

PACKAGING: Bag of 3

26" / 65 CM SYMMETRIC DOUBLE-END WIRE BRUSHES

ITEM #	COLOR	TOTAL LENGTH	BRUSH 1 (Ø)	BRUSH LENGTH	SHEATH (Ø)
40-65005	White	25.59" / 65 cm	0.19" / 5 mm	1.968" / 50 mm	0.118" / 3 mm
40-65006	Green	25.59" / 65 cm	0.23" / 6 mm	1.968" / 50 mm	0.118" / 3 mm
40-65008	Blue	25.59" / 65 cm	0.31" / 8 mm	1.968" / 50 mm	0.157" / 4 mm
40-65010	Red	25.59" / 65 cm	0.39" / 10 mm	1.968" / 50 mm	0.157" / 4 mm
40-65012	Yellow	25.59" / 65 cm	0.49" / 12.5 mm	1.968" / 50 mm	0.157" / 4 mm
40-65015	Green	25.59" / 65 cm	0.59" / 15 mm	1.968" / 50 mm	0.228" / 5.8 mm
40-65020	Blue	25.59" / 65 cm	0.78" / 20 mm	1.968" / 50 mm	0.228" / 5.8 mm

PACKAGING: Bag of 3

31" / 80 CM SYMMETRIC DOUBLE-END WIRE BRUSHES

ITEM #	COLOR	TOTAL LENGTH	BRUSH 1 (Ø)	BRUSH LENGTH	SHEATH (Ø)
40-80005	White	31.496" / 80 cm	0.19" / 5 mm	1.968" / 50 mm	0.118" / 3 mm
40-80006	Green	31.496" / 80 cm	0.23" / 6 mm	1.968" / 50 mm	0.118" / 3 mm
40-80008	Blue	31.496" / 80 cm	0.31" / 8 mm	1.968" / 50 mm	0.157" / 4 mm
40-80010	Red	31.496" / 80 cm	0.39" / 10 mm	1.968" / 50 mm	0.157" / 4 mm
40-80012	Yellow	31.496" / 80 cm	0.49" / 12.5 mm	1.968" / 50 mm	0.157" / 4 mm
40-80015	Green	31.496" / 80 cm	0.59" / 15 mm	1.968" / 50 mm	0.228" / 5.8 mm
40-80020	Blue	31.496" / 80 cm	0.78" / 20 mm	1.968" / 50 mm	0.228" / 5.8 mm

PACKAGING: Bag of 3

ASYMMETRIC BRUSHES

5 MM

6 MM



8 MM

10 MM



10 MM

12.5 MM



15 MM

20 MM




BATRIK REUSABLE IN-BRUSH™ ANTIMICROBIAL TWISTED WIRE BRUSHES

Our In-Brush™ Antimicrobial Reusable Twisted Wire Brushes are designed for use with Frazier/Baron suction tubes. The twisted wires have a fine/medium diameter yet are rigid enough to ensure excellent cleaning capabilities.

Our wide assortment of Nylon Twisted Wire Brushes are ideal for large bore, uniquely shaped, and endoscopic instruments.

GOOD TO KNOW: OUR TWISTED WIRE BRUSHES FEATURE A PROTECTIVE BLUE TIP FOR ADDED SAFETY.

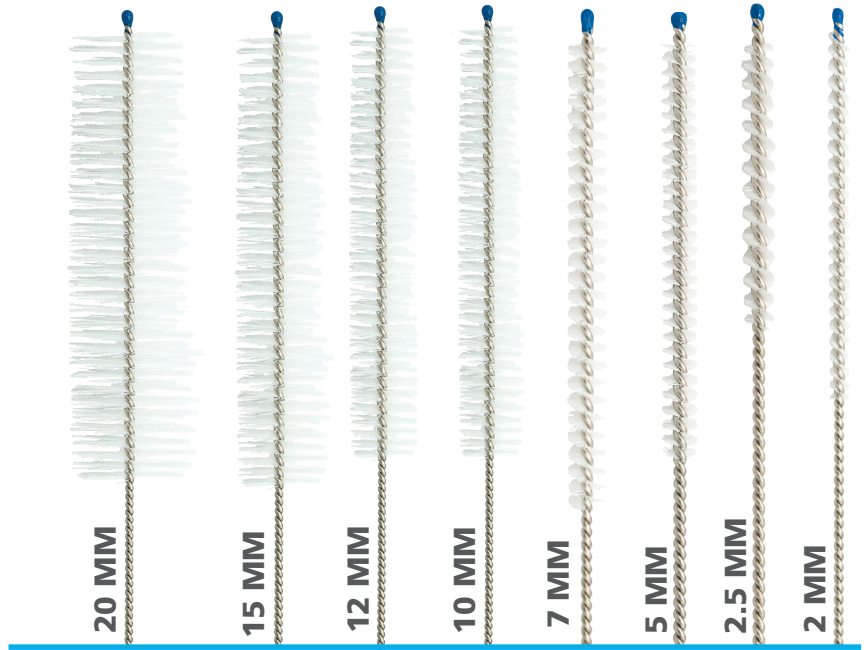

10" / 250 MM FRAZIER /BARON SUCTION TUBE BRUSHES

ITEM #	SIZE	BRUSH LENGTH	BRUSH (Ø)
21-10039	3 FR	0.5" / 12.7 mm	0.039" / 1 mm
21-10067	5 FR	0.5" / 12.7 mm	0.065" / 1.67 mm
20-10059	5 FR	1" / 25.4 mm	0.059" / 1.50 mm
20-10079	6 FR	1" / 25.4 mm	0.079" / 2 mm
20-10090	7 FR	1" / 25.4 mm	0.090" / 2.3 mm
20-10106	8 FR	1" / 25.4 mm	0.106" / 2.7 mm
20-10118	9 FR	1" / 25.4 mm	0.118" / 3 mm
20-10129	10 FR	1.5" / 38 mm	0.129" / 3.3 mm
20-10146	11 FR	1.5" / 38 mm	0.0146" / 3.7 mm
20-10157	12 FR	1.5" / 38 mm	0.157" / 4 mm

PACKAGING: Bag of 3.
Also available in bags of 50, add the suffix "-50".

12" / 305 MM NYLON TWISTED WIRE BRUSHES

ITEM #	BRUSH LENGTH	BRUSH (Ø)
22-12078	2" / 50.8 mm	0.078" / 2 mm
22-12118	2" / 50.8 mm	0.118" / 3 mm
22-12157	3" / 76 mm	0.157" / 4 mm
22-12197	3" / 76 mm	0.197" / 5 mm
22-12236	3" / 76 mm	0.236" / 6 mm
22-12393	3" / 76 mm	0.394" / 10 mm

PACKAGING: Bag of 3.
Also available in bags of 50, add the suffix "- 50".

16" / 406 MM NYLON TWISTED WIRE BRUSHES

ITEM #	BRUSH LENGTH	BRUSH (Ø)
23-16197	2.5" / 63.5 mm	0.197" / 5 mm
23-16236	2.5" / 63.5 mm	0.236" / 6 mm
23-16276	3" / 76 mm	0.276" / 7 mm
23-16394	3" / 76 mm	0.394" / 10 mm
23-16472	3" / 76 mm	0.472" / 12 mm
23-16591	3" / 76 mm	0.591" / 15 mm
23-16787	3" / 76 mm	0.787" / 20 mm

PACKAGING: Bag of 3.
 Also available in bags of 50, add the suffix "-50".

18" / 457 MM NYLON TWISTED WIRE BRUSHES

ITEM #	BRUSH LENGTH	BRUSH (Ø)
19-18010	0.47" / 12 mm	0.039" / 1 mm
19-18015	0.5" / 12 mm	0.06" / 1.6 mm

PACKAGING: Bag of 3.
 Also available in bags of 50, add the suffix "-50".

24" / 610 MM NYLON TWISTED WIRE BRUSHES

ITEM #	BRUSH LENGTH	BRUSH (Ø)
24-24098	2.5" / 63.5 mm	0.098" / 2.5 mm
24-24197	2.5" / 63.5 mm	0.197" / 5 mm
24-24276	3" / 76 mm	0.276" / 7 mm
24-24394	3" / 76 mm	0.394" / 10 mm
24-24472	3" / 76 mm	0.472" / 12 mm

PACKAGING: Bag of 3.
 Also available in bags of 50, add the suffix "-50".

TRACHEOSTOMY BRUSHES

ITEM #	DESCRIPTION	TOTAL LENGTH	BRUSH LENGTH	BRUSH (Ø)
25-08500	Large 6-12.3 3 / 4"	8" / 203 mm	3.75" / 95 mm	0.5" / 12.7 mm
25-37550	Small 000-5.1 7 / 8"	4" / 93 mm	1.88" / 48 mm	0.5" / 12.7 mm

PACKAGING: Bag of 3. Also available in bags of 50, add the suffix "-50".



BRONCHO BRUSHES

ITEM #	DESCRIPTION	BRUSH LENGTH	BRUSH (Ø) INCHES / MM
26-06175	Small	6" / 152 mm	1.75" / 44.5 mm
26-08325	Large	8" / 203 mm	3.25" / 82.6 mm

PACKAGING: Bag of 3.
 Also available in bags of 50, add the suffix "-50".

BATRIK STAINLESS STEEL ARC™ MICRO BRUSHES

Our Stainless Steel ARC™ Micro Brushes are some of the smallest brushes in the industry. They are designed to clean channels with an inner diameter ranging from 0.5 mm to 1.40 mm.

- Flexible stainless steel sheath will not bend like the regular metal twisted wire brushes
- Polyamide bristles prevents or slows bacterial growth between sterilizations
- Reusable*
- Available in two lengths

TIP: USING THE PULL-THROUGH
TECHNIQUE IS HIGHLY
RECOMMENDED WITH THESE
BRUSHES. SIMPLY INSERT THE
BLUNT, NON-BRUSH END FIRST.



PRODUCT SELECTION

ITEM #	TOTAL LENGTH	BRUSH SIZE (Ø)	BRUSH LENGTH	SHEATH (Ø)	PACKAGING
35-03015	12" / 30 cm	0.06" / 1.5 mm	0.6" / 15 mm	0.02" / 0.40 mm	6 / Tube
35-04515	18" / 45 cm	0.06" / 1.5 mm	0.6" / 15 mm	0.02" / 0.40 mm	6 / Tube

*Although all materials used in the fabrication of these brushes can be re-processed repeatedly, due to the delicate nature of these brushes re-use is conditional on their handling and application.

BATRIK INSTRUMENT CLEANING FLAT BRUSHES

Our Instrument Cleaning Flat Brushes are specifically designed for hospital use in the pre-clean and preparation of instruments.

- Long ergonomic handle permits ease-of-use
- Extra long length of brush head (60 [±5] mm in Single-Use and 75 [±5] mm in Reusable)
- Transparent brush fibers allow for quick verification of its cleanliness
- The plastic handle will not release any toxic components
- Available in 6 models (4 reusable/autoclavable and 2 single-use)



This brush is designed for large instruments that require a more gentle cleaning application such as optics, plastics, or hard to reach spot.



The very long, flexible bristles eliminate any stress during sanitation, and makes this the ideal brush for fragile instruments as well as both pliable and rigid endoscopes. The length of the bristles assist in cleaning areas that would otherwise be difficult to reach, such as the accessory channel or tip deflection of pliable endoscopes.



The rigidity of this brush, along with the length of the bristles and handle, improves and simplifies instrument cleaning.



Instrument safe polyamide fibers offer steel-like rigidity and allows deep cleaning of soiled instruments, particularly in grooved or uneven surfaces.



PRODUCT SELECTION

DELTRIN PLASTIC

ITEM #	BRUSH TYPE	DESCRIPTION	COLOR	BRUSH LENGTH	BRUSH SIZE	PACKAGING
50-40917	Extra Rigid	Reusable Autoclavable Flat Brush	Blue	220 mm	75 mm x 20 mm x 10 mm	Bag of 3
50-40918	Rigid	Reusable Autoclavable Flat Brush	Green	220 mm	75 mm x 20 mm x 15 mm	Bag of 3
50-40919	Flexible	Reusable Autoclavable Flat Brush	Yellow	220 mm	75 mm x 20 mm x 15 mm	Bag of 3
50-40920	Extra Flexible	Reusable Autoclavable Flat Brush	Red	220 mm	75 mm x 20 mm x 25 mm	Bag of 3
50-40921	Flexible	Single-Use Flat Brush	Orange	220 mm	60 mm x 20 mm x 15 mm	Bag of 20
50-40922	Rigid	Single-Use Flat Brush	Light Blue	220 mm	60 mm x 20 mm x 15 mm	Bag of 20

Reusable models are also available in bags of 50. Add the suffix "- 50".


BATRIK

REUSABLE PRECISION BRUSH

INSTRUMENT CLEANING

Our Reusable Precision Instrument Cleaning Brush is the ideal brush for the quick removal of difficult debris. It cleans externally as well as deep into instrument jaws. Although the firm and durable bristles feel like stainless steel, there is absolutely no risk of instrument damage or personal injury.

- Bristles have an extra firm density and superior durability
- Polyamide (TYNEX®) bristles
- Steamable thermoplastic ergonomic handle
- Steam / chemical sterilizable
- Latex free

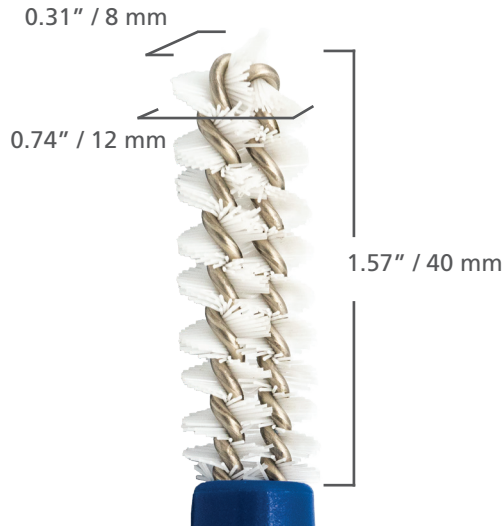

GOOD TO KNOW: ELIMINATES THE NEED FOR THE STAINLESS STEEL TOOTHBRUSH-STYLE BRUSH, WHICH CAN DAMAGE VALUABLE INSTRUMENTS!

PRODUCT SELECTION

BLUE DELRIN PLASTIC

ITEM #	DESCRIPTION	TOTAL LENGTH
50-40800	Reusable Precision External Cleaning Brush	6.29" / 160 mm

PACKAGING: Bag of 3

NEW!

TRIFECTA™ BRUSH

INSTRUMENT CLEANING

Introducing the **Trifecta™ Brush**, a versatile **3-in-1 tool** designed to meet all your instrument cleaning needs. This reusable instrument brush is available in two styles: one with **all-nylon** bristles for gentle cleaning and another with nylon bristles on two sides and **stainless steel bristles** on the third, ideal for tackling tougher buildup and stubborn cleaning tasks. This innovative design ensures effective, thorough maintenance for a wide range of instruments.

TRIFECTA™ BRUSH - NYLON

1 NYLON BRISTLES


Gentle, transparent nylon bristles provide an effective yet safe cleaning solution. This side allows you to verify cleanliness easily, without risking damage to delicate surfaces. Ideal for a light, non-invasive clean.

2 PRECISION TIP

Reach those hard-to-clean spots with the precision tip. This small, firm-bristle tip is designed for tight, hard-to-reach areas, such as instrument jaws and box joints. It ensures a deep clean where it's most needed, without the risk of damage or injury to delicate parts.

3 FIRM NYLON BRISTLES

For tougher residues, the firm nylon side delivers the cleaning power you need. With bristles that feel similar to stainless steel but are safe for your instruments, this side easily handles stubborn dirt and buildup.

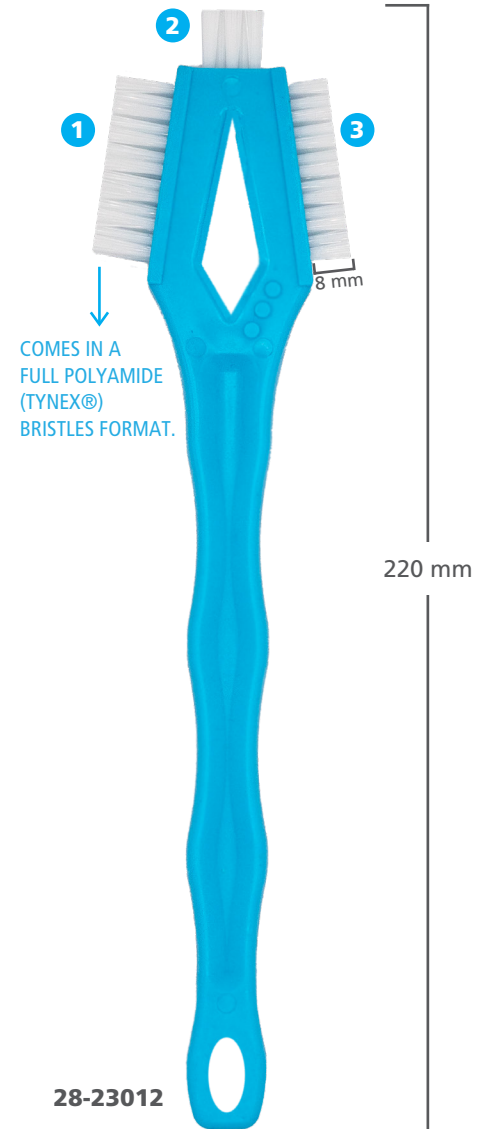

TRIFECTA™ BRUSH - NYLON & STAINLESS STEEL



29-23012

28-23012

PRODUCT SELECTION

ITEM #	COLOR	DESCRIPTION	LENGTH TOTAL	BRISTLE LENGTH	BRISTLES 1 Ø	BRISTLES 2 Ø	BRISTLES 3 Ø
28-23012	Blue	Triple Nylon Brush	8.66" / 220 mm	1.20" x 1.20" x 0.47" / 30.5 x 30.5 x 12.0 mm	0.47" / 12.0 mm	0.47" / 12.0 mm	0.315" / 8.0 mm
29-23012	Blue	Triple Nylon / Stainless Steel	8.66" / 220 mm	1.20" x 1.20" x 0.47" / 30.5 x 30.5 x 12.0 mm	0.47" / 12.0 mm	0.47" / 12.0 mm	0.47" / 12.0 mm

PACKAGING: Bag of 3 - For a bag of 50, add "-50" following the item code.

BATRIK ALL PURPOSE REUSABLE BRUSHES

Plastic Handle Brushes

We offer a wide selection of brushes for instrument cleaning applications. These brushes are ideal for large bore, uniquely shaped, and endoscopic instruments.

In-Brush™ Plastic Handle

Available in regular toothbrush style (nylon or stainless) and double-ended brushes.



28-72550



28-69394



29-72550



29-69394



29-70050



50-40911

PRODUCT SELECTION

ITEM #	DESCRIPTION	LENGTH TOTAL	BRUSH LENGTH	BRUSH HEIGHT
28-72550	Plastic Single End Nylon Reusable Handle Brush	7.25" / 184 mm	1.5" / 38 mm	0.5" / 12.7 mm
29-72550	Plastic Single End Stainless Steel Reusable Handle Brush	7.25" / 184 mm	1.5" / 38 mm	0.5" / 12.7 mm
29-70050	Plastic Single End Brass Reusable Handle Brush	7.25" / 184 mm	1.5" / 38 mm	0.5" / 12.7 mm
28-69394	Plastic Double End Nylon Reusable Handle Brush	6.9" / 175 mm	1.38" / 35 mm 1" / 25 mm	0.39" / 10 mm 0.2" / 5 mm
29-69394	Plastic Double End Stainless Steel Reusable Handle Brush	6.9" / 175 mm	1.38" / 35 mm 1" / 25 mm	0.39" / 10 mm 0.2" / 5 mm
50-40911	Hand Scrub Nylon Brush Reusable autoclavable	4.09" / 104 mm	3.74" / 95 mm	1.1" / 28 mm

PACKAGING: Bag of 3 - For a bag of 50, add "-50" following the item code.

BATRIK

ACETABULAR REAMER BRUSH

Our double-ended cleaning brush is manufactured exclusively for cleaning acetabular bone reamers. The "Pig Tail" rounded style circular head brush is designed to effectively clean inside the reamer, while the rigid straight brush on the reverse side is used for cleaning individual cutting holes. Special fusion of the brush heads to the hand piece along with its high quality antimicrobial medical grade bristles ensures a durable, reusable, and highly effective cleaning tool.



STRAIGHT BRUSH



1.18" / 30 mm length

ROUND BRUSH



2" / 50 mm length

PRODUCT SELECTION

GREEN DELRIN PLASTIC

ITEM #	DESCRIPTION	HANDLE LENGTH	STRAIGHT BRUSH	ROUND BRUSH HEAD
50-ACETAB	Acetabular Reamer Cleaning Brush	4.72" / 120 mm ø 0.47" / 12 mm	1.18" / 30 mm ø 0.3125" / 8 mm	2" / 50 mm ø 1.75" / 45 mm


PACKAGING: Bag of 3

NOTE: Specifications subject to small variances

BATRIK REUSABLE ORTHOPEDIC RASPS CLEANING BRUSH

Our Reusable Orthopedic Rasps Cleaning Brush is used to clean orthopedic rasps and other similar instruments used in orthopedic surgery.

- Long lasting and can withstand many sterilization cycles
- Ergonomic steamable thermoplastic handle
- Antibacterial polyamide bristles
- Reusable
- Latex free

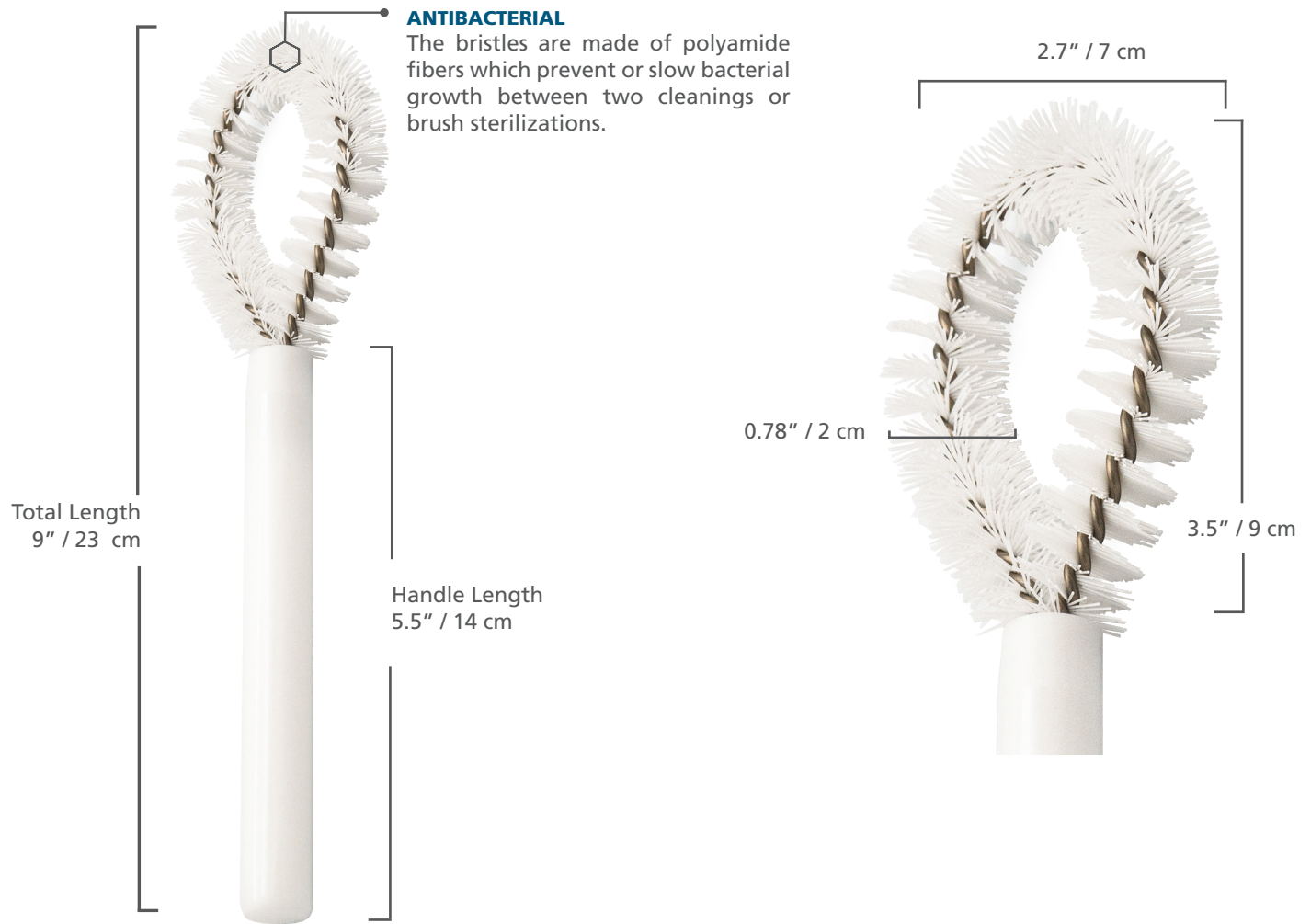

PRODUCT DESCRIPTION

ITEM #	TOTAL LENGTH	HANDLE LENGTH	BRUSH SIZE	PACKAGING
50-RASP	9" / 23 cm	5.5" / 14 cm	0.78" x 2.7" x 3.5" (2 cm x 7 cm x 9 cm)	3 / Bag

BATRIK BOTTLE & DRAIN CLEANING BRUSHES

Designed specifically for hospital applications, our Bottle and Drain Cleaning Brushes feature long and flexible bristles to maximize the effectiveness of large-surface cleaning.


- High quality polyamide clear bristles
- Easy to identify debris
- Autoclavable
- Latex free

REUSABLE JAR/BOTTLE BRUSH



Ergonomic
Handle

STERILIZER DRAIN BRUSH

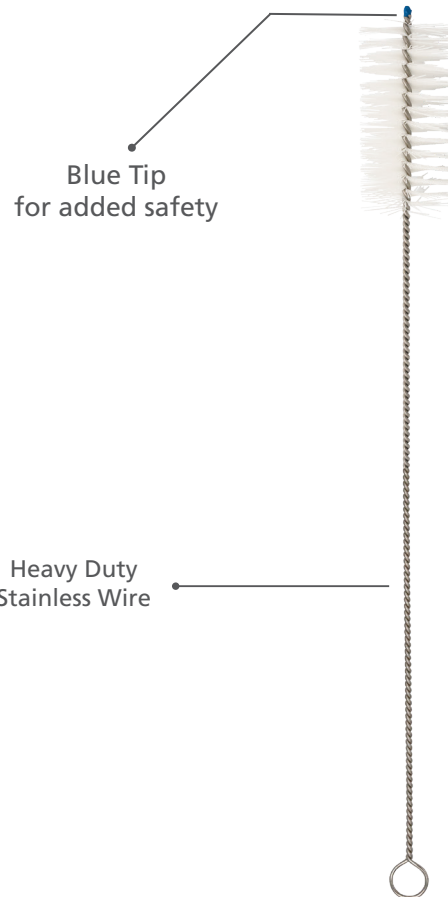

Blue Tip
for added safety

Heavy Duty
Stainless Wire

PRODUCT SELECTION

BLUE DELRIN PLASTIC

ITEM #	DESCRIPTION	HANDLE LENGTH	BRUSH LENGTH	BRUSH (Ø)	TOTAL LENGTH
50-10295	Reusable Jar / Bottle Brush	120 mm ± 2 mm	120 mm ± 5 mm	7 mm ± 2 mm	13.19" / 355 mm
50-10224	Sterilizer Drain Cleaning Brush	N/A	4" / 10.1 cm	2" / 5.08 cm	24" / 609 mm

PACKAGING: Bag of 3

BATRIK

AUTOCLAVE CLEANING BRUSH KITS



ROUND BRUSH

- Features 3 holes for adjustable pole placement
- Adjustable brush placement (parallel or perpendicular) and longer brush length provides easy access to hard-to-reach corners
- Made with FDA approved polyester bristles

HAND BRUSH

- Ergonomic handle for better hand comfort
- Longer, FDA approved polyester bristles
- Easy to clean (no cavity)

PRODUCT SELECTION

REGULAR AUTOCLAVE CLEANING BRUSH KIT

ITEM # (KIT)	DESCRIPTION	ITEM # (UNIT)	DESCRIPTION	DIMENSIONS
50-CLAVEREG	Autoclave REGULAR Brush Cleaning Kit	50-35900	1 Hand Brush	Length: 3.5" x 8" / 8.89 cm x 20.32 cm
		50-55800	1 Rounded Corner Brush (Head Brush)	Head size: 5.5 x 10" / 14 x 25.4 cm
		50-EXTPOLE	1 Anodized Aluminum Extension Pole	2-4 ft / 60-122 cm

SUPERIOR AUTOCLAVE CLEANING BRUSH KIT

ITEM # (KIT)	DESCRIPTION	ITEM # (UNIT)	DESCRIPTION	DIMENSIONS	CASE DIMENSIONS
50-CLAVESUP	Autoclave SUPERIOR Brush Cleaning Kit	50-35900	1 Hand Brush	Length: 3.5" x 8" / 11.5 cm x 20.32cm	
		50-55800	1 Rounded Corner Brush (Head Brush)	Head size: 5.5 x 10" / 14 x 25.4 cm	
		50-EXTPOLE	1 Anodized Aluminum Extension Pole	2-4 ft / 60-122 cm	
		SF-240	2 SAVEFACE™ Face Shields	9" / 23 cm	CASE OF 240 / 26.5" x 21" x 15"

REPLACEMENT COMPONENTS



HAND BRUSH
50-35900

Packaging: 1



BRUSH HEAD
50-55800

Packaging: 1



EXTENSION POLE
50-EXTPOLE

Packaging: 1



SAVEFACE™
FACE SHIELDS SF-240

Packaging: 100 / Box

SURGMED

GROUP 

Proudly Canadian, **The Surgmed™ Group of Companies** began operations in 1992 as a single distribution company known as **Advanced Surgi-Pharm Inc.** The company quickly became known as an industry leader and specialty distributor of value-added surgical products in the Canadian hospital market. Five years later, **Surgmed expanded upon its success with the launch of Batrik Medical Manufacturing Inc.** in 1997 and the subsequent acquisition of **Imperial Surgical Ltd.** one year later in 1998.




Together, our three unique companies serve the Canadian and global healthcare markets with high quality and competitively priced products, all backed by world-class customer care.



Founded in 1992, Advanced Surgi-Pharm Inc. is the distribution arm of Surgmed™ in Canada, representing leading US, international and domestic manufacturers for a wide range of specialty products with unique benefits for applications in the Operating Room, Women's Health, Labor & Delivery, ICU, and Sterile Processing areas. Equipped with a direct sales team, ASP provides effective nationwide coverage as the Canadian distributor of Batrik, Vygon, Rocket Medical, Cooper Surgical, SunMed, and other industry leading products including a wide-ranging portfolio of top-rated PPE products.

Established in 1997 and sold in 85 countries globally, Batrik Medical Manufacturing Inc. specializes in the manufacturing of endoscopy accessories, enzymatic solutions, infection control, central sterilization, and surgical devices. Batrik's in-house controlled environment production facility adheres to strictly enforced quality and regulatory control systems. They are registered as a device manufacturing establishment with the United States FDA and are in full compliance with CE and Health Canada Regulations. All Batrik products are rigorously tested, liability insured, and meet/exceed high international quality standards.

Established in 1935, Imperial Surgical Ltd. is one of Canada's oldest hospital capital equipment manufacturers, well-known for their stainless steel healthcare equipment for Sterile Processing, Operating Room, Examination and Treatment, Patient Handling, and Mortuary/ Pathology areas. Their growing portfolio of products also has applications in Infection Control & Prevention, an important area of focus for their future.

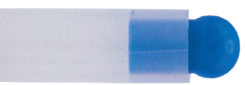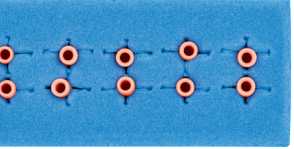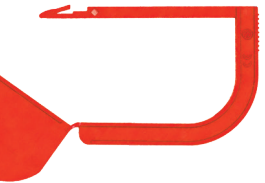

2165 46th Ave., Lachine, QC, Canada, H8T 2P1
Phone: +1.514.631.7988 | Toll Free: 800.661.5432 | Fax.: +1.514.631.9083
Email: info@surgmed.com | Website: www.batrik.com